# FOREO

ESPADA™

**FOREO.COM/SUPPORT**

---

English  
Lietuvių  
Latviešu  
Eesti  
Українська  
български  
ελληνικά  
Slovenčina  
Čeština  
Dansk  
Magyar  
Norsk  
Română

**EN GETTING STARTED:** This booklet contains all essential instructions; for additional product resources, visit [foreo.com/support](http://foreo.com/support).

**⚠ WARNING: No modification of this equipment is allowed.**

Please **READ ALL INSTRUCTIONS BEFORE USE**. Use this product only for its intended use as described in this manual.


**INTENDED USE:** The ESPADA™ has been developed to be used as an at-home device for the treatment of mild to moderate inflammatory acne on the facial area. The device is safe to use on all skin tones. The ESPADA™ quickly, noninvasively and effectively eliminates acne-causing bacteria to clear blemishes and breakouts for skin that looks clearer, healthier and smoother.

Suitable for:

- Mild acne (whiteheads, blackheads and papules)
- Moderate inflammatory acne (inflamed, sensitive bumps filled with pus)

#### ESPADA™ FEATURES

1. Treatment surface
2. Antibacterial silicone
3. Start/stop button
4. Charging port
5. Indication light

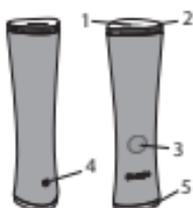- Blue (LED) light – For killing acne bacteria. Activates when device touches the skin. Flashes 3 times to indicate the start of each routine. Each time the ESPADA™ is switched on, it performs a self-check to ensure blue LED emissions remain at the appropriate level for use. If moved from the skin the light will automatically turn off.
- T-Sonic™ pulsations – when the device is turned on, its pulsations will be in a lower intensity and when it touches the skin, the intensity will increase so more blue light can penetrate the pores for increased effectiveness. After the device is moved from the skin, the pulsations will return to a lower intensity setting automatically.
- Targeting light (red cross-shaped light) – Targets the exact spot you wish to apply the treatment. Activates when device is away from the skin and turns off when the device touches the skin.
- Built-in timer – Every 30 seconds the device will pause its pulsations and flash to indicate the end of each routine. After indication, the device should be moved to the next treatment spot.
- If the device does not make contact with the skin when it is turned on, it will automatically switch off after 3 minutes.

**USING YOUR ESPADA™:** Cleanse and dry your face before using the device. Press the center button once to activate. Direct it over the blemish you wish to treat with the targeting light for ultimate precision. Gently place the device onto the blemish and hold it there for 30 seconds and wait for the built-in timer to indicate when the treatment is finished. Blemishes can be subjected to more than one 30 second routine, depending on severity. For mild acne, apply two 30 second routines. For moderate acne, 3 to 4 routines. Move the device to another blemish you wish to treat, if any, and repeat the procedure. There is no need to switch off the device between treatments.

**NOTE: SKIN WARMTH:** The skin may feel warm during treatment, which is normal and not cause for any concern. **TEMPORARY SKIN REDNESS:** The skin may develop a slight redness that fades within 24 hours. **TEMPORARY SKIN DISCOLORATION:** If the skin develops a slight discoloration that lasts more than 24 hours, stop using the ESPADA™.

**Before your 1st treatment perform a simple compatibility test.** Test the device on a less sensitive area of the skin, such as the inner arm, before beginning treatment on the face following the HOW TO USE instructions exactly. Wait 24 hours and check the treated area. If you don't notice any difference (redness, discoloration or irritation) between the treated and untreated skin, you may begin the Blue Light Acne Treatment.

#### WARNINGS:

- Consult your dermatologist prior to using the ESPADA™ if you have medical concerns.
- Do NOT use the device for severe, nodular or cystic acne.
- Do NOT use the device on or come in contact with injured skin or skin that has undergone surgery.

- Do NOT use the device on beauty spots/moles, warts, ingrown hairs or irritated skin.
- Do NOT use device in and around the eye and lip area.
- Do NOT use the device if you are particularly sensitive to light.
- Do NOT use the device if you have had sunburn within the past 2 weeks.
- Do NOT apply the treatment to a spot more than 3 times daily. Any more often than this could lead to skin burns or injuries. Leave a minimum time of 4 hours between each treatment.
- Do NOT use device longer than 2 minutes on one blemish. Overuse can lead to skin injury.
- Excessive use beyond the recommended treatment time may cause the device to become warm. The device's temperature will not go beyond 40°C if used at room temperature, but may cause discomfort or prolonged redness.
- If after using device you experience a worsening of the skin (skin feels tense, uncomfortable, irritated, red, burning, etc.), undesired side effects or prolonged redness, stop using the device and consult your doctor. However, a temporary, mild reddening after application is normal.
- Never look directly at the blue light as this could lead to eye injury.
- Do NOT let children use or play with the ESPADA™. The ESPADA™ has not been tested on children and may cause skin injury such as temporary redness or pigment changes (hypo- and hyperpigmentation) or affect their eyes.
- Proper supervision and instructions for safe use should be given to those with reduced physical and mental abilities.
- Do NOT share the ESPADA™ with other people. Sharing the device could transfer bacteria from one person to another.
- Do NOT use the device if it is visibly damaged, and never attempt to open or repair the device. Doing this could result in damage or malfunction to the device, serious eye or skin injury, fire, or possibly fatal electric shock.
- Pregnant or breastfeeding women should not use this device.
- Those who have had a migraine in the past two years should not use this device.
- Do NOT use over a thyroid gland as it may result in growth stimulation.
- If you are taking drugs or steroids for pain or skin conditions that cause light sensitivity, using this device may result in a photosensitive reaction.
- Do not use over or near a steroid injection received within the last six months. Doing so may result in a photosensitive reaction.
- In order to prevent fire or serious burns, immediately turn off device if it begins to overheat. Keep away from aerosol products.

#### **Contraindications**

- Do not use if you have a history of light triggered seizures.
- Do NOT use this device over a suspicious or cancerous lesion as it can delay proper medical treatment.

#### **Typical characteristics of potentially cancerous lesions include:**

- Asymmetry: one half of the abnormal skin area is different than the other half
- Borders: irregular borders
- Color: varies from one area to another with shades of tan, brown, or black (sometimes white, red, blue)
- Diameter: usually (but not always) larger than 6 mm in size (diameter of a pencil eraser)
- Bleeding: Any skin growth that bleeds or will not heal

**CAUTION:** Before charging, make sure the charging port and charger are free from water. Do NOT use your FOREO device while charging, and do NOT charge for more than 24 hours. Never use the device or charger if they are damaged, and only use the charger supplied by FOREO with the device. To reduce risk of electric shock, do not place the device where it can fall into water or other liquid. If device accidentally falls into water when charging, unplug immediately. The charging chord should be kept away from all heated surfaces. Never insert any object into any openings of the device. Never operate this device if its chord or plug has been damaged or is working improperly.

**MAINTENANCE AND SAFETY:** To keep your FOREO device in optimum condition, we recommend cleaning the ESPADA™'s treatment surface after each use. Clean the device with soap, rinse with warm water and dry with a soft, lint-free cloth. We recommend spraying the device with FOREO's Silicone Cleaning Spray for optimal results. NEVER use cleaning products that are silicone-based or containing alcohol, petrol or acetone as they may damage silicone.

#### **WARRANTY TERMS AND CONDITIONS**

**REGISTER WARRANTY:** To activate your 2-Year Limited Warranty, register the number provided on the FOREO Magnetic Card at [www.foreo.com/support](http://www.foreo.com/support).

**2-YEAR LIMITED WARRANTY:** FOREO warrants this device for a period of TWO (2) YEARS after the original date of purchase against defects due to faulty workmanship or materials arising from Normal Use of the device. The warranty covers working parts that affect the function of the device. It does NOT cover cosmetic deterioration caused by fair wear and tear, or damage caused by accident, misuse or neglect. Any attempt to open or take apart the device (or its accessories) will void the warranty. If you discover a defect and notify FOREO

during the warranty period, FOREO will, at its discretion, replace the device free of charge. Claims under warranty must be supported by reasonable evidence that the date of the claim is within the warranty period. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period. To claim your warranty, you must log in to your account at [www.foreo.com](http://www.foreo.com) then select the option to make a warranty claim. Shipping costs are nonrefundable. This undertaking is in addition to your statutory rights as a consumer and does not affect those rights in any way.

**DISPOSAL INFORMATION:** Disposal of old electronic equipment (applicable in the EU and other European countries with separate waste collection systems). The crossed-out dustbin symbol indicates that this device should not be treated as household waste, but rather be brought to the appropriate collection point for recycling of electrical and electronic equipment. Please see back page for instructions on how to remove the battery.

**BATTERY CARE:** It is recommended to regularly use the device to ensure proper performance. The ESPADA™ should be fully charged at least every 3 months to maintain the battery.

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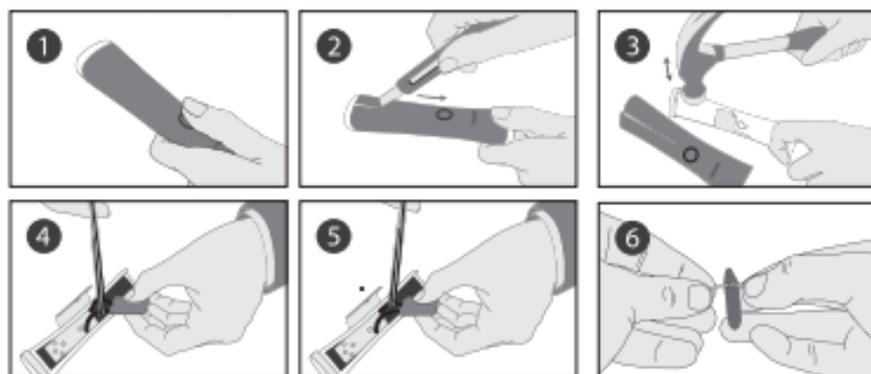
ESPADA™ complies with regulatory requirements of the following European Directive 93/42/EEC concerning medical devices and all the latest standards listed below.

EN60601-1/IEC60601-1, EN60601-1-2/IEC60601-1-2, IEC60601-1-11,  
EN60601-2-57/IEC60601-2-57, EN62471/IEC62471, EN ISO10993-5, EN ISO10993-10

Manufactured by Suzhou Sventech Co., Ltd. Address: 108 East Unit, No.77  
Suhong Middle Road, Suzhou Industry Park, Jiangsu Province, P.R. China  
Foreo AB, Birger Jarlsgatan 22,114 34, Stockholm, Sweden  
customer@foreo.com

MANUFACTURER

RESTRICTION OF HAZARDOUS  
SUBSTANCES COMPLIANT

AUTHORIZED REPRESENTATIVE  
IN THE EUROPEAN COMMUNITY

GARBAGE WASTE ELECTRICAL AND  
ELECTRONIC EQUIPMENT DIRECTIVE

TYPE BF APPLIED PART

IPX7 WATERPROOF

EUROPEAN CONFORMITY

CONSULT INSTRUCTIONS FOR USE

PACKAGING IS MADE OF  
RECYCLABLE MATERIAL

LS EQUIPMENT EMITS OPTICAL  
RADIATION WARNING

DC SOCKET POSITIVE AND NEGATIVE  
POLARITY /1.5VDC, 1.55W

SERIAL NUMBER

EU DISTR.: FOREO AB, BIRGER JARLSGATAN 22, 114 34 STOCKHOLM, SWEDEN.  
US DISTR.: FOREO INC., 3565 S. LAS VEGAS BLVD. #159, LAS VEGAS, NV 89109, USA.  
MANUFACTURED FOR / FABRIQUÉ PAR FOREO AB.

**FOREO.COM/SUPPORT**