
Body Cardio

Quick Installation Guide

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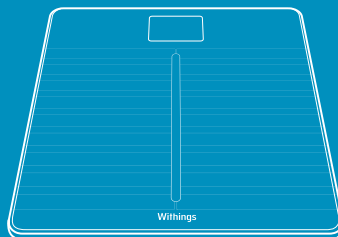
➔ <http://support.withings.com/bodycardio>

Thank you for choosing the Withings Body Cardio

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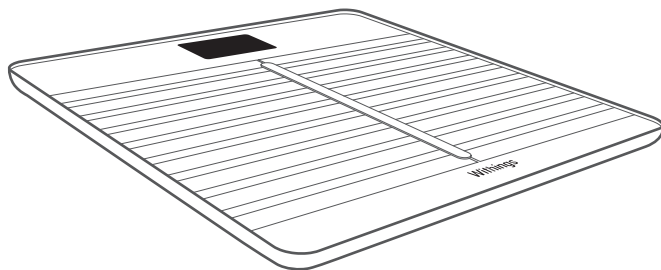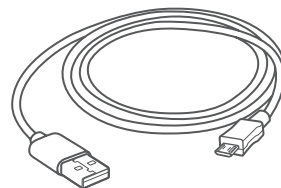
USB Cable

USB kabel | USB kábel
USB kábel

Connect the cable to a power source to charge your Body Cardio. Unplug the cable to weigh yourself again.

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

1

Press and hold the side button for a few seconds

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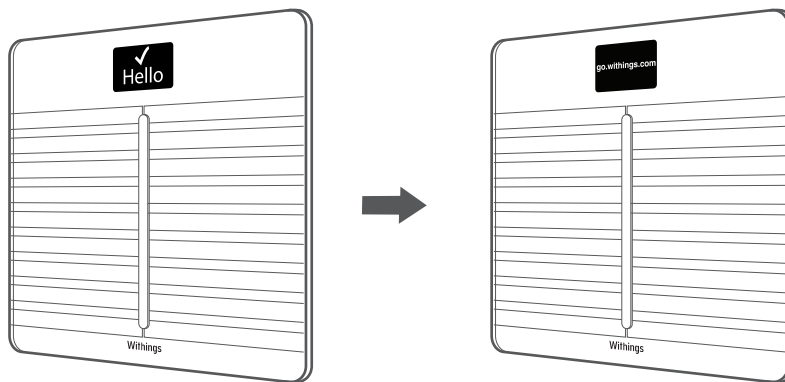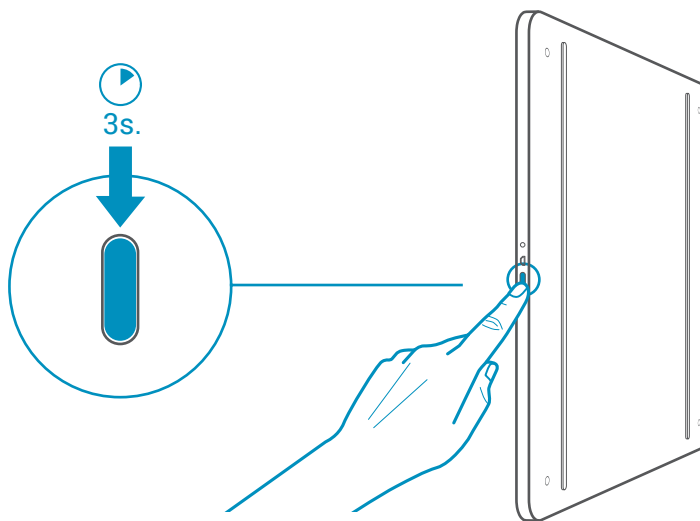2

Visit *go.withings.com*
to get the app and follow
the instructions

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
Analyze your Body Composition, Heart Rate, Pulse Wave Velocity and be notified of Weather Forecast for the day, while weighing

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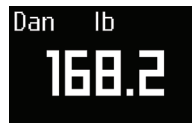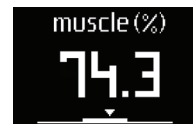
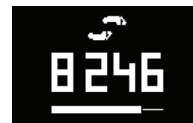
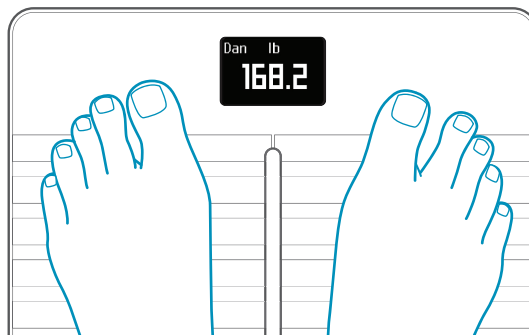
Change the order of the screens and choose what matters most to you

If you use a Withings activity tracker or track your steps with the Health Mate™ app, activate the Steps screen to know how much you walked the previous day!

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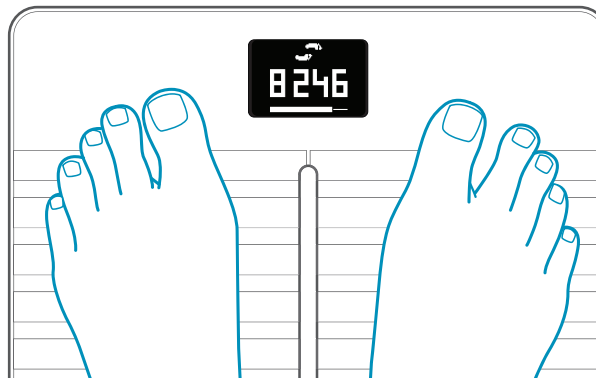
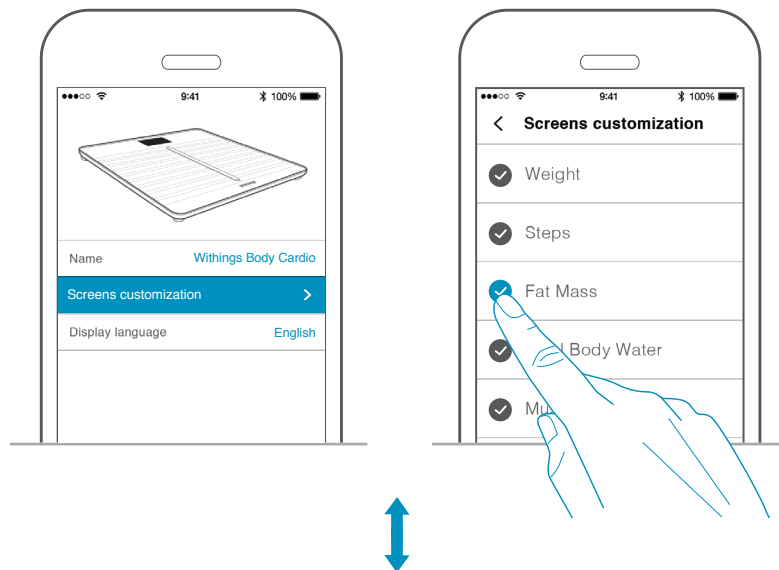
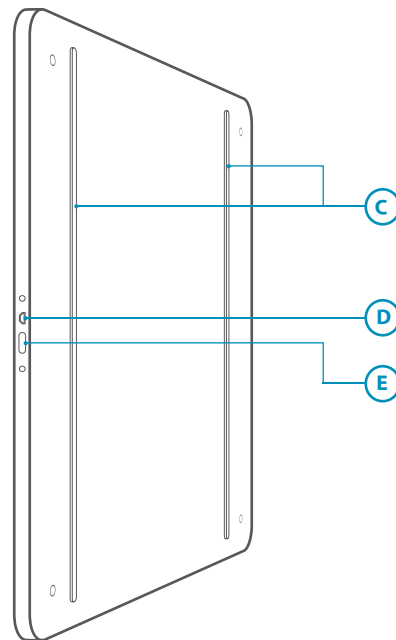
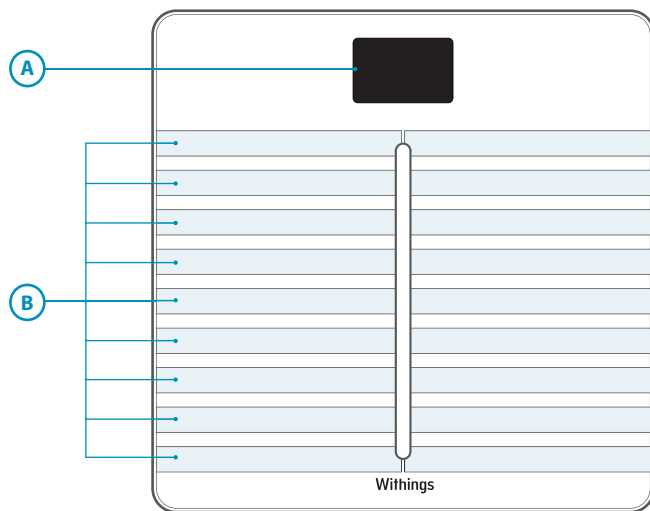
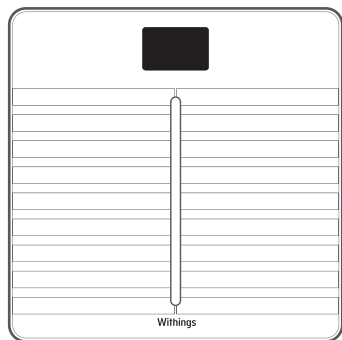
Use only a cloth dampened with water to clean your scale. The use of any detergents or cleaners may cause irreversible damage to the coating that protects the electrodes used to measure your body fat and heart rate.

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- This product is not suitable for people wearing a pacemaker or any other internal medical device.
- This product is not a medical device and should not be used to diagnose or treat any medical condition.

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US Withings™ One (1) Year Limited Warranty – Withings Body Cardio

Withings Inc. of 16192 Coastal Highway, Lewes, DE 19958 (“Withings”) warrants the Withings-branded hardware product (“Withings Product”) against defects in materials and workmanship when used normally in accordance with Withings’ published guidelines for a period of **ONE (1) YEAR** from the date of original retail purchase by the end-user purchaser («Warranty Period»). Withings’ published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product’s use.

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FCC ID: XNAWBS04 IC: 11411A-WBS04

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

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This device complies with the R&TTE Directive. A copy of the EU Declaration of Conformity is available online at www.withings.com/compliance

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2, rue Maurice Hartmann | 92130 Issy-les-Moulineaux | France

©2016 Withings. All rights reserved.

iPad, iPhone and iPod touch are trademarks of Apple Inc.

Android is a trademark of Google Inc.