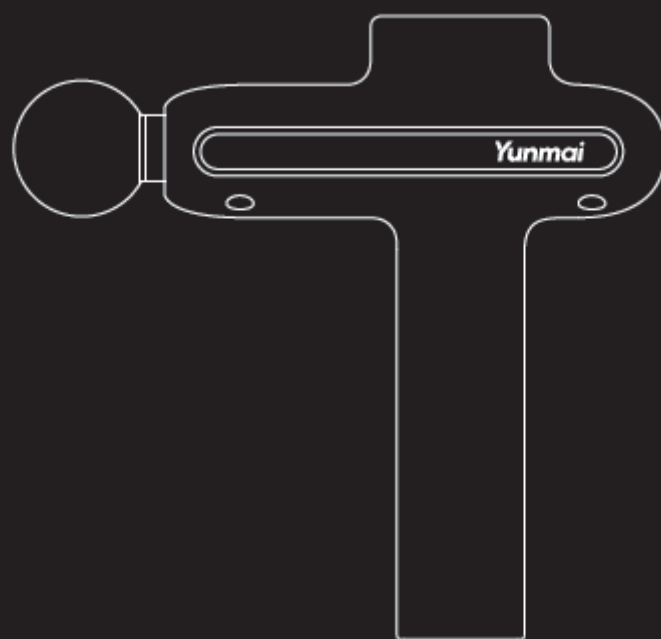
YUNMAI Fascia Massage Gun Pro Basic

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The product user guide contains product features, how to use, and the operating procedure. Read the user manual carefully to get the best experience and avoid unnecessary damage. Keep this manual for future reference. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Thank you for purchasing this product. Please read this manual carefully before use.
For future reference, please ensure this manual is properly kept.

Application scope:

The product can be used for physiotherapeutic massage before and after workout.

Key Information

To reduce risks of malfunction, electric shock, injury, fire, death and damage to equipment, please read all statements and warnings in this manual carefully before use.

Prohibitions

1. The use of this product by underaged and injured persons is strictly prohibited. In the event of muscle injury, please consult your physician before using this product.
2. Do not use this product on wet or dirty skin. This product can be used on clothed areas. Users can press this product on the skin gently, move it along different spots, letting it stay at each particular spot for about 60s.
3. Do not use this product on body parts above the neck, nor on body parts with less muscle or bone content. This product can be only used on soft tissues when the user does not experience pain or adverse reaction.
4. Do not use this product on body parts with bruises, pain, or adverse reaction.
5. Do not place foreign objects in the vent of the product.
6. Do not immerse the product in water. Do not allow water to enter the vent.
7. Do not place the product in or near sources of fire.
8. Do not throw or drop the product on purpose.
9. Do not recharge the product with chargers other than the one provided by Yunmai.
10. Do not disassemble the product.
11. Fingers, hair or other parts of the body may be trapped in the shaft and back side of the massage head. Do not come into close contact with these areas.
12. Not for children under 8 years old. Children are not allowed to use, clean, and maintain the product without supervision.

Do not touch this part to avoid pinching, when the machine is on without massage head installed.

Do not touch this part to avoid pinching, when the machine is on with massage head installed.

Statement

Those having any one of the following circumstances shall consult the doctor before using.

1. Pregnant women
2. Pacemaker wearer
3. Undergone surgery recently
4. Abnormalities in spines, joints or bones
5. Complications caused by other diseases (such as neuropathy or retinal damage)


Current research focuses on the influence of vibrational training on special medical conditions and may reduce the scope for the above-mentioned contraindications. Practical experience has shown that the combination of treatment and vibrational training also serves as a type of therapy in many circumstances. However, such an approach should be recommended by physicians, experts or therapists or be conducted in their presence.

Precautions


01. Unplug the charger after the completion of charging.
02. Do not use the product when it is charged.
03. Please use the product according to the instructions.
04. Before each use, carefully check whether the product can be used normally according to the operation methods stipulated in the instructions.

Product Introduction

Thank you for choosing Yunmai Fascial Gun Pro Basic

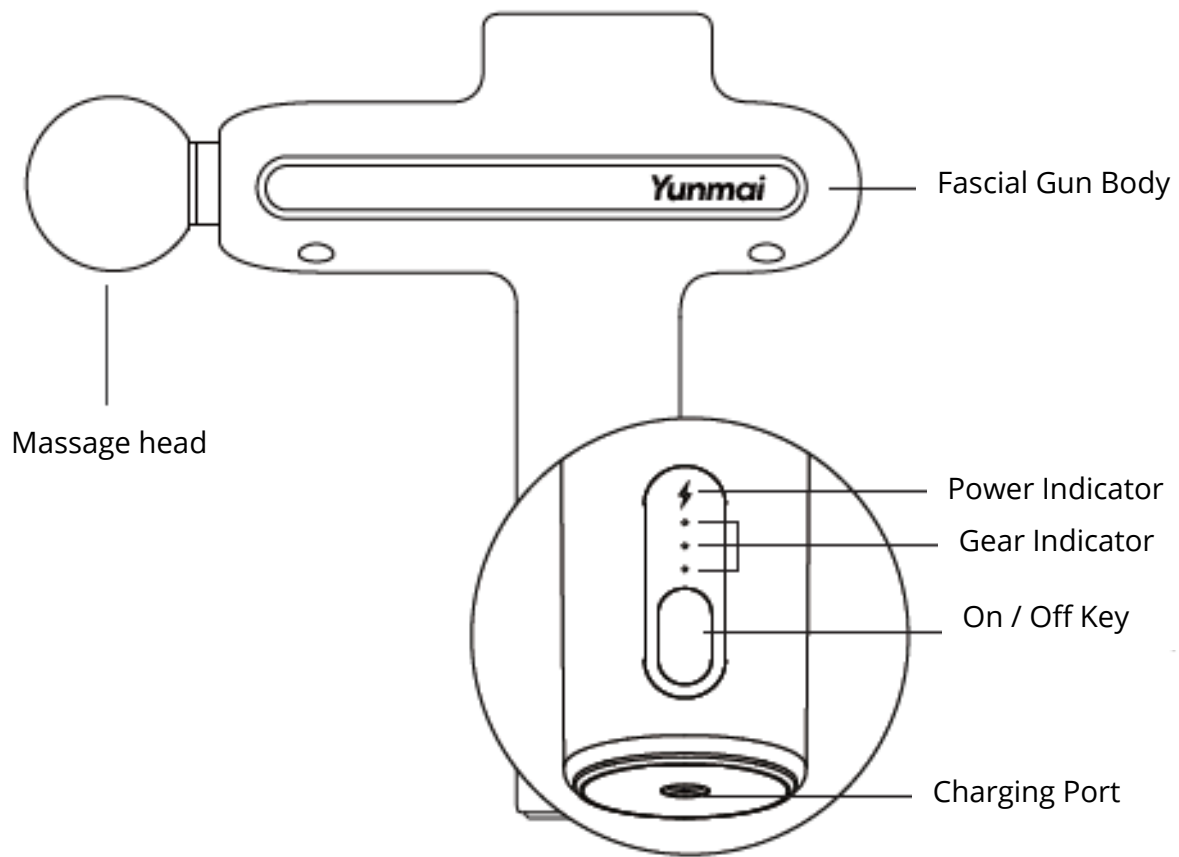Round Head
Mainly Used for large muscle groups

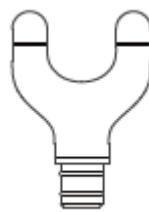U-shaped Head
Mainly used for the neck, both sides of spines and achilles tendons


Flat Head
Mainly used to relax and shape muscle parts

Bullet Head
Mainly used to impact deep tissues

Use Method

Operation

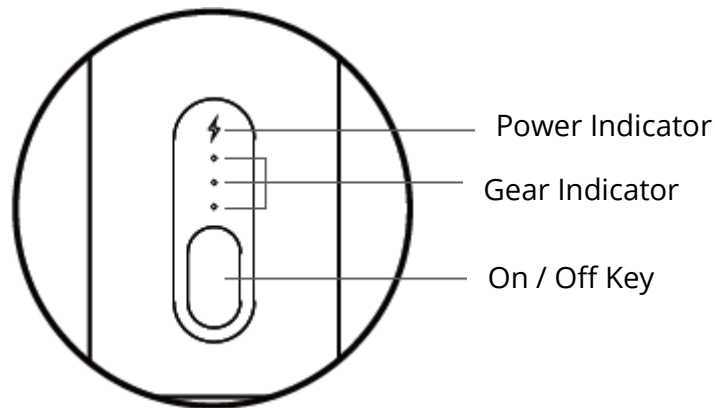1. Gear indicator: Gear 1 to Gear 3.


2. Power indicator:

- 1 When the indicator light turns white, the battery power is equal to or greater than 20%. to or greater than 20%.
2. When the indicator light turns red, the battery power is less than 20%, please charge in time. (In the process of charging, if the indicator light flashes white, it means the product is charging; if the white light is always, it means the product is charging; if the white light is always on, it means the product is fully charged)
3. When the indicator light flashes red, the battery power is less than 5%, and the product will stop working.


3. ON/OFF key:

1. After pressing and holding the ON/OFF key for 2s, the product will be turned on, and the power indicator and gear indicator will be on simultaneously to enter Gear 1 by default. Then press the ON/OFF key to enter Gear 2, followed by Gear 3, Neutral (Standby) and Gear 1 if pressed once again and again after the previous operation, and enter the recycle.
2. Press and hold for shutdown in any mode.
3. In neutral (standby) status, the product will automatically shut down in case of no operation for 15s. down in case of no operation for 15s.
4. The product will enter the intelligent power-off protection mode after continuous use of 10min. In this case, restart the product if you need to use it.

Massage Head Assembling

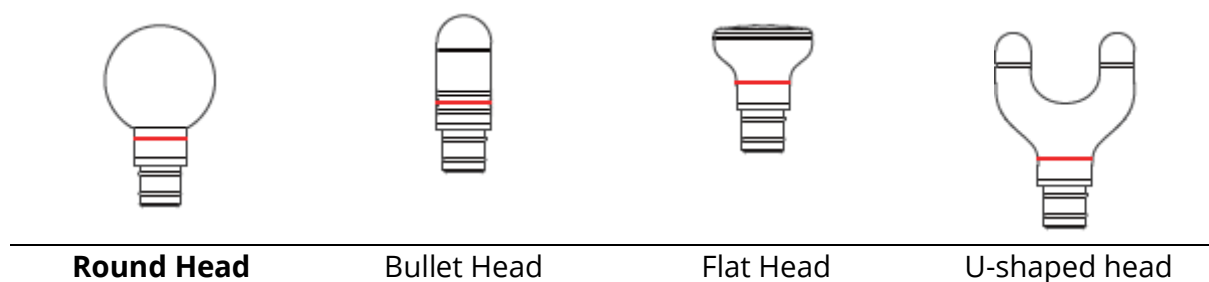
Notes: The red solid line in the figure is the line for alignment, which is bright black in the actual product. When assembling or replacing the massage head, please insert the massage head directly into the mouth of the fascial gun until the line for alignment coincides with the mouth (as shown in the figure); when replacing the massage head, pull out the massage head directly.

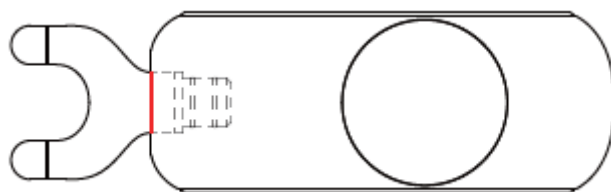Charging

The charging port is located at the bottom of the fascial gunner handle

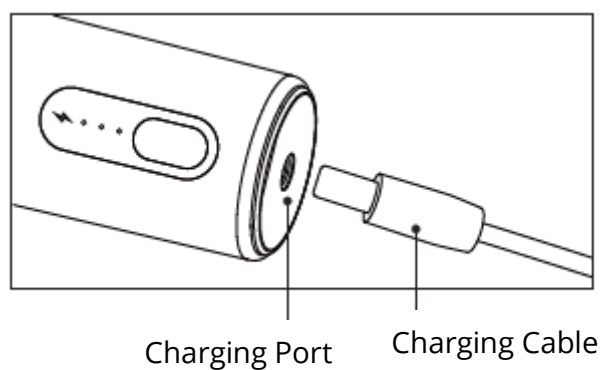1. Before the first use, please fully charge the battery of the product.
2. When recharging, connect the DC interface of 26V adapter to the recharging interface of the battery, and plug the AC cable in a 100-240V power socket.
3. When the power indicator on the product flashes regularly, it means the product is charging.
4. When the power indicator is always on, it means the product is fully charged.,
5. The battery can be charged at any time and at any power level, and the product cannot be started when charging.
6. If the power indicator turns red, please charge it.

Basic Parameters

Name: YUNMAJ Fascia Massage Gun Pro Basic	Model: YMJM-551S
Net weight: 1 kg	Input voltage: 100-240V-50/60Hz
Interface type: DC, 26V =1A	Battery capacity: 2600mAh
Intensity gears: 3 gears for adjustment	Charging time: About 4h
Power supply mode: Rechargeable Li-ion battery	Applicable environment: 0-40°C
Size: 183 X 230 X 71mm	Applicable age: Adults above the age of 18 (except pregnant women)

Frequently Asked Questions

Start-up failure

High temperature:

When used in the environment above 40°C, the fascia gun may trigger overheating protection, causing normal start-up failure of the fascia gun.

Solution:

Place the fascia gun in an environment below 40°C for 2h before starting it up for use.

2. Battery capacity is " 0 "

When the battery is out of power, if the user does not charge it for a long time, the battery will initiate low-battery protection, collectively referred to as over discharge protection.

The fascia gun cannot work normally if it is started at this time.

Solution:

Keep the fascia gun off, then connect the charger (power adapter), and the fascia gun will automatically enter the charging state. When it is fully charged, unplug the charger, and the fascia gun can work normally.

Charging failure

1. High fascia gun temperature:

When the temperature of the fascia gun body is too high, it cannot be charged when it is connected to the charger.

Solution:

Place the fascia gun in an environment below 40°C for 2h before charging.

2. High charger temperature:

When the fascia gun is charged in the environment above 40°C, the power adapter may trigger overheating protection, thus causing normal charging failure of the fascia gun.

Solution:

Place the fascia gun in an environment below 40°C for 2h before charging.

Troubleshooting of Fascia Gun Operation

Fault description	Troubleshooting
1.The massage head slides out and falls	Check whether the massage head is installed in place (see page 5 of the instructions for specific operation)
2.The product cannot be started	Check whether the power is sufficient (if the product has not been used for a long time, please charge it for more than 8 hours)
3. The product cannot be charged	Check whether the charging connector and power supply are plugged in place
4. The product stops suddenly after 10min of continuous operation	This phenomenon is normal as it is caused by the setting for long-term operation protection, which switches from operation to standby mode. And the product can be restarted through short pressing (see page 4 of the instructions for specific operation)
5. The product stops working after the massage head encounters the phenomenon of locked rotor	This phenomenon is normal as it is caused by the setting for abnormal operation protection, which switches from operation to standby mode. And the product can be restarted through short pressing (see page 4 of the instructions for specific operation)

Reasonable use of the product can effectively extend its service life

Suggestions:

- 1.Do not throw or cast the product during use.
- 2.Keep the maximum continuous service time lower than 30min.
- 3.Please charge the product every other month if it is not used for a long time.

Recycling

Do not discard this device together with general household waste when it is at the end of its service life.

Place it in the government-appointed recycling bin or send it back to after-sales service centre where workers will dispose of it in an environmentally friendly manner.

Maintenance, Cleaning, Storage, Transportation

1. Clean the device with a damp towel and wipe dry with a soft cloth.
2. Place the device inside the packaging box provided during storage or travelling.
3. This device should be stored in an area that is dry, well ventilated, and non-corrosive. Device should be placed no less than 250mm away from the floor, and the maximum number of layers on top of it should not be more than 5.
4. This device should be packed and stored properly if it is to be left unused for a long period of time. In such circumstances, the device should be recharged once every 6 months.
5. This device can be transported by car, train, sea, or air.
6. During transportation, this device should be gently handled and well-secured in an upright position to avoid damage due to strong vibration, collision, rolling or falling. Water-proofing measures should also be in place.

Name of parts	Hazardous and noxious substances					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Shell	0	0	0	0	0	0
Circuit board assembly	X	0	0	0	0	0
Metal parts	0	0	0	0	0	0
Rubber parts	0	0	0	0	0	0
Adapter	X	0	0	0	0	0
Connecting line	0	0	0	0	0	0
Handlebar grip	0	0	0	0	0	0

This table has been compiled in accordance with 5J/T 11364

✕: Circuit board assembly consists of parts of printed circuit board

machine, such as resistors, capacitors, integrated circuits, connectors.

0: It means that the content of hazardous substances contained in homogeneous materials of this device is within the limit specified in GB/T26572.

X : It means that the content of hazardous substances contained in some homogeneous material of this device is beyond the limit specified in GB/T26572.


Instructions on environmental service life

Hazardous substances or elements contained in this device will not leak if this device is used under the conditions stated in the manual. The device can be used for a period of 10 years without causing any environmental pollution or severe harm to the human body or property


Massage Position

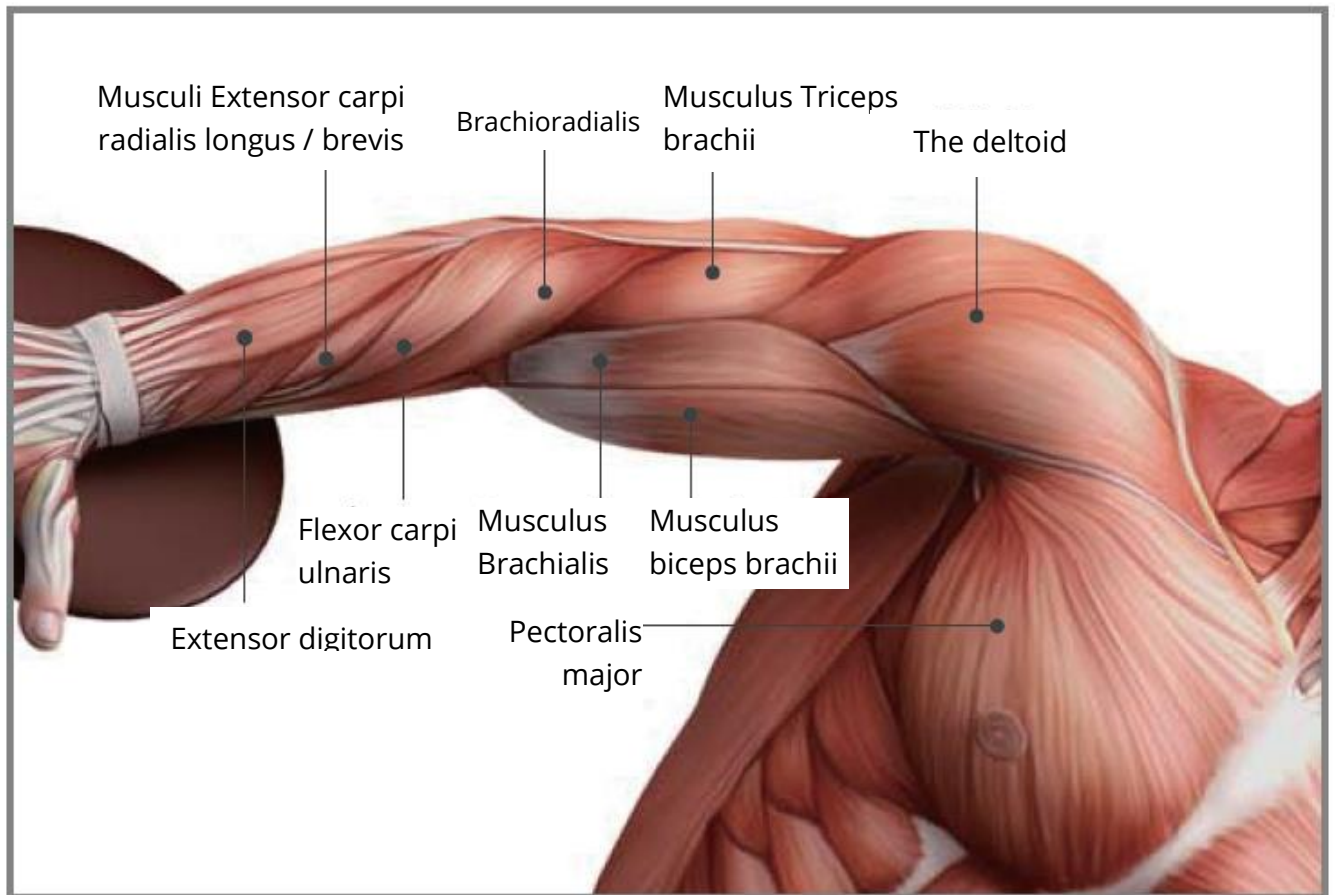Arms


Position	Duration (s)	Frequency
Brachioradialis muscle	10 – 15 s	3 times
Flexor carpi radialis	10 – 15 s	3 times
Palmaris longus	12 – 15 s	3 times
Flexor carpi ulnaris	10 – 15 s	3 times
Extensor digitorum	05 – 10 s	Twice
Extensor carpi longus radialis	05 – 10 s	Twice
Extensor carpi radialis brevis	05 – 10 s	Twice
The biceps brachii	12 – 20 s	3 times
Brachialis muscle	12 – 20 s	3 times
The deltoid	12 – 20 s	3 times
Pectoralis major	12 – 20 s	3 times

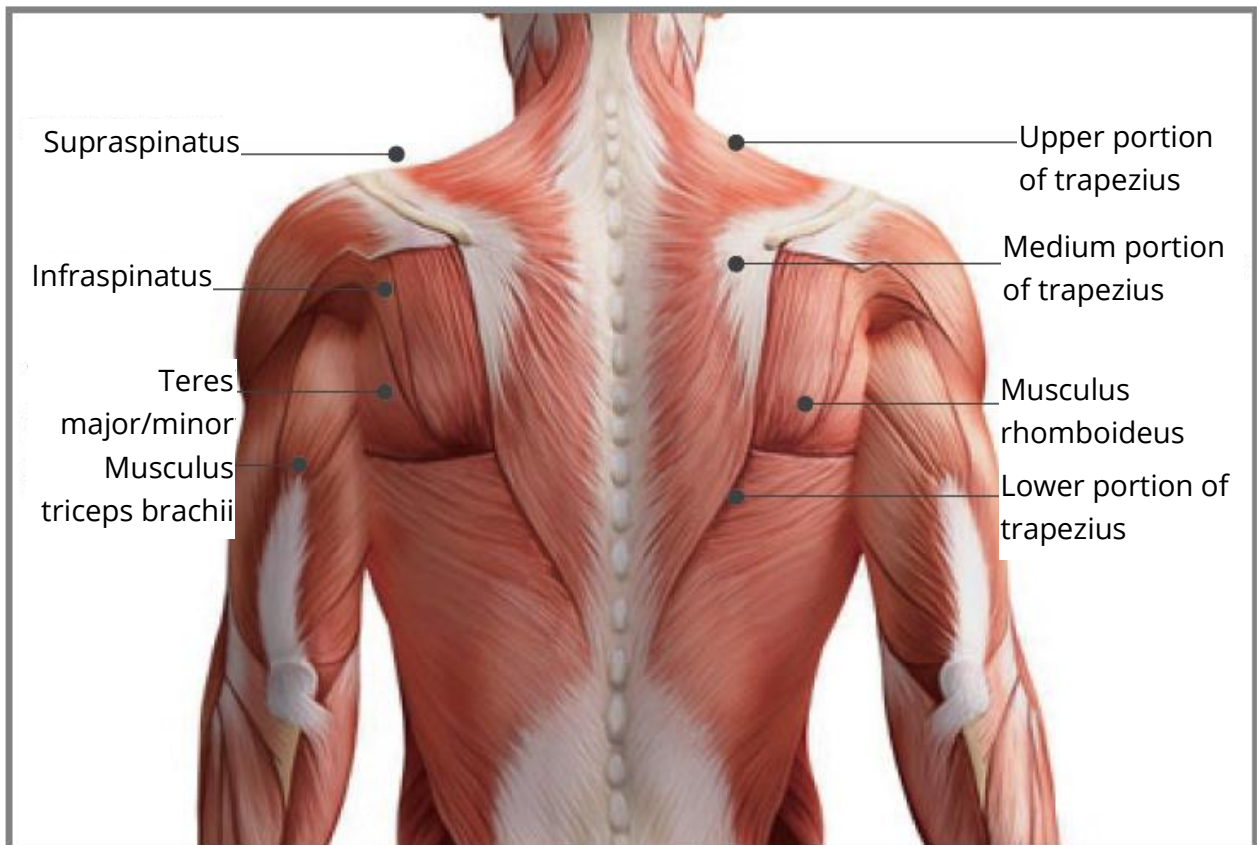Shoulder

Position	Duration (s)	Frequency
Triceps brachii	15 – 20 s	3 times
Musculus Rhomboideus	20 – 30 s	3 times
Upper portion of trapezius muscle	15 – 20 s	3 times
Middle portion of trapezius muscle	15 – 20 s	3 times
Lower portion of trapezius	15 – 20 s	3 times
Supraspinatus muscle	10 – 15 s	3 times
Elevator scapulae	10 – 12 s	Twice
Teres major muscle	10 – 12 s	Twice
Teres minor muscle	10 – 12 s	Twice
Infraspinatus muscle	10 – 12 s	Twice


Suggestions: Massage muscles from the middle to both sides and from the top to bottom.

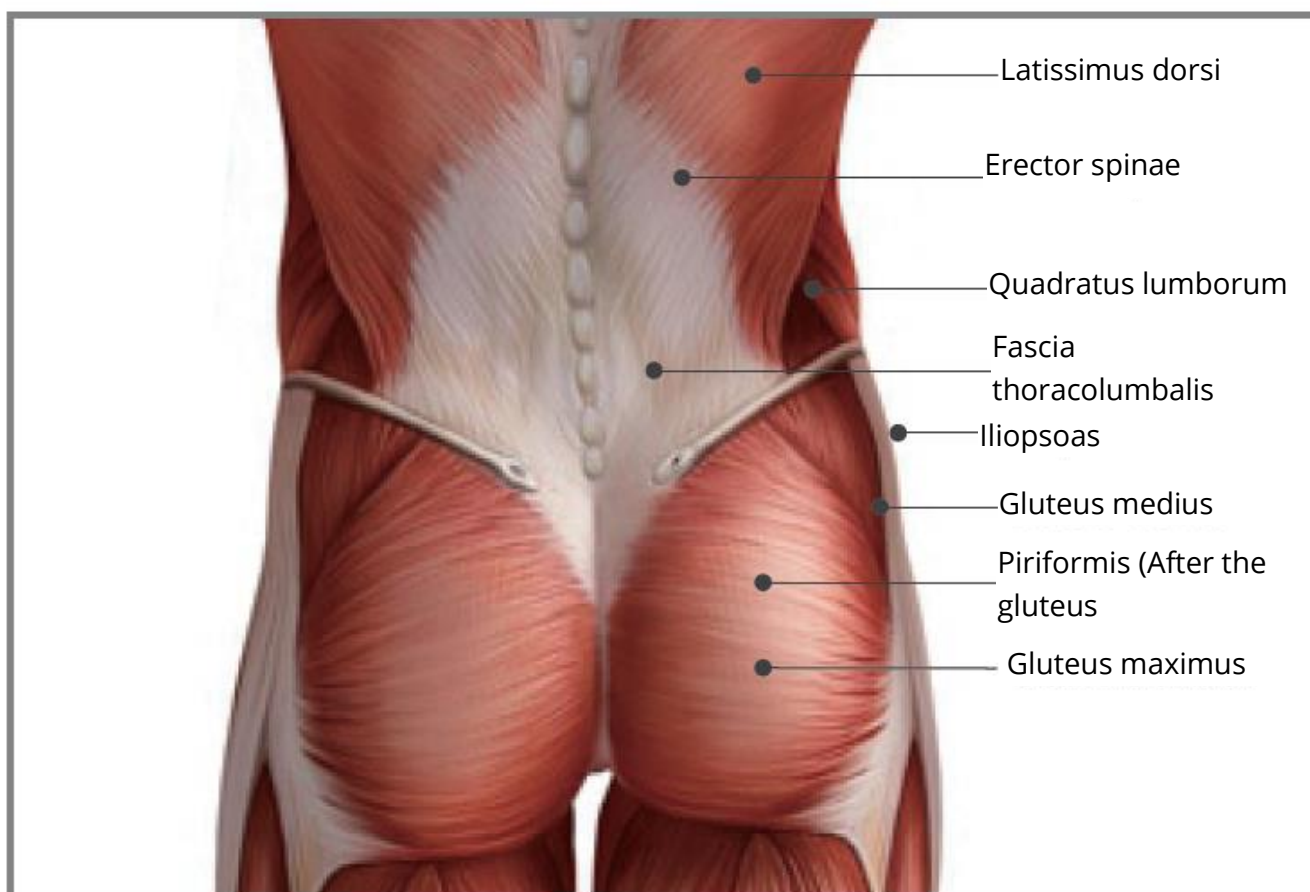Waist

Position	Duration (s)	Frequency
Erector spinae (lumbar segments)	30 – 60 s	4 times
Latissimus dorsi	30 – 60 s	4 times
Thoracolumbar fascia	40 – 60 s	3 times
Psoas major	30 – 60 s	3 times
Internal and external oblique muscle of abdomen	10 - 15 s	Twice
Erector spinalis (thoracic segments)	15 – 20 s	Twice
Quadratus psoas	20 – 30 s	Twice
Iliopsoas muscle	20 – 30 s	3 times
Gluteus maximus	30 – 60 s	4 times
Gluteus medius	30 – 60 s	4 times
Piriformis muscle	60 – 90 s	Once


Suggestion: Massage muscles from the middle to both sides and from the top to bottom. Move the device along the muscle when massaging large muscle groups.

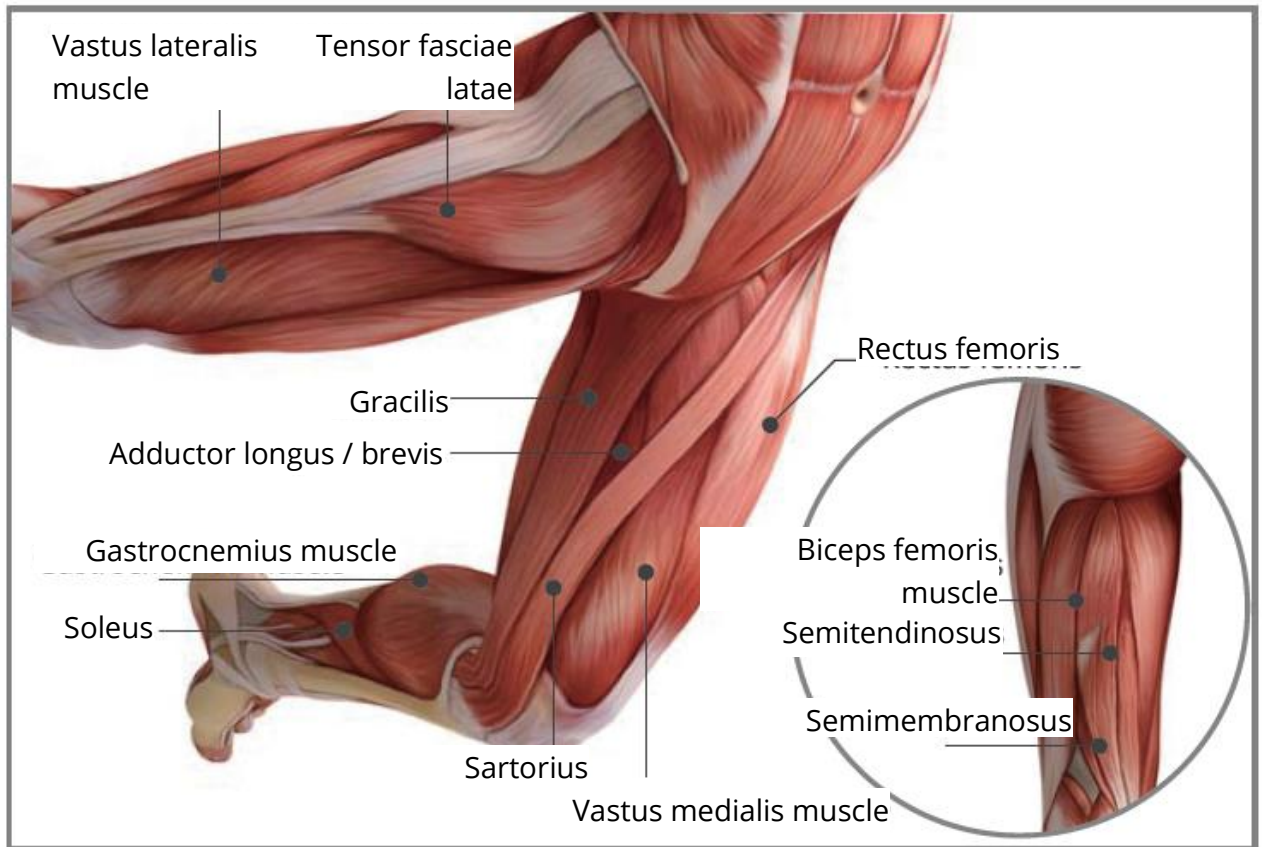Legs

Position	Duration (s)	Frequency
Biceps femoris muscle	20 – 30 s	3 times
Semitendinosus	20 – 30 s	3 times
Semi membranous muscle	20 – 30 s	3 times
Gracilis	10 – 15 s	Twice
Adductor brevis	10 – 15 s	Twice
Adductor longus	10 – 15 s	Twice
Rectus femoris	30 – 45 s	Twice
Vastus latarelis	30 – 45 s	Twice
Medial femoral muscle	30 – 45 s	Twice
Sartorius muscle	30 – 45 s	Once
Tensor fascia lata	15 – 20 s	Twice

Suggestions: Massage muscles from top to bottom.

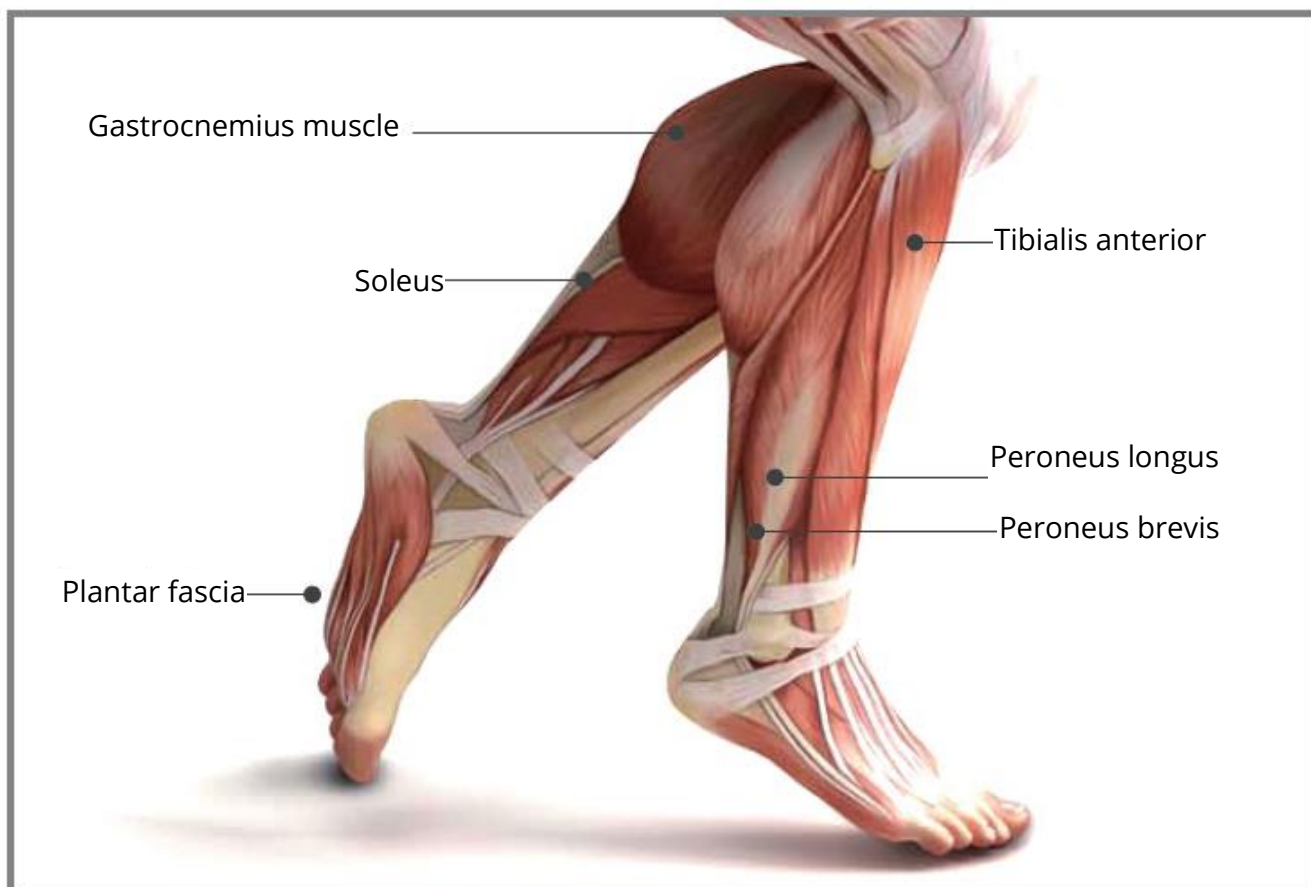Foot

Position	Duration (s)	Frequency
Gastrocnemius muscle	20 – 30 s	3 times
Soleus muscle	20 – 30 s	3 times
Tibialis anterior muscle	15 – 20 s	Twice
Peroneus longus	15 – 20 s	Twice
Peroneal brevis	15 – 20 s	Twice
Plantar fascia	20 – 30s	3 times

Suggestions: Massage muscles from top to bottom.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g. during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

IČO: 27082440

Subject of the declaration:

Title: Massage Gun

Model / Type: Yunmai Pro Basic

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. (EU) 2014/53/EU

Directive No. (EU) 2014/35/EU

Directive No. (EU) 2011/65/EU as amended 2015/863/EU

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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