X10

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important Safety Tips

- This product is recommended to use by person over 12 years of age.
- The product can only be used with Oclean replacement heads. Oclean bears no responsibility for any consequences resulting from the use of non-certified third-party brush heads.
- This product is for personal cleaning only and is not recommended for medical purposes and multi-user use.
- This product is for indoor use only. Do not place the toothbrush under direct sunlight for a prolonged period of time. It is recommended to be placed in the dry well-ventilated areas with normal temperature.
- Please keep the toothbrush body dry before charging and wipe it dry with a soft cloth. After charging, please install the bottom cover back to its original place and ensure that it is sealed and waterproof. Oclean is not responsible for any consequences caused by improper operation.
- It must be only supplied at safety extra low voltage corresponding to the marking on the product.
Instructions for products containing non-user replaceable batteries shall state the substance of the following:
- This product contains batteries that are only replaced by trained professionals.

Brush head

When replacing the brush head, please pull out and insert the replacement brush head vertically, and do not rotate left or right.

It is recommended by dentists to change the brush head every two to three months.

Instruction for use

- Up and down button control
The toothbrush can be set via buttons, tap the up button to switch the mode, and tap the down button to switch the strength level. The toothbrush can be selected from 5 modes: intense whitening, specialized cleaning, sensitive mode, polishing, and gum massaging. The toothbrush has 5 levels of strength to choose from.
- ON/OFF button control
Once the toothbrush turns on, tap the ON/OFF button to brush your teeth. Tap the button during its running state will cause it to pause; then tap again within 15 seconds will resume its operation; without operation for more than 15 seconds will turn off the display and hibernate the toothbrush. When expression is displayed, the toothbrush will hibernate directly with a tap on the ON/OFF button.
To reset language setting, please press the “ON” button for 8 seconds to revert to factory settings.
- Brushing program
The toothbrush will run for 2 minutes by default “intense whitening mode”, the toothbrush will remind the user of to change the brushing area every 30 seconds by vibrating and will automatically stop running after 2 minutes of operation.
- Smart brushing feedback
After brushing, the screen on the electric toothbrush will display the score and feedback on the results of the brushing, and the score is based on star scoring system. When the brushing section is scored 3 stars, the screen of the toothbrush will display an emoji, and ends the brushing. If the brushing section is scored less than 3 stars (Including 2 stars, 1 star and no star), the display of the toothbrush will then display a model of the teeth to remind the areas that need additional brushing and will automatically generate the brushing time the additional brushing it requires. After reminding for additional brushing, you will need to follow the tips and taps the up or down button to turn on the re-brushing mode. After the additional brushing is fully completed, the toothbrush will show and expression on the display, and the brushing process will end.

- Charging the toothbrush

This product adopts the Type-C charging mechanism. When charging is required, connect the charging cable to the power plug, remove the protection cover at the bottom, and insert the Type-C plug on one end of the charging cable into the USB port at the bottom of the toothbrush for charging. When fully charged, please unplug the charging cable as soon as possible and install the protection cover at the bottom. (As shown in the figure)

Leaving the fully charged toothbrush on the live charging cable for a prolonged period of time may lead to premature depletion of the built-in battery, thus affecting the usage experience.

Thanks for using Oclean Smart Sonic Electric Toothbrush.

Attention should be paid to those with the following conditions of oral sensitivities

- Those who have recently undergone periodontal surgery or are undergoing Oral-related treatments. Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy gum bleeding.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with above conditions and should proceed to consult a professional dentist before use.
- This product is not suitable for person with disabilities, mental illness or lack relevant experience and knowledge [including children (0-12 years old)]. Unless under the supervision of their guardians.

- New Message Icon
A new message icon will appear on the display when there is new activity message, software update notification or another alert information.

- **Low Battery Alert**

When the battery percentage is below 20%, an alert will pop up on the display of your toothbrush from time to time to indicate that the battery is running low, and you need to charge it soon. When the battery is seriously low, and the percentage is lower than 1%, your toothbrush will be off, and the toothbrush's motor will vibrate to alert you to charge the battery, and the screen will not be turned on either.

How to Charge

When charging is required, remove of the protection cover at the bottom of the toothbrush, connect the charging cable to the power plug, and insert the Type-C charging cable connected to the power plug into the charging port at the bottom of the toothbrush. During normal charging, the screen on the toothbrush will display a charging interface. Generally, and it will be fully charged in about 3 hours.

Recommended tooth brushing method

Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.


Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.

Please rinse with water after brushing.

Cleaning and Maintenance

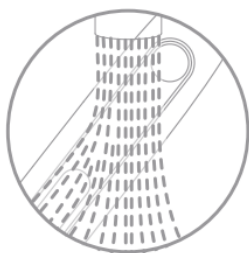Cleaning the toothbrush body and the brush head.
Rinse under the running water, put cover on the brush head bristles. After washing, air dry the brush head and do not pull the bristles.

Product Specifications

Product Name	Oclean XS Smart Sonic Electric Toothbrush
Model	R3100
IP Rating	IPX7
Input Power	2.5W
Charging	Type-C charging
Input	5V $\overline{=}$ 0.5A
Battery Life*	60 days
Battery Capacity	1200mAh
Charging Time	3 hours
Implementation standard	Q/SZYD 1-2021

* Calculated under whitening mode, 2 minutes/time, 2 times/day.

- Do NOT disassemble, strike, extrude or throw the device into a fire.
- Do NOT use the toothbrush if its battery is building or has expanded.
- Do NOT place in a high-temperature environment.
- Do NOT use the device if its battery has been subjected to immersion in water.
- Before the appliance is discarded, its battery must be removed from it.
- Before the battery is removed, the appliance must be powered off.
- The battery shall be disposed of safely.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: XS Smart Sonic Electric Toothbrush

Model / Type: R3100

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/30/EU


Directive No. 2014/35/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

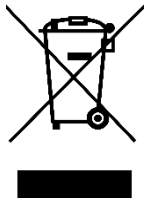This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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