Air 2

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important Safety Tips

- Due to the high-speed vibrations, it is normal to experience a tingling sensation while using the toothbrush for the first time. This feeling will gradually go away after 7 to 14 days.
- If the toothbrush is not used for a long time, the battery will run out due to the power consumption from standby mode. This can be solved by charging, but the charging indicator on the brush will not light up instantly. Please wait for 6 – 10 minutes for the toothbrush go back to the normal state.
- It is normal for the toothbrush to heat up when it is charging. If the toothbrush becomes extremely hot, immediately stop charging and/or turn it off and send it to Oclean's authorized customer service centre for checking.
- The toothbrush can only be used with the Oclean special heads. Oclean will not be able to take any consequences caused by non-official brush head.
- The toothbrush is for oral cleaning only. Please do not use it for other purposes.
- The toothbrush is a personal cleaning item, please avoid sharing it with others.
- Please store the toothbrush in a cool and dry place rather than anywhere exposes to direct sunlight.

Caution for Oral Sensitive Person

- The following 3 types of people are not recommended to use this toothbrush unless you have already consulted the professional dentists.
- The person who has undergone or is under treatment for periodontal surgery.
- The person with poor fitting dentures, continuously bleeding gums or any other oral related symptoms.
- The person whose oral is sensitive.
- Unless directed by a guardian, otherwise this toothbrush is not suitable for people with reduced mental and physical capabilities, and those who lack related experience, and the children under the age of 12 are included.

Notes:

- Please shut off the power when the charging base is needles or before cleaning.
- Please regularly clean the power plug to prevent the accumulation of dust and dirt which can reduce the charging effect.

Interaction Introduction

Power On

Start the toothbrush by short pressing the power button, the toothbrush will automatically shut down after a 2-mins brushing plan is finished.


Intensity Switch

Oclean Air2 has present 2 types of 2mins brushing plan: Standard Cleaning (Default) and Deep Cleaning.

To switch the intensity, just long press the power button for 2s under working condition.

Battery Indication


Please charge the toothbrush when the red indicator light is on.

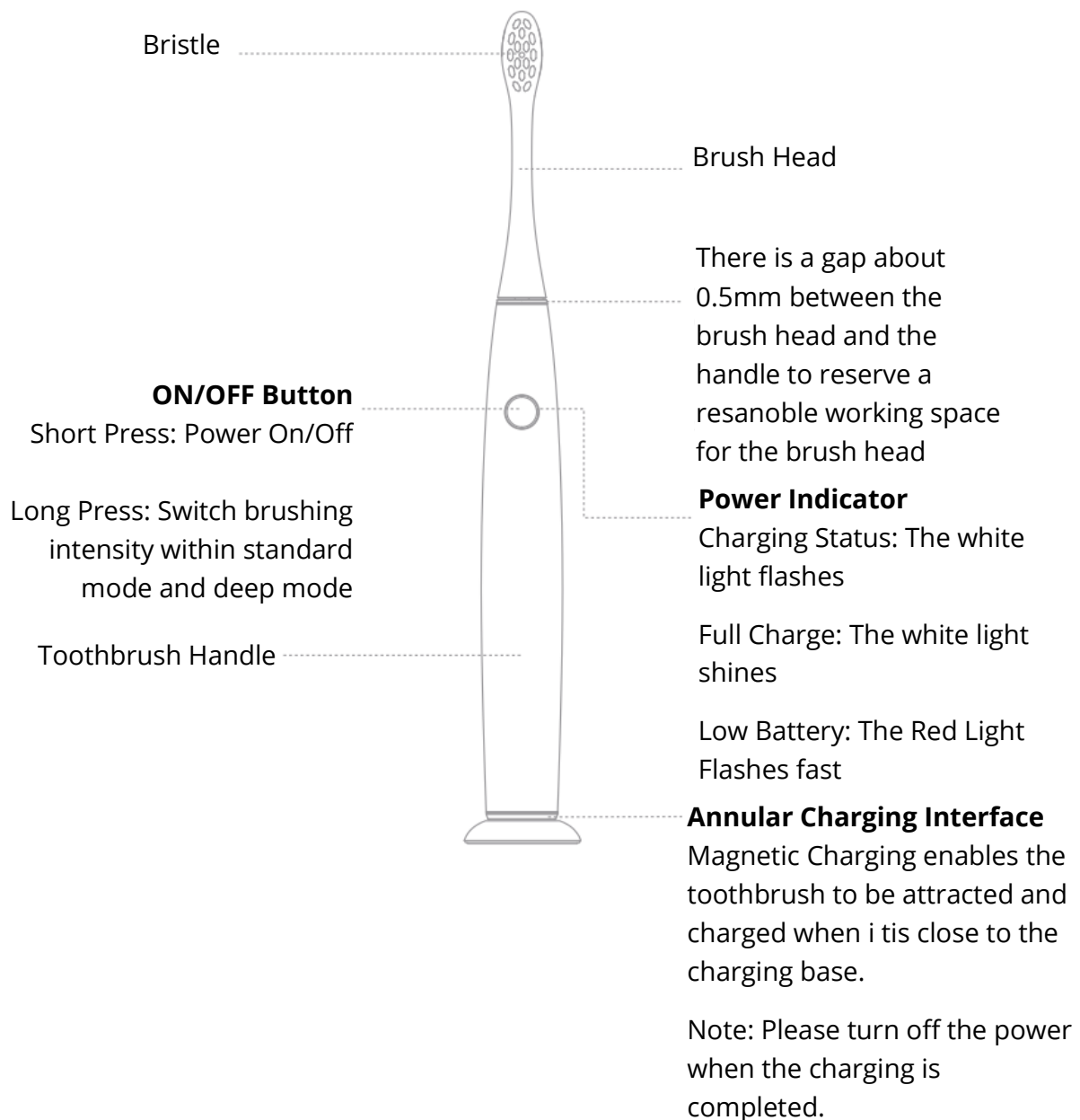
Recommended Tooth Brushing Method

Brushing Time

Oral experts (Dentists) recommended that each brushing should be 2 to 3 minutes long.

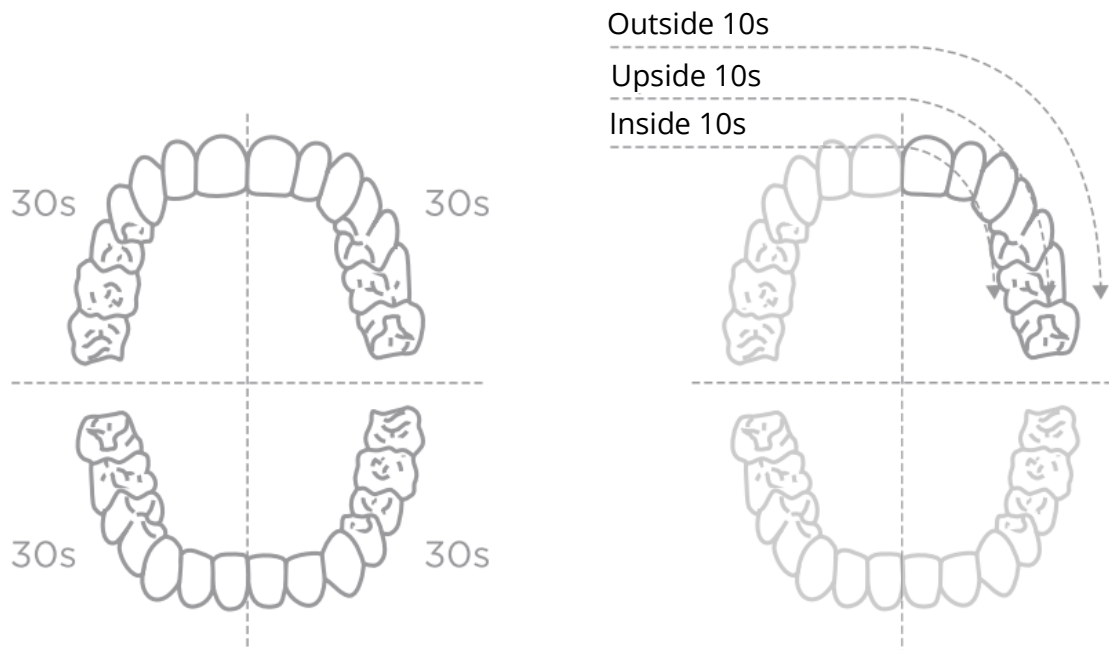
Recommended Gesture for Different Tooth Region

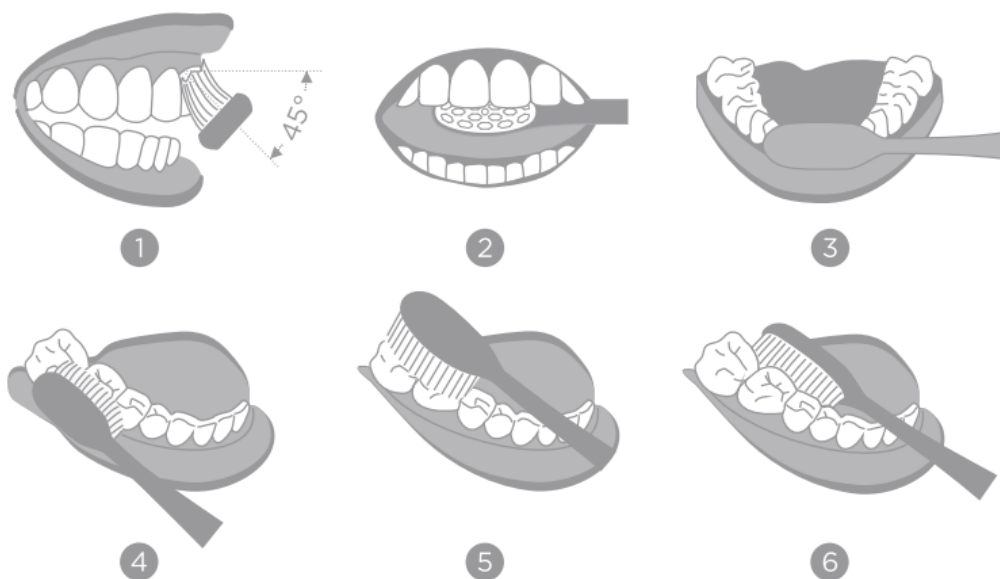Using Instructions

Firstly, please insert the brush head forwardly onto the toothbrush body well. A small gap of about 0.5mm should be left between the head and the main body so the brush head can function properly.

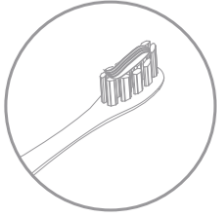Secondly, please turn on the toothbrush when it has been put into your mouth with a small amount of soaked toothpaste to prevent toothpaste from splashing.

Thirdly, please make sure the bristles at 45° from the gum line and move the toothbrush slowly.

Fourthly, hold the toothbrush handle and slowly move it in a circular motion. Do not press too hard.

Finally, please gargle with clean water after brushing.

Specifications

Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air 2
Input Power	2.5W
Charging	Oclean USB Charge
Input	5V $\overline{\text{---}}$ 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800mAh (**)
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V $\overline{\text{---}}$ 0.5A.

Charging Process

0%-5%, the power indicator light flashes.

95%-100%, the Power Indicator light is constantly on.

First Use

Charging time depends on the battery life demand of user.

Working Process

If battery power is less than 5%, the red indicator will flash as charging reminder.

Notes:

The battery cannot be replaced by users themselves and can only be replaced by the staff in Oclean Authorized Client Centre.


Unauthorized bodies are forbidden to disassemble this product. The USB cable cannot be replaced by unauthorized body. If the product was damaged, please contact the Oclean After-sale Centre for maintenance.

Please use the USB charger of 5V/1A or 5V/2A of mobile phone to charge. Such as the iPhone Charger, Tablet Charger etc.

If the charge current is less than 1A, it may be prolonging the charging time.


If the charger voltage is higher than 5V, it may damage the charging circuit.

Cleaning and Maintenance

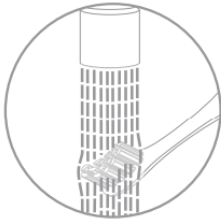
Cleaning The Brush Head

Clean the brush head after use and air dry it. Do not pull the bristles, because it may result in deformation of the brush head.

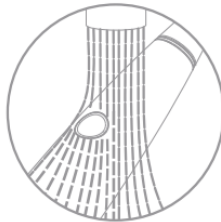Cleaning The Body

Use clean water to wash the surface of handed and wipe it with soft cloth and air dry it. Do not press the rubber sealing ring with sharp item or metal object, otherwise it will be damaged.

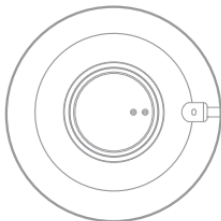Cleaning The Charger

Wipe the dirty parts of the charger with a soft cloth. Pull out the battery plug from the socket prior to cleaning charger. Wipe the surface of charger with wet cloth.

Charging The Head

It is suggested by dentist that the head should be changed every 2-3 months to maintain optimal brushing.

Common Problems and Solutions

The following are some common problems that users can face while using the toothbrush.

These problems can be easily solved by the users with the solutions given below. Please conduct the following trouble shooting techniques prior to contacting professional after-sale services.

1. The vibrations strength becomes weak

Possible Problems:

- The toothbrush is pressed too hard against the tooth.
- The gap between the brush head and body is too loose or too

Solutions:

Lightly place the brush head on tooth. Pull the brush head from the body or gently insert the brush head leaving enough space between the brush head and body.

2. The toothbrush functioning for short time

Possible Problems:

- The charging cradle is not in contacted with the charging ring of the body and may have resulted in improper charging which can also be detected by the LED light.
- Insufficient battery charge.

Solutions:

- Clean the charger cradle.
- Charge for longer period.

Battery has completed its service life. If the red light of the foundation is on, please contact the Oclean Official Authorized Centre for a replacement charger.

3. The product is unable to work

Possible Problem

- Insufficient battery power.

Solutions:

- Please charge the toothbrush
- Please contact Oclean Official Support Centre for a replacement battery if the first was doesn't help.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Sonic Electric Toothbrush

Model / Type: Oclean Air 2

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/30/EU


Directive No. 2014/35/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

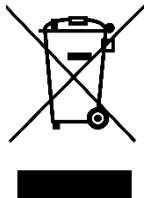This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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