Flow

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411


Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important Safety Notice

- When using the product for the first time, due to the high-speed swinging of brush, if it feels simulating or numb, please do not worry, it is normal, and the feeling will weaken and disappear gradually in around 7-14 days.
- If not using for long time (the product battery has no power due to the standby power consumption), when charging for the toothbrush again, the charging indicator doesn't show, please wait for 6-10 minutes patiently, it can recover normal charging.
- If the main part of toothbrush becomes slightly hot when charging, it is normal. If you feel hot with your hand, please immediately stop charging and submit the toothbrush to Oclean's After-sales Services Centre for testing.
- The product can only use the special brush head for Oclean Sonic electric toothbrush. Oclean will not be held liable for all the consequences caused using other unofficial or non-licensed brush head products.
- The product is only used for cleaning teeth and must not be used for other purposes.
- The product is only intended as a personal cleaning product. It is not recommended to use it by multiple persons simultaneously.
- The product is only used indoors, do not place the toothbrush under the blazing sun for quite a long time. It shall be placed in a dry and well-ventilated place with normal temperature.

Attention should be paid for those with the following conditions of oral sensitivities

- Those who have recently undergone periodontal surgery or are undergoing treatments.
- Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy bleeding from the gums.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with above group or consult a professional dentist before use.
- This product is not suitable for people with disabilities, mental illness or lack of relevant experience and knowledge (including children 0-12 years old). Unless used under the supervision of their guardians.

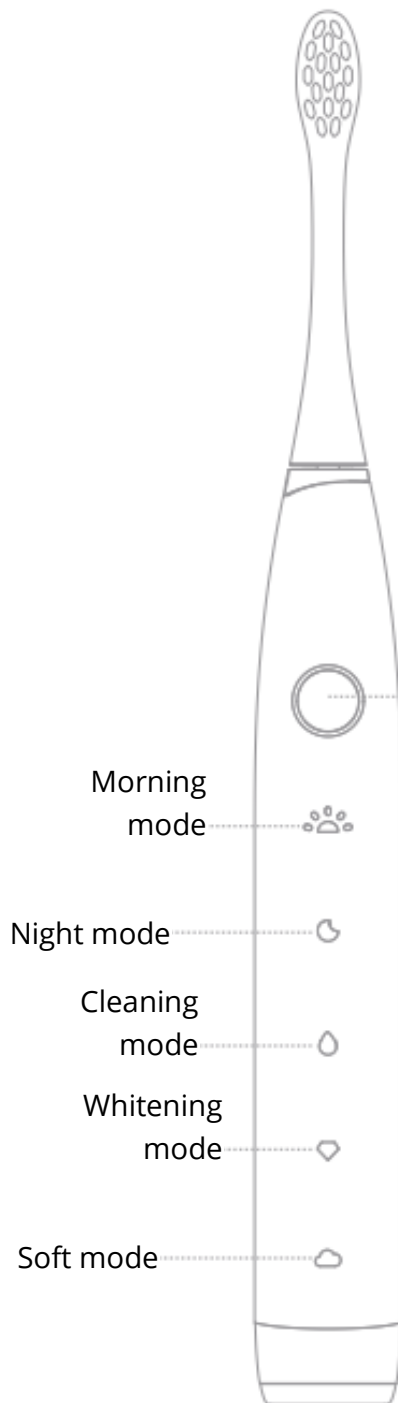ON/OFF button control

Enables:
short press the button

Switch mode:
Within 5 seconds after the button is short pressed, you can switch between modes by rapidly short pressing the button when the mode indicator is in the state of "breathing light"

Turn off:

1. You can turn this product off by short pressing the button again after circularly short pressing the button until this product operates in the "soft mode".

2. You can turn off this product by short pressing the button again when the mode indicator is in the "normally on" state within 5 seconds after short pressing the button to start up the toothbrush.

3. This product is automatically turned off after it is in operation for 2 minutes after you short press the button to start up it.

Low Battery reminder: When the battery level is lower than 10% the button circular lamp is in the state of "flickering in red" after you short press the button to turn on the toothbrush.

Instruction for use

Default Brushing program

Five toothbrushing modes are built in this product: Morning mode, night mode, cleaning mode, whitening mode and soft mode.

- Morning mode: compared with “night mode”, the pressure is low. As a cleaning and care tool, it can effectively remove dental plaque caused during the sleep.
- Night mode: deeply remove food residues and dental plaque on the teeth with high pressure.
- Cleaning mode: meet requirements for routine and ordinary cleaning. Compared with “night mode”, the pressure is relatively tender.
- Whitening mode: effectively clean bite marks in those who have habits of smoking or frequently drink coffee, tea, etc.
- Soft mode: protect the health of gingiva while fulfilling the fundamental function or removing dental plaque in those who have oral allergy.

Type-C direct connection charging

Open the bottom cover of the electric toothbrush, connect Type-C interface into the port at the bottom of the toothbrush, connect one end of the USB interface to the power supply, the button circular lamp flickers in white and this indicated that the toothbrush is being charged.

Charging process

Charging status: the button circular lamp is on in white and flickers.

Use for the first time

The user may charge the toothbrush depending on his or her need for battery life to ensure normal use of this product.

Application process

When the battery level is lower than 10%, the button circular lamp is in the state of “flickering in red” after you short press the button to turn on the toothbrush.

Precautions

Battery cannot be replaced by users themselves, it can only be replaced by personnel in the customer centre authorized by Oclean, non-specialized personnel cannot dismantle the product by themselves, or it will be regarded as that the user give up the warranty automatically.

Please charge the Oclean toothbrush with 5V/1A or 5V/2A mobile phone charger.


Such mobile phone charger includes iPhone charger, tablet charger, etc. If the charging current is lower than 1A, this may cause charging time to be extended, if you charge the toothbrush with a mobile phone charger with a voltage of higher than 5V, this may cause the charger circuit to be damaged.

- Do NOT disassemble, strike, extrude or throw the device into a fire.
- Do NOT use the toothbrush if its battery is bulging or has battery to a high-temperature environment.
- Do NOT subject the battery to a high-temperature environment.
- Do NOT use the battery if it has been subjected to moisture.

Before the appliance is discarded, its battery must be removed from it. Before the battery is removed, the appliance must be powered off. The battery shall be disposed of safely.

Recommended tooth brushing method

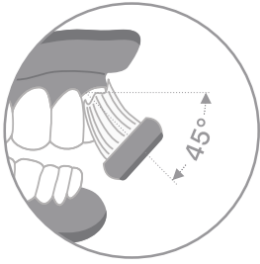Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.


Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.

Please gargle with clear water after brushing.

Cleaning and Maintenance

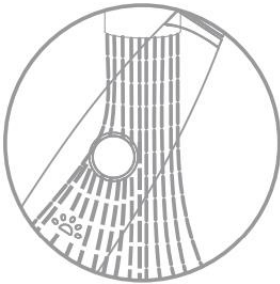Rinse with tap water and air dry, and do not pull the bristles of the brush head

Product specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	F5002
IP Rating	IPX7
Power of Subject	1.85W
Charger	Type-C direct Recharge
Input	5V $\equiv$ 1A
Service Time*	$\leq$ 180 Days (*)
Battery Capacity	2500mAh (**)
Charging time	$\leq$ 5 hours (***)
Implementation standard	GB 4706.1-2005 GB 4706.59-2008

* Validated by 2 brushing processes per day and 2 minutes of each process. (Morning mode)

** The Battery capacity theoretically has a capacity tolerance of $\pm 50\text{mAh}$.

*** It is suggested to adopt the USB charger or 5 $\equiv$ 1A in the phone industry.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Sonic Electric Toothbrush

Model / Type: F5002

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/30/EU


Directive No. 2014/35/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

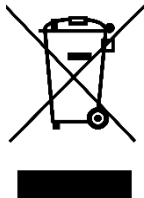This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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