X-Pro

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
Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

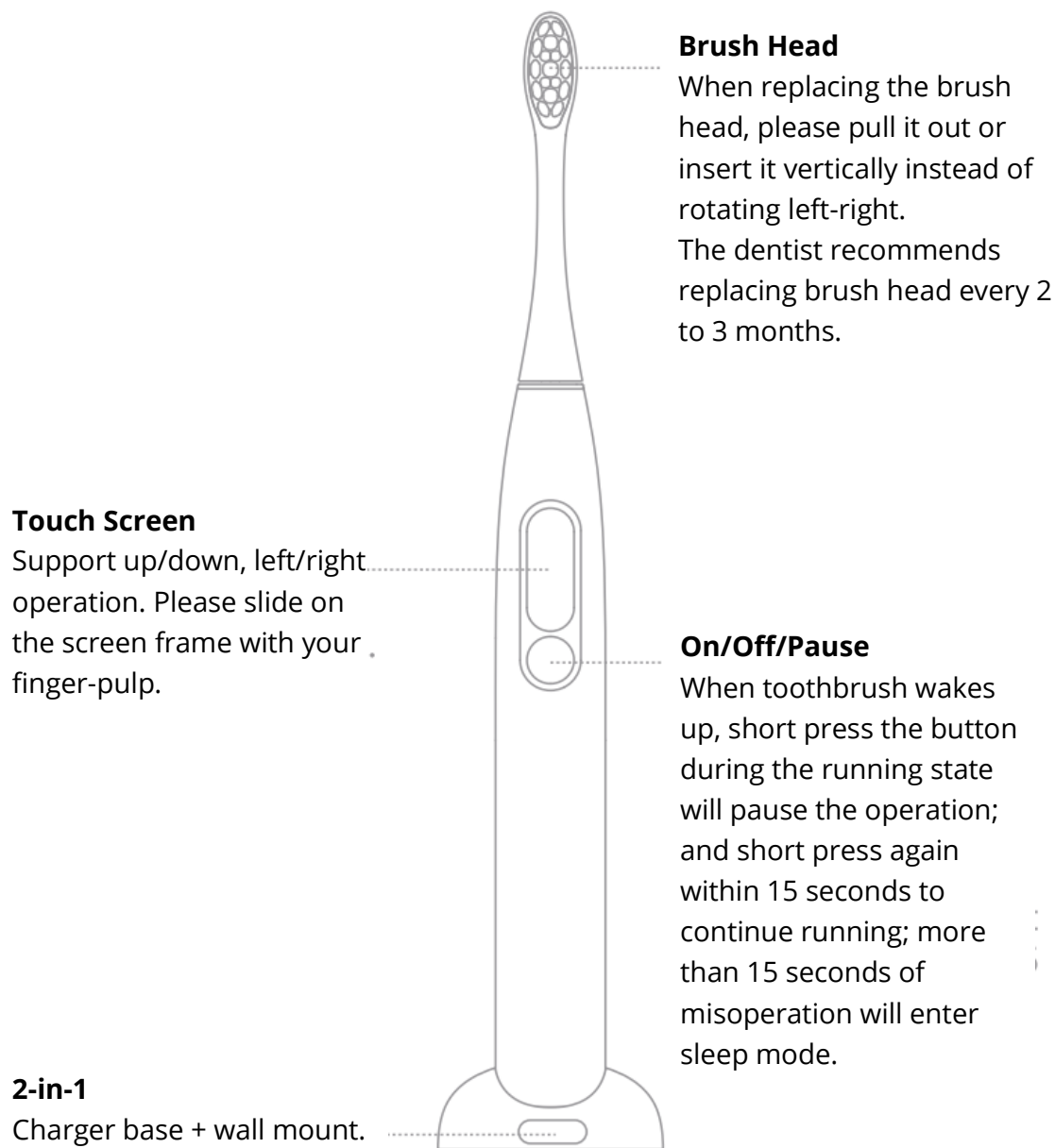First time use


Short press on ON/OFF button to activate toothbrush and Bluetooth, after being activated, the device will be constantly in stand-by status.

Brushing Plan

Toothbrushing will start default mode with "Standard Cleaning Plan" with 2 mins if no new customized brushing plan setting in APP. There will be a reminder with vibration to remind user to change the brushing area in every 30s. Toothbrush will stop operating automatically after 2 mins.

Charging

Please wipe the toothbrush to dry before charging. Please keep an eye on the battery percentage, only charge it when battery is low, do NOT overcharge or repeat charging frequently.

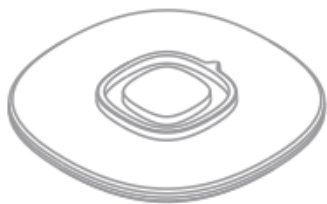Wall Magnet Sticker
(Paste on the wall)

+


=

Bottom Base
Put in a flat position when
using it for charging.
Keep upright and insert to
wall magnet sticker when
using it for a holder

Hanging Base
2-in-1 feature

Notes

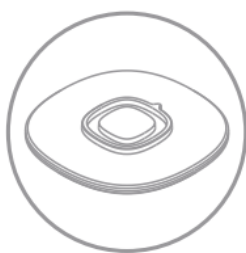Wall mount

When installing the wall hanging, be sure to place the base surface. The small arrow on the board is facing up.

Important Safety Tips

- This product is recommended for people over 12 years old.
- This product can only use Oclean brush heads and will not bear all the consequences of non-certified third-party brush heads.
- This product is only a personal cleaning product and is not recommended for medical use or multi-user use.
- This product is for indoor use only. Do not place it in direct sunlight. It is recommended to be place in a dry and ventilated place.
- Before charging, please keep the body dry, and regularly clean the charging base with a dry cloth to prevent the deposition of stains.

Not recommended for use

- Individuals who have recently undergone oral surgery. It is recommended that consumers with oral disease, such as dental cavities, periodontal (gum) disease, oral cancer, oral infectious disease, consult with a dentist before using this toothbrush.
- Individuals with mental and/or physical disabilities should only be used under supervision.
- People who have mouth diseases, loose false teeth, or gums with haemorrhage.
- People with sensitive and susceptible oral cavity, they should consult a professional dentist before using this product.
- The toothbrush is not suitable for people with reduced mental and physical capabilities. Not suitable for infant and children under the age of 12. Children who use this product should be strictly supervised by an adult.

Instructions

- **First Time Use**

The toothbrush can be turned on by pressing the ON/OFF button for the first time on a brand-new toothbrush. The Bluetooth function will be activated and remain on. After installing the Oclean App, please complete the prompts to complete the connection and settings for the toothbrush.

- **Install App & Toothbrushing Settings**

You can use the app to customize your brushing programs – automating tooth cleaning, tooth whitening and gum massaging in one brushing session. There is no need to open the Oclean App each time you brush. Simply open the app every few days to sync the brushing data.

- **Direct Use without App Settings**

You can also use the toothbrush without connecting to the app. It has touchscreen display and an ON/OFF button. It offers touch control via the touch-sensitive screen.

- **Brushing Performance**

After brushing, the display will show the effective brushing time, the brushing scores, and the tips of the brushing area.

- **Check Brushing Records in App**

Only after the toothbrush is paired with the app and connected to the Bluetooth, can the app get access to the brushing data stored in the toothbrush. Due to limited storage space, the brushing data can usually be stored for no more than 7 days.


- **Mode Setting**

Enter the mode-setting section by sliding up and down on the display. Each of the modes, such as Cleaning, Whitening and Massage, has 32 levels of cleaning intensity. After the new setting is completed, it will automatically register in the memory.

Brushing intensity: Sliding up/down on the screen, the intensity will be changed by 1 level (+1 or -1); Side-sliding up/down on the screen side, the intensity will be changed by 5 levels (+5 or -5).

Slide up/down/right/left

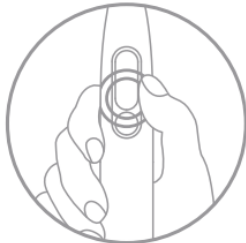Long press on screen or short press on power button to confirm the settings

Time Adjustment

Connect with the app and pair with the toothbrush to synchronize. After the message, the toothbrush will be updated automatically for local area time.

- No disassembly, impact, extrusion or throwing into fire.
- Please stop using the toothbrush if the battery has expanded.
- Please don't put the battery in high temperature environment.
- Do not use the battery if water entered the device.

Brushing Plan & App Settings

You can either choose or customize a brushing plan through the supporting Oclean App. You can also upgrade the firmware to continuously update its performance.

Low Battery Prompt

When battery percentage is lower than 20%, you'll see low battery prompt appears on your toothbrush display frequently, you will need to charge the toothbrush. When the battery percentage is lower than 1%, you will be prompted by motor vibration, and the screen will not be on.


How to Charge

Remove the charging base from the wall mount and connect it to the mobile phone charger. When the bottom of the toothbrush is placed gently above the charging base, it will automatically position itself. The display will show the battery percentage when the toothbrush is charging. A full charge usually takes about two hours.

Important Charging Tips

The toothbrush runs with a fast-charge mechanism. It is recommended that the toothbrush be placed on the charging base only when it is necessary. Please do not put the toothbrush on the charging base for a long time, as it may cause the built-in battery to decay prematurely so that the user experience is affected.

Recommended Tooth Brushing Method

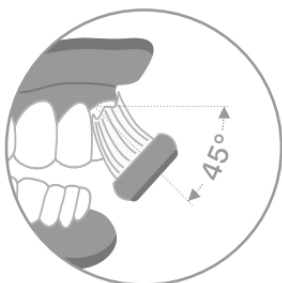Slowly move the bristles at 45° to the gum line to brush teeth properly. No need to move the bristles up and down manually.


Hold the toothbrush handle and slowly move it in circular motion. Do not press too hard.


Please rinse with clean water after brushing.

Cleaning and Maintenance

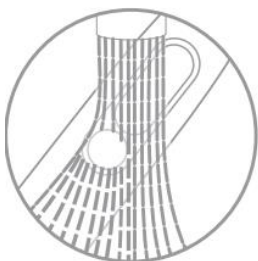
Cleaning the body and brush head

Rinse with water and air dry, do not pull the bristles.

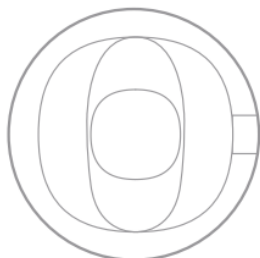Cleaning Charger

Clean the dirty parts of the charger with a soft cloth. Be sure to remove the battery plug from the socket prior to cleaning the charger and use it again only when it is fully dry.

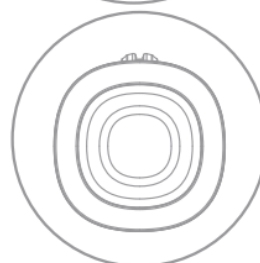Cleaning the Handle Base

Always keep the bottom clean and dry before charging. Clean the dirty parts of the handle with a soft cloth.

Specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Oclean X
Edition	Pro
Max Output	0dBm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V --- 1A
Battery Life*	20 days (Level 1), 40 days (Level 32)
Battery Capacity	800 mAh
Charging Time	2 hours
Bluetooth	4.2 above

* Calculated by 2 minutes/time, twice/day

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Smart Sonic Electric Toothbrush

Model / Type: X Pro

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/30/EU


Directive No. 2014/35/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

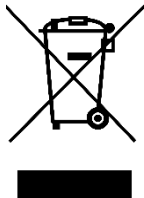This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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