W1

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.


✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important Safety Tips


- This product can only be used with Oclean special nozzle. Oclean will not bear the responsibility caused using any third-party accessories.
- To prevent water from spraying toward eyes, nose, ears, or throat to cause accidental injury, do NOT use for purpose other than oral cleaning.
- People who have just undergone oral surgery for treatment. It is recommended those who have mouth disease such as seriously haemorrhaged gums to consult a dentist before using it.
- It's normal to feel a bit discomfort or bleeding slightly when use it for the first time, situation will be much better as gum gets healthier. If there is a large amount of bleeding or bleeding after 3-4 weeks still, it's recommended to consult a dentist.
- To avoid injury caused by accidental inhalation by children or infants, this product is only recommended for people over 12 years old.

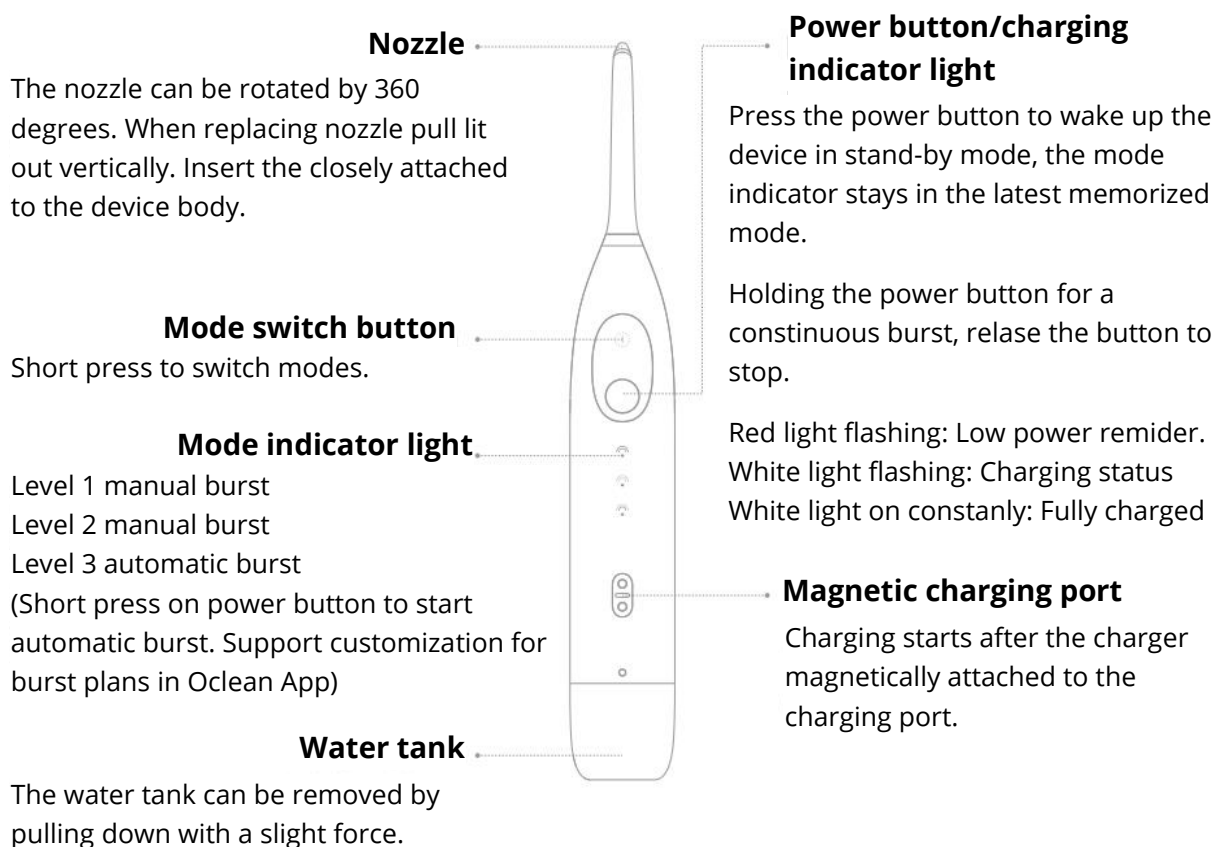Instructions


1. Refill the water tank

Pull down the water tank with a little force to detach it. Fill water or mouthwash into the water tank (The liquid temperature should not be over 40°C), then insert the water back into the device.

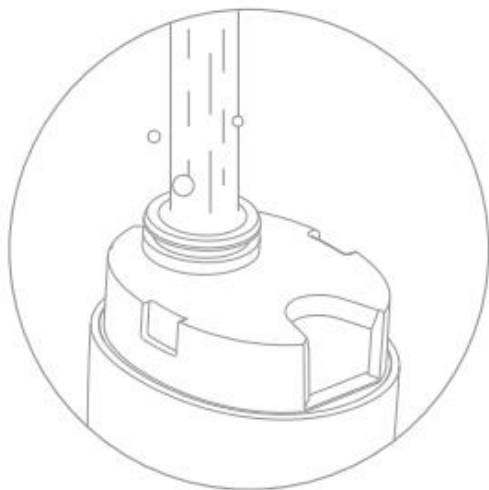
2. Installation of nozzle

1. Fix the silicone cap to the nozzle.
2. Insert the nozzle to the device body.

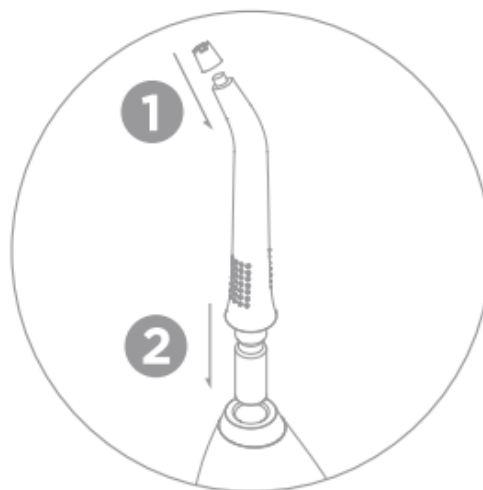
3. Mode switch

Short press on modes switch choose the best fit mode, the device will memorize the latest mode.

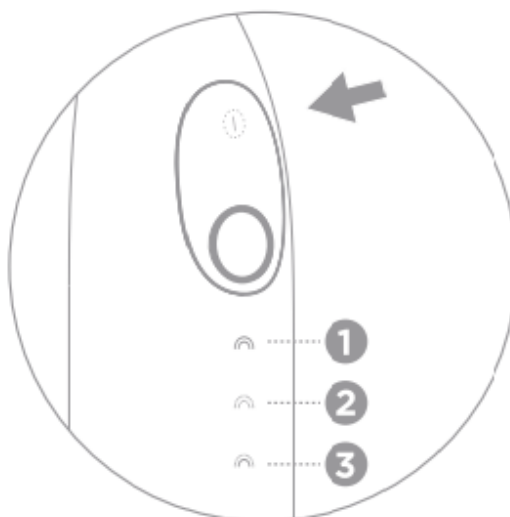
4. Bursting Instructions

Put the nozzle intermittently close to the slit between teeth. Press and hold the power button to start cleaning. Slowly move along the gum line to the next zone, repeat the cleaning for all slit.

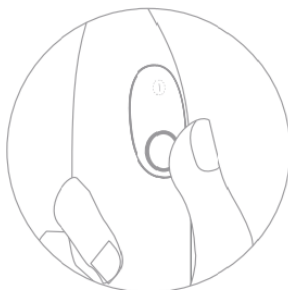Device Cleaning and Maintenance

1. Remove the nozzle, clean the nozzle with running water after each use.

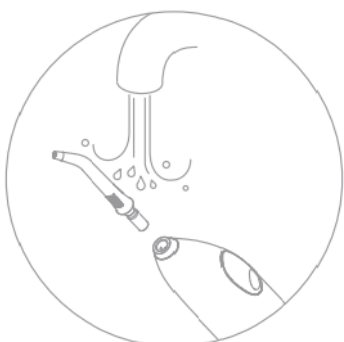2. Detach the water tank, get rid of the remaining liquid and clean the water tank with running water.

3. Wipe the surface of the device with a soft wet cloth to keep it clean.

Important Tips

- It's normal that the first few bursts may be only gas, it will become a mixture of air and water burst when water reaches the nozzle.
- Do NOT soak the device body in water, use dry cloth to clean it after use.
- If the device will be long time shelved, please empty the water tank to dry. It's recommended to use spray function to clear the remaining water. Clean and store the product in a cool and dry place. Do NOT expose it in direct sunlight.
- It's normal that the device heats up a little during cleaning process.
- It's normal there's a little plastic smell in the first few tanks of cleaning, the smell will disappear after several times of use.
- Please empty the water tank and wipe it dry before putting into the travel case.

Specifications

Product Name	Oclean W1 Smart Oral Irrigator
Model	Oclean W1
Water Tank Capacity	30 ml
Waterproof level	IPX5
Input Power	2.5 W
Charging	Oclean magnetic charging
Input	5V  0.5A
Battery Life*	30 days
Battery Capacity	1200 mAh
Charging time	3 hours
Standard	GB 4706.1-2005 GB 4706.59-2008

* Remark: Battery life is calculated based on default mode, working once/day, 1 tank/time

Other Tips

- It's recommended to change the nozzle every 4-6 months. If the nozzle is seriously blocked or cannot be installed on the device body normally and tightly, it also needs to be replaced immediately.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Smart Oral Irrigator

Model / Type: W1

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):


Directive No. 2014/30/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

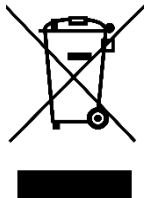This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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