FITTOP

Super hit Lite Version

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important information

Statement

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings in Super Hit massage gun user manual before using.

Prohibition

- In case of muscle injury please consult your doctor before using Super Hit message gun. Minors are advised to use Super Hit massage gun under the adult supervision.
- Do not use wet and dirty skin. It can be used over the clothes, just gently press, and move on the skin for about 60 seconds on each part.
- Do not use over the neck, bones, and body parts with less muscle coverage. It's only applicable on the body soft tissue without causing pain or discomfort.
- Do not use when the body has bruise, pain, and discomfort.
- Do not put any objects in the vent holes of the product.
- Do not immerse the product in water and please keep water out of the product vents.
- Do not put the product into the fire or close to the fire source.
- Do not drop this product intentionally.
- Do not use other chargers to charge Super Hit message gun. It must be charged with the charger provided by Fittop Company.
- Do not use this product while charging.
- Do not disassemble this product by any means.
- Do not close to the axis and back of the message head by fingers, hair or other parts of the body, as it may be clipped.

Do not touch this part when starting the Super Hit massage gun without message head installed


Do not touch this part when starting the Super Hit massage gun with message head installed.

Warning and reminders

- Keep away from liquids or heat sources and keep the vents dust free.
- After the charging is completed, pull out the charger.
- Please use this product according to the instructions.
- Carefully check whether the product can be used normally according to the instructions in this book before each use.

Medical note:

Super Hit massage gun cannot be used in any of the following cases without the prior consent of the doctor. Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, spondylitis, vertebral detachment, or spinal joint stiffness, recent joint replacement surgery or IUD, metal needle or steel plate, and any other serious problems related to your health. These contradictions do not mean that you cannot use this product. We recommend that you consult a doctor before use.

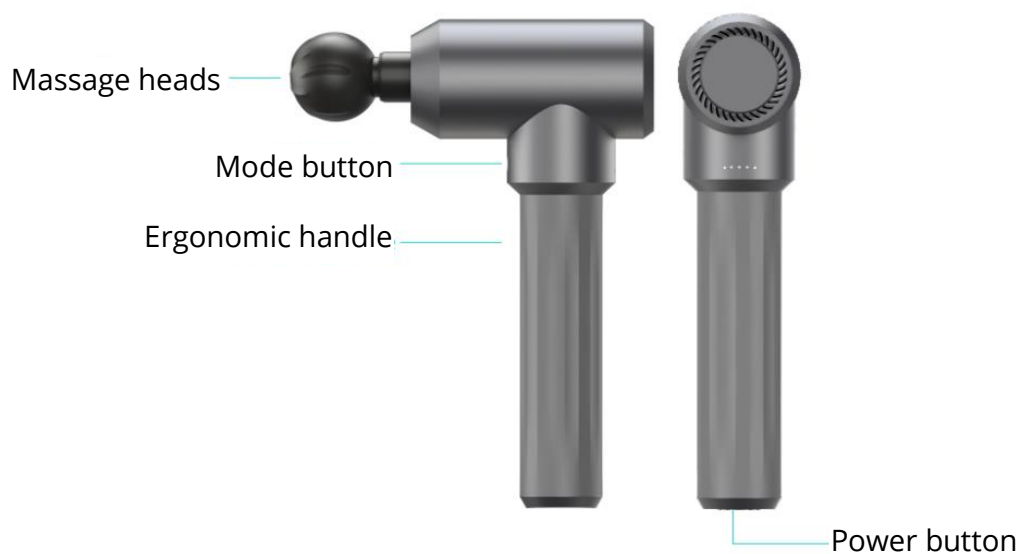- Mute technology use adopting dynamic force balancing for silence work.
- Skin-friendly silicone massage heads provides a comfortable massage experience.
- High quality rechargeable lithium-ion battery.
- Compact handle for comfortable gripping.
- Wireless, lightweight, and portable.
- Anti-sweat and moisture-proof humanized design.
- Battery and the body integrated design gives ergonomic handle with no risk of loose.

Super Hit massage gun 5 changeable massage heads equipped

Achilles tendon
massage head

Triger point
massage head

General
massage head

Joint massage
head

Warm-up
activation
massage head

Massage heads installation method:

While pulling/installing, slightly rotate the massage head to pull out/install. Do not pull/install directly. Choose massage head position for most comfortable massage.

Basic parameter

Product Name	Percussion Gun
Product Model	FSL961
Size	145*62*253.5mm
Product Weight	About 0.8kg
Rated Voltage	18V
Power adapter specifications	DC, 18 $\equiv$ 1A
Rated current	1000mA
Massage heads	5 different massage heads
Mode	Five speeds optional, manual speed control
RPM	1800-3200rpm

Charging

1. Charge the battery before initial use. (If the product is sufficiently charged, there's no need for charging)
2. When charging, connect the Fittop provided adapter plug to the general household power supply, and connect the other end to the battery charging port on the bottom of the handle.
3. Do not start Super Hit massage gun when charging.

Battery display:

- The battery indicator light is always on when is fully charged, and flashes when is not.
- There are three indicated battery power levels, low, medium, and high.
When battery discharging specific level indicator light on:
Fully charged: all three indicators,
Medium: medium and low-level indicator,
Low: only low-level indicator is on or flashes.

Instructions

1. Turn on the power button  ON position, the device is standby, Speed 3 indicator light breathes.
2. Mode: Short press the mode button  to start up and start with Speed 1. The Speed 1 indicator light is always on. Short press the mode button  to level up one speed each time (the current speed indicator light is always on). When the device reaches Speed 5, short press the mode button  again to enter standby status.
3. Manual speed control mode: long press the model button  for 3 seconds to enter manual speed control mode, and the speed will increase with the button pressing time (the same goes for the speed indicator lights) until the maximum speed. Short press the mode button  to enter standby status.
4. When the product is not in use, please make sure that the power button  is switched to the OFF position, otherwise the power will be consumed.

Maintenance, Cleaning, Storage, Transport

1. Clean the device with a slightly wet towel and dry with a soft cloth.
2. Place the product in the package provided while storing or travelling.
3. Stored in a dry, well-ventilated place with non-corrosive gas and no less than 250mm height from the ground. The number of stacking layers does not exceed 5 layers.
4. If not use for a long time, store in the package after cleaning and charge every six months.
5. This product is suitable for car, train, ship, and aircraft transport.
6. This product must be fixed during the transit to prevent severe vibration, collision, roll and fall and place it upward, lightly to somewhere rainproof.

General Massage Head

Unique wave-curve structure design for better fascia stimulation. Applicable for multiple muscle groups.

Warm-up activation message head

Hitting your muscles with hand-imitation structures. Quickly penetrate stimulation and activate warm-up status.

Joint massage head

Distinctive inward concave design, for better joints and coccyx protection during percussion.

Trigger points massage head

Characteristic spiral texture design for deeper percussion into pain points and resolving connected pain points one by one!

Achilles tendon message head

Special bionic V-shape design fit adequately for Achilles tendon, hind, neck and back.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Percussion Gun

Model / Type: FSL961

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):


Directive No. 2014/30/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

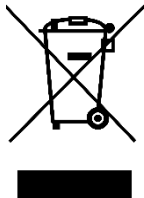This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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