FITTOP

Super hit Pro Version

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important information

Statement

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings in Super Hit massage gun user manual before using.

Prohibition

- In case of muscle injury please consult your doctor before using Super Hit massage gun. Minors are advised to use Super Hit massage gun under the adult supervision.
- Do not use on wet and dirty skin. It can be used over the clothes, just gently press, and move on the skin for about 60 seconds on each part.
- Do not use over the neck, bones, and body parts with less muscle coverage. It's only applicable on the body soft tissue without causing pain or discomfort.
- Do not use when the body has bruise, pain, and discomfort.
- Do not put any objects in the vent holes of the product.
- Do not immerse the product into the fire or close to the fire source.
- Do not drop this product intentionally.
- Do not use other charges to charge Super Hit massage gun. It must be charged with the charger provided by Fittop Company.
- Do not disassemble this product by any means.
- Do not close to the axis and back of the massage head by fingers, hair, or other parts of the body, as it may be clipped.

Do not touch this part when starting the Super Hit massage gun without massage head installed.

Do not touch this part when starting the Super Hit massage gun with massage head installed.

Warning and reminders

- Keep away from liquids or heat sources and keep the vents dust free.
- After charging is completed, pull out the charger.
- The Super Hit massage gun cannot be used during charging.
- Please use this product according to the instructions.
- Carefully check whether the product can be used normally according to the instructions in this book before each use.

Medical note:

Super Hit massage gun cannot be used in any of the following cases without the prior consent of the doctor. Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, spondylitis, vertebral detachment, or spinal joint stiffness, recent joint replacement surgery or IUD, metal needle or steel plate, and any other serious problems related to your health. These contraindications do not mean that you cannot use this product. We recommend that you consult a doctor before use.

- Mute technology use adopting dynamic force balancing for silent work.
- Skin-friendly silicone massage heads provides a comfortable massage experience.
- High quality rechargeable lithium-ion battery.

- Compact handle for comfortable gripping.
- Wireless, lightweight, and portable.
- Anti-sweat and moisture-proof humanized design.
- Battery and the body integrated design gives ergonomic handle with no risk of loose.

Super Hit massage gun 5 changeable massage heads equipped

Massage heads installation method:

While pulling/installing, slightly rotate the massage head to pull out/install. Do not pull/install directly. Choose massage head position for most comfortable massage.

Basic parameter

Product Name	Intelligent Chip Percussion Gun
Product Model	FSP951
Size	153.5*63*253.5mm
Product Weight	About 1kg
Rated Voltage	29V
Power adapter specifications	DC, 29V $\approx$ 1A
Rated current	About 1000mA
Timing time	Automatic shutdown after 15-minute working
Massage heads	5 different massage heads
Working mode	Intelligent mode, manual mode

Charging

- Charging the battery before initial use. (If the product is sufficiently charged, there's no need for charging).
- When charging, connect the Fittop provided adapter plug to the general household power supply, and connect the other end to the battery charging port on the bottom of the handle.
- Charging process battery indicator: battery power <25% (red light flashing), battery power <50% (yellow light flashing), battery power <75% (cyan light flashing), battery power >75% (green light flashing), battery power 100% (while charging is finished and battery saturated, green light is always on).
- Do not start Super Hit massage gun when charging.

Instructions

1. Short press the power button  to turn on the power, then short press the mode button  to start percussion.
2. In any intelligent mode, long press (3 seconds) the mode button  to switch to manual mode. Then you can manually control speed by holding the mode button.
3. The system will enter into intelligent protection mode after 15 minutes of each use, and Super Hit massage gun will automatically stop. If you want to use it again, you only need to press mode button  again.

Voice broadcasting function: When the device is on, the voice broadcasting opens by default. When the voice is on, long press power button  for 5 seconds to turn off the voice broadcasting function, and “beep, beep” indicates voice broadcasting closes. When the voice is off, long press the power button  for 5 seconds to restart the function. And “beep” indicates the voice broadcasting opens. When the battery is no more than 25% during the voice mode, it would broadcast “low battery, please charge in time” ever once in a while. (Note: when the voice is turned off, it can still indicate low power when the power is low)

Mode switch:

1. Intelligent mode: short press the mode button  to start the device, and the product defaults to the intelligent mode. The proximity sensor, massage head recognition function works in the intelligent mode. After starting up, the device automatically identifies massage head and enters into the intelligent slow speed. When the massage head is close to the object, it turns into the rated reservation speed of the current massage head. When the massage head is always from the body, it enters into the slow speed, along with the corresponding massage head and icon flashing in display screen. If there is no massage head detected, Super Hit will shut down in 5 seconds.

2. Manual mode: long press the mode button ○ in any modes after turning in the device, the product will automatically switch to the manual mode. The proximity sensor, and massage head recognition function would be off by default.
In this mode you can control speed manually. Press and hold mode button “○” Super Hit will accelerate infinitely until you release the button, or until reach the highest possible speed. Press the mode button short to shut down.

Battery status indicator:

Battery at 75% ~ 100% charged level (shows 4 lines, and power indicator shows green)

Battery at 50% ~ 75% charged level (shows 3 lines, and power indicator shows cyan)

Battery at 25% ~ 50% charged level (shows 2 lines, and power indicator shows yellow)

Battery at 5% ~ 25% charged level (shows 1 line, and power indicator shows red)

Battery at <5% charged level (shows empty battery symbol, and the power indicator shows flashing red).

Maintenance, Cleaning, Storage, Transport

- Clean the device with a slightly wet towel and dry with a soft cloth.
- Place the product in the package provided while storing or travelling.
- Stored in a dry, well-ventilated place with non-corrosive gas and no less than 250mm height from the ground. The number of stacking layers does not exceed 5 layers.
- If not in use for a long time, store in the package after cleaning and charge every six months.
- This product is suitable for car, train, ship, and aircraft transport.
- This product must be fixed during the transit to prevent severe vibration, collision, roll and fall and place it upward, lightly to somewhere rainproof.

4 Black Techs

Intelligent Control Black Technology 1

Intelligent switching massage head percussion modes

Other than traditional indistinguishable percussion, Super Hit is more intelligent and attentive for muscle needs.

When inserting different massage heads, the system will automatically identify and configure the corresponding percussion mode, intensity and amplitude. Precise configuration of percussion for different muscle area is naturally out of necessity.

Intelligent Control Black Technology 2

Intelligent voice broadcast

Instead of the traditional “mute” operation, Super Hit provides convenient and instant voice guidance.

When inserting different massage heads, the system will automatically switch to the corresponding percussion mode and voice broadcasting.

For example, when inserting the general massage head, the system will automatically switch to multi zone relaxation modes and voice broadcast: “General mode suitable for all muscle fascia parts.” You can turn on and off the intelligent voice broadcast any time you want. Press and hold power button for three seconds when device is turned on. You will hear “beep” sound for turning on and “beep, beep” sound for turning off.

Intelligent Control Black Technology 3

Visual Intelligent Control Interface

Different from the traditional one-way operation, Super Hit massage gun comes with unique user-friendly interface with clear visible icons:

- Percussion modes (Intelligent, Manual)
- Massage head icons (Multi part, Joint massage, Run & Jump, Warm up, Trigger point)
- Massage areas (Different body parts flashing)
- Time countdown (15 minutes protection mode countdown digital clock)
- Battery reminder (Battery shape with 4 lines in)
- Speed indicator circle (Green and red lines in circle)

Fittop engineers developed user-friendly interface which helps a lot and allows beginners to easy understand how Super Hit massage gun PRO works!

Intelligent Control Black Technology 4

Intelligent idling system for power saving

Distinguish from the traditional non-power-saving percussion, Super Hit possesses intelligent idling technology for power saving.

When the massage gun hovering over or away from the body, it will automatically switch to idling to or non-operating mode. Intelligent power saving not only for battery life extending, but also for longer production lifespan.

General massage head

Maximum RPM: 3500rpm

Unique wave-curve structure design for better fascia stimulation. Applicable for multiple muscle groups.

Voice broadcast entry: General mode suitable for all muscle fascia parts.

Warm-up activation massage head

Maximum RPM: 3150rpm

Hitting your muscles with hand-imitation structures. Quickly penetrate stimulation and activate warm-up status.

Voice broadcast entry: Warm-up mode for quick entering into motion state.

Trigger point massage Head

Maximum RPM: 2975rpm

Characteristic spiral texture design for deeper percussion into pain points and resolving connected pain points one by one!

Voice broadcast entry: Dedicated Trigger points release mode for whole body trigger points relieving.

Achilles tendon massage head

Maximum RPM: 2625rpm

Special bionic v-shape design fit adequately for Achilles tendon, hind neck and back.

Widely very well received by exercisers and practitioners.

Voice broadcast entry: Achilles tendon and ankle joint relaxation mode for relaxation before and after intense jumping and running exercise.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Intelligent Chip Percussion Gun

Model / Type: FSP951

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):


Directive No. 2014/30/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

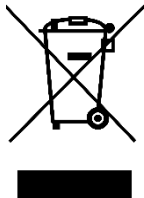This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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