# SIGURO


**Vacuum  
blackhead  
remover**

SGR-SK-V640W

<b>ENGLISH</b>	<b>5</b>
<b>DEUTSCH</b>	<b>8</b>
<b>ČESKY</b>	<b>12</b>
<b>SLOVENSKY</b>	<b>15</b>
<b>MAGYAR</b>	<b>18</b>
<b>FRANÇAIS</b>	<b>21</b>
<b>ITALIANO</b>	<b>24</b>
<b>ROMÂNĂ</b>	<b>27</b>
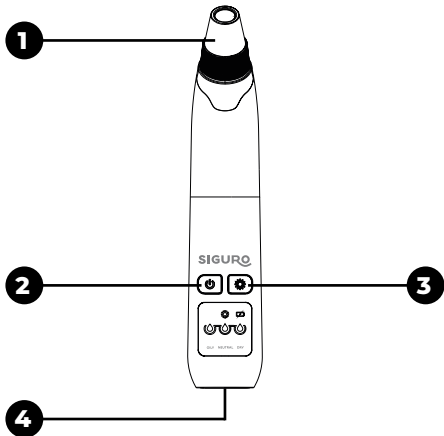
## **SAFETY INFORMATION**

---

1. This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, even if they have been given supervision or instruction concerning use of the appliance in a safe manner and understand the potential hazards. Children should be supervised to ensure that they do not play with the appliance.
2. Never immerse the appliance in water!
3. The appliance is intended for household and similar use only. It is not intended for use in a beauty salon or other commercial use. Use of this device is not intended as a substitute for medical care.
4. Do not use the appliance if you have hypersensitive, damaged or skin with eczema (e.g. fresh scars, various wounds, burns, redness, swollen cheeks or skin diseases). To prevent the spread of micro-inflammation, do not use it even if you have acne on your skin.
5. If you feel pain or discomfort during cleaning, or if you experience any other health problems, switch off the appliance immediately and consult your doctor or dermatologist about its further use.
6. Do not use the appliance on children, legally incompetent persons or people who are unable to express their feelings.
7. Do not use the appliance with exfoliating or otherwise aggressive gels (e.g. with alcohol)!
8. Do not use the appliance close to sensitive body parts, not intended for this device.
9. Do not use the appliance in one time longer than specified in the operating instructions! Prolonged use may cause skin injuries.

- 1 **Vacuuming head**
- 2 **Power / Intensity button**
- 3 **Heating function button**
- 4 **Heating pad**

**Included in package and not shown:** interchangeable heads and charging cable.

Dimensions (H x W x D)	172 × 35 × 35 mm
Weight	120 g
Power	5 V / 1 A
Power consumption	5 W
Battery	Li-ion, 500 mAh / 3.7 V

## USE

### Charging

- You can charge the device using the included USB cable.
- Fully charge the device before using it for the first time.
- The charging status and remaining battery capacity are shown on the display.

### Method of use

The device is used to remove blackheads, comedones and

clean the pores by suction. The device is equipped with three suction intensities shown by indicators **OILY** (lowest), **NEUTRAL** (medoim) and **DRY** (highest). To facilitate cleansing, the device is also equipped with a heating pad, with which you can loosen the facial skin and expand the pores.

### **Controls**

- Attach the appropriate head to the device as needed.
- Press the **2** button to start suction. Press repeatedly to switch the intensity in the order OILY → NEUTRAL → DRY → Off. The current mode will be indicated by a symbol on the display.
- Press the **3** button to turn on the heating pad. **Attention:** you cannot start the heating and suction functions in the same time.

## **CLEANING AND MAINTENANCE**

---

- Always make sure that the appliance is switched off before cleaning or storing it.
- Rinse the device under running water. Do not immerse in water or any other liquid. Do not use any cleaning agents, lubricants or polishes.
- Do not dry the device in direct sunlight.

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## ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on [www.alzashop.com/doc](http://www.alzashop.com/doc). This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic

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© 2022 Alza.cz, a. s. All Rights Reserved.