Infračervená lampa

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Obsah:

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User Manual

Table of Contents

1. Product Introduction	1
2. The effects of the Infrared Lamp	1
3. Specification	1
4. Components	2
5. Functions and Features	2
6. Safety notices	3
7. How to use	5
8. Product Cleaning and Storage	7
9. Problem and Solution	7

1.Product Introduction

Infrared Lamp gives effective pain relief for muscles and joints. Its 300 Watt infrared warmth is very comfortable and penetrates deep into the skin, stimulates the blood circulation.


2.The effects of the Infra red Lamp

The Infrared Lamp, can penetrate deeply into the skin and warms the tissues. The warmth modulates the neuronal signals to the brain and thus reduces pain. Warmth also stimulates the blood circulation and dilates the blood vessels, which accelerates the transport of substances necessary to rebuild and nourish the body's tissues. Furthermore, the infrared lamp speeds up the metabolic process and the removal of waste substances from the body. Because it makes tissue more pliable, warmth also reduces stiffness and makes joints more flexible. Because of these effects, treatment with the Infrared Lamp can temporarily relieve pain due to muscle and joint problems.

3. Specification

Power:	300W
Voltage:	220V
Frequency:	50Hz
Lifetime of lamps:	500 hours
Product dimensions:	22.5 x 22.5 x 33.5 (L x W x H)cm
Product weight:	1.65kg
Box dimensions:	27.5 x 27.5 x 38.5 cm

3. Components

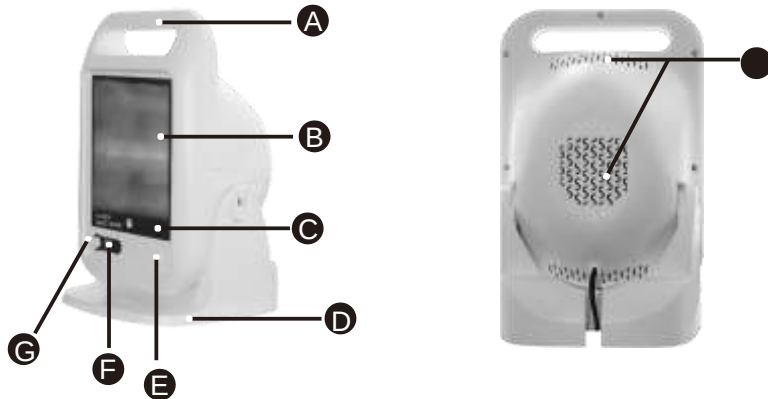A.Handgrip	E.On/Off
B. Lamp	F.Timer
C. Ceramic Glass	G.Time Selecting Button
D.Base	H. Cooling Hole


5.Functions and Features

1.0-35 degree adjustable angle.

2.An automatic protection against overheating.

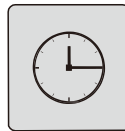3.15 minutes automatic cut off.

4.Ceramic glass, 100% UV proof.

5.Automatic protection against falling down.

◆When the appliance was falling down, angle $\geq 90^\circ$, the lamp stop working, while the cooling fan still on work, cool down the temperature of the appliance.

7.Safety notices

(1) Danger

- ◆Water and electricity are a dangerous combination!
- ◆Do not use this appliance in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).
- ◆Do not let water run into the appliance or spill water on to the appliance.

(2) Warning

- ◆Check if the voltage indicated on the appliance corresponds to the local mains voltage before you connect the appliance.
- ◆If the mains cord is damaged, you must have it replaced by a

service centre authorised by the manufacturer or similarly qualified persons in order to avoid a hazard.

- ◆Do not leave the appliance unattended when it is switched on.
- ◆This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- ◆This appliance is not intended for use by children. Children should be supervised to ensure that they do not play with the appliance.
- ◆This appliance is intended for the treatment of muscle and joint pain and muscle and joint stiffness, not for the treatment of abdominal, facial and other sorts of pain.
- ◆To prevent overheating of the skin, avoid very short treatment distances. The warmth should be comfortable throughout the treatment.
- ◆To prevent eye damage, do not look into the light of the lamp when it is on.
- ◆Do not fall asleep during treatment.
- ◆Prevent the treatment area from cooling down too fast immediately after treatment.

(3) Contraindications for the Infrared Lamp treatment:

- ◆ Do not use the appliance without consulting your doctor first if you suffer from severe diseases such as heart disease, acute inflammatory diseases, thrombosis, blood coagulation disorders, adrenal suppression, systemic lupus erythematosus or malignant diseases.
- ◆ Do not use the appliance if you are oversensitive to infrared lamp or if you have oedema. If in doubt, consult your doctor.
- ◆ Do not use the appliance if you are insensitive to infrared radiation and heat due to a disease (e.g. advanced diabetes) or the use of painkillers. If in doubt, consult your doctor.
- ◆ Do not use the appliance if you suffer from conditions associated with nerve demyelination such as carpal tunnel syndrome, multiple sclerosis and ulnar neuropathy. If in doubt, consult your doctor.
- ◆ Do not use the appliance if you have an implanted device, particularly not if it is located in the vicinity of the treatment area. If in doubt, consult your doctor.
- ◆ If you are pregnant, do not use the appliance to treat areas close to your abdomen or your lower back area, to avoid overheating of these areas. If in doubt, consult your doctor.
- ◆ Do not use the appliance if you have a swelling or an inflammation, as heat could aggravate the complaints. If in doubt, consult your doctor.
- ◆ Do not use the appliance on damaged or irritated skin or on wounds. If in doubt, consult your doctor.
- ◆ Do not use the appliance to treat areas on which you have recently applied topical counter irritants (creams or ointments containing e.g. menthol), as this can cause a mild inflammatory reaction.

(4) Cautions

- ◆ If the wall socket used to power the appliance has poor connections, the plug of the appliance becomes hot. Make sure you plug the appliance into a properly connected wall socket.
- ◆ Do not subject the appliance to heavy shocks.
- ◆ Always unplug the appliance after use and in case of a power failure.
- ◆ Place the appliance on a stable, level surface. Make sure there is at least 25cm free space around it to prevent overheating.

- ◆ Do not cover the appliance with, for example, a cloth or a piece of clothing when it is switched on.
- ◆ Let the appliance cool down for approx. 15 minutes before you clean or store it.
- ◆ If you have just been swimming or have just taken a shower, make sure you dry your skin properly before you use this appliance.

7. How to use

(1) Using the Infrared Lamp

- A: Place the appliance on a stable, level surface and make sure there is at least 25cm free space around it to prevent overheating.
- B: Put the mains plug in the socket. Press the on/off button to switch on the appliance (Fig. 1). The red power-on light goes on.
- C: Press the time setting button to select the treatment time (15-0 minutes). The appliance starts to work after 5 seconds. (Fig. 2)

(Fig. 1)

(Fig. 2)

- ◆ Make sure the distance between your body and the appliance is large enough to prevent Accident.
- ◆ Choose a distance at which the heat is comfortable and increase the distance if the heat becomes too intense.
- D: To adjust the direction of the light beam, tilt the lamp housing forwards or backwards.
- E: Sit or lie down at the right distance from the appliance.
- ◆ At a distance of 30-40 cm to the body, the Infrared Lamp treats an area of about 40 cm by 30 cm (w x h). At a longer distance, the area treated becomes larger, therefore the infrared lamp heat less intense. At a shorter distance, the area treated is smaller and the infrared lamp heat more intense.

F: The appliance can be used in many different ways. Below you find some examples.

- ◆To treat a sore shoulder, a sore neck or part of your upper back area, place the appliance on a table and tilt the lamp housing backwards. Sit down on a chair or stool with your back towards the lamp.
- ◆To treat a lower-back ache, place the appliance on a table beside a chair. Sit down on a chair or stool with your back towards the lamp.
- ◆To treat a sore elbow, place the appliance on a table and sit down beside it.
- ◆To treat a sore knee, place the appliance on a bedside table and lie down on the bed with one leg pulled up towards the lamp.

G: After use, press the on/off button to switch off the appliance.

Shoulder

Neck

Back

Arm

Leg

Stomach

(2) Duration of treatment

The duration of the treatment depends on the type of treatment and the person to be treated. Treatment of a single muscle or joint should last approximately 15 minutes. If necessary, the treatment can be repeated several times a day. A longer treatment time than 15 minutes does not lead to better results. You can repeat the course of several treatments a day on several consecutive days to get a result. For chronic problems, the appliance can be used for longer periods of time.

8. Product Cleaning and Storage

Cleaning:

- A: Unplug the appliance and let it cool down for about 15 minutes before you clean it.
- B: Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.
- C: Clean the lamp housing and the base with a moist cloth. Do not let any water run into the appliance.
- D: You can clean the surface of the appliance with a soft cloth sprinkled with some methylated spirits.

Storage:

- A: Do not let the mains cord come into contact with the surface of the appliance when it is hot.
- B: Unplug the appliance and let it cool down for about 15 minutes before you store it. Always lift or carry the appliance by its handgrip. Store the Infrared Lamp in a dry place.

9. Problem and Solution

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, contact the selling center of the Infrared Lamp. Don't try to disassemble the appliance by yourself.

Problem	Possible Cause	Solution
The appliance does not go on.	The plug is not inserted into the socket properly	Put the plug in the socket properly
	The protective tube inside the appliance was broken	Change the protective tube
No light	The temperature controller, Cooling fan or the lamp was broken	Change the temperature controller, cooling fan or the lamp
When the appliance was falling down, angle $\geq 90^\circ$ C, the lamp still working	The inversion switch was broken	Change the inversion switch

QualityWarrantyCardofAK-2012-RInfraredLamp

- a. We hereby guarantee that for customers purchased product, if it is used normally as per this user manual, in case of any failure, all the parts of the device(except the lamp) will be free of charge within one year from the date of purchase.
- b. When it's above one year less than two years, if the part replacement is needed, we only charge for the part cost and the repair will be done free of charge.
- c. Please mail the failure goods to us or to the local dealer, and we will handle immediately.

Model: AK-2012-R

One year of warranty period from the date of purchase (main unit) Purchase date:	
User's address:	
Name:	Telephone:
Dealer:	
Address:	
Telephone:	
Name:	Stamp:

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