HAVIT

MG1505

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

The deep tissue message gun is a portable rechargeable vibration device with large torque and long adjustment. It can effectively reduce the pain caused by the accumulation of lactic acid after intense exercise, and relieve the discomfort of muscles in the sedentary and long-lying population.

Packing List

- Massage Gun x 1
- Massage head x 6
- Carrying bag x 1
- AC Adapter x 1
- User's Manual x 1
- After-sales card x 1

Production Structure

Charging

- To charge, connect the DC end of the supplied 17V adapter to the battery's charging port and plug AC side cable into the wall outlet.
- LED bank on battery lights pulsate different colour to show battery level and indicate active charging. Battery level 10-20 is red, 30-50 is yellow, 60-100 is green.
- Full charge is indicated when the LED lights show solid green.
- It is not recommended to fully deplete the battery to the red LED level.
- 2500mAh Long Working Time Battery-Muscle massage gun, equipped with a rechargeable Lithium battery, lasts 6-11 hours massage after a full charge (depending on different massage intensity). Auto shut-off after 15 minutes of continuous use to prevent overheating. If you do not use it for 3 minutes, it will enter to sleep mode automatically. You need to touch the power button to restart it.

Operating

1. Power device on by moving the power switch on the bottom of the battery to the ON position. Power device off by moving the power switch to the OFF position.
2. Press "+" and "-" button to adjust machine speed. Total have 20 levels.
3. If the use process force is too large, when the machine is stopped and the speed is displayed as E2, the machine does not run. This is to press the middle pause key to reset and restart.
4. When the battery power is too low, the speedometer will display "E3" prompt, at this time, you need to charge the machine for a long time.
5. When the battery is full, the light line at the bottom of the machine changes from red to green.

Massage head instructions:

	Ball head For large, muscle groups, such as arms, legs, calves, back, hip.
	Flat head For all muscle groups, relaxation and body shaping.
	Rubber head For gentle massage.
	U-shape head For massaging neck, spine and achilles tendons.
	Bullet head For meridians, palms, plantars, joints, and trigger points.
	Mushroom head For limbs and shoulders.

Technical Specs

Work voltage:	13.2-16.8V
Rated power:	60W
Battery capacity:	2500mAh
Product weight:	750g
Product size:	16.2*24.1*6.8cm

Safety Instructions

DANGER: To reduce risks of electric shock, fire, and personal injury, this product must be used in accordance with the following instructions.

- For adult use only.
- Use only on a dry clean surface on the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use only on the body's soft tissue and desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Use only the applicators that provide the best desired effect.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes on the Massage Gun.
- Do not immerse in water or allow water to enter the ventilation holes of the Massage Gun.
- Do not drop or misuse the Massage Gun.
- Only recharge with supplied 17V DC Massage Gun charger.
- Carefully examine device and battery before each use.
- Do not tamper or alter Massage Gun in anyway. Never leave the Massage Gun operating or charging unattended.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Massage Gun

Model / Type: MG1505

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/30/EU


Directive No. 2014/35/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

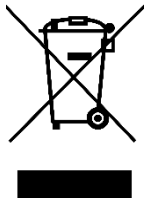This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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