medivon®

User manual

COLLAR MASSAGE

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1. Setup and operation

1. Connect the adapter cable to the cable on the massager.
2. Plug the home adapter into an outlet.
3. Turn on the massager and select your desired settings through below controller.

Controller support:

1. Power Button

Press to turn on the massager and activate shiatsu clockwise rotation at high level, press again to turn off the massager.

2. Direction Button


When the massager is turned on, press to shift Shiatsu clockwise rotation or Shiatsu anti-clockwise rotation

3. Intensity Button

When the massager is turned on, press to set the shiatsu intensity at high or low level


4. Heat Button

When the massager is turned on. Press to turn on or off the heat.

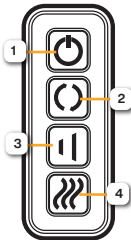Note: Press the Power button to turn on the massager and activate the default setting massage mode i.e. the massage heads will rotate clockwise at high level of intensity and with heating, until it auto shuts off or press the corresponding button to change it. **This device will switch off automatically after 15 minutes timer runs out.**


2. How to use your massager

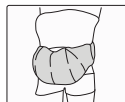
Abdomen


Neck & Shoulder

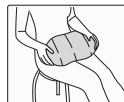
Thigh

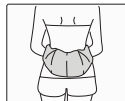Waist

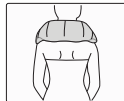Neck

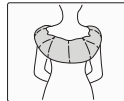Back

- Remove the adapter from the outlet when not in use.
- Do not setup or use the massager in a bathroom or similar wet / damp areas.

• For your own safety we equip the massager with a protection against over-heating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk and after a runtime of 15 Minutes. The massager will switch off automatically and will only allow switching on again after an appropriate downtime. Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 Minutes.

3. Safety instruction

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the IEC safety regulations.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.

- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel. In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.

• If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.

• The packaging material supplied should not be used as a toy.

4. Care and Cleaning Instructions

Always unplug the unit before cleaning it. Never touch the power plug with wet hands.

Before cleaning, allow unit to cool.

Wipe the unit with a soft, dry cloth. Do not use cloths containing any kind of chemical or alcohol and solvent liquids.

Never submerge any part of the unit in liquid.

If the unit is particularly dirty, use the following instructions:

1. Dampen a cloth in water or a 3%–5 % mild detergent solution.
2. Wipe the dirty areas with the wet cloth.
3. Clean the cloth in water and wring it out, then wipe away any remaining detergent off the surface of the unit.
4. Wipe the areas with a soft, dry cloth to dry.
5. Wait for the unit to become fully dry before use.

Place massager in a safe, dry, and cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the surface.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health damage or environmental pollution during further disposal or recycling of the equipment.

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