NURSENSE IPL Epilátor

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Laser Hair Removal

Instruction Manual

Please read and keep this manual properly before use. The pictures in this manual are for illustration only and the actual product shall prevail

CATALOG

<u>Safety precautions</u>	28
<u>Product structure</u>	33
<u>Accessories & packing list</u>	35
<u>Product parameter</u>	35
<u>Prepare before use</u>	36
<u>Instructions</u>	38
<u>skin care</u>	42
<u>cleansing and maintenance</u>	43
<u>Hair removal essentials</u>	44
<u>Possible effects</u>	46
<u>Problems</u>	47

Product features

Our company has militarized and homeized the medical beauty freezing point painless hair removal technology. The innovative combination of semiconductor refrigeration technology, integrated refrigeration technology, get rid of the burning sensation in the traditional hair removal process. The unique cooling technology can quickly soothe the skin, reduce skin irritation and bring a comfortable and ice-free painless hair removal experience.

- Refrigeration technology makes the whole hair removal feed ice comfortable and get rid of the burning sensation

The chiller refrigeration technology can achieve a safe and comfortable ice-feeling painless hair removal experience. In the process of hair removal, the cooling sheet can reach 10OC when it comes into contact with the skin, which can effectively cool the epidermis and reduce the burning injury of the skin, making the whole process comfortable and painless.

- Intense pulsed light technology is perfectly applied in the field of household hair removal instrument

Our company uses intense pulsed light (IPL) technology, which is the most effective way to restrain continuous hair regeneration on the market. Through close cooperation with skin experts over the years, we have improved this optical technology, which was originally and effectively at home. The light of the hair removal instrument is soft, the irradiation is convenient and effective, and the light intensity makes people feel very comfortable. After many times of use, the excess hair will eventually cease to exits, allowing you to enjoy the pleasure of getting rid of the hair and glow with charming brilliance every day.

- Long-life one-piece lamp holder

This product has a long-life integrated lamp holder that does not need to be replaced, and is suitable for multiple hair removal areas, such as face, lip hair, armpit, bikini line, arm, leg etc., saving you expensive lamp holder replacement costs and no need to purchase accessories frequently.

SAFETY PRECAUTIONS

In order to treat your skin better, please use it closely

- Do not apply any gel or care products (aloe gel, etc.) before use
- Please be sure to use the product in accordance with the regulations contained in the Instruction Manual
- Please read the Instruction Manual carefully and keep it for inquiry
- If there is any question found in use, please contact the online customer service in the desired language

DANGER

Misoperation may lead to personal injury

NOTICE

Misoperation may lead to danger

BE CAREFUL

Misoperation may lead to personal safety or property damage.

DANGEROUS

- ⊘ Warning: do not use this product near bathtubs, shower basins, washbasins or other utensils that hold water. Please keep this product and power adapter dry.
- ⊘ If the product is damaged, do not touch any internal parts to avoid electric shock.
- ⊘ Please do not insert any objects into this product, and please do not place this product at will.
- ⊘ Do not use this product around the eyes, wearing protective glasses in a common operation to protect the user, but if you operate around the eyes, the protective glasses may not be able to completely block the energy, and the light will be absorbed by the iris melanin in your eyes at close range.
- ⊘ We recommend that you wear protective glasses every time you use it. Although the light does not directly affect your eyes, high-frequency flash light may affect your eyesight.
- ⊘ When you need a back hair removal or ice compress operation, you may need an assistant or partner to avoid the risk of using it alone
- ⊘ If your hair is longer than 1mm, please shave first. We recommend that you use a hair knife. Excessive hair absorbing heat on the skin surface will cause skin redness and swelling, and hair burning will produce a strange smell, which will affect your experience of using product.

- ⊗ If the skin of the irradiated area is not dry and clean, it will affect the effect of use, and it may also cause the product not to flash or divination.
- ⊗ Please pay attention to prevent children under the age of 15 from coming into contact with this product, and please stay away from children under the age of 15 when operating this product.
- ⊗ Please use different light intensity stalls for different body parts as needed, because the skin color of different body parts is not the same.
- ⊗ If the power cord is damaged, it must be replaced by a professional from the manufacturer, product maintenance department, or similar department in order to avoid danger.

! NOTICE

- ! Do not use antiperspirant immediately after exposure to the armpit.
- ! Please confirm and clean the dirt of the irradiation mouth regularly to avoid scalding and failure.
- ! Do not put the irradiation mouth near items that are easy to degauss, so as not to cause failure and damage.
- ! Do not use it when the vent is blocked, so as not to cause any smoke and fire
- ! Do not point the irradiation port toward the black object, so as to prevent the machine from working properly.
- ! Do not irradiate clothes, hair and other items, so as to avoid scorching, discoloration and damage.
- ! Do not use it on pets to avoid injury.
- ! Do not use this product around the eyes to avoid eye injury.
- ! Please wipe the hair and skin care products attached to the depilator after use to avoid failure and skin damage.
- ! Please do not use it around children, the hair removal instrument will give off a strong light and easily hurt your eyes.
- ! Please set the irradiation intensity according to your skin condition, and be sure to use it from gear 1 for the first time
- ! Make the skin shout slightly got and red, which is a normal phenomenon, which may turn on the ice compress function and relive skin discomfort.
- ! The product cannot be left unattended after it is turned on. Please turn it off in time after use.
- ! Clean the skin and make sure there is no residual hair, completely dry, and no oily substances or stains on the surface of the skin.

- ❗ If the ambient temperature difference is more than 20°C (from low temperature environment to high temperature environment or vice versa), please turn off the instrument and rest for 2 hours and then use it again.
- ❗ Do not use in the vicinity of flammable and explosive materials (including sprayers, gas stations) which can easily lead to fire, electric shock injury and other dangers.
- ❗ The ambient temperature of this depilator is 0mm 40X. Please confirm the room temperature before using it.
- ❗ Please do not touch the light outlet after use, because the light outlet produces high temperature, it is easy to cause burns.
- ❗ Please avoid direct ultraviolet radiation within 24 hours after irradiation.
- ❗ Do not place the product in a high-temperature and high-humidity environment, which will cause product deterioration and failure.
- ❗ If the built-in fan does not work after turning on the power supply (you can not hear the bussing sound), do not use it, so as not to cause smoke and fire failure.

■ Unsuitable People and Related Situations

- Do not use it for pregnant women, breast-feeding women, menstruation, diabetes, fragile skin and light-sensitive skin.
- Severe heart disease, moderate and severe hypertension, malignant tumor, folliculitis, infectious dermatosis, burn skin, anticoagulant disorders are prohibited.
- Do not use for those with abnormal collagen metabolism, including history of scar and poor wound healing.
- If you have vascular disease,
-
- do not use it if you have varicose veins and vasodilation at the irradiated site.
- Do not use if the skin is sensitive to light and is prone to rashes or allergies.
- If the irradiation site infection, herpes simplex, skin lesions, hematoma, etc., do not use.
- Do not use to suppress immune response diseases (including HIV infection, lupus erythematosus, dysphagia, etc.,)
- It is forbidden for any patients with hemorrhagic disease, for topical use of isotretinoin, (AHA), (salicylic acid, (BHA) or azelaic acid.
- For those who use chemical hair removal within 30 days, use skin whitening products or cosmetics containing vitamin A in a short period of time.

- Do not use essential oils at the same time. Please make sure the skin in the illuminated area is clean, dry and does not use creams, lubricants and other cosmetics.
- If the skin surface is not clean and dry, it may affect the experience of using the product, and may even damage the product.
- It is prohibited during the period of taking aspirin drugs, using exemption preparations, taking isotretinoin and other drugs.
- Recent use of light-sensitive drugs (such as tetracycline, sulfonamide, etc.) is not recommended.
- Light damage (sunburn) Please do not use sensitive parts such as mucous membrane.
- It is forbidden to use this product if any type of wound is not healed.
- Skin hirsutism caused by hirsutism and other diseases is not recommended.
- If you take painkillers, it will reduce the skin's sensitivity to heat, so it is recommended to use it with caution before the effect is removed.
- If the color of you skin or hair is not suitable, do not use this product.

■ Do not use the following body parts or situations

- | | |
|---|------------------------------|
| 13. Inside the mouth, lips, genitals,
mucous membranes, etc. | 20. Eyebrows, hair, forehead |
| 14. Plastic surgery, surgical site | 21. Ears |
| 15. Adam's apple | 22. Nipple and areola |
| 16. Eyeballs, eyelids, around the eyes | 23. Inside the belly button |
| 17. Head | 24. Large veins |
| 18. Chicken skin, acne area | 25. Tattoos, body painting |
| 19. Spots and moles | |


(Note: this list is not exhaustive. If you are not sure whether you can use this product, we recommend that you consult a doctor.)

BE CAREFUL

Do not use this product for naturally dark skin

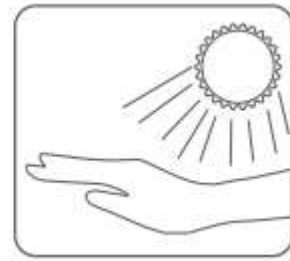Do not expose to the sun before and after hair removal

It is forbidden to use in a dark environment

Please wear protective glasses when using

Before hair removal, please shave your hair first (do not use chemical hair removal methods)

When using for the first time, be sure to start with the first gear

Do not use in direct sunlight to avoid machine failure


Do not violently hit or shake this product.

PRODUCT STRUCTURE

Part Name

- 10) Flash window
- 11) Switch button
- 12) Button (flash/auto flash, cooling switch)
- 13) Display
- 14) DC line interface

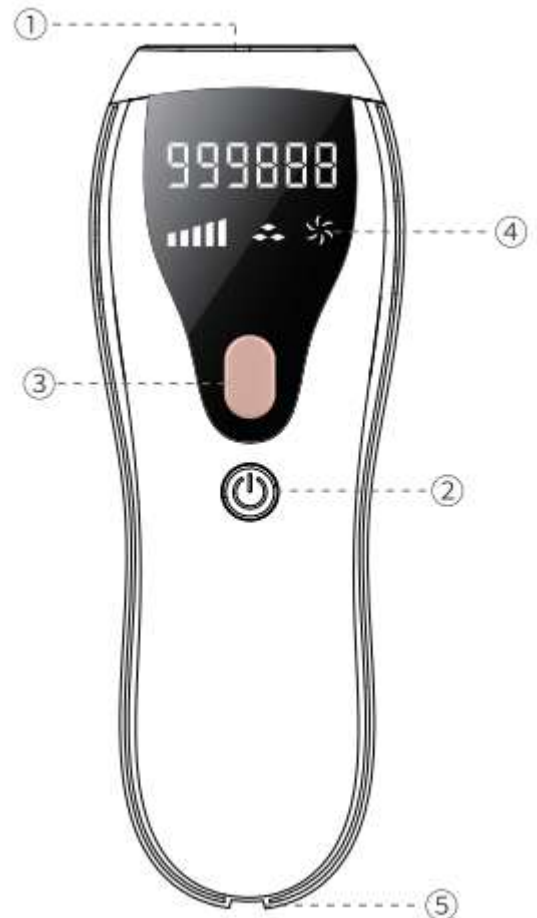
Display screen

- 15) Number of flashes
 - 16) Refrigeration symbol
 - 17) Continuous flashing symbol
 - 18) Energy mark
- (1-5 gears, the higher the gear, the

greater the energy)


It is strictly forbidden to block the air inlet and outlet with sundries, otherwise it will cause serious damage to the product.

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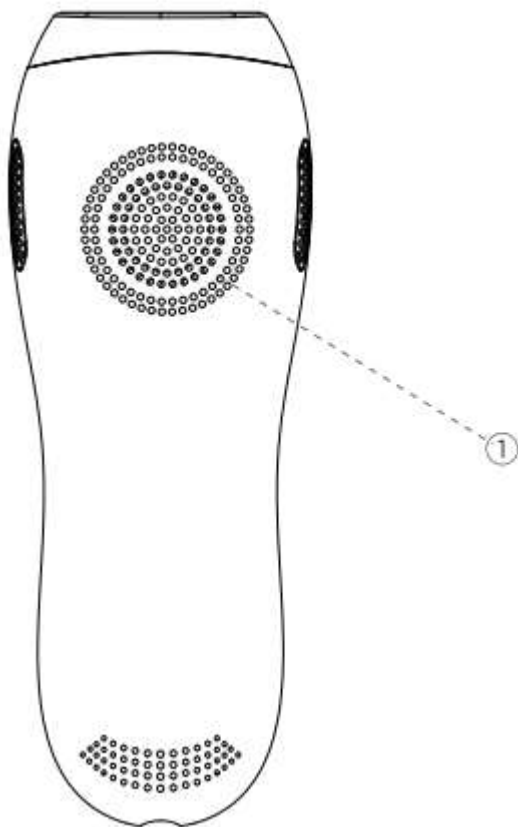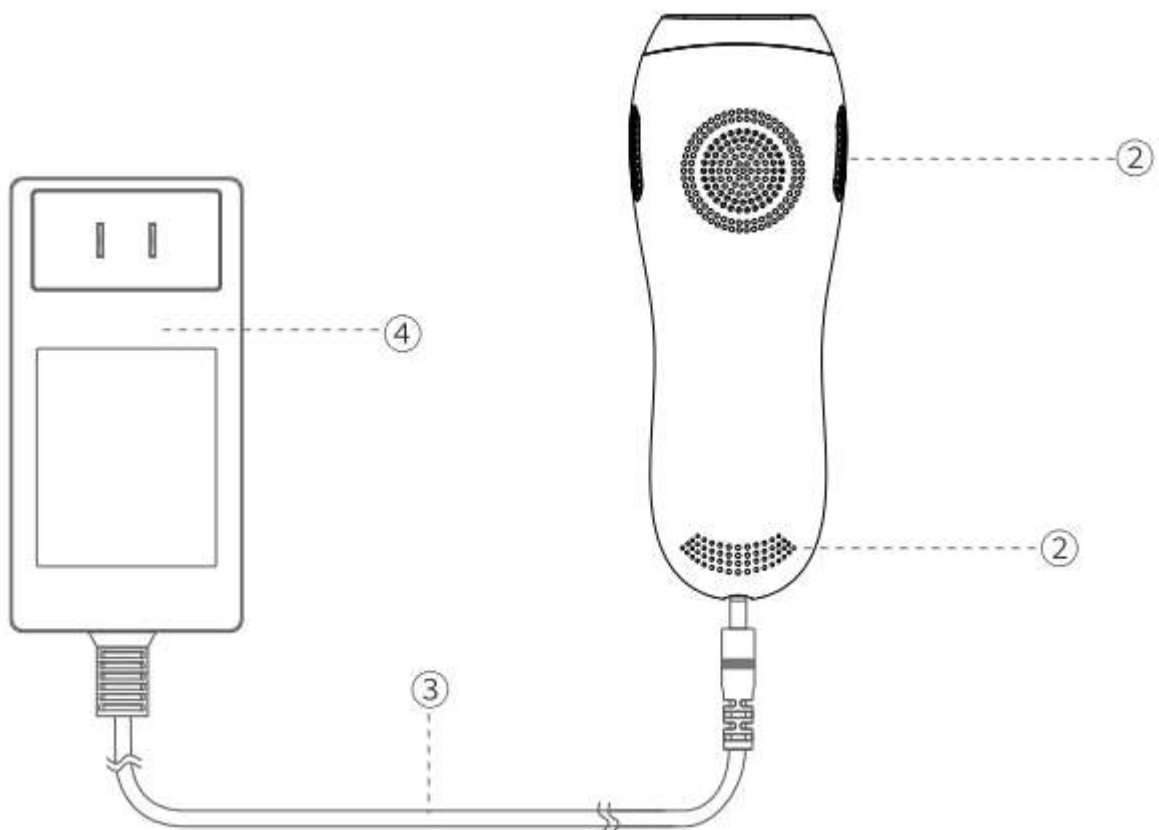DC line x 1

PRODUCT PARAMETER

Product Name	IPL freezing point hair removal
Product Size	185 X 69 X 58 mm
Weight	217.3 g
Enter	12 V = 3 A
Power	36 W
Light Emitting Area	3.36 cm ²
Gear	5 levels adjustable

PREPARE BEFORE USE

- Thoroughly clean the skin to ensure that the skin is hairless, completely dry and free of oily substances (perfume, milk, body fluid, sunscreen, etc.).
 - Make sure that there is no dirt residue such as dust and hair at the irradiation mouth, and there is no damage or abnormality. If there is dirt, please keep it clean so as not to cause product failure, deterioration and skin problems.
- Please do not apply any gel or care products (aloe gel, etc.) before hair removal. This product has ice removal technology, which can take care of the skin in time, making it more comfortable and painless to use

Irradia3tion intensity setting

This product is specially developed for family hair removal, so its light energy is softer, providing 5 different kinds of light energy. You can adjust the light energy to suit your skin and body hair color and make yourself feel comfortable. The light energy meter recommended below will help you determine the appropriate light energy.

3. Please refer to the following table to determine which light intensity is appropriate for your skin and hair color.

Tip: you can put the skin color card on your skin and compare it with light intensity is appropriate for your skin and hair color

4. If you feel any pain or discomfort, please lower the energy gear. People with darker skin are more likely to develop skin allergies because the skin absorbs more energy. Therefore, people with darker skin are advised to use lower stalls. Pulsed light may feel warm or scorching on the skin, but it should not be painful.

Color Hair Color						
White						
Cream-Colored						
Light Brown						
Medium Brown						
Dark Brown						
Brown-Black/Darker						
White Grey	×	×	×	×	×	×
Red	×	×	×	×	×	×
Light Gold	×	×	×	×	×	×
Dark Gold/Beige Brown	3/4/5	3/4/5	3/4/5	1/2	×	×
Brown	3/4/5	3/4/5	3/4/5	1/2	×	×
Dark brown	3/4/5	3/4/5	3/4/5	1/2	×	×
Black	3/4/5	3/4/5	3/4/5	1/2	×	×

(Note: x indicates not applicable)

Note: Please refer to this table to select the appropriate gear

Note: Due to various reasons, your skin may react differently due to different exposure time/situations.

Note: Regardless of the skin color, it is recommended to start with the first gear for the first use and gradually increase the gear.

■ Skin Light Intensity Test

If you are using this product for the first time or have recently been tanned, you should perform a skin test on each position to be irradiated. The skin test is necessary to check whether the skin responds to radiation and to determine the correct light intensity setting for each body part.

8. Select a part near the part to be irradiated.

Note: Do not try to use this product on difficult or sensitive areas (ankles or bones).

9. Press and hold the power button for 3 seconds to turn on, and make sure to select the first gear.

10. Place the device vertically on the skin so that the flashing port is in contact with the skin.

11. Press the flash button to emit a flash.

12. Slide the device against the skin to the next part to be irradiated.

13. Increase the setting step by step, emit a flash, and then slide the device to the next position.

Repeat this operation for all levels within the recommended range of your skin type.

14. After the skin test, confirm the skin condition after 2 hours. If there are no abnormalities such as redness or blisters on the irradiated area, you can use this product with confidence.

INSTRUCTIONS

Steps To Use Intermittent Flash Mode

01 / Shaving

Use a razor to shave off the hair in the area that need care.

- ▲ Skin that has not been shaved or is not cleaned may catch bad things during hair removal.

02 / Moist skin


After shaving, you can moisten the skin with toner, and then start to irradiate hair removal, the experience will be more comfortable.


03 / Disinfection of irradiated surface

Please disinfect the irradiated surface with a damp cotton pad or cotton swab.

- ▲ Cleaning is carried out without power supply

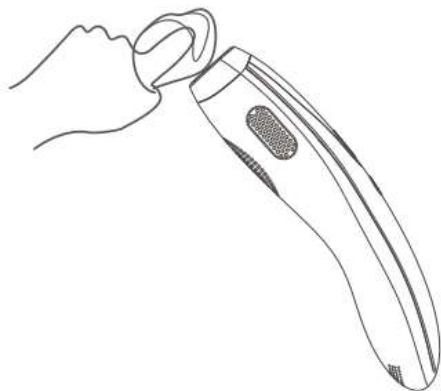04 / Connect to the power supply

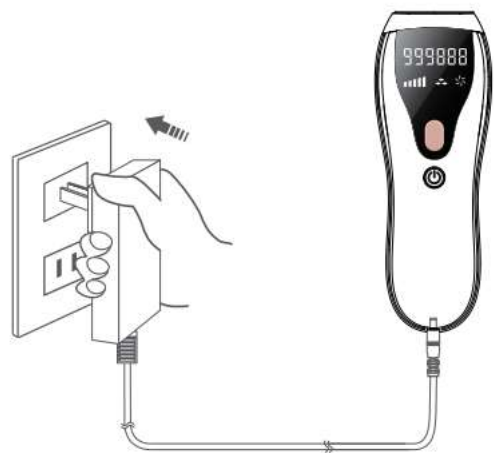05 / Power on

Press and hold the witch key "⏻" for 3 seconds to light up the operation window

06 / Power On Indication

After power-on, it will display the gear (one gear by default) and intermittent flash mode.


Long press "⏻" to turn on/off the continuous flash mode, and short press twice to turn on/off the refrigeration.

07 / Compare skin color card, select gear

Use the skin color card to compare your skin and choose the most suitable gear. If your skin color is the same as the color x on the color card, please do not continue.

- ▲ Short press the energy plus/minus key "⏻" to adjust the energy level. The darker the skin tone, the lower the level. It is recommended to use level 1 for the first use.
- ▲ Under the condition of not damaging the skin, the higher the gear, the better the hair removal effect.

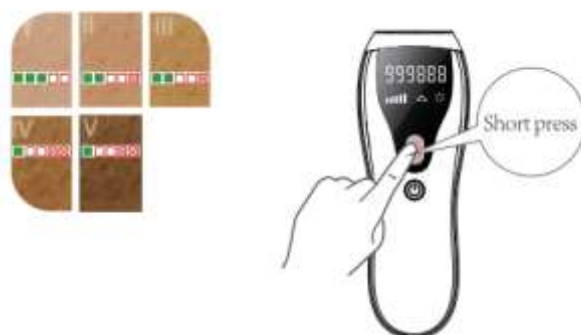
08 / Wear protective goggles

Wear professional protective goggles during operation to avoid eye injury from flashes.

09 / Hair removal while applying ice

At this time, place the product vertically on the part that needs to be depilated, and the lamp head of the body fits perfectly with the skin, press the flash button to emit a flash, and start hair removal while applying ice.

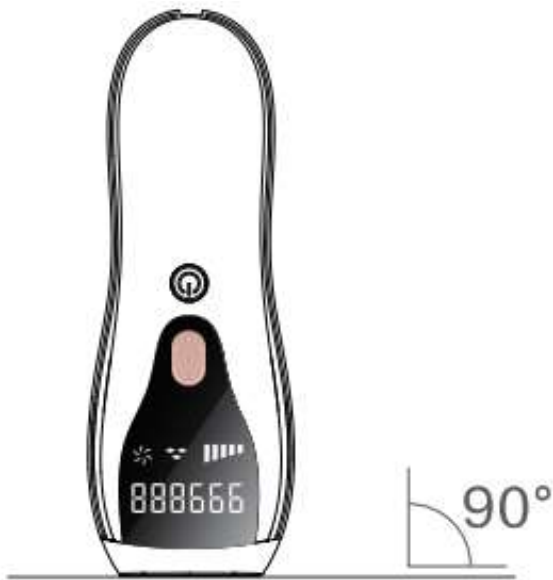10 / Skin care after hair removal

After hair removal, apply aloe essence to moisturize the skin and make it more skin-care.

Steps To Use Continuous Flash Mode

01 / Enter continuous flash mode

After powering on, press and hold the flash button for 3 seconds to enter the continuous flash mode.


02 / Continuous flashing hair removal while applying ice

The ice pack mode is entered by default when the machine is turned on. Press and hold the light button for 5 seconds to turn off the freezing point mode, but it is recommended to turn on the freezing point to achieve painless hair removal.

03 / Start of hair removal

Hold the machine tightly against the skin while slowly sliding to change the position, it can automatically light up without pulling the trigger button. Long press the flash button again to exit the continuous flash and enter the intermittent flash mode.

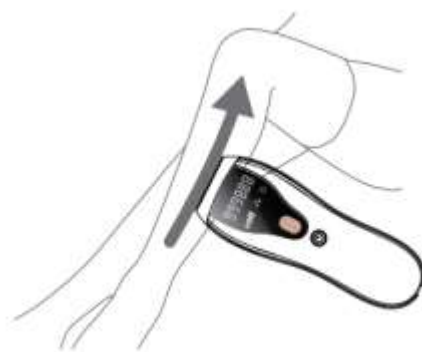
▲ The hair removal will automatically shut down if no one uses it for 5 minutes in all modes.

SKIN CARE

- If the skin turns red after hair removal by light irradiation, please use the nursing mode for nursing. It is recommended to avoid touching water within 6 hours after hair removal.
- Do not use perfume, mask and cosmetics containing scrub beads after skin hair removal.
- Please avoid ultraviolet exposure within 24 hours after hair removal. Apply sunscreen if necessary.
- If you have any discomfort, please consult a surgeon.
- Sunbathing is likely to cause skin blackening or sunburn on the surface of the skin. When irradiated by intense pulsed light at this time, it will cause redness, swelling, inflammation and so on. So during the period when you still need light irradiation to depilate, please be careful not to expose yourself to the sun. Sunburn repair time is not less than 1 week.
- It is forbidden to expose the hair removal area within 48 hours of the hair removal operation, even if you need to be exposed to the sun after 48 hours, you should always pay attention to whether there is a bad condition in the skin. Please use sunscreen after 72 hours.
- The use of artificial light tanning skin (changing the physical or chemical formula of the skin to bronzing) is not allowed during hair removal.
- If you need plastic surgery after hair removal (only refers to the hair removal area), please consult your doctor first.

Note: even if you follow all the instructions, you may still have a skin reaction. In this case, please stop using this equipment and contact the customer service center in your area.

Recommended Exposure Period

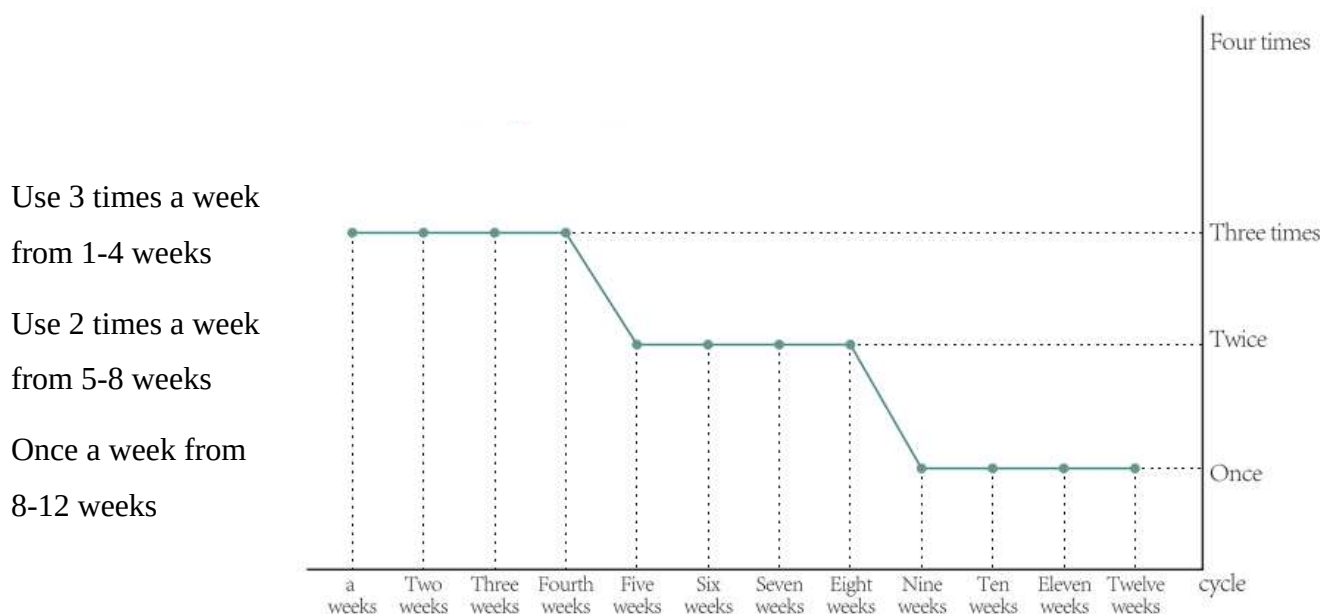- ▲ Please use it for multiple times according to the growth of the hair. For heavy hair, the gear and frequency can be appropriately increased on this basis. The specific situation varies from person to person.

CLEANSING AND MAINTENANCE

Storage And Transportation

In order to ensure that the product achieves the best results and the longest service life, please clean the product before and after each hair removal, and clean the product during hair removal if necessary. If you do not clean normally, this product will not get its due effect.

- Check this product, especially the light-transmitting window and body parts. Make sure they are all clean and tidy.
- Do not clean the product and its accessories under the tap or in the dishwasher.
- Do not use steel wool, abrasive cleaners (toothpaste,) or corrosive liquids (gasoline, engine oil, acetone, etc.) to clean the product.

Note: If you can no longer clean the light-transmitting window, please stop using this product.

- Turn off the product, unplug the power plug, and let it cool down.
- When using this product, the flash window will heat up, make sure it cools down before cleaning.
- If necessary, clean the outside of the product with a clean soft cloth.
- Damp the soft cloth with a few drops of water, and then use it to clean the flash window.

Note: if the water cleaning effect is not good, you can use a few drops of high-concentration alcohol (not more than 96%) to clean the above parts.

Environmental Protection

Then disposing of the product, please do not throw it together with general household garbage, but should give it to the designated recycling center. This will be beneficial to the environment.

Electromagnetic Field (EMF)

This product complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

HAIR REMOVAL ESSENTIALS

Before using the machine, shave off the obvious hair on all parts, and then use a toner to moisten the depilated area with a towel to avoid discomfort due to dry skin.

Hair Removal Essentials

When using a strong pulse light depilator, please follow the following points

- Underarm hair removal
 - ❖ Before depilating, please clean the antiperspirant, perfume, cream and other cosmetics applied to the armpits.
 - ❖ Where the skin is loose, stretch the skin with your fingers first, and then depilate. It is normal for the skin to have a slight burning sensation when irradiating the mouth.
- Bikini line hair removal
 - ❖ For areas with accumulated pigments and dark colors, please start using the low gear first, and use the first gear before adapting.
- Hands and feet hair removal
 - ❖ Bones and pits and bumps are more prone to redness and pain than other parts, but this is normal.
- ▲ Do not use the following parts

Areas with burn marks and dark colors, eyeballs, eyelids, around the eyes, forehead and moles, eyebrows, forehead, hair, sun burns, etc.

 - ❖ After hair removal, avoid the place with ultraviolet rays for 24 hours to prevent the sun from burning the skin during hair removal.
 - ❖ Be careful not to look directly at the flash.

POSSIBLE EFFECTS

Adverse Reactions

- Skin slightly red or slightly itchy, warm, these are normal phenomena, are harmless to the body, and will soon disappear.
- After hair removal, the reaction of the hair removal laser for a short time is similar to that of exposure to the sun, If it is not eliminated after 48 hours, it is recommended that you consult a doctor.
- Shaving and hair removal work together in a similar period of time, which will cause dryness and itching of the skin. This phenomenon will gradually reduce the discomfort of the normal reaction. You can apply ice to cool this one. After 24 hours, you can apply a mild and non irritating moisturizer.

Rare Serious Adverse Reactions

5. Burns, severe redness, and bloating pain: these reactions are rare and are usually caused by the use of light that is too high for the color of the skin. If these reactions do not disappear within three days, you are advised to consult your doctor. Wait until the skin has healed before a second exposure, and be sure to use a lower light intensity.
6. Skin discoloration: this phenomenon rarely occurs. Skin discoloration usually appears darker or brighter than the surrounding area. This is due to the fact that the light intensity used is too high for skin color. In the case of auxiliaries usually due to the use of too high not disappears and the skin turns gray back to normal complexion.
7. Skin infection and inflammation: this rarely occurs and may be caused by the use of this product on shaving wounds or incisions, existing wounds, or inward-growing hair.
8. Extreme pain: if you use this product on unshaved skin: the light intensity you use is too high for your skin color; you flash many times in the same position; and if you use this product on open wounds, inflammation, infection, tattoos or burns, this reaction may occur during or after hair removal.

PROBLEMS

This chapter outlines the problems encountered when using this product. If you cannot use the following information to solve the problems encountered, please contact the customer service center of your country/region.

Fault	Possible Reasons	Solution
Unable to power on.	Power is not connected	Check if the power supply is plugged in.
Can't get the light.	Whether the temperature of the host head and body is too high. Is the host head close to the skin.	Wait for the machine to cool down naturally before using it. Place the host close to the ferret and wait for the light indicator to turn on.
The work indicator is not on.	Bad machine.	Check whether the power supply is working properly, repair or return the product immediately.
There is a strange smell in the flash, and the part of the hair turns red.	The filter is dirty and contaminated, and the body surface hair remains.	Wipe the filter to shave off the hair and reduce the light intensity.
When depilating the pain is unbearable.	Did not shave the hair in the area to be depilated.	Before using this product, shave the hair on the part to be irradiated. If shaving causes skin irritation, please cut the hair as short as possible and use it in a comfortable setting.
	The part being depilated is not suitable or using this product.	Do not use this product on the inner labia, vagina, anus, nipples and areola. Men should not use this product on the face or scrotum.
	The light intensity is too high.	Reduce the light intensity to a comfortable level.
	Damaged light-transmitting window.	Repair or return the product immediately.
After hair removal, the depilated area turns red.	Slight redness is normal, it is harmless to the body, and it will subside quickly.	Use aloe vera gel to soothe the skin or use a lamp cap for ice soothing.
The hair removal effect Is not satisfactory.	You used a light intensity setting that is too low for you.	Please choose a higher setting next time.
	You didn't shine a flash of light on the area adjacent to the previously irradiated area.	You should emit flashes in areas close to each other, and there should be some overlap between flashing areas.
	This device has no effect on your hair color.	If your hair is light blonde, gray, red or white, the irradiation will not be effective.

	Insufficient frequency of using this device.	If your hair is light blonde, gray, red or white, the recommended irradiation plan.
	Your response to light exposure is slow.	This product should be used continuously for at least 6 months, because it can still inhibit hair regrowth during this period.
The product becomes hot and humid during use.	May be due to long working hours of the instrument.	Turn off the instrument and let the instrument rest and fully dissipate heat before using it again.

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