
STORMRED

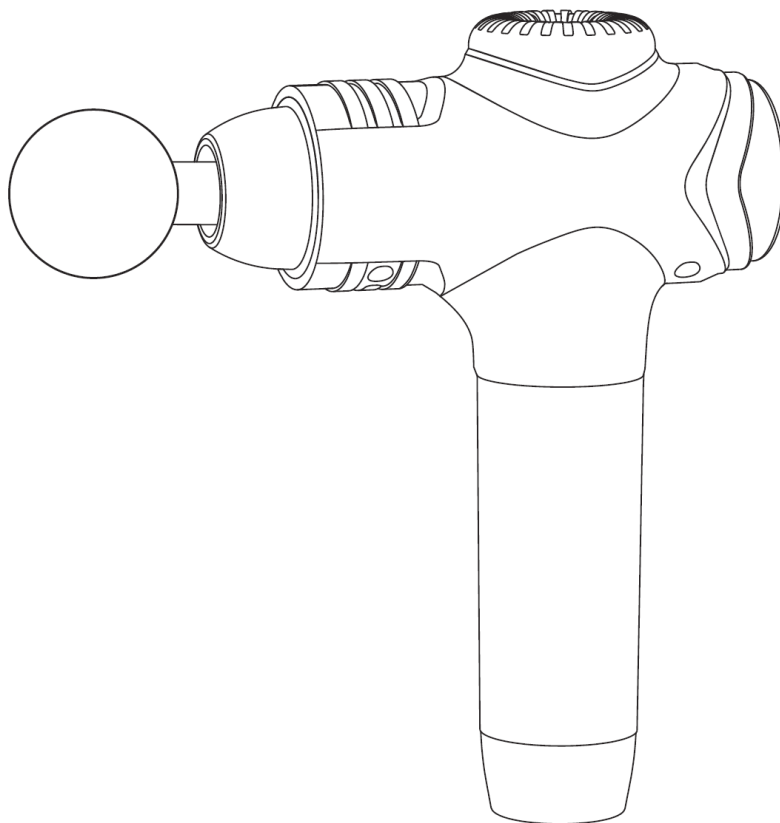PRO REHAB

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Important Safety Instructions

To reduce the risk of fire, electric shock, and personal injury, please follow these instructions carefully when using this product.

- For adult use only.
- Use only on a dry, clean surface of the body or clothing by lightly pressing and moving across the body for approximately 60 seconds per area.
- Use the massage gun only on the body's soft tissue as desired without producing pain or discomfort.
- Do not use on head or bone areas of the body.
- Use only the attachment that provides the best desired effect.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers and hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the massage gun.
- Do not immerse in water or allow water to enter the ventilation holes of the massage gun.
- Do not drop or misuse the device.
- Carefully examine device and battery level before each use.
- Do not tamper or alter massage gun in any way.
- Never leave the massage gun operating or charging unattended.

Please do not use the massage gun or any vibration device without first obtaining approval from your doctor if any of the following apply:

- Pregnancy, diabetes with complications such as neuropathy or retinal damage, use of pacemakers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis or spondylosis, recent joint replacement or IUDs, metal pins or plates or if there are any concerns about your physical health.
- Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device, but we advise you to consult a doctor first.
- Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above.
- Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in number of cases. This must be done on the advice or and in accompaniment with a doctor, specialist, or physiotherapist.

Note: You should understand that when participating in any exercise or therapy, there is the possibility of physical injury. Any exercise program including this one may result in injury. Consult with a medical professional before beginning this or any exercise or therapy program.

Features

Relieve muscle soreness and stiffness with the Stormred PRO REHAB muscle massage gun. Ideal for relieving shoulder, neck, waist, and back pain, this massager is your answer for soothed and relaxed muscles. It includes 10 different massage heads for a comfortable massage experience. The muscle massage gun is compact, portable, and ultra-quiet. Choose from 10 massage speeds which allow you to get the right and appropriate amount of pressure to relieve your pain and improve muscle strength, flexibility, and coordination.

By the rhythmic hammering of the massage head against your muscles, it helps to accelerate athlete recovery from tension workout, and make it easy for those who need to relieve stress, feelings of anxiety, soreness, pain, or muscle tension.

Massage guns help with the following:

- Treats muscle bruises, sprains, strains, and pains.
- Increases circulation and lymphatic flow.
- Relieves muscle spasms and stiffness.
- Reduces the accumulation of lactic acid in muscles.
- Activates the nerve system and muscles.
- Gentle muscle stretch.
- Increases joint mobility.
- Relaxation and thickening of connective tissue and fascia.



Specifications

Model	PRO REHAB
Amplitude	12mm
Massage head	10× adjustable
Working noise	≤40dB
Charging time	~ 3-4 hours
Motor run speed	1200-3200 RPM
Charging voltage	100-240V
Running time	2-5 hours (depending on speed level)
Battery	Undetachable
Adjustable levels	10 levels
Weight	839g

Package Contents

- Stormred PRO REHAB massage gun
- 10 massage heads
- USB cable
- Portable case

Product Diagram and Functions

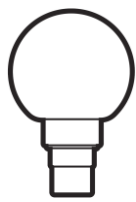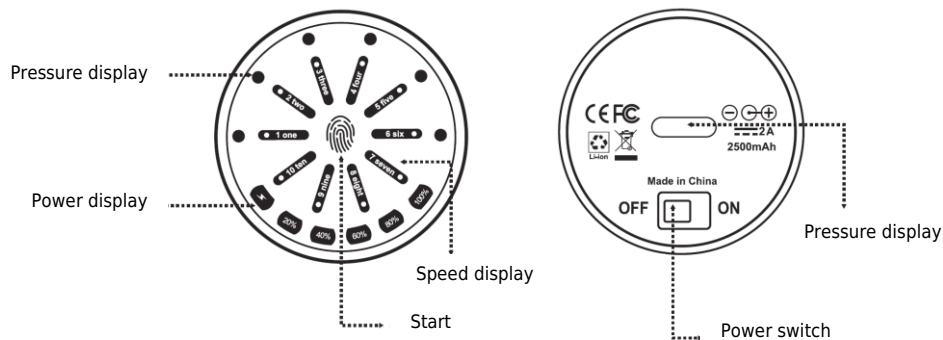Round head

Suitable for thighs, calves, back.

Flat head

Suitable for all body parts.

Bullet head

Suitable for local pain points.

Spade head

Relaxes hamstring, glutes, improves flexibility.

Shovel head

Suitable for abdominal muscles and low back.

Arrow head

Suitable for joint deep tissue, palm, and foot massage.

Fork head

Suitable for neck, spine, and Achilles tendon.


Elbow head

Suitable for palm and sole massage.

Convex head

Suitable for quads, gluteal.

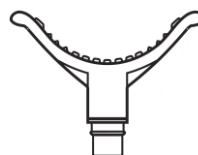Y-shaped head

Suitable for arm and leg muscles.

To **attach** the massage head, choose a suitable massage head according to your needs and rotate the massage head firmly into the massage gun.

To **detach** the massage head, rotate the massage head to loosen it and pull it out of the massage gun.

Charging

- It is recommended to use a 5V/2A charging plug for charging.
- Fully charge the massage gun before initial use.
- Connect the USB charging cable to the charging port.
- Battery can be recharged at any time/any battery level.
- Battery should not be fully depleted prior to recharging.
- Battery charge level is indicated on the bottom of the LED screen.
- For your safety, you should never turn on the massage gun when the battery is charging.
- For security reasons, the massage gun cannot be turned on while charging. It only works properly when charging is completed and charging adapter is disconnected.

Note: When using electrical appliance, always observe safety regulations to reduce the risk of fire, electric shock, and personal injury. Always check that power source corresponds to the voltage on the rating label of the device.

Cleaning and Maintenance

- Before cleaning, be sure to turn off the massage gun.
- Gently clean the massage gun with slightly damp towel and dry it with a soft cloth.
- Do not attempt to clean internal surfaces – cleaning should not require disassembling the device. Disassembling the device to clean internal mechanisms could result in device failure.
- Turn off the power and store it in the carrying case when not using it.
- Keep the massage gun away from heat sources, liquids, and children.

Troubleshooting Guide

If you encounter any issues with your massage gun, try to solve them using the troubleshooting guide. Should the problem prevail, please contact Alza customer support by means mentioned at the beginning of this user guide.

Why can't I turn on my massage gun?

- Check whether the power switch is in the "ON" position.
- Check the battery status. Before the initial use, charge the massage gun for about six hours.

Why does the massage gun stop working in the course of a massage?

- The massage gun automatically stops after 10 minutes in order to remind you not to massage one body area for too long, as it may cause a bruise. Restart the massage gun by pressing the Start button.

How hard should I press the massage gun against my muscles?

- The best practice is to start the device away from your body and on lowest settings, then lightly contact your body with the massage head and allow the device to float across your tissue. Once you have a feel for the pressure you desire, you can increase the pressure slowly until you reach a comfortable balance.

Why the massage gun stops working after applying pressure to the body?

- Restart the massage gun. Should the problem prevail, please contact Alza customer support by means mentioned at the beginning of this user guide.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

This equipment is in compliance with the essential requirements and other relevant provisions of EU directives.

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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