

SENCOR®

SOC 0910BL / SOC 0911RD
SOC 0912GR



EN ■ Kids sonic electric toothbrush

Translation of the original manual

READ CAREFULLY AND STORE FOR FUTURE USE.**General warnings**

- Toothbrushes may be used by children and persons with reduced physical or mental capabilities or with a lack of experience and knowledge, if they are supervised or have been instructed to use the appliance in a safe manner and understand the potential dangers.
- Children must not play with the appliance. Cleaning and maintenance performed by the user must not be performed by unsupervised children.
- In the event that the power cord is damaged, have it repaired at a professional service centre to prevent causing a dangerous situation. It is forbidden to use the appliance if it has a damaged power cord.
- The parts of the appliance that need to be mounted or stood upright, must be mounted or stood upright in such a way that they cannot fall into water.
- The appliance must only be powered using a small safe voltage corresponding to the markings on the appliance.

**Attention:**

The toothbrush must only be used with the original charging station that is supplied with it. In order to prevent a dangerous situation from arising, never use a different type of charging station.

Electrical safety

- Before connecting the charging station to a power socket, make sure the voltage stated on its rating label corresponds to the voltage in your power socket.
- Connect the only to a properly grounded socket. Do not use an extension cord.

- Do not disconnect the charging station from the power socket by pulling on the power cord. This could damage the charging station or the power socket. Disconnect the charging station from the power socket by gently pulling on the plug of the power cord.
- Do not place the power cord over sharp objects. Make sure that the power cord does not hang over the edge of a table and that it is not touching a heated surface.
- Use only the original charging station with the power cord to charge the toothbrush. Never use non-original accessories.
- Never use the charging station supplied with this toothbrush for any other purposes.
- Keep the charging station and its power cord dry.
- Always locate the charging station on an even and stable surface.
- Do not wind the power cord around the handle of the toothbrush.
- Do not connect and disconnect the power cord of the charging station to/from the power source with wet hands.
- Store, charge and use the toothbrush at an ambient temperature of 15 to 35 °C. Do not charge the toothbrush for longer than necessary for it to fully charge. Once it is charged, disconnect it from the power source.
- The toothbrush is equipped with a built-in long-life battery that can only be removed from the toothbrush at a professional service centre. Never attempt to remove it from the toothbrush yourself.



Warning:

Never disassemble the toothbrush for the purpose of removing the battery. Never short circuit the toothbrush or throw it into a fire. Incorrect use, risks releasing toxic substances or causing a fire or explosion.

- Do not expose the toothbrush to temperatures higher than 50 °C. Otherwise, the built-in battery could be damaged.
- Do not charge, use or store the toothbrush or leave it lying in a bathroom or in locations where it could fall into a bath, wash basin or other containers filled with water.
- In the event that the toothbrush or the charging station fall into a bath, under no condition reach for it, rather immediately disconnect the power cord from the power socket. Never use a toothbrush or a charging station that has been submerged in water. Hand it over to an authorised service centre for inspection.
- In the event that the toothbrush, charging station or its power cord are damaged in any way, do not use them. The power cord of the charging stand cannot be replaced. In the event that it is damaged, the charging stand must be disposed of


and replaced with an original type. For replacement of damaged accessories and for all repairs of the toothbrush, contact an authorised service centre.

Safety during use


- This toothbrush is intended for domestic use only and is used for cleaning teeth and gums. This toothbrush is not intended for commercial purposes or for treating animals.
- For hygiene reasons, the toothbrush head should be used by only a single user. Regularly replace the toothbrush head every 3 months or sooner if it shows signs of wear. Use only original replacement toothbrush heads.
- When the toothbrush is used by children, supervise them while they are brushing their teeth until they acquire the proper tooth brushing technique.
- Use the toothbrush and toothbrush head in accordance with the instructions contained in this user's manual.
- Use only the original accessories supplied with this toothbrush.
- This toothbrush and its charging stand must only be used in indoor areas. Do not use it outdoors.
- Do not use the toothbrush in rooms where aerosol sprays are used or where oxygen is administered.
- Do not expose the toothbrush and its accessories to the outdoor environment, extreme temperatures, direct sunlight, excessive humidity or excessively dusty environments.
- Do not place the toothbrush and its accessories near heaters, open fires or other sources of heat.
- Before every use, check the toothbrush and the toothbrush head for damage. Never use a toothbrush with a damaged toothbrush head.
- Keep the toothbrush at a safe distance from eyes.
- Keep the toothbrush clean. Perform cleaning and maintenance according to the instructions in chapter Cleaning and maintenance.
- Do not perform any maintenance on the toothbrush other than cleaning the individual parts as described in chapter Cleaning and maintenance. Do not wash the toothbrush or its accessories in a dishwasher.
- Do not use the toothbrush if it has fallen on the floor, if it is not functioning correctly or if it is damaged in any way.
- In order to prevent a hazardous situation from arising, do not repair the toothbrush or any of its accessories yourself or modify it in any way. For all repairs, contact an authorised service centre. By tampering with the product you risk voiding your legal rights arising from unsatisfactory performance or quality warranty.
- When using a whitening toothpaste containing peroxide, baking soda or sodium bicarbonate, thoroughly clean the toothbrush head using soap and water. This will prevent the plastic from cracking.

- In the event that dental work or oral surgery was performed in the past two months, consult the use of a sonic toothbrush with your dentist.
- In the event that gum bleeding occurs when the toothbrush is used and if bleeding does not stop even after a week of its use, contact your dentist.
- In the event that using the toothbrush causes unpleasant or painful feelings, stop using the toothbrush and consult this issue with your dentist.
- In the event that you have a pacemaker or another implanted device, consult the use of the toothbrush with your doctor or manufacturer of the implant.


A

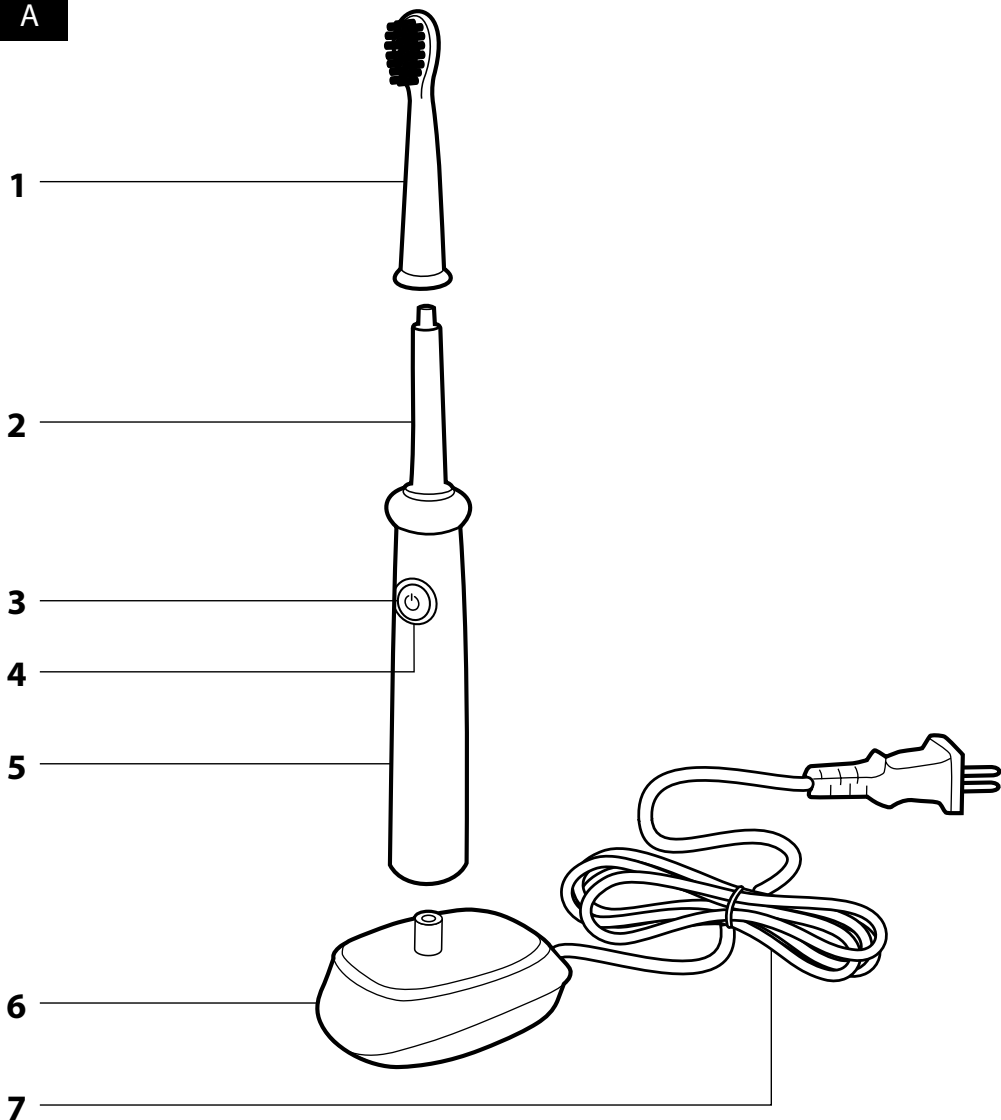

B

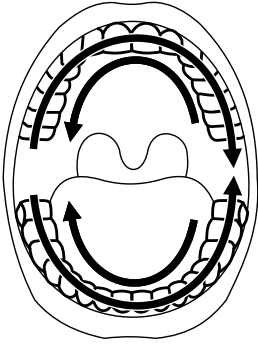

C

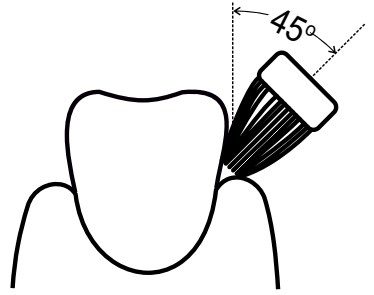

- Prior to using this appliance, please read the user's manual thoroughly, even in cases, when one has already familiarised themselves with previous use of similar types of appliances. Only use the appliance in the manner described in this user's manual. Keep this user's manual in a safe place where it can be easily retrieved for future use. In the event that you hand this appliance over to somebody else, make sure to also include this user's manual.
- Carefully unpack the appliance and take care not to throw away any part of the packaging before you find all its parts. We recommend saving the original cardboard box, packaging material, purchase receipt and responsibility statement of the vendor or warranty card for at least the duration of the legal liability for unsatisfactory performance or quality. In the event of transportation, we recommend that you pack the appliance in the original box from the manufacturer.

DESCRIPTION OF THE APPLIANCE AND ITS ACCESSORIES

A1 Toothbrush head with a soft rubber rear side

A2 Vibrating shaft

A3 Mode / power toggle switch

A4 Two-colour LED cleaning / charging mode indicator

A5 Ergonomically-shaped handle

A6 Charging station

A7 Power cable

PURPOSE OF USE

- This sonic toothbrush is to only be used for personal dental hygiene of children 6–12 years of age. It is intended for home use and for private purposes. It is not intended for business activities or other professional areas of use nor for the treatment of animals.
- This sonic toothbrush is able to remove dental plaque several times better than when teeth are brushed using a standard toothbrush. Using the two preset programs, care for your gums using specially set speed of 30000 vibrations per minute, whilst consuming an extremely low amount of power.

CHARGING

- The toothbrush is equipped with a built-in battery that serves as an internal power source and enables cordless operation. Before first use or after not using it for a long time, first allow the toothbrush battery to fully charge. Only use the original charging station **A7** to charge the toothbrush.
- Locate the charging station **A7** on an even and stable surface and plug the power cord into a power socket.
- Check that the toothbrush is turned off and put it down on to the pin of the charging station. The built-in battery will start to charge, which is indicated by the white LED indicator **A4** turning on and off. The first charge takes approximately 3 hours. Each subsequent charge will take approximately 2 hours to fully charge the battery. If you are charging the toothbrush for the first time or if you have not used it for a long time, it may take up to 10 minutes before the LED indicator **A4** lights up.
- As soon as the battery is fully charged, the indicator **A4** will remain lit in blue colour. Pull the power cord out of the power socket. When the toothbrush is disconnected from the power socket, the LED indicator will turn off.



Attention:

Do not leave the toothbrush connected to the power source longer than is necessary to charge it.

- A fully charged toothbrush will provide up to 15 days of operation when teeth are brushed for 2 minutes twice a day.
- When the indicator **A4** starts flashing white or blue during operation (depending on the mode), it means that the battery is almost flat and that it needs to be recharged.
- Recommendations for achieving the longest possible battery lifetime:
 - In order for the built-in battery to serve you as long as possible, charge, use and store the toothbrush in a dry environment at an ambient temperature of 15–35 °C. Lower or higher temperatures have a negative impact on the lifespan of the battery.
 - Do not wait for the battery to become completely flat. It is better to charge the battery before it is completely flat.
 - When the toothbrush is not used for an extended period of time, charge it at least once every 6 months.

USING THE TOOTHBRUSH

ATTACHING AND REMOVING THE TOOTHBRUSH HEADS

- Slide the toothbrush head **A1** on to the vibrating shaft **A2** from above. Prior to this, turn the toothbrush head so that the grooves inside the brush head slide on to the tabs on the shaft so that the brush head is completely seated on the toothbrush handle.



Attention:

Do not attempt to push the toothbrush head down on to the toothbrush in a different direction than permitted by the tabs on the shaft. This could damage the toothbrush head and the shaft.

- Remove the toothbrush head by pulling it upwards.

SETTING THE BRUSHING MODE

- Repeatedly press switch **A3** to switch between the CLEAN and SENSITIVE modes. The CLEAN mode is intended for standard tooth brushing and is indicated by the white colour of the indicator **A4**. The SENSITIVE mode is intended for brushing sensitive teeth and is indicated by the blue colour of the indicator.



Tip:

If you leave the toothbrush running until it turns off automatically, the last selected mode will be saved to memory and when the toothbrush is turned on the next time it will run in the mode that was set previously.

CLEANING

- Moisten the toothbrush bristles and apply toothpaste on to them.
- Place the toothbrush into your mouth so that the toothbrush bristles are at an angle of approximately 45° to the gum line, see figure **B**. Then turn on the toothbrush.
- Start cleaning the first dental quadrant and slowly move the toothbrush from one tooth to the next. Do not brush your teeth as you would with a standard toothbrush. To achieve maximum efficiency, only push gently on the toothbrush and allow it to clean your teeth for you.

QUADPACER FUNCTION

- The toothbrush has an intelligent timer function that momentarily pauses the toothbrush every 30 seconds to inform you that it is time to start brushing the next dental quadrant. This will ensure that all the dental quadrants are brushed for the same amount of time. To correctly utilise this function, it is recommended to divide the dental area into four quadrants, see fig. **C** and then after every 30 seconds of brushing move from one quadrant to the next until all your teeth have been brushed.

AUTOMATIC SHUT-OFF TIMER

- The automatic shut-off timer automatically turns the toothbrush off after 2 minutes of operation. This is the tooth brushing time recommended by dentists for the sonic toothbrush.

CLEANING AND MAINTENANCE

- For hygiene reasons, regularly clean the toothbrush after every use.
- Make sure that the toothbrush turned off before cleaning it.
- Thoroughly rinse the toothbrush head **A1** under lukewarm running water to remove any toothpaste or other remains. If whitening pastes are used when brushing teeth, it is recommended to apply a drop of soap on to the toothbrush bristles and then thoroughly wash it under running lukewarm water. At least once a week, remove the toothbrush head from the toothbrush and rinse out its inner part under running water.
- Rinse the toothbrush handle **A5** under running lukewarm water and then dry it.



This toothbrush can be washed under a water tap.

- In the event that the charging station **A6** is soiled, wipe it with a clean dampened cloth and then wipe dry. Ensure that the power cord does not come into contact with water or moisture.
- For cleaning, do not use abrasive cleaning products, thinners, etc. Otherwise the surface of the toothbrush may be damaged.
- Before you start cleaning the toothbrush, turn it off and disconnect it from the power source.

REPLACING THE TOOTHBRUSH HEAD

- Regularly replace the toothbrush head **A1** every 3 months or sooner if it shows signs of wear. Use only original toothbrush heads.

STORAGE

- When not using the toothbrush for an extended period of time, disconnect the charging stand **A6** from the power source and put the toothbrush in a dry, clean place. Keep the stored toothbrush and its accessories away from extreme temperatures, direct sunlight and excessive humidity and do not store it in an overly dusty environment.

TECHNICAL SPECIFICATIONS

Toothbrush	
Built-in battery.....	Ni-MH, 1.2 V, 600 mAh
Operating time when fully charged	up to 15 days
Frequency	up to 30000 / min.
Noise level	65 dB(A)
Electrical shock protection class.....	II
Protection class level against damaging entry of water	IPX7

Charging station	
Rated voltage range	100–240 V AC
Rated frequency.....	50 / 60 Hz
Power input	1.4 W

The declared noise emission level of the appliance is 65 dB(A), which represents a level A of acoustic power with respect to a reference acoustic power of 1 pW.

EXPLANATION OF TECHNICAL TERMINOLOGY

Safety class for electrical shock protection:
Class II – Electrical shock protection is provided by double or heavy-duty insulation.

Protection class level against damaging entry of water:
IPX7 – The appliance is protected against damaging entry of water when submerged at a depth of 1 metre for 30 minutes.

Changes to text and technical parameters are reserved.

DISPOSAL OF USED BATTERIES

Prior to handing over the toothbrush at a specialised collection facility, it is necessary to ensure that the built-in battery is removed. Battery removal will be performed at the authorised service centre. In the event that the product has been disassembled for the purpose of ecological disposal then it must not be reassembled.



This symbol indicates that used batteries do not belong in ordinary municipal waste. Hand over used batteries to designated collection points for proper disposal. Correct disposal of batteries will help prevent negative impacts on the environment and human health.

Disposal in other countries outside the European Union.

This symbol is valid in the European Union. If you wish to dispose of used batteries, request the necessary information about the correct disposal method from the local authorities or from your vendor.

INSTRUCTIONS AND INFORMATION REGARDING THE DISPOSAL OF USED PACKAGING MATERIALS

Dispose of used packaging material at a site designated for waste in your municipality.

DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT



This symbol on products or original documents means that used electric or electronic products must not be added to ordinary municipal waste. For proper disposal, renewal and recycling hand over these appliances to determined collection points. Alternatively, in some European Union states or other European countries you may return your appliances to the local retailer when buying an equivalent new appliance. Correct disposal of this product helps save valuable natural resources and prevents potential negative effects on the environment and human health, which could result from improper waste disposal. Ask your local authorities or collection facility for more details. In accordance with national regulations penalties may be imposed for the incorrect disposal of this type of waste.

For business entities in European Union states

If you want to dispose of electric or electronic appliances, ask your retailer or supplier for the necessary information.

Disposal in other countries outside the European Union.

This symbol is valid in the European Union. If you wish to dispose of this product, request the necessary information about the correct disposal method from the local council or from your retailer.



This product meets all the basic requirements of EU directives related to it.