RENPHO

Body Composition Scale

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Important Safety Instructions

Please read all instructions carefully before using this scale. Following these guidelines will ensure proper use and prevent damage or injury.

- **Household Use Only:**
 - This scale is intended for household use only. It is not designed for medical or commercial purposes. The measurements provided should not be used for diagnosing medical conditions.
- **Special Considerations for Children and Certain Individuals:**
 - Children under the age of 8, pregnant women, and individuals with certain medical conditions—such as osteoporosis, pacemakers, artificial joints, or portable ECG devices—should use the scale only for weighing purposes and avoid using the body composition features. These individuals should wear shoes to avoid potential interference or electric current risks when stepping onto the scale.
- **Body Composition Data is Reference Only:**
 - All body composition information provided by the scale is for reference purposes only. It should not replace professional advice from a licensed healthcare provider.
- **Avoid Overloading the Scale:**
 - Do not exceed the maximum capacity of 180 kg. Overloading the scale can lead to damage and inaccurate measurements.
- **Inspect for Damage Before Each Use:**
 - Always inspect the scale for visible signs of damage before using it. If the scale appears cracked, damaged, or unstable, do not use it as this may lead to inaccurate readings or potential injury.
- **Check the Batteries for Malfunction:**
 - If the scale malfunctions, check the battery level. Replace the batteries if needed, following the manufacturer's instructions.
- **Keep Batteries Out of Reach of Children:**
 - Batteries pose a significant risk if ingested. Ensure that batteries are stored in a safe place away from children. If a battery is swallowed, seek medical attention immediately.


- **Ensure Proper Surface for Use:**
 - To get accurate results, always place the scale on a flat, hard, and dry surface such as tile or wood. Avoid using the scale on carpets or rugs, as this may affect accuracy.
- **Avoid Electromagnetic Interference:**
 - Do not use the scale in environments with strong electromagnetic fields, such as near electronic devices or appliances that emit radiation. This could interfere with the scale's sensors and lead to inaccurate results or device malfunction.
- **Use the Scale Safely:**
 - Do not stand on the edges of the scale, jump on it, or stomp on it, as this can cause injury or damage the scale's internal mechanisms.
- **Handle with Care:**
 - The scale is a delicate and precise device. Handle it carefully to avoid breakage or malfunction.
- **Keep the Scale Dry:**
 - The scale is **not waterproof**. Never submerge it in water or expose it to excessive moisture. Clean the surface with a damp cloth or glass cleaner, but avoid using soap, chemicals, or abrasive cleaning products.
- **Avoid Exposure to Extreme Conditions:**
 - Do not expose the scale to direct sunlight, high temperatures, or high humidity. Store it in a cool, dry place when not in use to maintain its performance and durability.
- **Do Not Disassemble:**
 - Do not attempt to disassemble or repair the scale yourself. Tampering with the internal components may cause permanent damage and void the warranty.
- **Dispose of Batteries Properly:**
 - When batteries are depleted, remove them from the scale immediately and dispose of them following your local recycling guidelines. Do not dispose of batteries in fire. If old batteries are left inside the scale for too long, they may leak electrolytes, causing damage to the scale and potentially leading to personal injury.

- **Battery Safety:**
 - Do not attempt to heat, short-circuit, force discharge, or otherwise modify the batteries. This could lead to leakage, explosion, or other hazards.
- **Do Not Solder or Weld Batteries:**
 - Never attempt to solder or weld connections directly onto the batteries, as this can cause them to overheat or explode.
- **Store Batteries Safely:**
 - Store unused batteries in their original packaging and away from metal objects to prevent accidental short-circuiting. Do not mix unpacked batteries or store different battery types together.
- **Remove Batteries for Long-Term Storage:**
 - If the scale will not be in use for an extended period, remove the batteries to prevent leakage or corrosion.
- **Avoid Extreme Temperatures:**
 - Do not expose the scale or its batteries to extreme temperatures, whether hot or cold, as this can affect performance and battery life.

Installing the Batteries

Before using the scale, you need to install the provided batteries. Follow these steps for proper installation:

- Insert the batteries into the battery compartment, ensuring that you align them according to the correct polarity. The positive (+) and negative (-) ends should correspond to the markings on both the battery and the scale.

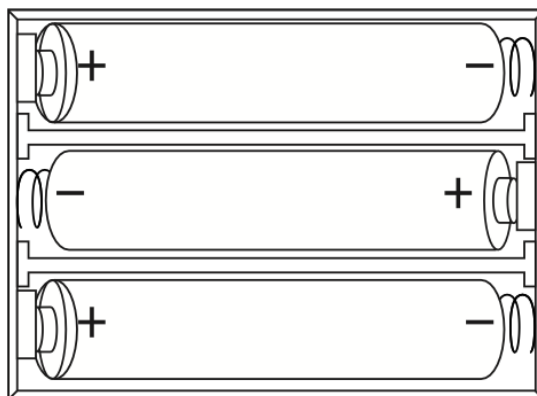
Note:

- Always use the appropriate size and type of batteries for the scale as recommended by the manufacturer. Never mix old and new batteries or combine batteries of different types or brands, as this could lead to performance issues or potential damage.
- Do not attempt to charge non-rechargeable batteries. Charging non-rechargeable batteries can cause leakage, overheating, or even explosion.
- Before inserting the batteries, clean the contacts on both the batteries and the scale to ensure a proper connection.
- If the display shows "Lo," it indicates low battery power. In this case, replace all batteries at once to avoid uneven power distribution, ensuring a longer battery life.

Preparing the Scale for Use

To ensure accurate measurements, always place the scale on a flat, hard surface. Avoid placing the scale on soft surfaces like carpets, rugs, or uneven flooring, as this can affect the accuracy of the results.

- **Correct Placement:** Use the scale on a hard, flat surface like tile or wood flooring to ensure the best results.
- **Incorrect Placement:** Avoid using the scale on soft or uneven surfaces, such as carpets or rugs, which can interfere with the scale's sensors and lead to inaccurate readings.

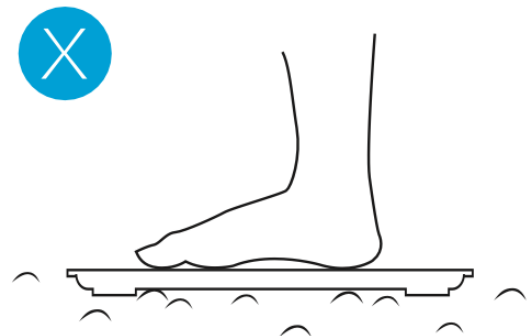App Guides

Download the RENPHO Health App

To start tracking your health metrics with your scale, you first need to download the RENPHO Health App. You can do this by following one of the methods below:

- Search for "RENPHO Health" in either the Apple App Store or Google Play Store.
- Alternatively, you can scan the QR code provided in your manual or packaging for direct access to the app.

System Requirements:

- For **iOS**: Requires iOS 12.0 or later.
- For **Android**: Requires Android 7.0 or later.

Note: Due to continuous updates and improvements, the RENPHO Health App interface might appear slightly different from the screenshots shown.

Log In or Sign Up

Once you've downloaded the app, follow these steps to log in or create an account:

- **Log in:**
 - If you already have a RENPHO Health account, simply enter your credentials to log in.
- **Sign up:**
 - If you don't have an account yet, tap Sign Up to register a new account using your email address. Complete your profile information to start tracking your health data.

- **Guest Mode:**

- You can also choose Guest Mode if you prefer not to create an account. This allows you to use the scale without registering. However, keep in mind that data will not be saved if you log in from another device.

Note: The Guest Mode does not support Wi-Fi pairing, and your measurements will not be saved across devices.

Add Device (Pairing the Scale)


Method 1: Quick Pairing

- **Enable Bluetooth** on your smartphone.
- **Step onto the Scale:** The scale will automatically detect your presence and a pop-up window will appear on your smartphone.
- **Tap 'Pair'** to connect your smartphone with the Body Fat Scale.

Method 2: Manual Pairing

- Tap the “+” icon in the top-right corner of the Home Page.
- Tap Body Fat Scale from the list of available devices.
- Follow the In-App Instructions: Step on the scale as prompted to initiate the connection process.
- Tap “+” to connect the scale to your device.
 - Once the pairing is complete, the Bluetooth icon on the scale’s display will turn solid white, confirming the successful connection.

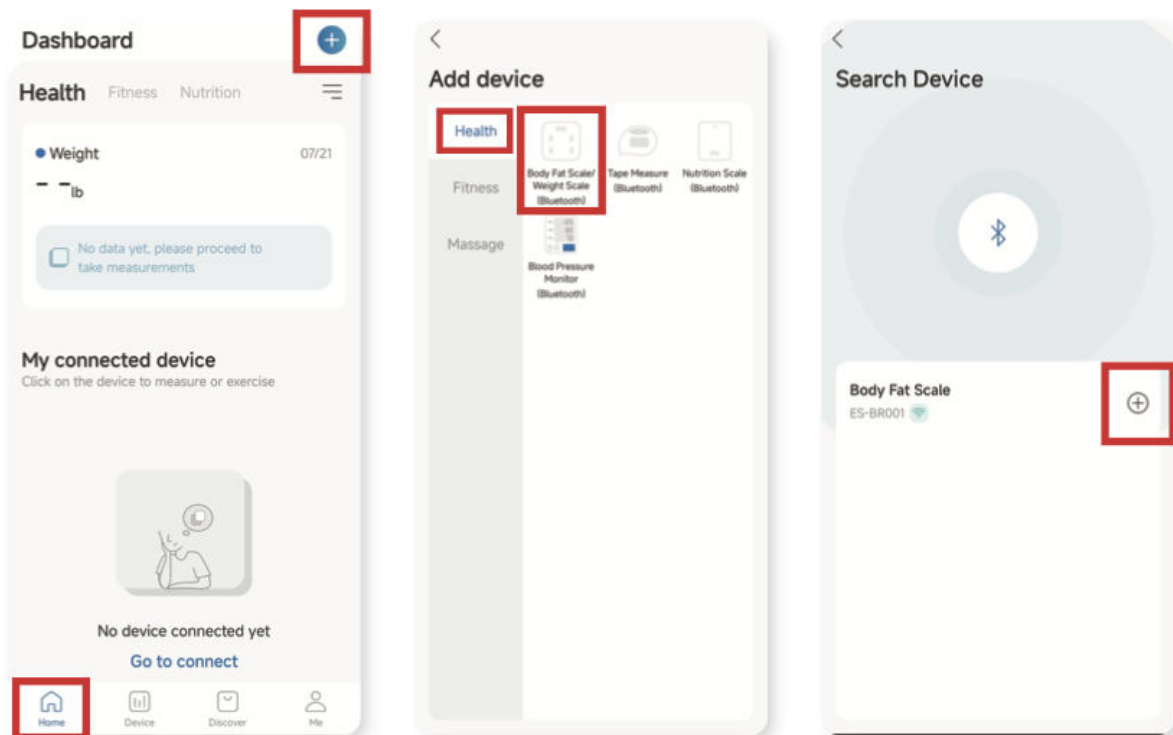
Note on First-Time Pairing: When pairing the scale with your smartphone for the first time, you must use the Bluetooth connection to complete the initial measurement. After successfully completing the first measurement, you can then set up the scale to use Wi-Fi for future measurements.

Pairing the Scale with Your Phone via Wi-Fi

Once Bluetooth pairing is complete, you can connect your scale via Wi-Fi to enhance convenience, allowing measurements without needing the app open.

1. Go to the “Device Page” on the app.
2. Tap the  Wi-Fi icon to enter the “Select Wi-Fi Page”.
3. Choose the appropriate 2.4 GHz Wi-Fi network and enter the correct Wi-Fi password.
 - (Note: 5 GHz Wi-Fi is not supported by the scale).
4. Step onto the Scale and remain standing until the Wi-Fi icon on the scale’s display turns solid white, indicating a successful connection.

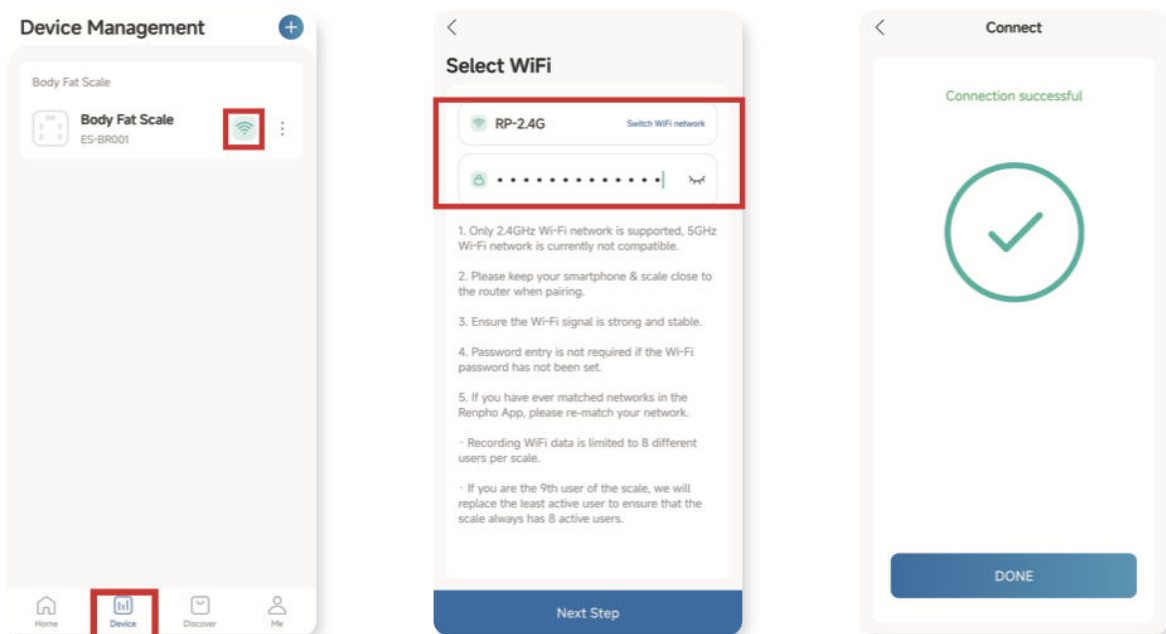
Note: Wi-Fi mode enables the scale to display your weight, BMI, and body fat percentage directly on the scale’s screen, without needing to run the app.

Calibration Instructions

If the scale has been moved, turned upside down, or hasn't been used for an extended period, it must be recalibrated to ensure accurate measurements. Follow these steps to calibrate your scale:

- **Place the Scale on a Stable Surface:**
 - Ensure that the scale is on a flat and hard surface such as tile or wooden flooring. Avoid placing it on soft surfaces like carpets or rugs, as this can lead to inaccurate readings.
- **Wake the Scale:**
 - Gently place one foot on the scale's surface to wake it up from sleep mode.
- **Remove Your Foot:**
 - Once the display shows digits or symbols, step off the scale. This step allows the calibration process to initiate.
- **Wait for Calibration to Complete:**
 - The display will first show "CAL," then it will shift to either "0.0" or "0.00." This indicates that the calibration process has been successfully completed.

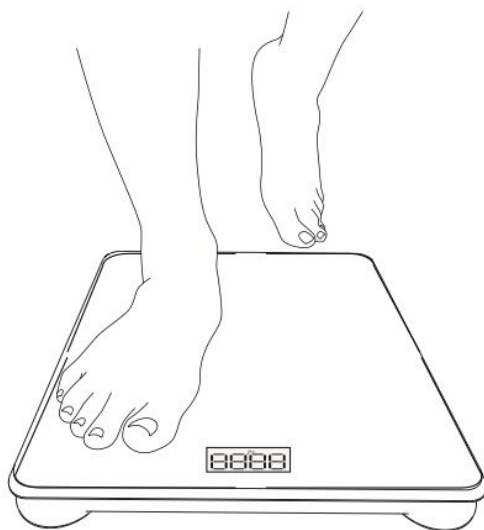CAL → 0.0

or

CAL → 0.00

For more detailed information, you can scan the QR code on your scale's packaging or visit our website for additional support and detailed instructions.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

This product complies with the legal requirements of the European Union directives.

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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