
RENPHO

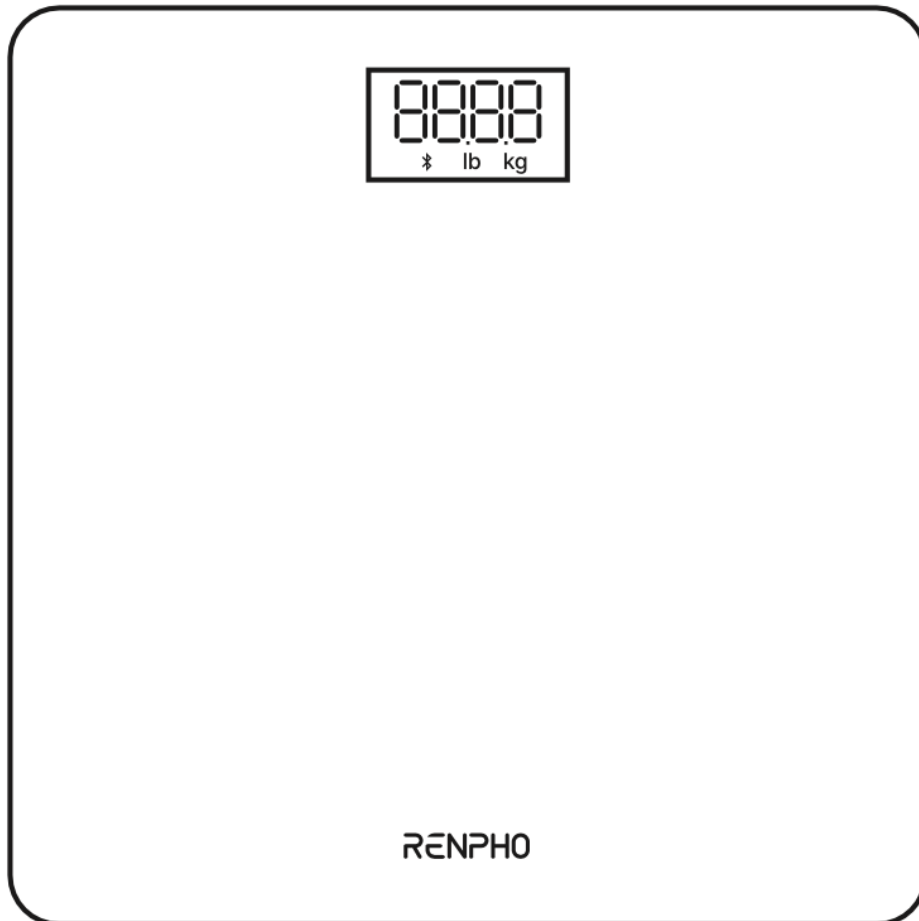Body Composition Scale

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Important Instructions

Read all instructions carefully before using this appliance.

- **Household Use Only:** This scale is designed exclusively for household use and is not intended for medical or commercial applications.
- **Restricted Use:** The scale should not be used by individuals (including children) with reduced physical, sensory, or mental capabilities, or by those who lack the necessary experience and knowledge, unless they are supervised or have received instructions from a responsible person regarding the safe use of the appliance.
- **Supervise Children:** Children should be monitored to ensure they do not play with the scale, as improper use could lead to injury or damage.
- **Reference Only:** The body composition data provided by the scale is intended for general reference and should not replace professional medical advice. Consult a licensed healthcare provider for accurate health guidance.
- **Maximum Capacity:** Do not exceed the maximum load limit of 400 lb (180 kg). Overloading the scale may cause permanent damage and result in inaccurate measurements.
- **Check for Damage:** Before each use, inspect the scale for any signs of damage. Do not use the scale if it is cracked or damaged, as it may provide incorrect readings or cause injury.
- **Battery Check:** If the scale malfunctions, check the battery power first. Replace the batteries if necessary to restore proper function.
- **Proper Surface for Use:** For accurate results, always place the scale on a hard, flat, and dry surface. Using the scale on uneven or soft surfaces can lead to inaccurate measurements.
- **Avoid Interference:** To prevent interference or potential danger, do not use the scale in environments with strong electromagnetic fields, such as near electronic equipment or appliances.
- **Use with Care:** To avoid injury, do not stand on the edges of the scale, jump on it, or apply uneven pressure. Always use the scale as intended to prevent damage.

- **Handle with Care:** The scale is a sensitive and precise measuring device. Handle it with care to avoid damage or breakage.
- **Water and Humidity:** The scale is not waterproof. Keep it dry and never immerse it in water. To clean the surface, use a damp cloth or glass cleaner. Avoid harsh chemicals such as soap or abrasive cleaners.
- **Storage Conditions:** Do not expose the scale to direct sunlight, extreme heat, or high humidity. Store it in a cool, dry place to prevent damage to the electronic components and batteries.
- **No Disassembly:** Do not attempt to disassemble the scale or tamper with its internal components, as doing so could void the warranty and potentially cause malfunction.

About the Battery


- **Temperature Considerations:** Do not expose the scale or its batteries to extreme temperatures. Excessive heat or cold can cause battery failure or malfunction.
- **Battery Replacement:** Replace all batteries at the same time with new, matching batteries to ensure consistent performance.
- **Battery Mixing:** Never mix different types of batteries or old and new batteries, as this can affect performance and cause leakage.
- **Proper Disposal:** Dispose of exhausted batteries according to local regulations. Do not dispose of them in regular trash as they are considered hazardous waste.
- **Storage:** If you plan to store the scale for an extended period without use, remove the batteries to prevent leakage or damage to the internal components.
- **Battery Condition:** Do not use modified, damaged, or compromised batteries, as they can present a risk of short-circuiting or failure.
- **Short Circuit Risk:** Be cautious of the risk of short-circuiting the battery terminals with metal objects such as jewelry or keys, which could cause the battery to overheat or fail.
- **Non-Rechargeable Batteries:** The batteries provided with the scale are non-rechargeable. Do not attempt to recharge them.
- **Battery Installation:** Always insert batteries according to the correct polarity (+ and -), as marked on both the battery and the scale.

- **Battery Safety:** Avoid exposing batteries to excessive heat, such as fire or direct sunlight, and do not attempt to deform or dismantle them.
- **Child Safety:** Keep batteries out of reach of children. Do not allow children to handle or replace batteries without adult supervision.
- **No Charging:** Do not attempt to charge non-rechargeable batteries.
- **Battery Disposal:** Never dispose of batteries in fire, as they can explode.
- **Supervision:** Ensure that children do not attempt to replace or tamper with the batteries.
- **Storage of Unused Batteries:** Store unused batteries in their original packaging and away from metal objects. If the packaging is opened, do not mix different battery types together.

Installing the Batteries

Before using the scale, you need to install the provided batteries. Follow these steps for proper installation:

- Insert the batteries into the battery compartment, ensuring that you align them according to the correct polarity. The positive (+) and negative (-) ends should correspond to the markings on both the battery and the scale.

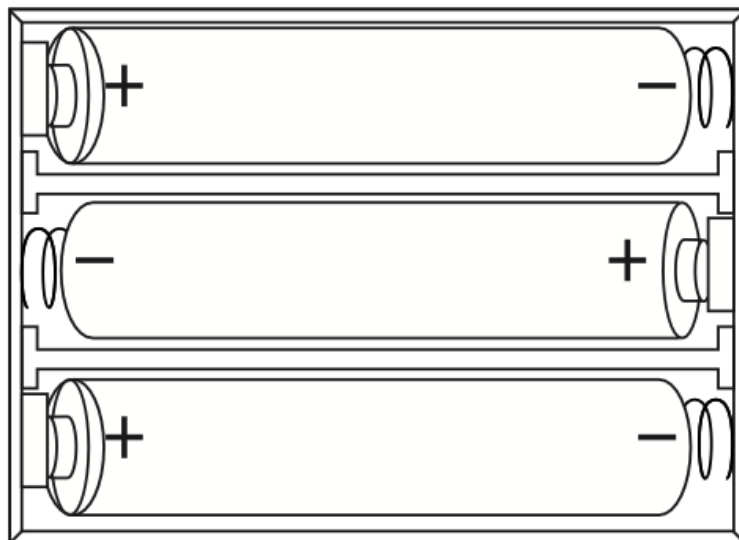
Note:

- Always use the appropriate size and type of batteries for the scale as recommended by the manufacturer. Never mix old and new batteries or combine batteries of different types or brands, as this could lead to performance issues or potential damage.
- Do not attempt to charge non-rechargeable batteries. Charging non-rechargeable batteries can cause leakage, overheating, or even explosion.
- Before inserting the batteries, clean the contacts on both the batteries and the scale to ensure a proper connection.
- If the display shows "Lo," it indicates low battery power. In this case, replace all batteries at once to avoid uneven power distribution, ensuring a longer battery life.

Preparing the Scale for Use

To ensure accurate measurements, always place the scale on a flat, hard surface. Avoid placing the scale on soft surfaces like carpets, rugs, or uneven flooring, as this can affect the accuracy of the results.

- **Correct Placement:** Use the scale on a hard, flat surface like tile or wood flooring to ensure the best results.
- **Incorrect Placement:** Avoid using the scale on soft or uneven surfaces, such as carpets or rugs, which can interfere with the scale's sensors and lead to inaccurate readings.

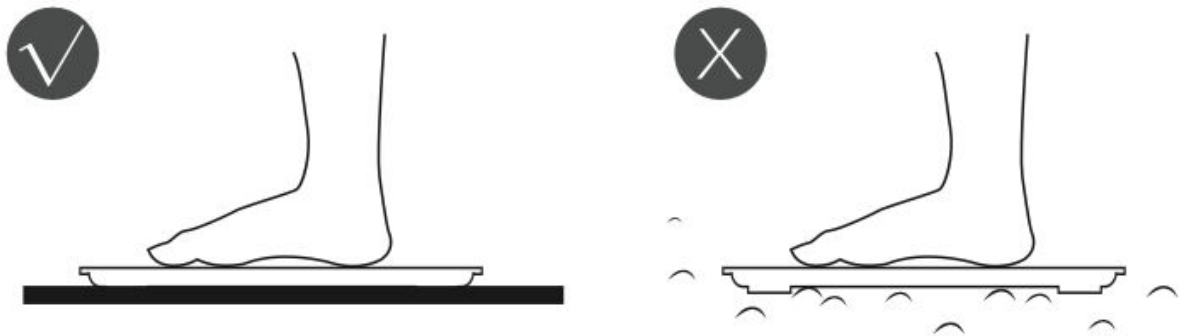Downloading the RENPHO Health App

To fully utilize your RENPHO scale, download the "RENPHO Health" app on your smartphone. This app allows you to track your body composition data, such as weight, body fat percentage, and more.

- You can download the app by scanning the provided QR code with your smartphone's camera, or by searching for "RENPHO Health" in your app store (available for both iOS and Android devices).
- The app may update periodically, so the interface could appear slightly different from the images shown.

Pairing the Scale with the RENPHO Health App

Once the app is installed, follow these instructions to pair the scale with the app:

- **Register Your Account:** Open the RENPHO Health app and register using your email address. If you prefer, you can log in using "Guest Mode," which allows you to skip the registration process, though you won't be able to save data across multiple sessions.
- **Enable Bluetooth on Your Smartphone:** Ensure that Bluetooth is turned on for a successful connection between the scale and the app.
 - For iOS users (iOS 12.0 or later): Go to your phone's settings, locate the RENPHO Health app, and make sure Bluetooth permissions are enabled.
 - For Android users (Android 7.0 or later): In addition to Bluetooth, you must enable location services for the scale to connect with your phone. You can do this by turning on location permissions in your device settings.

Note: Your phone's Bluetooth settings cannot directly connect with the scale. Always go through the RENPHO Health app to pair the device successfully.

Methods to Pair the RENPHO Scale with the App

Method 1: Quick Pairing

- **Step onto the Scale:** With bare and dry feet, step onto the scale and wait for the display to illuminate. This indicates that the scale is ready to pair.


- **Pair with the App:** Once the scale is active, a prompt should appear within the app. Tap "Pair" to complete the connection.

Method 2: Manual Pairing


- **Add Device:** On the main screen of the RENPHO Health app, tap the "Add Device" option located in the top right corner of the "Healthy" page.

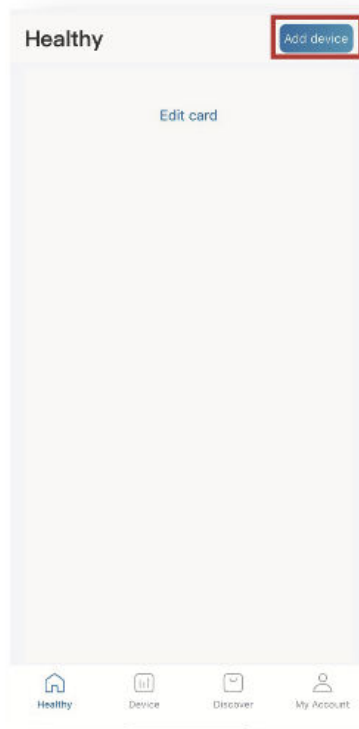
- **Select Scale Type:** In the list of available devices, select "Body Fat Scale" to initiate the pairing process.


- **Stand on the Scale Again:** With bare and dry feet, step onto the scale to ensure the app can fully detect your weight and initiate the connection. This will trigger the scale to communicate with the app for accurate pairing.

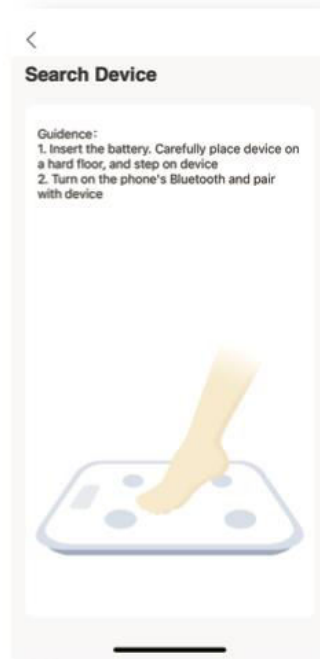- **Add the Scale:** Once the scale is detected by the app, its name will appear on your phone's screen (e.g., "body fat scale"). Tap the "+" button to confirm and add the device to your list of connected devices.

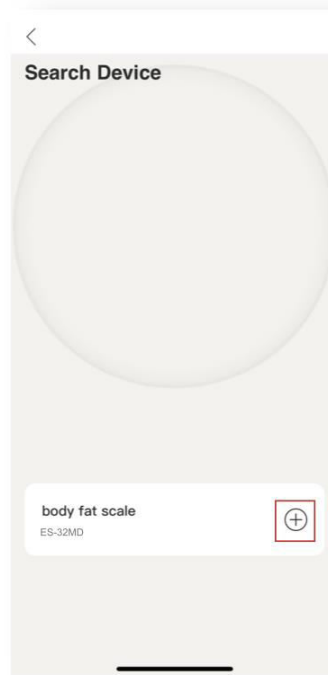
- **Wait for Pairing to Complete:** Remain standing on the scale for approximately 5 to 10 seconds. The pairing is complete when your weight and body metrics appear on the RENPHO Health App's display. The app will automatically show various metrics, including weight, body fat percentage, BMI, and other health data. These readings will now sync every time you step on the scale.

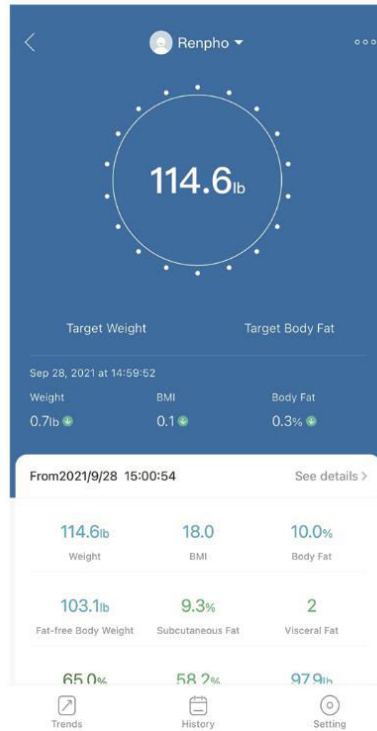Calibration Instructions

If the scale has been moved, turned upside down, or hasn't been used for an extended period, it must be recalibrated to ensure accurate measurements. Follow these steps to calibrate your scale:

- **Place the Scale on a Stable Surface:**
 - Ensure that the scale is on a flat and hard surface such as tile or wooden flooring. Avoid placing it on soft surfaces like carpets or rugs, as this can lead to inaccurate readings.
- **Wake the Scale:**
 - Gently place one foot on the scale's surface to wake it up from sleep mode.
- **Remove Your Foot:**
 - Once the display shows digits or symbols, step off the scale. This step allows the calibration process to initiate.
- **Wait for Calibration to Complete:**
 - The display will first show "CAL," then it will shift to either "0.0" or "0.00." This indicates that the calibration process has been successfully completed.

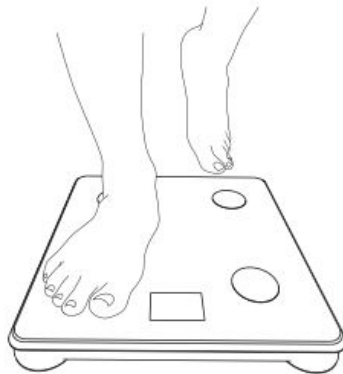CAL → 00

or

CAL → 000

For more detailed information, you can scan the QR code on your scale's packaging or visit our website for additional support and detailed instructions.

Frequently Asked Questions

Why isn't the scale detecting body fat and/or other metrics?

If the scale isn't measuring body composition data, such as body fat percentage, follow these troubleshooting steps:

- **Ensure Dry, Bare Feet:** Make sure your feet are both dry and bare for proper conductivity between the scale's sensors and your body.
- **Correct Stance:** Stand firmly on the scale with feet slightly apart, and ensure both feet are making full contact with the scale's sensors. Face the digital LED display to view the readings clearly.
- **Wait for Measurement:** Allow the scale about 10 seconds to complete the weighing and measurement process. Be patient as the sensors work to gather data.
- **Reapply Moisturizer:** If you have particularly dry skin, applying a small amount of moisturizer to your feet may help improve the conductivity and sensor contact. Try stepping on the scale again after applying moisturizer.
- **Try Another Person:** Have another person step onto the scale to see if the body fat measurement works for them. This can help confirm if the issue lies with the sensors or with individual variations.

I have registered an account on the RENPHO website, but I can't log into the RENPHO Health App. Why?

The RENPHO website and the RENPHO Health App operate on two separate systems. While you may have created an account on the website, you will need to create a separate account for the app. If you encounter Bluetooth connectivity issues, follow these steps:

- **Reset the Scale:**
 - Remove one battery from the scale.
 - Wait for 5 seconds, then reinstall the battery.
- **Reset Bluetooth Connection:**
 - Ensure that both Bluetooth and GPS (for Android devices) are enabled on your smartphone.
 - **Open the RENPHO Health App** and follow these steps:
 - Step one foot on the scale to wake it up.
 - After a few seconds, the scale's model number will appear on the app's screen.
 - Tap the scale's model number to complete the connection.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.


EU Declaration of Conformity

This product complies with the legal requirements of the European Union directives.

WEEE

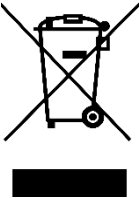This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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