RENPHO

Handheld Massager

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt


+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Diagram

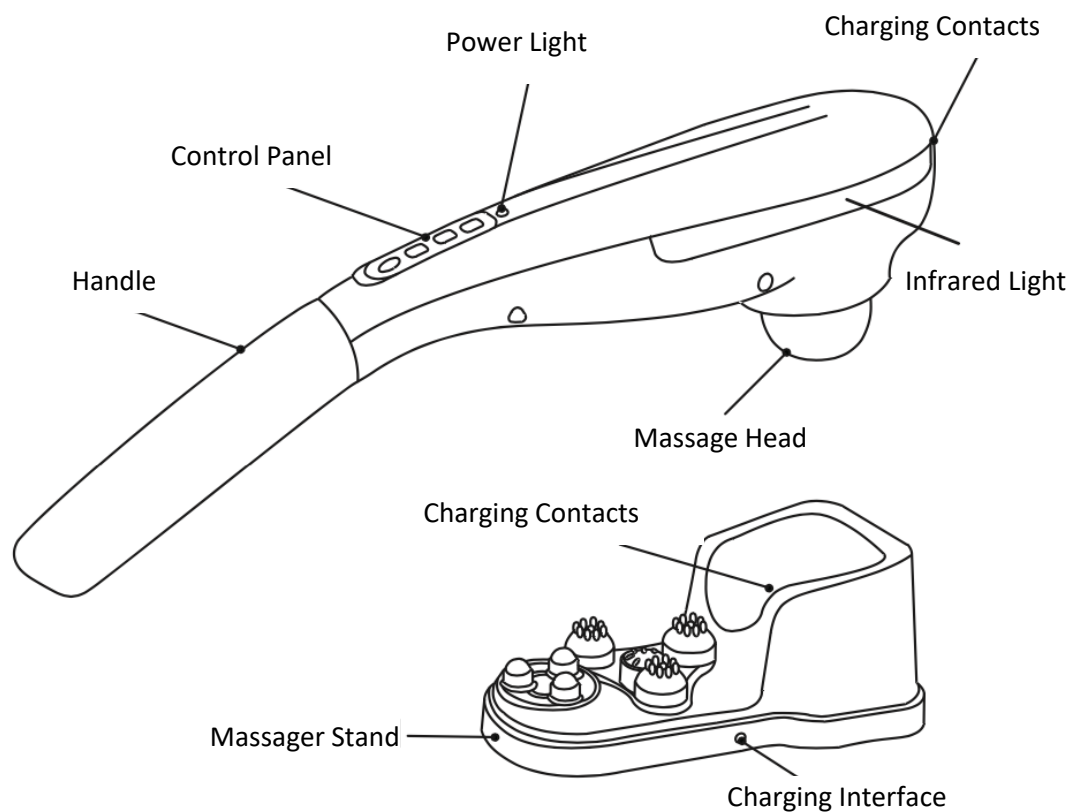Charging and Battery Life

Charging Instructions

Before using the massager for the first time, it is essential to charge the device fully to ensure optimal performance.

- **Important Safety Note:** To prevent any potential hazards, the massager is designed to not function while it is connected to the charger. Therefore, always unplug the charger before using the massager. Additionally, the massager is equipped with an automatic safety feature that will shut the device down after 20 minutes of continuous use to prevent overheating and prolong the life of the unit.

- **Charger Specifications:**


Specification	Details
Charger Input Voltage	100 V – 240 V, 50 / 60 Hz, 0.5 A
Charger Output Voltage	9 V = 1A
Rated Power	9 W
Battery Capacity	2600 mAh (19.24 Wh)
Battery Voltage	7.4 V

Please ensure you use the included charger or a compatible one with the same specifications to avoid damage to the device.

Steps for Charging the Massager

Follow these steps carefully to charge your massager:


- Step 1: Connect the charger to an available power outlet. Once plugged in, the indicator light on the charger will turn green, signaling that the power supply is connected correctly.

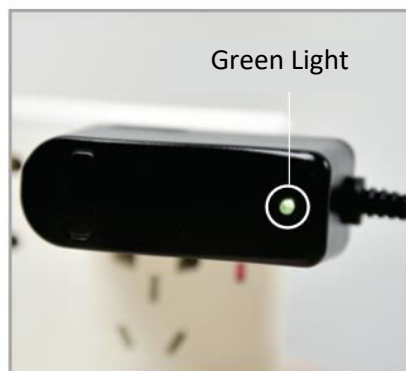- Step 2: Insert the charger plug into the massager's designated charging port. Ensure the plug is firmly seated for proper connection.

- Step 3: Place the massager on the provided stand so that the charging contacts are aligned and touching. When the contacts are properly aligned, the indicator light on the charger will turn red, confirming that the device is now charging.

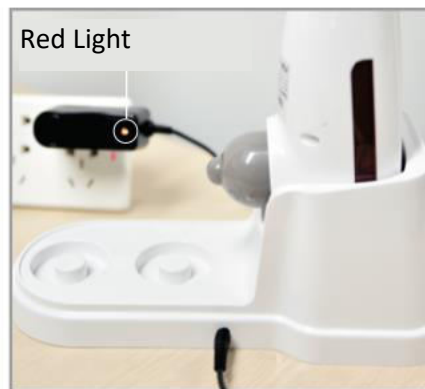- Step 4: Once charging is complete, the indicator light on the charger will switch back to green, signaling that the battery is fully charged.

- Step 5: To prolong battery life and maintain performance, it is recommended that you unplug the charger once charging is finished.

Battery Life and Maintenance

The massager is powered by two 18650 lithium-ion battery cells with a total capacity of 2600mAh. A full charge typically requires about 3 hours. After fully charging the unit, you can expect up to 140 minutes of continuous operation, depending on usage.

Battery life can vary based on the settings you use and the conditions under which the massager is operated. To ensure the longevity of the battery:

- If you do not plan to use the massager for an extended period, turn off the device and store it in a cool, dry place.
- Charge the massager at least once every month, even if it is not in use, to keep the batteries in good condition.

Important Usage Tips:

- Keep the massager away from high temperatures and avoid using it in environments where the surrounding temperature is excessively warm. Extreme heat can reduce battery efficiency and may cause damage to the device.
- Avoid overcharging the device by unplugging it once it reaches full charge, as this can help preserve the battery's shelf life.

Operating Instructions

To ensure the best experience with your massager, please follow these step-by-step instructions carefully:


- **Charge the Device:** Before using the massager, ensure that it is fully charged. This will guarantee optimal performance and battery life.
- **Disconnect from Charger:** Once the indicator light on the charger turns green, indicating a full charge, unplug the charger from both the device and the power outlet. This ensures that the massager is ready for operation and prevents unnecessary power consumption.
- **Attach a Massage Head (Optional):** If desired, you can select a massage stick or attachment suited to your massage needs. To remove an existing massage stick, twist it counterclockwise. To attach a new one, twist clockwise until it is securely in place. Make sure the attachment is fastened properly to avoid any interruptions during use.
- **Power On the Device:** Press the Power On / Off button to activate the massager. Once powered on, the device's power indicator light will illuminate, signaling that the unit is operational and ready to use.

- **Select a Mode:**

- Press the Mode button to choose between different massage modes, each offering a unique type of massage rhythm and intensity.
- To adjust the intensity of the massage, use the up (Λ) or down (V) buttons to either increase or decrease the power of the percussion.
- There are four available percussion modes, each offering five different speed levels. Experiment with different modes and speed levels to find the combination that best suits your preferences or therapy needs.

The control panel features the following buttons:

- **Mode:** Allows you to switch between massage modes.
- **Up Arrow (Λ):** Increases the intensity of the massage.
- **Down Arrow (V):** Reduces the intensity.
- **Power On/Off:** Turns the device on or off.

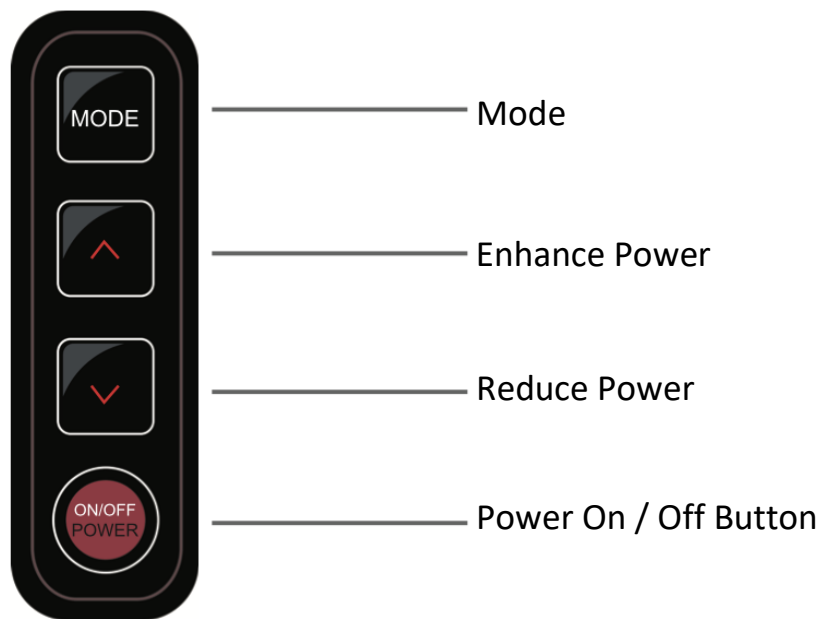
Targeted Massage Areas

This massager is designed for versatile use and can be applied to various areas of the body, providing relief and relaxation. Below are some of the recommended target areas and their corresponding illustrations:


- **Shoulder Massage:** Relieves tension in the shoulders, often a result of poor posture or stress.

- **Back Massage:** Eases muscle tightness in the upper or lower back.

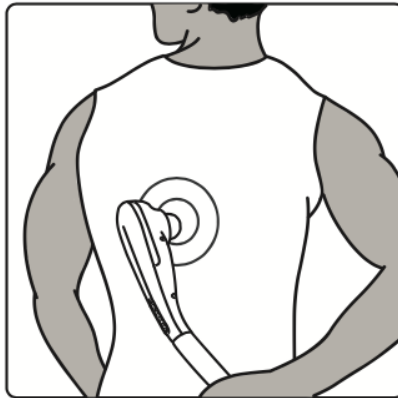
- **Foot Massage:** Provides soothing relief after a long day on your feet.

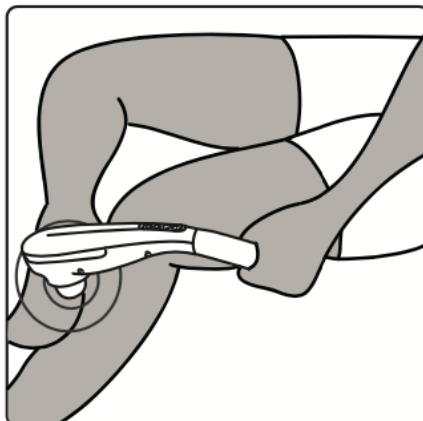
- **Arm Massage:** Targets muscle strain or soreness in the arms.

- **Waist Massage:** Helps to alleviate stiffness or discomfort in the waist region.

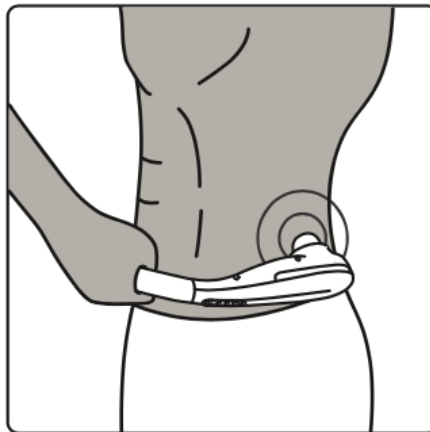- **Leg Massage:** Relaxes sore leg muscles, perfect after workouts or long periods of standing.

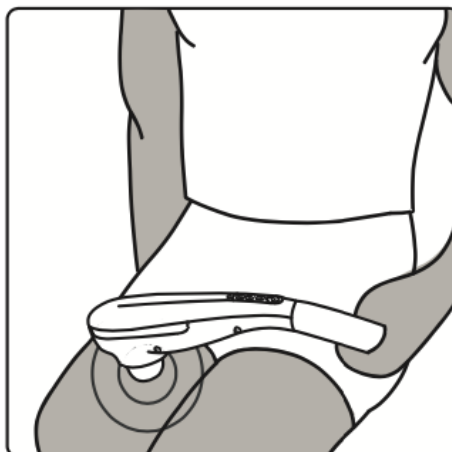Safety Feature: Automatic Power Off

For your safety, the massager is equipped with a protective feature that automatically shuts the device down if too much pressure is applied. This prevents damage to both the massager and your body.

If the device powers off due to excessive pressure, follow these steps to reset it:

- Turn the device off if it has not already powered down.
- Place the massager into the charging stand and connect the charger.
- Wait for approximately 3 seconds.
- Unplug the charger from the massager and the power outlet simultaneously.
- Turn the device back on and resume your massage session.

Massage Sticks

Your full-body massager comes equipped with five unique massage sticks, each designed to provide a different style of massage. Depending on whether you prefer a firmer, deeper massage or a softer, more gentle approach, you can select the attachment that best suits your needs.

Below is a guide to help you choose the most appropriate massage stick for different parts of the body and massage techniques:

- **Round Stick:**
 - The Round Stick is designed to deliver a powerful, deep-tissue massage. It is ideal for targeting larger muscle groups where you may need stronger pressure to relieve tension and soreness. This stick provides broad, firm contact, making it effective for intense massage sessions.

- **3-Points Stick (or 3-Head Stick):**

- The 3-Points Stick is perfect for improving blood circulation and relieving muscle stiffness. Its three rounded contact points allow for a focused massage, making it ideal for areas where there is a significant amount of muscle over bone, such as the shoulders and upper back. The design helps to evenly distribute pressure, enhancing relaxation and promoting muscle recovery.

- **Curved Head:**

- The Curved Head is shaped to comfortably contour to the arms and legs, providing a more targeted massage for those regions. It's great for releasing muscle tension in the limbs and can be particularly useful for athletes or those with frequent muscle strain. Its curved shape makes it perfect for applying even pressure along the curves of your body.

- **One-Point Head:**

- The One-Point Head is intended for more precise application, such as massaging acupuncture points. This stick delivers focused pressure to small, specific areas, making it ideal for relieving tension in tight spots. Use this attachment carefully, especially when applying it to sensitive areas.

Important Safety Note:

None of the massage sticks are designed for use on the head. Please avoid using any of the attachments on your head or face to prevent injury or discomfort.

Warnings & Safety Precautions

Important – Please Read All Instructions Before Use

- **Intended Use:** This appliance is not designed for use by individuals (including children) with reduced physical, sensory, or mental capabilities, or those who lack experience and knowledge, unless they have been given proper supervision or instruction by a person responsible for their safety.
- **Supervision of Children:** Children should be closely supervised to ensure they do not play with the massager. It is not a toy and can pose a risk if mishandled.
- **Maintenance:** Cleaning and user maintenance should not be performed by children without supervision to ensure proper care and avoid accidents.

Safety Notice: To prevent serious injury, death, or damage to property, read this manual thoroughly and use the massager only as instructed. Improper use can result in hazardous situations.

Usage Restrictions

- **Age Recommendation:** Due to the high power of this device, it is recommended for use only by individuals aged 16 and above. Keep the massager out of the reach of children and pets to avoid accidental injury.
- **DO NOT Use the Massager On:**
 - **Sensitive or Damaged Skin:** Avoid using the massager on areas with sensitive, swollen, or inflamed skin. Also, refrain from using it where skin infections, rashes, or eruptions are present, or where you experience unexplained abdominal pain, as this could worsen the condition.
 - **Head:** The massager is not designed for head use. Applying pressure to the head can be dangerous.
 - **Varicose Veins:** Do not use the massager on areas with varicose veins, as this could exacerbate the condition.
 - **Throat Area:** The throat is a sensitive area, and using the massager here could pose a risk to your airway and vital structures.
 - **Genital Area:** This is not a personal or sexual massager. Using it on the genital area can cause severe bodily harm due to its strong power.
 - **Unconscious or Sleeping Individuals:** Never use the massager on someone who is unconscious or sleeping, as they cannot provide feedback on pressure or discomfort.
 - **Pregnant Women:** The powerful vibrations and intensity of this massager may not be safe for use during pregnancy. Always consult a healthcare provider before using the device if you are pregnant.
- **Use Only Original Accessories:** Always use the original charger and massage sticks provided by RENPHO. Using third-party chargers, adapters, or massage sticks can result in overheating, explosion, electric shock, or fire. Additionally, non-compatible attachments may cause damage to the massager and void the warranty.
- **Avoid Wet Conditions:** Never operate the massager with wet hands, and do not allow moisture to come into contact with the charger or the device itself. Water can cause severe damage to the massager and may increase the risk of electric shock. Do not insert foreign objects into the charging port, as this can lead to electric hazards.

- **Keep Away from Heat Sources:** Do not place the massager near heat-emitting appliances, such as ovens, microwaves, or high-pressure cookers. Avoid storing it in overly hot environments, such as inside a parked car on a hot day or near windows with direct sunlight. Excessive heat can damage the internal batteries, degrade the product's performance, and cause malfunction.
- **No Disassembly:** Under no circumstances should you attempt to disassemble or modify the massager or its batteries. Drilling into the device or tampering with its components can lead to electric shock, fire, or irreparable damage. Altering the product will also void your warranty, and you may be held responsible for any resulting damages.
- **Temperature Guidelines:** For optimal performance and safety, the massager should only be used in environments where the temperature ranges between 5°C (41°F) and 35°C (95°F). Avoid exposing the device to extreme temperature fluctuations or high humidity, as this may cause damage to the electrical circuits or batteries.
- **Handle with Care:** Do not drop or apply excessive force to the massager. Applying heavy weight or impact can cause internal damage, malfunctions, and reduce the lifespan of the device. Always store and handle the massager with care to avoid unnecessary wear and tear.
- **Intended Use Only:** The massager is intended solely for personal, non-commercial use. Using it for any purpose other than that described in the manual could result in malfunction and premature product failure. This product is not designed for continuous industrial or commercial operation.
- **Magnetic Fields:** Keep the massager away from areas with strong magnetic fields. Prolonged exposure to magnetic fields can affect the device's functionality and cause battery discharge or other issues.
- **Warning Against Overheating:** Do not operate the massager under blankets, pillows, or any enclosed surface where heat cannot dissipate properly. Excessive heat buildup may result in fire, electric shock, or injury to the user.

Additional Warnings

- **Supervision Required:** Close supervision is necessary when the massager is being used near children, individuals with disabilities, or those who require assistance. The massager should not be used by anyone unable to fully understand its operation or risks.
- **No Wet Environments:** Do not use the massager while bathing, showering, or near any water source. The risk of electric shock is significantly increased in wet environments. Never allow the massager to come into contact with water or any other liquid.
- **Ventilation & Placement:** Do not operate the massager on soft surfaces (such as beds or couches) where the ventilation openings could become blocked. This could cause the device to overheat and result in malfunction or fire.
- **Original Power Supply Only:** Always use the power supply unit that came with the massager. Using a non-recommended power source can lead to electrical issues, short-circuiting, or potential injury.

Battery Safety Information

- **Battery Handling:** Be mindful of the risk of short-circuiting the terminals of the battery-operated device. Metal objects such as keys, coins, or jewelry can cause a short circuit if they come into contact with the battery terminals.
- **Non-Replaceable Battery:** This massager contains non-replaceable batteries. If the batteries reach the end of their life, the device should be properly disposed of following local guidelines for electronic waste. The batteries should only be removed or replaced by qualified personnel.
- **Temperature Considerations:** Avoid exposing the massager and its batteries to extreme temperatures. Excessive heat or cold can damage the batteries, reducing their efficiency or causing them to fail altogether.

Product Care & Maintenance

- **Storage:** Always store the massager in a safe, cool, and dry place when it is not in use. This will help preserve its functionality and prolong its lifespan.
- **Unplug Before Cleaning:** Before performing any cleaning or maintenance, make sure the massager is unplugged from the charger to avoid electrical hazards.
- **Cleaning:** After each use, gently wipe the massager with a soft, damp cloth. It is important to use only non-abrasive cleaners to avoid scratching or damaging the surface of the device. If necessary, you may also sterilize the massager with appropriate sterilizing products to ensure cleanliness, but avoid harsh chemicals.
- **Avoid Water Submersion:** Never submerge the massager in water or any other liquid. Doing so may cause internal damage to the device's electrical components, leading to malfunction or safety hazards.
- **Keep Away from Harsh Chemicals:** Avoid exposing the massager to solvents, bleach, or harsh detergents. These substances can degrade the materials and affect the massager's performance.
- **No User Repairs:** Do not attempt to repair the massager yourself. There are no user-serviceable parts inside the device, and opening it could void your warranty. If the product malfunctions, please contact customer service for assistance or repairs.
- **Proper Disposal:** The massager contains a lithium-ion battery and other electronic components that are classified as hazardous E-waste. Do not dispose of the device in a landfill. Instead, follow local recycling programs or E-waste disposal guidelines to ensure the unit is properly and safely discarded.

Troubleshooting

Issue 1: The Massager Won't Turn On

If your massager is not powering on, try the following steps:

- **Check the Charge:** Ensure that the massager is fully charged by connecting it to the charger and verifying that the indicator light turns on.
- **Reset the Massager:**
 - Turn off the massager completely.
 - Place it in the charging stand and connect the charger.
 - After waiting approximately 3 seconds, unplug both ends of the charger simultaneously (from the massager and the wall socket).
 - Press the Power button to turn the massager back on.
- **Battery Life:** If the massager still does not turn on after charging, it is possible that the lithium-ion battery has reached the end of its life cycle. In this case, you may need to replace the unit or consult customer support for further assistance.

Issue 2: The Massager Doesn't Start Immediately

- **Start-Up Delay:** The massager may take up to 3 seconds to start after pressing the power button. During this time, the motor is initializing.
- **Adjusting Speed/Intensity:** If the device is slow to start, press the Mode button to increase the speed or intensity settings. The Power Light on the massager should turn solid green when the device is ready for use.

Frequently Asked Questions (FAQ)

If you encounter issues with your massager, consult the following troubleshooting guide to resolve common problems.

Problem	Status	Solutions
No power on charger's indicator light.	The charger's indicator light does not turn on.	Try using a different power outlet. If the problem persists, the outlet may be faulty, or the charger itself may need replacement.
Charger's indicator light still not lighting up.	The indicator light on the charger fails to light up, even after trying different outlets.	The charger may be defective. Please contact the seller or support team for assistance with a replacement or repair.
Charger light turns from red to green too quickly or never turns green.	The charger's indicator light behaves unusually, changing too quickly or failing to change to green when fully charged.	This may indicate an issue with the battery. The massager could be defective. Please contact the seller or support for troubleshooting and possible replacement.
Massager shuts off automatically.	The massager powers down on its own during use.	The device's built-in overheat protection may have been triggered. Allow the massager to cool down for a few minutes, then try turning it back on. If the issue recurs, consider checking the environment or usage duration to ensure the device isn't being overheated.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

This product complies with the legal requirements of the European Union directives.

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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