RENPHO

Massage Gun

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

www.alza.co.uk/kontakt

+44 (0)203 514 4411

Importer

Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Important Safety Instructions

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions carefully before using.

Before using the device, it is crucial to consult your doctor if you fall into any of the following categories:

- **Pregnant or postpartum women**
Pregnancy or postpartum recovery may pose certain risks when using the device, so consulting with a healthcare provider is necessary.
- **Persons with blood stasis, fractures, or severe inflammatory injuries**
Conditions involving circulatory issues or severe inflammation can exacerbate with the device's stimulation.
- **Persons with hypertension (high blood pressure)**
Individuals with severe varicose veins or other blood vessel conditions should avoid using the device without medical consultation.
- **Persons with heart disease**
If you have a pacemaker, defibrillator, or metal implants, avoid using the device near the heart or chest area due to potential interference with the implanted devices.
- **Persons recovering from joint surgery**
After surgery, certain areas might be sensitive to pressure and vibrations, so the device should not be used on or near recent surgical incisions.

If you experience any discomfort or unusual sensations while using the device, discontinue use immediately and seek medical advice.

Caution

- The device is intended for household use only and is not designed for medical purposes. Do not use it to treat or diagnose any conditions.
- Children aged 8 years and above, as well as those with reduced physical, sensory, or cognitive capabilities, should only use the device under supervision and guidance from someone who understands the hazards.
- Children should not play with the device. Cleaning and maintenance should not be performed by children without adult supervision.

Danger – To reduce the risk of electric shock:

- DO NOT use the device while bathing or showering.
- DO NOT place or store the device in places where it may fall into water, such as bathtubs, sinks, or pools. Avoid dropping the device into water or other liquids.

Warning – To reduce the risk of fire, electric shock, or injury to persons:

- Always use the device for its intended purpose as described in this manual. Do not use attachments that are not recommended by the manufacturer, as these could cause damage or injury.
- Supervision is required when the device is used by, or near, children, individuals with disabilities, or persons with reduced cognitive abilities.
- Keep the air openings of the device clear. DO NOT block the air openings with hair, lint, or other objects. Never drop or insert items into the openings.
- DO NOT operate the device under blankets or pillows. Excessive heating can occur, leading to a fire hazard or risk of electric shock.
- DO NOT use the device in areas where aerosol spray products are in use or where oxygen is being administered.
- DO NOT keep the device near magnetic fields as it can interfere with the device's functionality or cause battery discharge.
- Always store the device away from extreme heat. Avoid using it in environments where the temperature is too high.


About the Battery

- Chemical Exposure: If the battery pack leaks, do not allow the liquid to come into contact with your skin or eyes. If exposed, wash the area immediately with water and seek medical advice.
- Battery life varies based on usage. If the device is not in use for an extended period, power it off and charge it once a month to maintain the battery's lifespan.
- DO NOT expose the battery or device to extreme temperatures, as this may cause malfunction.
- Short-circuit risks: Be cautious of metal objects that may come into contact with the battery's terminals, potentially causing a short circuit.
- This device contains non-replaceable batteries. When the battery reaches the end of its life, it must be properly disposed of. Do not throw the device into a fire or expose it to high heat. Place the device in a sealed bag and dispose of it according to local environmental regulations.

Packaging

- 1 × Massage Gun
- 5 × Massage Heads
- 1 × Type-C Charging Cable
- 2 × Rubber Sleeves
- 1 × Travel Case
- 1 × User Manual

Diagram

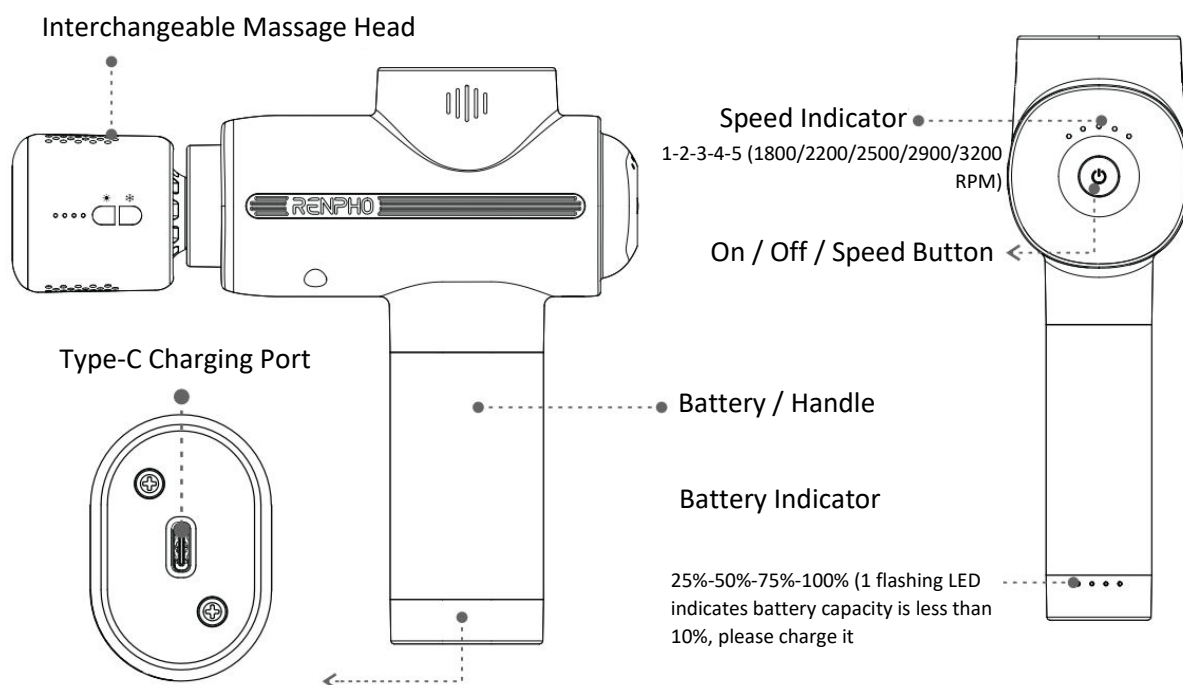Massage Head Attachments

- **Heating & Cooling Massage Head (Rechargeable):**

- Suitable for use on all body parts, providing both heating and cooling therapy for enhanced muscle recovery and relaxation.
- Note: This head has its own charging system.

- **Ball Head:**

- Ideal for large muscle groups such as the quadriceps, glutes, and muscles prone to extreme soreness.

- **Bullet Head:**

- Designed for precise deep tissue massages, targeting joints, trigger points, and smaller muscles like those in the wrist or feet.

- **Flat Head:**

- Best for working on the shoulder blades and IT bands, offering a relaxing, full-body massage experience.

- **U-shaped Head:**

- Specially designed for the spine, neck, and Achilles tendon, allowing safe and effective massage around these sensitive areas.

- **Rubber Sleeve:**

- Each massage head is pre-installed with a rubber sleeve that ensures a secure fit. If the sleeve becomes damaged, replace it with the included spare.

- **Warning:**

- Do NOT use massage heads on your head.


Charging Instructions

- **Before First Use:**

- Make sure to fully charge the massage gun and the heating & cooling massage head before using them for the first time.

- **How to Charge:**

- Plug the provided Type-C charging cable into the charging port on the device.
- Connect the other end to a power source such as a laptop, mobile charger, or power bank.
- The device's battery indicators will flash while charging, and the respective indicator lights will show charging progress.

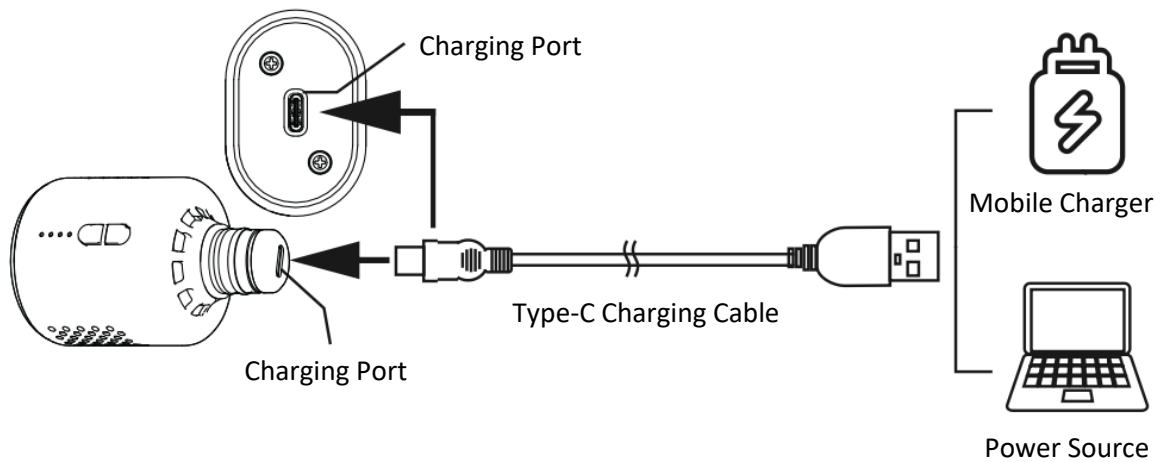
- **Charging Indicators:**
 - **Massage Gun:**
 - During charging: 4 battery indicators will flash BLUE in sequence.
 - When fully charged: 4 indicators will remain solid BLUE.
 - **Heating & Cooling Massage Head:**
 - During charging: 4 indicators will flash RED in sequence.
 - When fully charged: 4 indicators will remain solid BLUE.
- **Important Note:**
 - Do not use the device while it is charging to prevent any damage or safety hazards.

Attaching the Massage Head

- **Select the Desired Massage Head:**
 - Choose the appropriate massage head based on the type of massage you want (refer to the previous section on massage head types for guidance).
- **Insert the Massage Head:**
 - Align the joint of the massage head with the slot on the device.
 - Firmly press and fully insert the head until it clicks into place.

- **Ensure Safety:**


- Make sure the device is turned off before switching or attaching any massage heads to avoid injury or malfunction.

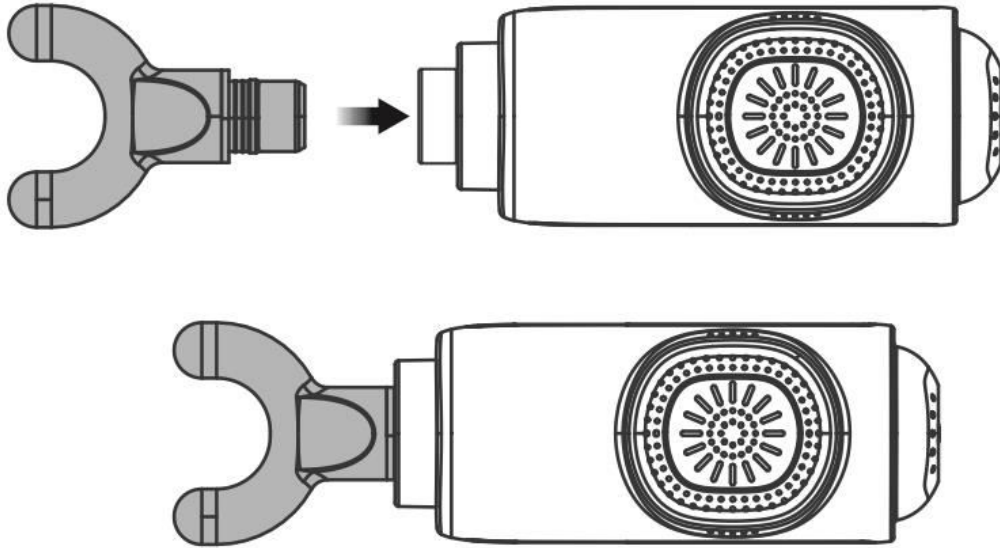Using the Product

Massage Gun

-  **Power On/Off:**
 - Long press to turn on the device. The speed indicator will activate.
 - Short press to cycle through the speed levels: 1-2-3-4-5.
 - Long press to turn off the device.

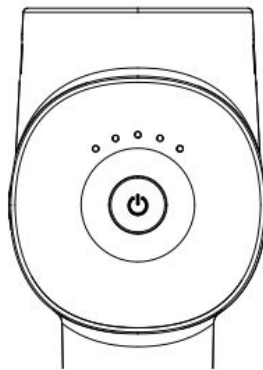
Heating & Cooling Massage Head

Heating Function:

- Long press to enable the heating function. The red indicator will turn on.
- Short press to cycle through the heating levels: 1-2-3-4.
- Long press to disable the heating function.
 - Heating Temperature: 99 – 113 °F / 37 – 45 °C.
 - Auto-off Timer: 10 minutes.

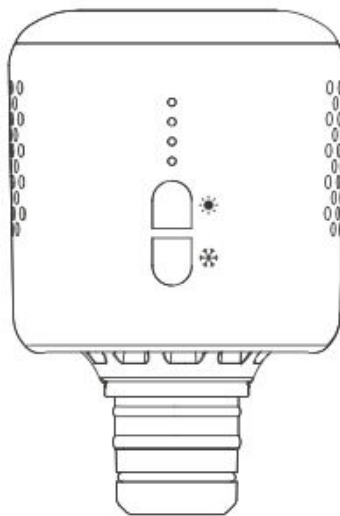- **Cooling Function:**

- Long press to enable the cooling function. The blue indicator will turn on.
- Short press to cycle through the cooling levels: 1-2-3-4.
- Long press to disable the cooling function.
 - Cooling Temperature: 47 – 61 °F / 8 – 16 °C.
 - Auto-off Timer: 5 minutes.

Notes:

- **Auto Shut-off Function:**

- The device and the heating function will automatically shut off after 10 minutes of continuous use.
- The cooling function will shut off automatically after 5 minutes of continuous use.

- **Locked-motor Protection:**
 - If excessive pressure is applied during a session, the device will stall to prevent injury or bruising.
- **Recommendation:**
 - Start each massage session with minimal intensity before increasing it. Do not use for more than 30 minutes at a time.

Cleaning and Maintenance

- **Step 1:** Before cleaning, ensure the device is turned off. This will help prevent any electrical issues during the cleaning process.
- **Step 2:** Use a soft towel or cloth to gently wipe the outer surface of the device. Avoid using any harsh chemicals or abrasive materials that could damage the exterior. If needed, a damp cloth with mild soap can be used, but make sure no water enters any openings.
- **Step 3:** Avoid disassembling or attempting to repair the device yourself, as it contains no user-serviceable parts. If any part of the device seems damaged or malfunctioning, contact the manufacturer for proper support.
- **Step 4:** When not in use, store the device in a cool, dry place. High heat or humid environments can degrade the device's materials and battery life over time.

Troubleshooting Guide

Problem	Cause	Solution
The device cannot be turned on	The battery is empty.	Charge the device using the provided Type-C charging cable.
The device stops working suddenly	Auto-off timer: The device automatically shuts off after 10 minutes of continuous use.	Simply restart the device and resume use if necessary.
	Built-in locked-motor protection: When the massage head is pressed too firmly, the motor protection activates to prevent damage.	Release pressure on the massage head and restart the device.
	Built-in overheat protection: The device will shut off if it becomes too hot during use.	Wait for the device to cool down before turning it back on.
	Built-in overload protection: When the input voltage exceeds the safe limit (6.5 -30 V), this protection will activate to prevent damage.	Make sure that the input voltage does not exceed 5 V.

Specifications

Specification	Massage Gun	Heating & Cooling Massage Head
Product Name	Massage Gun	
Model	R-C003H	
Input	5 V == 3 A	5 V == 1 A
Amplitude	10 mm	
RPM	1800 / 2200 / 2500 / 2900 / 3200 r / min	
Auto-off Timer	10mins	10mins (heating); 5mins (cooling)
Battery Voltage	DC 7.4 V	
Battery Capacity	2000 mAh 14.8 Wh	900 mAh
Charging Time	2 H	1 H
Working Time	2-3 H	1 H
Product Dimensions	7.3 × 5.7 × 2.3 inch / 186 × 147 × 60 mm	
Product Weight	1.5 lbs / 680 g	

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

This product complies with the legal requirements of the European Union directives.

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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