SIGURO

Baby Scale

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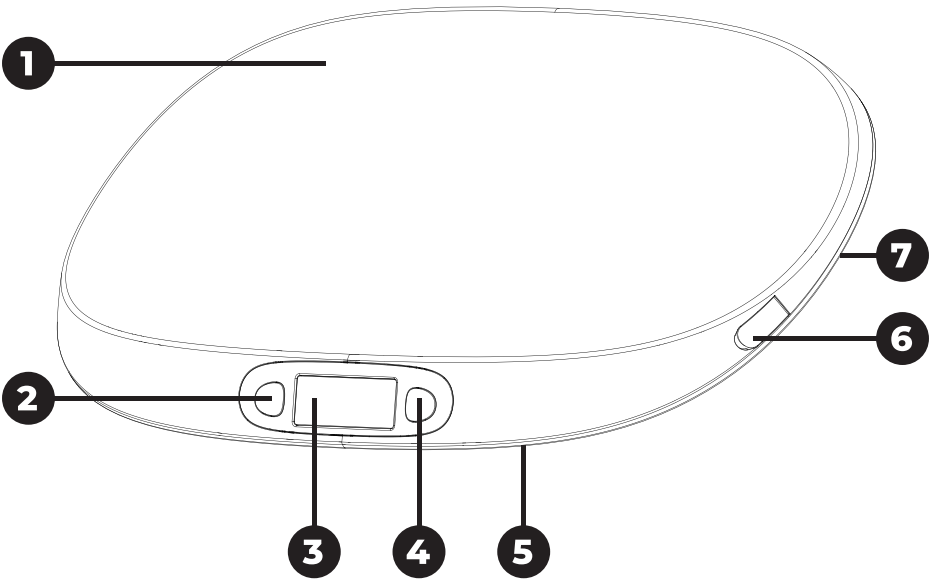
Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are pleased to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, exploitation prevention and fair working conditions.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's helpline. For your convenience in solving potential problems with the product, the above contact information can also be used to lodge complaints or inquire about post-warranty servicing.

The current version of this manual is available at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. This product is intended for household use. Do not use it for medical or other purposes.
2. This product may be used by children and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, when supervised.
3. Children should be supervised to ensure that they do not play with the product.
4. Do not use the product outdoors.
5. Always place it on a stable, level, and firm surface. Do not place it on unstable surfaces, carpets, etc.
6. Do not place the product on or near an electric or gas stove, open flame or other heat source.
7. Do not expose the product to direct sunlight.
8. Do not use the product in highly humid or dusty environments.
9. Do not use this product within strong electromagnetic field.
10. Make sure the platform is clean and dry before weighing. We recommend placing a soft nappy or towel on the surface to lay the baby on. The surface can be cold, and this can be uncomfortable for babies.
11. Take extra care during weighing to keep the baby from falling off. Never leave the scale unattended with a baby on it.
12. The measured values are for personal use only. No medical conclusions can be drawn from the measurement results. If you suspect a possible worsening from the measured values, contact your general practitioner.

13. We recommend weighing at the same time every day. This allows you to achieve more relevant values.
14. Do not place any objects on the product. Store it in a horizontal position.
15. Remove the batteries during longer periods of disuse. There is a risk of damage by electrolyte leakage.
16. Always use new batteries of the same type. Never mix old and new batteries or batteries of different types. Never put batteries in fire, short circuit or damage them in any way. After the end of their lifetime, discard them in environmentally friendly manner. Never dispose of batteries in household waste.
17. Do not use rechargeable batteries to power the appliance.
18. Keep batteries out of the reach of children, as they can be dangerous to them. If a child swallows a battery, seek medical attention immediately.
19. Keep the product clean. Follow the instructions in this manual.
20. Do not immerse the product in water or any other liquids or rinse it under running water.
21. The product does not contain parts which can be serviced by a user. Have all repairs or adjustments carried out by an authorized service centre.

DISPOSAL OF USED BATTERIES

Batteries contain materials harmful to the environment. Therefore, do not dispose of batteries in regular municipal waste. Dispose of used batteries properly at collection points.

OVERVIEW

- 1 Weighing surface
- 2 **HOLD UNIT** button saves measured weight / switches units
- 3 Display
- 4  **TARE** button turns scale on / off and resets
- 5 Battery compartment (located at the bottom of the scale)
- 6 Silicone cover for measuring tape
- 7 Measuring tape rewind switch (at bottom)

TECHNICAL SPECIFICATION

Power	4x 1.5 V battery type AAA
Weight limit	Max 25 kg / 55 lb
Weighing accuracy	± 5 g / 0.01 lb
Selection of weighing units	kg / lb
Dimensions	4.5 × 55.2 × 33 cm (H × W × D)
Weight	1.34 kg

USE

Before first use

Remove the scale from packaging. Before disposing of the packaging, be sure to remove all components. Check the scale for damage.

Wipe the scale with a soft sponge slightly dampened in warm water. Wipe dry with a clean cloth.

Inserting and replacing batteries

The product comes with batteries already inserted. However, they are disconnected by an insulating foil for transport. This must be removed before first use.

1. Turn the scale upside down and remove the battery compartment cover.
2. Remove the insulating foil and put the cover back on.

When replacing batteries, always replace all batteries at the same time. Do not use alkaline and rechargeable batteries at the same time.

1. Turn the scale upside down and remove the battery compartment cover.
2. Remove dead batteries and insert new ones. Observe the polarity markings on the bottom of the battery bay.

3. Replace the cover.

Warning:

Never put batteries in fire, short circuit or damage them in any way. After the end of their lifetime, discard them in environmentally friendly manner. Keep batteries out of the reach of children.

Note:

If you will not be using the scale for a long time, remove the batteries to prevent electrolyte leakage and damage to the scale.

Weighing

1. Place the scale on a solid and level surface. Do not place the scale on carpets or other soft surfaces, as the measurement result may be skewed.
2. Switch the scale on with the button $\frac{\text{ON}}{\text{TARE}}$ button and wait until the display shows "0.000".
3. Carefully place the baby on the plate.
4. The display will show their current weight.
5. Press $\frac{\text{HOLD}}{\text{UNIT}}$ to store the measured weight in memory.

Reset

To reset the displayed value, e.g. to discount the weight of a nappy or towel on which you will set the baby, press $\frac{\text{ON}}{\text{TARE}}$. The value on the display will reset to zero. Proceed to weigh the baby.

Saving Measured Weight

To save the measured weight in the scale's memory, press $\frac{\text{HOLD}}{\text{UNIT}}$. The value will be saved. When the next measurement is saved, the difference from the previous one is indicated by up or down arrows. The new saved value then overwrites the old one.

Note:

Replacing batteries will cause saved values to be lost.

Weight unit change

While the scale is on, long-press $\frac{\text{HOLD}}{\text{UNIT}}$ to switch weight units between pounds and kilograms.

Turning off the scale

To turn off the scale manually, long-press $\frac{\text{ON}}{\text{TARE}}$. The scale automatically switches off after 2 minutes of inactivity.

Tape Measure

The scale is equipped with an extendable tape measure, which can be used to measure the baby's height, head size, etc.

1. Carefully remove the silicone cap on the side of the scale.
2. Pull out the tape measure and take the measurement.
3. When the measurement is complete, press the switch at the bottom of the scale to rewind the tape.

Display reports

Err Overload

The scale has been overloaded. The baby weighs more than the scale's maximum load. Take the baby off the scale to avoid damaging it.

L0 Low battery

The batteries are very low and need to be replaced.

Tips

1. To ensure the measurement accuracy, we recommend taking measurements at the same time each day. Do not use the scale near an unstable and strong electromagnetic field.
2. Do not use the scale in a humid and overly hot environment. The recommended operating temperature is 5-35 ° C.
3. When the batteries are low, replace them. Also, remove the batteries if you will not be using the scale for an extended period of time. There is a risk of electrolyte leakage and damage to the scale.
4. Protect the scale from falls, shocks or jolts.

CLEANING AND MAINTENANCE

Keep the scale clean. Wipe the scale with a soft sponge slightly dampened in warm water. Wipe dry with a clean cloth.

Warning:

Do not immerse the appliance in water or other liquids.

Do not use plastic or metal scouring pads, harsh cleaners, chemicals, petrol, solvents or other similar substances to clean any part of the scale. The surface finish could be seriously damaged.

Storage

Make sure that the scale is clean before you store it. Keep it in a dry and well-ventilated place out of the reach of children or pets. Do not place any objects on the scale.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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CE

SIGURO

Alza.cz a.s., Jankovcova 1522/53, 170 00 Prague 7, Czech Republic
© 2024 Alza.cz a.s. All Rights Reserved.

www.siguro.net