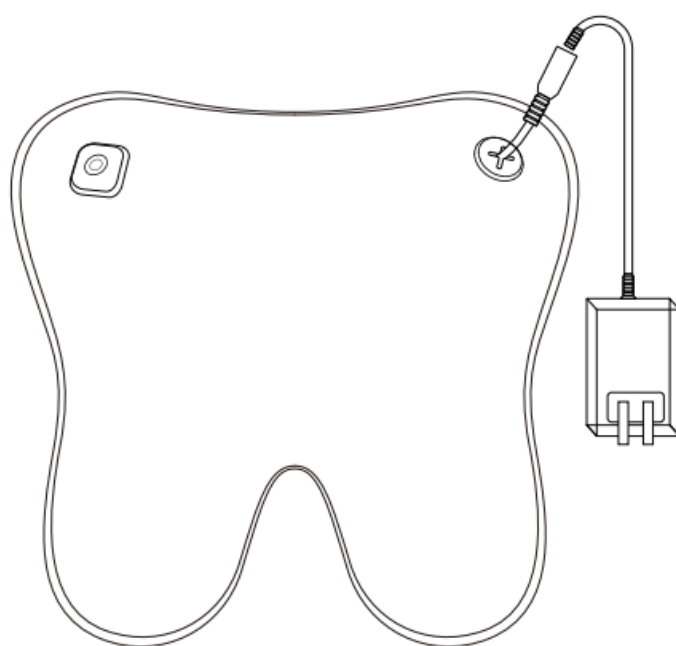
EZback

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


Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please visit our customer support page: www.alza.cz/EN/kontakt.

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Instructions for Users

- Please read this manual thoroughly before using the product and **SAVE THESE INSTRUCTIONS** for future reference.
- This product is designed with an ergonomic build and high-quality components. It offers a relaxing massage experience for the back and neck and can also be used on other body parts such as shoulders, calves, thighs, abdomen, and feet.

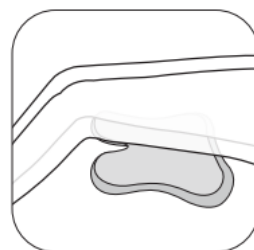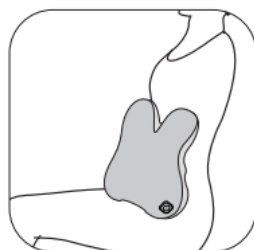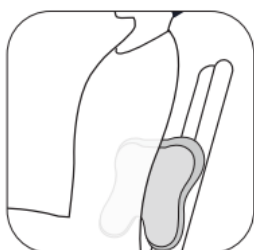Note:

- **Do not use the product on your head or chest.**
- **Do not use the product on sensitive areas, such as eyes, teeth, breasts, implants, genital organs, artificial limbs, and similar areas.**

Important Safety Instructions

DANGER – To Reduce the Risk of Electric Shock:

- Always unplug the appliance immediately after use and before cleaning.
- Do not reach for the appliance if it falls into water. Unplug it immediately.
- Do not place or store the appliance where it may fall into water or a sink.
- Do not immerse the appliance in water or any other liquid.

WARNING – To Reduce the Risk of Burns, Fire, Electric Shock, or Injury:

- Never leave the appliance unattended while plugged in. Unplug it when not in use and before attaching or removing parts.
- Do not operate the appliance under a blanket or pillow, as excessive heating may cause fire, electric shock, or injury.
- Always supervise the use of the appliance when children, invalids, or individuals with disabilities are present.
- Use the appliance only as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never use the appliance if it has been damaged, dropped, or submerged in water. Contact a service center for examination or repair.
- Do not carry the appliance by its cord or use the cord as a handle.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- Do not handle the plug or appliance with wet hands.
- Do not operate the appliance near aerosol (spray) products or where oxygen is being administered.
- When disconnecting, turn all controls to the "OFF" position before unplugging the appliance.

Precautions for Use

- Use only the original adapter intended for this product.
- Confirm the power supply voltage matches the adapter rating before connecting.
- Inspect the product before each use. Do not use if damaged.

- Remove any dust from the plug regularly, as dust accumulation can cause fire risks.
- If the product malfunctions or feels abnormal, stop use immediately and unplug it.
- In case of a power outage, unplug the appliance.
- The surface of the appliance may become hot. Individuals sensitive to heat should use caution.
- Do not stand on or misuse the product as a toy.
- Never disassemble, modify, or attempt repairs on your own. This voids the warranty.
- Do not wear heavy clothing or let fabric block the direct contact of the product with your skin.
- Do not use the product on sensitive, inflamed, or wounded areas.
- Do not use the product without its covering fabric.
- Do not forcefully pull or stab the product with sharp objects.
- Do not attempt to fold the product, as it is not designed for folding.
- Do not plug the appliance into an overloaded outlet.
- This product is not a medical device and should not replace medical treatments.

Disclaimer: This product is designed for relaxation and massage only. It is not intended for medical purposes. Always consult a physician if you have any concerns or experience discomfort while using the product.

Precautions for Maintenance

- Before and during maintenance, always disconnect the product from the power supply.
- **Do not drop the product or its charger.** Avoid hitting the product or charger with strong force.
- **Do not spray water on the product or immerse it in water or other liquids.**
- **Do not use the product in bathrooms or areas with high humidity.**
- If the product has not been used for an extended period, ensure the power button is functioning normally before reuse.
- During maintenance, **do not use volatile oils, diluents, alcohol, strong alkaline detergents, bleach, or abrasive cleaning agents.**
- Keep the product out of the reach of children.
- Ensure the product has fully cooled before storing it.

Groups That Should Exercise Caution When Using the Product:

Individuals undergoing medical treatment or those with the following conditions **must consult a doctor or physiotherapist before using this product:**

- Individuals with a pacemaker or other implanted medical electrical devices.
- People diagnosed with malignant tumors or heart disease.
- Pregnant women.
- Individuals with nerve disorders caused by diabetes.
- People suffering from skin injuries or trauma.
- Individuals with a body temperature above 38°C (100.4°F).
- People prone to numbness in their hands or feet.
- Individuals requiring a respirator or breathing assistance.
- Those with high blood pressure (hypertension).
- Individuals with acute symptoms of visceral diseases (e.g., gastritis, hepatitis, gastroenteritis).
- People suffering from lumbar disc herniation, spondylitis, or similar spinal conditions.
- Individuals with degenerative arthritis, arthritis, or gout.
- Those with spinal fractures or swollen bones.
- Individuals prone to low-temperature burns.
- People with abnormal spinal curvature or conditions affecting the back.
- Individuals with vascular diseases or conditions related to poor blood circulation.
- People using medical restorative devices.
- Those experiencing abnormal body sensations, such as unexplained pain, tingling, or discomfort.

Groups Prohibited from Using the Product:

- Individuals prohibited from exercising or stretching by their doctor, such as those with:
 - Thrombosis
 - Severe aneurysms
 - Acute venous aneurysms
 - Skin infections or dermatitis
- Individuals unable to operate the product independently or express their needs.
- Children aged 14 or younger.
- Individuals with impaired skin sensitivity or reduced heat perception.

Cleaning Instructions

- **Preparation:** Ensure the product is turned off before cleaning. Do not remove the fabric covering.
- Use a vacuum cleaner to remove dust from the fabric instead of patting it. Avoid applying force.
- Wipe dust off the surface with a clean, dry cloth.
- For deep cleaning, dampen a cloth with a neutral detergent and gently wipe the surface.
- Do not submerge or spray water on the product.
- Allow the product to air-dry naturally.
- Avoid storing the product in direct sunlight or humid environments. If not in use for long periods, clean and cover the product to prevent dust accumulation.

Note:

- Disconnect the power cord and ensure it is dry before cleaning.
- The comfort cover is removable and washable. Hand wash it in cold water. **Do not bleach.** Air-dry or hang it to dry.

About the Product


What's in the Box

- 1 × Back and Neck Massager
- 1 × Comfort Cover
- 1 × User Manual
- 1 × Power Adapter

Specifications

Specification	Details
Product Name	Back and Neck Massager
Product Model	EZBack
Adapter	Input: AC 100-240V 50/60Hz
	Output: 12V
Rated Voltage	12V
Rated Power	24W
Net Weight	1.2 kg / 2.6 lb
Dimensions	290 mm × 290 mm × 110 mm (11.4" × 11.4" × 4.3")
Auto-Shutoff Timer	15 minutes

Product Introduction

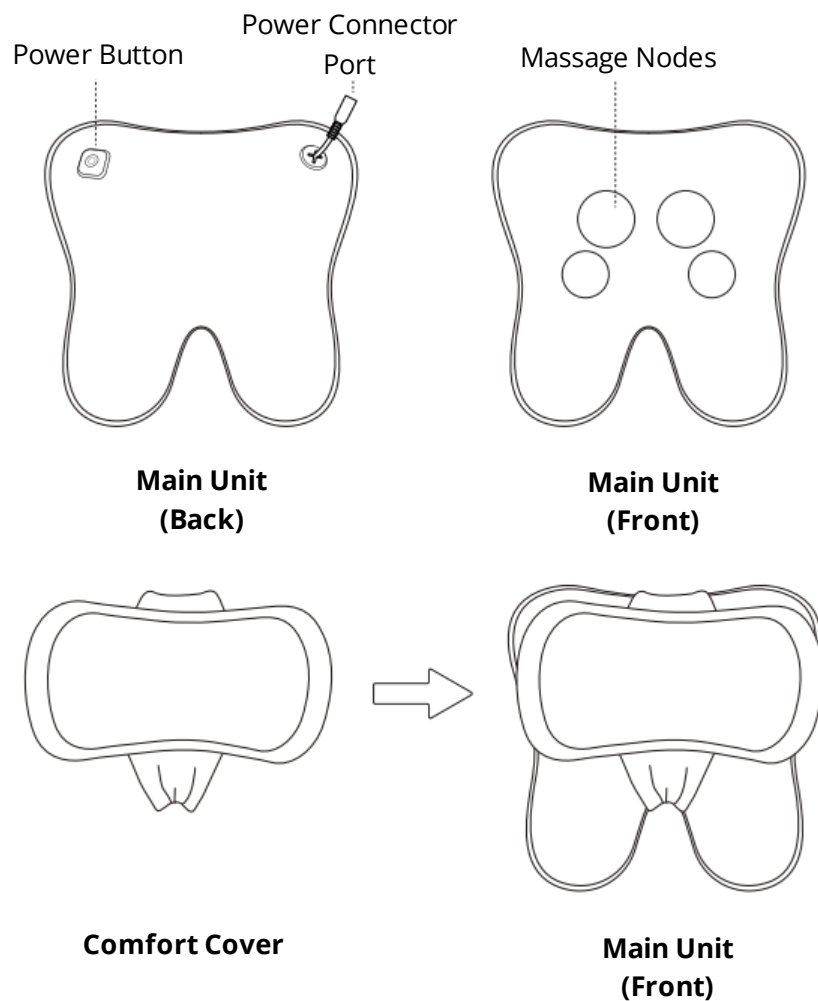
Preparations for Use


- Inspect the covering fabric to ensure it is not damaged before switching on the product.
- Verify that the product and its accessories are complete and in good condition.
- Remove any pendants, jewellery, or other accessories before using the product to avoid interference with the massage.

How to Use

- Plug the wall adapter into an electrical outlet.
- Ensure the adapter plug and power connector port securely "snap" into place when connected.

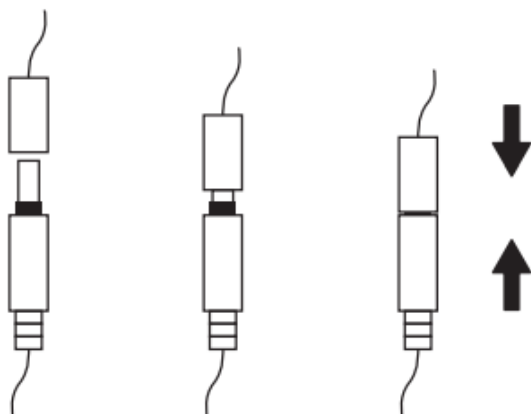Basic Functions

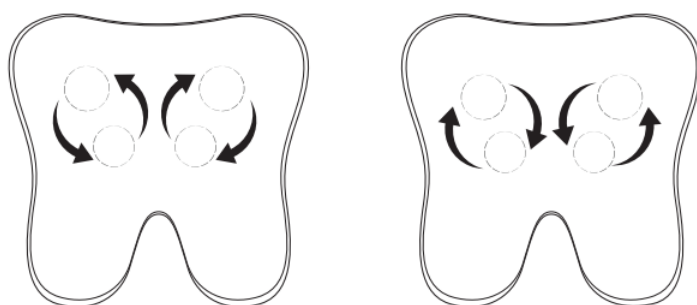Feature	Details
Massage Modes	- A1: Massaging in Direction 1 + Heating - A2: Massaging in Direction 2 + Heating - A3: Massaging in Direction 1, No Heating - A4: Massaging in Direction 2, No Heating
Heating	Heating is enabled in A1 and A2 modes but not available in A3 and A4 modes.

Power Button	- Press and hold for about 3 seconds to turn the product ON/OFF. - Short press to change the mode.
LED Indicator Light	- A1 Mode: Orange Light
	- A2 Mode: Orange Light - A3 Mode: Blue Light - A4 Mode: Blue Light
Auto Shut-Off	The product automatically shuts off after 15 minutes of continuous use.

Troubleshooting

Problem	Cause	Solution
You hear a mechanical sound when the product is in use.	Normal operation: The massage heads make a mechanical sound.	This is normal operation; no action is needed.
The product cannot be turned on.	Overheating protection is activated.	Stop using the product for 30 minutes to cool it down.
	The power adapter is not properly plugged into the electrical socket.	Ensure the electrical cord is securely connected to the socket.
	The power connector is not properly connected to the adapter plug.	Confirm the power connector and adapter plug are fully clicked together.
The product stops working suddenly during use.	Overheating protection is activated.	Stop using the product for 30 minutes to allow it to cool down.
	The built-in 15-minute auto-off timer is activated.	Restart the product.

Regular Checks

If the product has been left unused for a long time, inspect its condition before use to ensure safety and functionality:

Check the product for the following hazards:	If hazards are found:
Whether the power adapter is overheating.	Stop using the product immediately and disconnect the power supply.
Whether the power cord is scratched, deformed, or has a burnt smell.	
Whether there are abnormal sounds or vibrations during operation.	
Whether the product operates properly after being turned on.	
Any other abnormalities or malfunctions.	

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.

- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.


EU Declaration of Conformity

This equipment is in compliance with the essential requirements and other relevant provisions of Directive 2014/35/EU, 2014/30/EU and Directive 2011/65/EU amended (EU) 2015/863.

WEEE

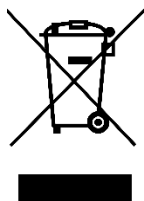This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19/EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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