
SIGURO

**Ultrasonic
skin scrubber**

SGR-SK-U560W

ENGLISH	5
DEUTSCH	10
ČESKY	16
SLOVENSKY	21
MAGYAR	26

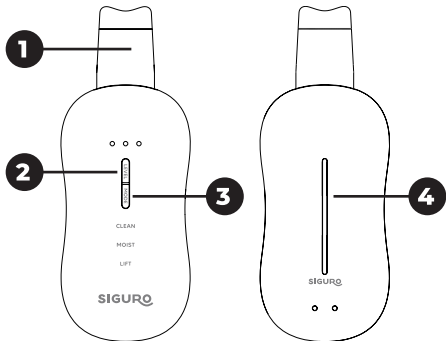SAFETY INFORMATION

1. This device must not be used by children or persons with reduced physical, sensory or mental abilities or lack of experience and knowledge, even if they are supervised or have been instructed in the safe use of the device and understand the possible dangers. Children should be supervised to ensure that they do not play with the device.
2. Never immerse the device in water!
3. The device is intended for household and similar use only. It is not intended for use in a beauty salon or other commercial use. Use of this device is not intended as a substitute for medical care.
4. Do not use the device if you have hypersensitive, damaged or skin with eczema (e.g. fresh scars, various wounds, burns, redness, swollen cheeks or skin diseases). To prevent the spread of micro-inflammation, do not use it even if you have acne on your skin.
5. If you feel pain or discomfort during cleaning, or if you experience any other health problems, switch off the device immediately and consult your doctor or dermatologist about its further use.
6. Do not use the device on children, legally incompetent persons or people who are unable to express their feelings.
7. Do not use the device with exfoliating or otherwise aggressive gels (e.g. with alcohol)!
8. Do not use the device close to sensitive body parts, not intended for this device.
9. Do not use the device in one time longer than specified in the operating instructions! Prolonged use may cause skin injuries.
10. Only use the supplied charging cable. If it is damaged, it must be replaced with the same type by the manufacturer or its authorized service.

- 1 Ultrasonic skin spatula
- 2 Intensity setting button
- 3 Power and mode switch
- 4 Sensor ring

Package parts not shown: charging cable, cleaning cloth and cosmetic bag.

TECHNICAL SPECIFICATIONS

Dimensions (H x W x D)	158 x 58 x 28 mm
Weight	116 g
Power	5 V/1.0 A
Power input	2.4 W
Battery	Rechargeable 3.7 V, 500 mAh

USE

Charging

- You can charge the device using the included charging cable.
- Fully charge the device before using it for the first time, which takes about 2 hours. When fully charged, the battery lasts approximately 35 minutes of continuous operation.

- The charge status will be indicated by the intensity indicators above the **2** button. When charging is complete, all 3 indicators will be lit continuously.

Method of use

The device is used to treat the face using ultrasonic waves. Always remove all make-up and moisten your face with water or tone your skin before using the device.

You can use the ultrasonic skin spatula several times a week. However, the recommended maximum time of use is 5 minutes.

The device offers 3 modes of operation:

- **CLEAN** (cleaning mode): serves for deep pore cleansing and cell regeneration. We recommend using this mode 1 to 2 times a week. Use the mode selection button – **3** to select the **CLEAN** mode. Use the **LEVEL** button to set the desired intensity level. Apply the ultrasonic spatula to the face at an angle of approx. 30 ° and move it slowly along the contours of the face.
- **MOIST** (moistening mode): improves the absorption of facial cosmetics. We recommend using this mode 3 to 5 times a week. Apply facial cosmetics to the skin and gently spread over the face with your fingers. Use the mode selection button – **3** to select the **MOIST** mode. Use the **LEVEL** button to set the desired intensity level. Apply the ultrasonic spatula to the face at an angle of approx. 30 ° and move it slowly along the contours of the face.
- **LIFT** (firming and lifting mode): firms the skin and smoothes wrinkles. We recommend using this mode 3 to 5 times a week after completing the moisturising mode. Use the mode selection button – **3** to select

the **LIFT** mode. Use the **LEVEL** button to set the desired intensity level. Apply the ultrasonic spatula to the face at an angle of approx. 30 ° and move it slowly along the contours of the face.

Recommendation:

- If you are using the ultrasonic skin spatula for the first time, we recommend setting the lowest intensity level and then adjusting as necessary.

Caution:

Stop using the ultrasonic skin spatula if you are uncomfortable using it, feel uncomfortable, have sudden redness on your skin, feel a pinching or uncomfortable pain during use, etc.

Controls

- Turn on the device by pressing the **3** button. The default mode is **CLEAN**.
- Press briefly the **3** button and you can switch between the individual modes. The active mode will be shown by the indicator.
- Use the **2** button to set the required intensity. The set intensity will be shown by indicators.
- You can turn off the device by pressing and holding the **3** button. It switches off automatically after 5 minutes of use and after a short press of button 3 at the end of the mode switching cycle (that is after the **LIFT** mode).

CLEANING AND MAINTENANCE

- Always make sure that the device is switched off before cleaning or storing it.
- You can wash the ultrasonic skin spatula under a weak stream of warm water with the device pointing downwards. Otherwise, water could leak in and cause irreversible damage. Wipe the ultrasonic skin spatula with a dry cloth.
- Wipe the outer body with a dry cloth.
- Do not immerse the device in water or any other liquid. Do not use any cleaning agents, lubricants or polishes.
- Do not dry the device in direct sunlight.

Storage

Store the device in the cosmetic bag provided. Store the device in a dry, well ventilated place. Even if you do not use the device for a long time, charge it regularly.

DISPOSAL OF THE BUILT-IN BATTERY

This symbol indicates that used batteries must not be disposed of in general waste. The product is equipped with a built-in battery. Remove the battery before environmentally friendly disposal of the product or its delivery to an electronics disposal point. The battery must only be removed by a qualified person at an authorised workshop. The battery contains environmentally harmful substances. Take the battery to an appropriate collection point for environmentally friendly disposal.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc. This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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CE

SIGURO

Alza.cz a.s., Jankovcova 1522/53, 170 00 Prague 7, Czech Republic
© 2025 Alza.cz a.s. All Rights Reserved.

www.siguro.net
www.alza.cz/contact