SIGURO

Oral Irrigator

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SGR-OC-I730W

ENGLISH

Safety information	4
Overview	9
Technical specifications	9
Use	10
Cleaning and maintenance	15
Troubleshooting	16

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


51

66

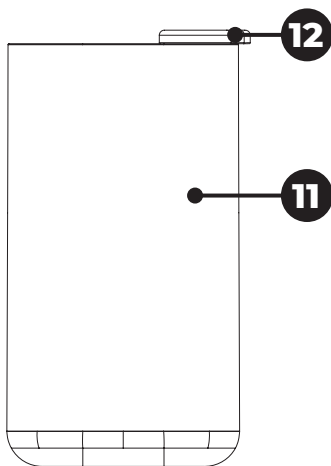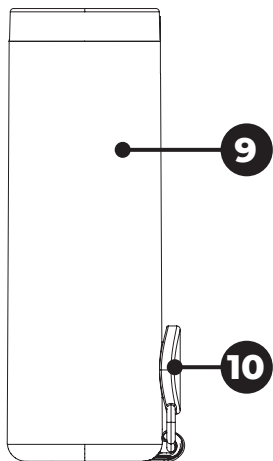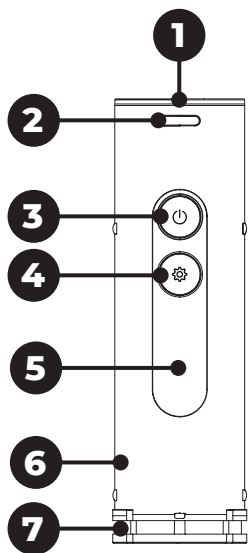Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are pleased to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, exploitation prevention and fair working conditions.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's help-line. For your convenience in solving any problems with the product, we have unified these contact points and the above contacts can be used in case of any complaints or post-warranty service.

The latest version of the operating instructions can be found at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. Children aged at least 8, persons with physical, sensory or mental impairment, and the untrained or inexperienced may use this product with instruction or supervision in its safe use, provided they understand the potential dangers.
2. Children are not allowed to play with the product.
3. Cleaning and user maintenance may not be

performed by unsupervised children.

4. If the charging cable is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a dangerous situation.
5. Before use, make sure that the oral irrigator is correctly assembled.
6. The oral irrigator must only be used with cold or warm water (up to 40 °C) or with special solutions that are intended for use in the oral irrigator only and according to the manufacturer's instructions.
7. Only use the oral irrigator with the charging cable that comes with it. Do not use the charging cable to charge another device. Do not charge the oral irrigator with another charging cable either.
8. Use the product as described in these operating instructions.
9. The oral irrigator is designed for oral hygiene. Do not use it for any other purpose. It is not intended for use in a commercial environment or for laboratory purposes.
10. Do not use the oral irrigator without a head.
11. Check the heads regularly, especially when they are fitted or after use. Do not use a damaged or worn head. Doing so risks severe injury.
12. Do not press the button to release the nozzle during use. The head could come loose and cause injury due to the high pressure.
13. Do not use the oral irrigator if you have an injury in the oral cavity or on the tongue, after a medical procedure, etc. In such case, consult your doctor.
14. Users with periodontal disease or those who have had oral surgery in the last six months should use it under the supervision of a dentist.

15. Do not use the head with heavy pressure on the gums or teeth to prevent damage to the gums.
16. Do not disassemble or modify the oral irrigator. This could cause damage, a dangerous situation (e.g. fire) or personal injury. If it is necessary to repair the product, contact an authorized service centre.
17. Do not place the oral irrigator in an environment with high temperature or humidity and avoid direct sunlight.
18. After use, pour the remaining water out of the water tank. Clean the oral irrigator as described in "Cleaning and maintenance". Water remaining in the tank could cause formation of bacteria, odours and malfunctions.
19. For charging, use a safe and compatible power adapter with a recommended DC 5 V 1 A rated power. Adapter is not included. Make sure that you connect a power supply that matches the rated voltage, otherwise electric shock or fire may occur.
20. Ensure that the charging port and charging cable are dry before charging, otherwise a short circuit may occur.
21. Do not insert foreign objects into the charging port. This could cause a short circuit or electric shock.
22. Only touch the charging cable and adapter with dry hands.
23. Hold the power adapter, not the charging cable, when plugging and unplugging the power supply.
24. Do not damage, modify, bend or pull the charging cable. Make sure that the charging cable does not touch sharp or hot objects. Do not let the charging cable hang over the edge of the table. Route the charging cable so that it is not tripped over or stepped on.
25. Disconnect the charging cable before use, cleaning

and maintenance.

26. To prevent the product from tipping over while charging, place it on a horizontal surface.
27. If you do not use the oral irrigator for a long time, please disconnect the power adapter and charging cable, otherwise the insulation performance may be reduced and electric shock may occur.
28. If you will not be using the oral irrigator for a long time, it should be fully charged before storage and again every 6 months thereafter.

Special instructions regarding the rechargeable Li-ion battery

1. The oral irrigator is equipped with a built-in rechargeable Li-ion battery.
2. Do not expose the oral irrigator to direct sunlight, frost or high temperatures (e.g. sunlight in the car).
3. Keep the oral irrigator in a dry place. Avoid moisture and contact with water.
4. Check the condition regularly. If you notice a bulge, odour, unusual heat, or other suspicious signs, stop using the oral irrigator immediately.
5. Avoid mechanical damage. Do not puncture, squeeze or deform the battery.
6. Do not leave the oral irrigator unattended while charging.
7. Do not use damaged or old batteries.
8. Do not disassemble the battery.
9. Do not expose the battery to fire or high pressure. Do not throw batteries in the fire.
10. Charging and discharging outside the recommended limits (e.g. full zero discharge) can irreparably damage the battery.
11. Do not use or store the oral irrigator where

flammables, volatiles or explosives are used or stored.

12. Avoid contact with leaked electrolyte. Wear protective gloves, goggles, and possibly a protective breathing mask. Place the oral irrigator on a non-flammable surface, ideally outdoors or in a well-ventilated area. Do not handle the battery directly with your hands. Contact a professional service or hazardous waste collection facility. Never dispose of a damaged battery like this in normal waste. In case of skin or eye contact, rinse with plenty of water and seek medical attention. In case of ingestion, seek professional medical attention immediately.

OVERVIEW

- 1 Hole for attaching the head
 - 2 Button to release the head
 - 3  button to turn on / off
 - 4  button to select the mode
 - 5 Display
 - 6 Oral irrigator base
 - 7 Connector for attaching to the tank
 - 8 Charging socket (rear)
 - 9 Small water tank
 - 10 Small water tank lid
 - 11 Large water tank
 - 12 Large water tank lid
- Not pictured: cloth bag, charging cable

List of supplied heads

- 3× standard head
- 1× orthodontic head
- 1× periodontal head
- 1× plaque removal head

TECHNICAL SPECIFICATIONS

Rated voltage	5 V = 1 A
Rated power	5 W
Battery	Li-ion 3.7 V 1650 mAh
Protection class	IPX7

Water tank capacity	250 ml (large) 150 ml (small)
Dimensions (large tank)	150 (retracted oral irrigator body, without head) / 245 (extended oral irrigator body, without head) / 335 (extended oral irrigator body, with head) × 45 × 65 mm (h × w × d)
Dimensions (small tank)	150 (retracted oral irrigator body, without head) / 245 (extended oral irrigator body, without head) / 335 (extended oral irrigator body, with head) × 45 × 45 mm (h × w × d)
Weight	170 g

USE

Before First Use

Unpack the oral irrigator and accessories. Before disposing of packaging, make sure no components are left inside. Check the oral irrigator and accessories for damage.

Wipe the outer surface of the oral irrigator with a soft and dry sponge.

We recommend flushing the water tank with a small amount of water.

Charging

- The oral irrigator must be charged before first use.
- Plug the end of the charging cable into the charging socket at the back of the oral irrigator and plug the other end of the charging cable into the connector of a suitable mains adapter (not included, recommended voltage 5 V / 1 A).
- The battery symbol on the display will gradually fill in.
- Charging time is usually 2.5 hours.

- When the battery symbol appears in full, the oral irrigator is fully charged.
- Disconnect the end of the charging cable from the charging socket and the AC adapter.
- When the battery symbol flashes on the display, the oral irrigator needs to be charged.

**Warning:**

The oral irrigator cannot be used while it is charging. Only use the supplied charging cable.

Notes on charging

If you are using the oral irrigator for the first time or have not used it for more than 6 months, the charging time will vary and the battery symbol may not light up immediately. After a few minutes of charging, the normal charging status is displayed. Charge your oral irrigator even if you don't use it regularly. We recommend that you charge it at least once every 6 months.

Checking Charge Level

Press and hold the  button to check the battery level. The battery symbol appears on the display and its filler indicates the charge status. For example, if the filler is half empty, the battery is approximately 50% charged.

Filling the water tank

Before using the oral irrigator, you must fill the tank with water.

1. Grasp the oral irrigator with one hand and turn it to a horizontal position.
2. Release the tank lid and fill the tank with clean tepid water. The maximum recommended water temperature is 40 °C so using the oral irrigator is comfortable.
3. Close the lid. Make sure it is tightly fitted to prevent water leakage during use.

**Note:**

You can also fill the tank with diluted mouthwash. The dilution ratio depends on personal preference. If you are using a mouthwash for the first time, we recommend a dilution of 1 : 4 (1 part mouthwash, 4 parts tepid water). Adjust the ratio as necessary.

**Warning:**

Do not put toothpaste, salt, mouthwash tablets, or mouthwash containing isopropyl myristate or other chemicals in the tank as they may clog or damage the internal piping and head.

Heads and their use

Standard head is designed to remove food debris, plaque; prevents cavities and tartar in the gaps and grooves of the teeth; suitable for daily cleaning or deep cleaning.

Orthodontic head is designed for cleaning around braces and other orthodontic appliances where a regular toothbrush won't easily reach.

Periodontal head is gentle and suitable for cleaning deeper gum pockets in clients with periodontal disease or other gum problems.

Plaque remover head effectively removes plaque from between the teeth and along the gum line to help prevent cavities and inflammation.

Attaching and removing the head

- Insert the end of the head into the opening at the top of the oral irrigator. Push down carefully until you hear a click.
- If you need to remove the head, first press the button to release the head, then pull the head upwards away from the oral irrigator.
- Check the head for damage or wear before fitting or after removal.

Warning:

- Do not use the oral irrigator without a head.
- Do not use a damaged or worn head. Doing so risks severe injury.
- Do not press the button to release the nozzle during use. The head could come loose and cause injury due to the high pressure.

Assembling the oral irrigator

- The oral irrigator is supplied with two water tanks.
- **Large tank:** attach the base of the oral irrigator to the connector at the top of the tank. Push down gently until you hear a click.
- **Small tank:** remove the top cover of the tank and attach the oral irrigator base to the connector at the top of the tank. Push down gently until you hear a click.
- If you need to remove the oral irrigator base, grasp the base with one hand and hold the tank with the other hand. Pull the base upwards to release it from the tank.

Using the oral irrigator

- Hold the oral irrigator in an upright position, lean over the sink and insert the head into your mouth. Align the head with the teeth or gum line, slightly clamp the lips to prevent water from splashing, and leave a gap to allow water to flow from the mouth into the basin.
- Press the  button repeatedly to select one of the operating modes:
 - **gentle mode:** suitable for sensitive teeth or if you are using the oral irrigator for the first time;
 - **normal mode:** suitable for everyday use;
 - **pulse mode:** suitable for cleaning the interdental spaces and massaging the gums;
 - **manual mode:** press the  button and the oral irrigator sprays water once. If you hold the button down, the oral irrigator will spray water continuously. This mode is suitable if you have a cavity in your tooth, or to clean hard-

to-reach areas such as around your wisdom teeth;

- **DIY (do-it-yourself) mode:** In the DIY mode, you can switch between three levels of cleaning intensity by long pressing the  button.
- Press the  button briefly and the oral irrigator will switch on in the selected mode.
- Adjust the angle at which the water is directed towards the teeth and gums. The recommended effective angle is 90°.
- We recommend starting from the back teeth, slowly moving the head towards the front teeth along the gum line and stopping briefly in the gap between the teeth. Do not spray water directly on the gums.
- When switched on, a 2-minute timer starts and the oral irrigator automatically switches off when the timer expires.
- If you need to switch off the oral irrigator earlier, press the  button briefly.
- After use, clean the oral irrigator as described in the “Cleaning and Maintenance” section.

Note:

- Do not apply force to the teeth or gums with the head, as the anti-rotation function may be activated. The oral irrigator will automatically pause and the three mode lights will flash for 10 seconds. After they go off, switch the shower back on to resume operation.
- The oral irrigator will start in the mode in which it was switched off when switched on again.

Caution when using

Distilled water from factory testing may remain in the water tank or in the tube. It is a normal occurrence. Rinse the water tank with a small amount of water before first use.

Even if you have healthy gums, you may feel itching or bleeding when you first use a oral irrigator. This is because the gums are not yet used to this stimulation. This condition usually resolves within 1–2 weeks.

The first time you use the oral irrigator, we recommend choosing the gentle mode for one week before choosing another mode.

**Warning:**

If discomfort persists for more than 2 weeks, stop using the oral irrigator and consult your dentist.

CLEANING AND MAINTENANCE

Turn off the oral irrigator before any cleaning or maintenance. For hygiene reasons, clean the oral irrigator and the head after each use.

Warning:

- Never immerse the oral irrigator base or the power cable in water or any other liquid.
- Do not use water at temperatures above 40 °C, corrosive liquids, thinners, alcohol, etc. to clean the oral irrigator.

Daily cleaning

- Remove the oral irrigator base from the tank.
- Open the tank lid and pour out the remaining water. Rinse the tank with clean water and leave the lid open to allow the tank to dry.
- Turn on the oral irrigator to remove any remaining water from the tube.
- Remove the head, rinse with clean water and wipe dry.
- Wipe the oral irrigator base with a soft, dry sponge.

Descaling

Regular descaling of your oral irrigator is important for its proper function and longevity. Hard water can gradually clog its internal parts, which reduces performance and can lead to malfunctions.

Descale with a small brush (such as a toothbrush). Rinse the tank thoroughly before next use.

Storage

Clean the oral irrigator according to the instructions above. Make sure the tanks are completely dry before closing the lid.

The body of the oral irrigator can be inserted into a large or small tank.

Place the oral irrigator and heads in a cloth bag. Keep the cloth bag in a safe place away from children and pets.

The top cover of the small tank is equipped with brackets on the bottom to which you can attach the heads. These are then placed in a small tank.

Recharge your oral irrigator regularly (at least once every six months), even if you don't use it for a long time. This will extend the life of the built-in battery.

TROUBLESHOOTING

Problem	Possible cause	Solution
The oral irrigator does not work properly and no button responds to being pressed.	Weak battery.	Charge the oral irrigator.
	Internal malfunction.	Contact authorized service centre.
The battery symbol did not light up during charging.	Unusual connection to power supply.	Connect to the correct power source.
	Internal malfunction.	Contact authorized service centre.
The operating time after charging is too short.	The charging time is too short.	Let the oral irrigator fully charge.
The oral irrigator cannot be switched on.	Weak battery.	Charge the oral irrigator.

Problem	Possible cause	Solution
Water pressure is too low.	Weak battery.	Charge the oral irrigator.
	Clogged or damaged nozzle.	Replace the nozzle.
	Clogged filter nozzle.	Clean the filter nozzle.

Removing the battery

The oral irrigator uses a lithium battery for its operation. If the oral irrigator reaches the end of its useful life, the batteries should be removed and recycled or disposed of in accordance with local waste management regulations before disposal. If you have any questions, please contact your local waste management authority regarding the removal of rechargeable batteries. If the battery cannot be removed, it is recommended to dispose of it strictly in accordance with local regulations for the separation and management of household waste.

Warning:

We recommend to have the battery removed by an authorized service technician or other experienced person.

When removing the battery, proceed as follows:

1. First remove the head, water tank and suction hose.
2. Using a screwdriver, loosen the two screw caps on the bottom cover, remove the screws, and then remove the bottom cover by pulling downwards.
3. Push the mechanism upwards from the handle.
4. Remove the battery from the battery holder.
5. For proper recycling or disposal of the battery, please contact your local municipality and waste management department to dispose of the battery in an environmentally friendly manner.

6. Be careful not to short-circuit the positive and negative poles of the battery.

Warning:

- Do not throw the battery into a fire, it could explode or release toxic and harmful substances. Do not short-circuit the battery, it could burn and cause a fire.
- If the product is not to be disposed of, do not dismantle it yourself, as premature dismantling will damage it.
- Basic safety precautions must be observed during operation and the power must be disconnected before removing the battery.
- Before removing the battery, make sure that the battery is discharged. If it needs to be discharged, you can briefly press the switch button to switch on the appliance and let it run until it stops automatically. Repeat this procedure several times until the battery is discharged.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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© 2025 Alza.cz a.s. All Rights Reserved.