SIGURO

Personal Scale

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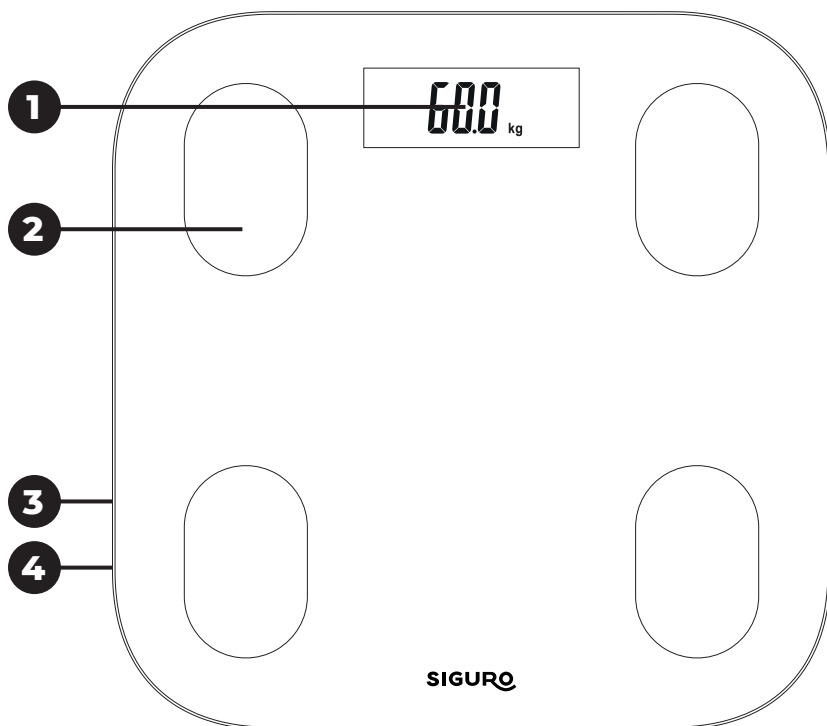
Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are pleased to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, exploitation prevention and fair working conditions.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's helpline. For your convenience in solving potential problems with the product, the above contact information can also be used to lodge complaints or inquire about post-warranty servicing.

The current version of this manual is available at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. This product is intended for household use. Do not use it for medical or other purposes.
2. This product may be used by children and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, when supervised.
3. Children should be supervised to ensure that they do not play with the product.
4. Do not use the product outdoors.
5. Always place it on a stable, level, and firm surface. Do not place it on unstable surfaces, carpets, etc.
6. Do not place the product on or near an electric or gas stove, open flame or other heat source.
7. Do not expose the product to direct sunlight.
8. Do not use the product in highly humid or dusty environments.
9. Do not use this product within strong electromagnetic field.
10. Make sure the platform is clean and dry before weighing.
11. The measured values are for personal use only. No medical conclusions can be drawn from the measurement results. If you suspect a possible worsening from the measured values, contact your general practitioner.
12. We recommend weighing at the same time every day. This allows you to achieve more relevant values.
13. Do not place any objects on the product. Store it in a horizontal position.
14. Remove the batteries during longer periods of disuse. There is a risk of

damage by electrolyte leakage.

15. Always use new batteries of the same type. Never mix old and new batteries or batteries of different types. Never put batteries in fire, short circuit or damage them in any way. After the end of their lifetime, discard them in environmentally friendly manner. Never dispose of batteries in household waste.
16. Do not use rechargeable batteries to power the appliance.
17. Keep batteries out of the reach of children, as they can be dangerous to them. If a child swallows a battery, seek medical attention immediately.
18. Keep the product clean. Follow the instructions in this manual.
19. Do not immerse the product in water or any other liquids or rinse it under running water.
20. The product does not contain parts which can be serviced by a user. Have all repairs or adjustments carried out by an authorized service centre.

DISPOSAL OF USED BATTERIES

Batteries contain materials harmful to the environment. Therefore, do not dispose of batteries in regular municipal waste. Dispose of used batteries properly at collection points.

OVERVIEW

- 1 Display
- 2 Weigh plate with electrodes
- 3 Battery compartment (located at the bottom of the scale)
- 4 Button for switching weight units (located at the bottom of the scale)

TECHNICAL SPECIFICATION

Power supply	4 × 1.5V battery type AAA
Weight limit	Max 180 kg / 396 lb
Weighing accuracy	± 10 g / 0.02 lb (0.1 – 10 kg) ± 50 g / 0.1 lb (10 – 100 kg) ± 100 g / 0.2 lb (100 – 180 kg)
Selection of weighing units	kg / lb
Dimensions	2.9 × 28 × 28 cm (H × W × D)
Weight	1.54 kg

USE

This scale transfers data to your smartphone app via Bluetooth. The measured values are stored in the app, where you can e.g. put them in a chart for a better grasp of your physical health.

Download the app to your smartphone and with the help of a scale you can take the first step to control your weight.

You can also use this scale to measure and evaluate body fat (in %) and hydration rate (in %). In order to always achieve accurate measurement results, it is important to enter the correct information about your gender, height or age in the first place.

The measured values are for personal use only. This scale is not designed for medical use and should not be used to weigh patients in healthcare facilities. This scale should not be used by children, pregnant women, people with symptoms of oedema or people with heart bypass.

Before First Use

Remove the scale from packaging. Before disposing of the packaging, be sure to remove all components. Check the scale for damage.

Wipe the scale with a soft sponge slightly dampened in warm water. Wipe dry with a clean cloth.

Inserting and replacing batteries

The product comes with batteries already inserted. However, they are disconnected by an insulating foil for transport. This must be removed before first use.

1. Turn the scale upside down and remove the battery compartment cover.
2. Remove the insulating foil and put the cover back on.

When replacing batteries, always replace all batteries at the same time. Do not use alkaline and rechargeable batteries at the same time.

1. Turn the scale upside down and remove the battery compartment cover.
2. Remove dead batteries and insert new ones. Observe the polarity markings on the bottom of the battery bay.
3. Replace the cover.

Warning:

Never put batteries in fire, short circuit or damage them in any way. After the end of their lifetime, discard them in environmentally friendly manner. Keep batteries out of the reach of children.

Note:

If you will not be using the scale for a long time, remove the batteries to prevent electrolyte leakage and damage to the scale.

Weighing

Place the scale on a solid and level surface. Do not place the scale on carpets or other soft surfaces, as the measurement result may be skewed.

Stand barefoot with your feet so that your feet are on the electrodes. Stand still and straight as the measurement takes place.

Once the values on the display stabilize and flash, they are transferred to the application on your smartphone. After a successful transfer, the scale displays the weight and measured values, then switches off automatically. The received data is evaluated in the application and stored in the measurement history.

You can search and compare the measured values in the application. At the same time, you will receive a notification asking if these are your personal results.

If the measured values appear to be in order, click on the request to confirm it.

Once the request is complete, the application will display the weighing screen.

Click on the weight value to see more information about your body's measured values.

The detailed record includes various physical evaluations and results.

Use with application

The scale allows you to manage and store weighing and measurement results in the Tuya mobile app. The measured values are transferred to your smartphone via Bluetooth technology.

Downloading the App

Download the Tuya app to your smartphone using the QR code at the top of this manual.

Device registration and connection

1. Open the application and create your account or log in to an already created account. Follow the instructions on the screen of your smartphone.
2. On the app homepage, tap "Add Device" or "+" to find and add a new device.
3. Start an automatic search for a new device. Once the scale is found, click on the scale icon to start the connection and pairing. The scale display must be lit during

connection and pairing. If the application did not find the personal scale, press and hold the change units button for 3-5 seconds and check if the personal scale icon is displayed.

4. Click "Finish" after successfully adding a new device.
5. Return to the home screen and select the personal scale you want to use.
6. Select the user ID.
7. Enter personal data and save them.

Note:

To ensure comprehensive and methodical measurements, enter true and real personal data.

Display reports

	Bluetooth	Kg lb	Weight unit
0000	Weighing in progress	Err	Overload
ErrL	Very low body fat percentage	Lo	Low battery
ErrH	Very high body fat percentage		

Weighing Items

The scale can also be used to weigh smaller items.

1. Step lightly on the scale to turn it on.
2. Place the item in the centre of the weighing area and wait for the display to show its weight.
3. Remove the item from the scale.

Note:

Items lighter than 100 g will not be registered by the scale, and their weight will not be displayed.

Weight unit change

You can change the weight units either directly in the application or on the scale. However, the units in the application are not synchronized with the units on the scale. To change the weight units on the scale, turn the scale upside down and switch to the desired weight units.

Tips and advice

1. To ensure the measurement accuracy, we recommend taking measurements at the same time each day. Do not use the scale near an unstable and strong electromagnetic field.
2. If you will be using the scale after a long time, first step on it lightly, wait until the display lights up and "0.0" is displayed. Then wait until the scale switches off automatically. Only then can you perform the weighing.
3. Stand barefoot with your feet so that your feet are on the electrodes. Stand still and

straight as the measurement takes place. If you stand on the edge of the scale, you could lose your balance and fall.

4. Do not use the scale in a humid and overly hot environment. The recommended operating temperature is 5-35 ° C.
5. When the batteries are low, replace them. Also, remove the batteries if you will not be using the scale for an extended period of time. There is a risk of electrolyte leakage and damage to the scale.
6. Protect the scale from falls, shocks or jolts.

CLEANING AND MAINTENANCE

Keep the scale clean. Wipe the scale with a soft sponge slightly dampened in warm water. Wipe dry with a clean cloth.

Warning:

Do not immerse the appliance in water or other liquids.

Do not use plastic or metal scouring pads, harsh cleaners, chemicals, petrol, solvents or other similar substances to clean any part of the scale. The surface finish could be seriously damaged.

Storage

Make sure that the scale is clean before you store it. Keep it in a dry and well-ventilated place out of the reach of children or pets. Do not place any objects on the scale.

TROUBLESHOOTING

Problem	Solution
No values were transferred to the app.	Make sure the scale is added to the device list and is successfully paired with your smartphone.
	Make sure your smartphone's Bluetooth is switched on and the scale is also on and within Bluetooth range of your phone.
	Check that the batteries are not empty. Replace them if necessary.
	Close and restart the app.
	Turn the Bluetooth function on your phone off and on again and check the pairing with the scale.

Problem	Solution
It is not possible to pair the scale with a smartphone.	The scale can only be paired with one smartphone at a time.
	Check the Bluetooth version of your smartphone. The scale supports Bluetooth version 4.0 and higher.
The display is still off.	Check the condition of the batteries and replace them if necessary.

If the problem with the scale is not listed in the table above or the problem cannot be solved, switch off the scale and contact an authorized service centre.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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R20250828

SIGURO

Alza.cz a.s., Jankovcova 1522/53, 170 00 Prague 7, Czech Republic
© 2025 Alza.cz a.s. All Rights Reserved.

www.siguro.net
www.alza.cz/contact