SIGURO

Massage Pillow

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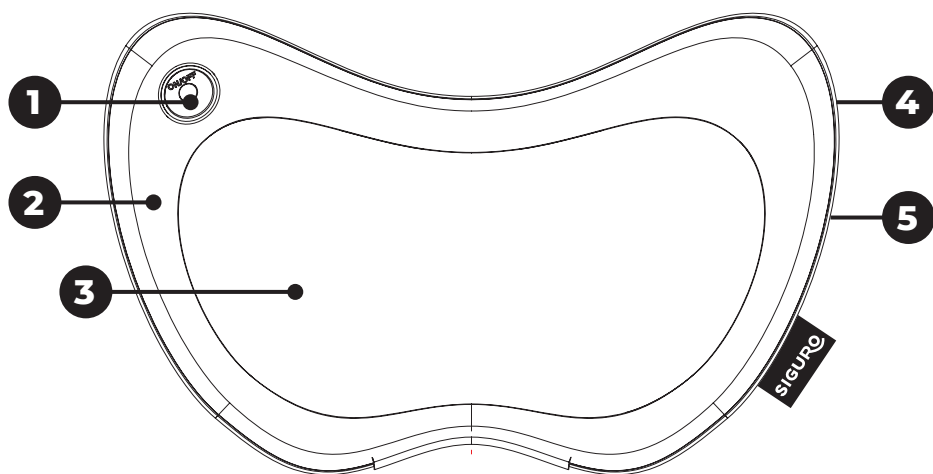
Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are pleased to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, abuse prevention, and workplace fairness.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's helpline. For your convenience in solving potential issues with the product, the unified contact points above can also be used for any complaints or post-warranty servicing.

The current version of this manual is available at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. Children aged at least 8, persons with physical, sensory or mental impairment, and the untrained or inexperienced may use this appliance with instruction or supervision in its safe use, provided they understand the potential dangers.
2. Children are not allowed to play with the product.
3. Cleaning and user maintenance may not be performed by unsupervised children.
4. If the power adapter is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons to avoid danger.
5. This product has a heated surface. Persons with impaired heat sensitivity should take care when using this appliance.
6. Use only the adapter supplied with the product. Do not use different adapters. There is a risk of electric shock or property damage.
7. Do not use the supplied adapter to power another product.
8. Use the product as described in this manual.
9. This product is designed as a massage pillow to massage the cervical spine, shoulders and the top of the shoulder blades. Do not use it for any other purpose. It is not intended for use in a commercial environment or for laboratory purposes.
10. The product is designed for indoor use. Do not use it outdoors.
11. Do not place the product in an environment with high temperature or humidity or in direct sunlight. Do not use the product at temperatures above 40 °C.

12. Avoid using sharp objects.
13. Do not use it if damaged.
14. Use only the accessories supplied with the appliance.
15. Use the product only on dry parts of the body. Do not use massage cream.
16. Turn off and unplug the product after use and before cleaning, storage or moving.
17. Do not cover the product while in use. Otherwise, overheating and damage may occur.
18. Do not disassemble or modify the product. It may become damaged, create a dangerous situation (e.g. fire) or cause personal injury. If necessary, contact an authorized service centre.
19. Stop using the product if it makes you feel unwell, dizzy or uncomfortable.
20. Do not touch the product with wet or damp hands, especially when plugging or unplugging it.
21. Consult your physician before using the product if:
 - your legs are injured, post-surgery, or you have another leg condition;
 - you have high blood pressure, heart disease, cerebral arteriosclerosis, cerebral thrombosis, malignant tumour, or osteoporosis;
 - you are pregnant;
 - you use a pacemaker or other medical aids.
22. Do not use the appliance continuously for more than 20 minutes.
23. Do not use the product while sleeping.
24. Use only dry cloth for cleaning. Follow the instructions in this manual.

Protection class II equipment

The "interior use only" symbol means the appliance is not designed for outdoor use.

OVERVIEW

- 1 ON/OFF button
 - 2 Massage pillow
 - 3 Massage head (under cloth protector)
 - 4 Adapter port (at back)
 - 5 Elastic strap (at back)
- Power adapter (not pictured)

TECHNICAL SPECIFICATION

Appliance

Nominal power	24 W
Input voltage	12 V 
Adapter model	SGR-MM-X002
Dimensions	310 × 100 × 190 mm
Weight	1.3 kg

Adapter

Manufacturer name or trademark, trade registration number and address	Alza.cz a.s. Jankovcova 1522/ 53,170 00 Prague 7, Czech Republic
Product identifier	PN: SGR-MM-X002 Model: ZL-024WL1202000EU01
Input voltage	100–240 V~
Input frequency	50/60 Hz
Output voltage	12.0 V
Output current	2.0 A
Output power	24 W
Mean efficiency in active mode	89.03 %
Efficiency at low load (10%)	86.21 %
Standby power consumption	0.039 W

BEFORE FIRST USE

Remove the massage pillow and its accessories from the packaging. Before disposing of packaging, make sure no components are left inside. Check the massage pillow and accessories for damage.

Wipe the outer surface of the massage pillow with a soft and dry sponge.

USE

1. Plug the device into a properly grounded mains outlet.
2. Insert the adapter plug into the input terminal at the back of the pillow.
3. Press **ON/OFF** briefly to switch on the massage pillow. The massage heads rotate clockwise or counterclockwise, reversing direction every minute.
4. Press **ON/OFF** again to turn on the massage pillow's heating feature.
5. Another press of **ON/OFF** will turn off the heating.
6. Press and hold **ON/OFF** for 3 seconds to turn the massage heads off.

Clean the appliance as described in the "Cleaning and maintenance" section.

Caution:

We recommend using the massage cushion for no more than 20 minutes at a time. Do not use it for longer periods.

CLEANING AND MAINTENANCE

Always switch off and unplug the appliance before cleaning or any maintenance.

Warning:

Do not immerse the appliance, adapter, or power cord in water or other liquids.

Caution:

Do not use plastic or metal scrapers, harsh cleaners, chemicals, petrol, solvents, or other similar substances to clean any parts of the appliance. The surface finish could be seriously damaged.

Wipe the outer surface of the massage pillow with a soft and dry sponge.

Storage

Make sure the cushion is cool and clean before storing.

Store in a dry, airy place out of reach of children and pets. Do not place any objects on top of the product.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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R20250611

SIGURO

Alza.cz a.s., Jankovcova 1522/53, 170 00 Prague 7, Czech Republic
© 2025 Alza.cz a.s. All Rights Reserved.

www.siguro.net
www.alza.cz/contact