SIGURO

Electric toothbrush

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

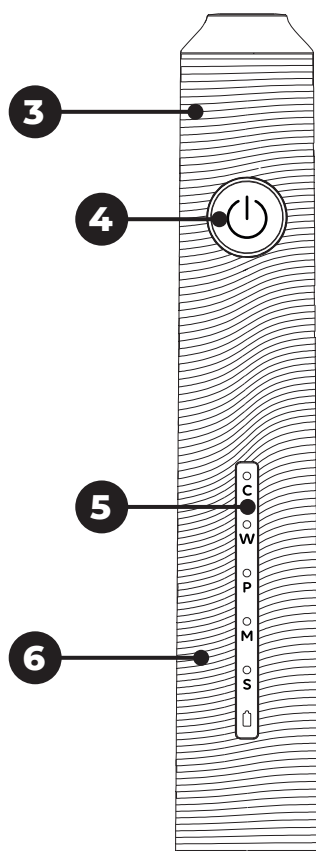
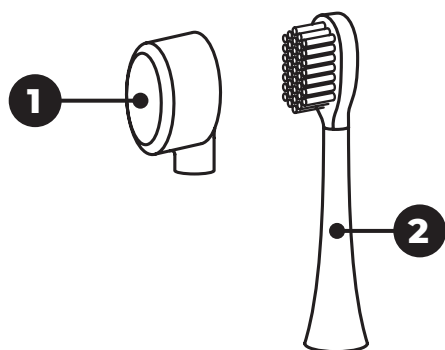
A

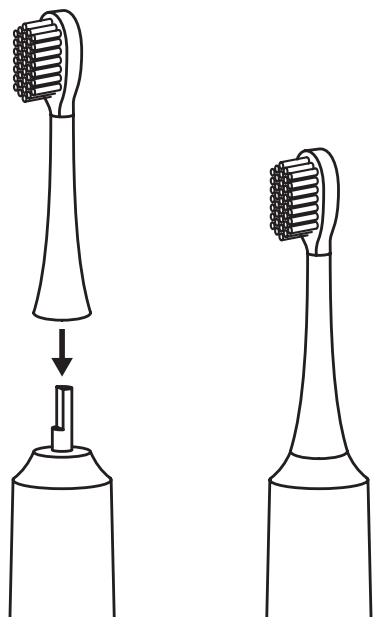B

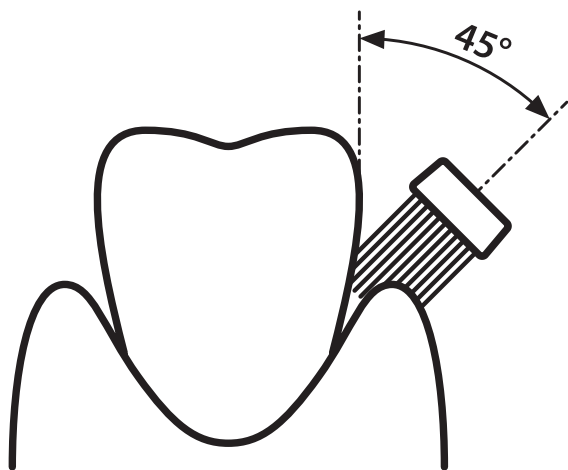Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are pleased to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, exploitation prevention and fair working conditions.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's help-line. For your convenience in solving potential issues with the product, the unified contact points above can also be used for any complaints or post-warranty service.

The current version of this manual is available at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. This appliance is intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, provided they have been given supervision or instruction concerning use of the appliance in a safe manner and understand the potential hazards.
2. Children are not allowed to play with the product.
3. Cleaning and user maintenance may not be performed by unsupervised children.
4. If the charging cable is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a dangerous situation.
5. Before use, make sure the toothbrush is correctly assembled.
6. Only charge the toothbrush with the supplied charging cable. Do not use the charging cable to charge other devices. Do not charge the toothbrush with another charging cable.
7. Use the product as described in this manual.
8. This toothbrush is designed for oral hygiene. Do not use it for any other purpose. It is not intended for use in a commercial environment or for laboratory purposes.
9. Do not use the toothbrush without the brush head on.
10. Check the heads regularly, especially when attaching them and after use. Do not use a damaged or worn head. Doing so risks severe injury.
11. Do not apply pressure to the brush head during use. High pressure will cause it to stop.

12. Do not use the toothbrush when your oral cavity or tongue is injured, after a medical procedure, etc. In such cases, consult your physician.
13. Do not use the head with heavy pressure on the gums or teeth to prevent gum injury.
14. Do not disassemble or modify the toothbrush. It may become damaged, create a dangerous situation (e.g. fire) or cause personal injury. If necessary, contact an authorized service centre.
15. Do not place the product in hot or humid environments or in direct sunlight.
16. Clean the toothbrush after use as described in "Cleaning and Maintenance."
17. Charge with a safe and compatible power adapter. We recommend DC 5 V, 0.5 A. Adapter not included. Make sure your power source matches the rated voltage, otherwise electric shock or fire may result.
18. Ensure the charging port and cable are dry before charging, otherwise a short circuit may occur.
19. Do not insert foreign objects into the charging port. This may cause a short circuit or electric shock.
20. Only touch the charging cable and adapter with dry hands.
21. Hold the power adapter, not the charging cable, when plugging and unplugging the power supply.
22. Do not damage, modify, bend or pull the charging cable. Make sure the charging cable does not touch sharp or hot objects.
23. Do not let the charging cable hang over the edge of the table. Route the charging cable so that it is not tripped over or stepped on.
24. Disconnect the charging cable before use, cleaning and maintenance.
25. Place the product on a horizontal surface to prevent it from tipping over during charging.
26. When not using the toothbrush for longer periods, disconnect the power adapter and charging cable, otherwise the insulation performance may degrade and electric shock may occur.
27. When not using the toothbrush for longer periods, it should be fully charged before storage and again every 3 months thereafter.

Special instructions for rechargeable Li-ion battery

1. The toothbrush is equipped with a rechargeable Li-ion battery.
2. Do not expose the toothbrush to direct sunlight, frost or high temperatures (e.g. car parked in the sun).
3. Keep the toothbrush in a dry place.
4. Check appliance condition regularly. If you notice bulging, odour, unusual heat or other suspicious signs, stop using the appliance immediately.
5. Avoid mechanical damage. Do not puncture, squeeze or deform the battery.
6. Do not leave the toothbrush unattended when charging.
7. Do not use damaged or old batteries.
8. Do not disassemble the battery.
9. Do not expose the battery to fire or high pressure. Do not throw batteries in the fire.

10. Charging and discharging outside the recommended limits (e.g. full discharge) may irreversibly damage the battery.
11. Do not use or store the toothbrush where flammables, volatiles or explosives are used or stored.
12. Avoid contact with leaked electrolyte. Wear protective gloves, goggles, and possibly a breather mask. Place the toothbrush on a non-flammable surface, ideally outdoors or in a well-ventilated area. Do not handle the battery directly with bare hands. Contact a professional battery disposal service or hazardous waste collection facility. Never dispose of damaged batteries in household waste. In case of skin or eye contact, rinse with plenty of water and seek medical attention. In case of ingestion, seek medical attention immediately.

BATTERY DISPOSAL

This symbol indicates that used batteries must not be disposed of in general waste. The product is equipped with a built-in battery. This must be removed before the product is disposed of in an environmentally friendly manner or taken to a collection point. Only qualified personnel at authorized service centres may remove the battery pack. The battery contains environmentally harmful substances. Take it to a collection point for environmentally friendly disposal.

This symbol indicates the device is protected against temporary immersion in water (up to 1 m, for at most 30 minutes). The device can be cleaned under running water, but is not designed for use in wet environments (e.g. bath or shower).

OVERVIEW

- 1 Brush head cover
- 2 Replacement brush head (3 pcs included)
- 3 Toothbrush motor unit
- 4  button switches device on/off
- 5 Operating lights
- 6 Handle
- 7 Charging stand

Not pictured: charging cable, protective case

TECHNICAL SPECIFICATION

Input voltage	5 V  1 A
Output	2 W
Battery	Li-Ion, 3.7 V 800 mAh
Protection class	IPX7
Dimensions	28 × 28 × 252 mm
Weight	109 g (including head)
Power consumption in off mode / in standby mode / in standby mode when connected to the mains	Not available
Time required for the appliance to reach the corresponding mode or state with low power consumption	Not available
Power consumption (other conditions)	0.41 W (measured according to sampling method point 5.3.2 EN 50564:2011 + (EU) 2023/826; the product is connected to the mains.)

USE

Before First Use

Remove the toothbrush and its accessories from the packaging. Before disposing of packaging, make sure no components are left inside. Check the toothbrush and accessories for damage.

Wipe the outer surface of the toothbrush with a soft and dry sponge.

Charging

- The toothbrush must be charged before first use. Always switch off the toothbrush before charging.
- Plug the end of the charging cable into the socket at the back of the charging stand. Plug the other end of the cable into a suitable mains adapter (not included, we recommend 5 V / 1 A).
- Place the toothbrush in the charging stand. Check that the electrical connection is good.
- The battery light will flash red as the toothbrush charges.
- Charging time is usually ≤ 12 hours.
- When the battery light turns blue, the toothbrush is fully charged.
- Disconnect the charging cable from the stand and the AC adapter.
- When the battery light flashes red, the toothbrush needs charging.

Warning:

Do not use the toothbrush while it is charging. Only use the supplied charging cable.

Note:

Charge the toothbrush even when not using it regularly. We recommend charging at least once every 3 months.

Inserting and Removing the Brush Head

Place the brush head on the shaft at the top of the toothbrush (Fig. A). Push down gently until you hear a click.

When removing, pull the brush head upwards away from the toothbrush.

Three brush heads are included. Each has different bristle shape and is suitable for a different mode:

- blue – Clean mode
- green – Sensitive mode
- grey – White mode

Toothbrush modes

Mode	Indicator light	Use
Clean	C	Standard mode suitable for everyday brushing
White	W	Suitable for thorough cleaning, helps remove surface stains (e.g. from coffee or tea)
Polish	P	Helps make your teeth appear smoother and shinier
Massage	M	Suitable for gum massage
Sensitive	S	Gentle mode suitable for sensitive teeth, exposed roots, periodontal disease, etc.

Using the Toothbrush

1. Remove the toothbrush from the charging stand, then remove the protective cap from the brush head.
2. Apply toothpaste to brush head bristles and wet under a weak stream of water.
3. Place the brush head in your mouth.

4. Press  to start the toothbrush. The indicator light for the most recently used mode will flash. While it is flashing, press  repeatedly to cycle modes. When the indicator light stops flashing, pressing  will instead switch off the toothbrush.
5. During cleaning, the brush stops about every 30 seconds. At this time, you should move it to a new position. After 4 such stops (2 in Polish mode), the toothbrush will automatically switch off.
6. Rinse your mouth and wash the brush head under running water to remove any residual paste and clean the head. Wipe the handle with a dry cloth and place the toothbrush in the storage case.

Note:

- To manually stop the toothbrush at any time, press .
- The toothbrush is equipped with a memory chip. When switched on, it will start in the most recently used mode.

When using an electric toothbrush, it is important to place the brush head against the teeth at the correct angle – about 45° to the gum line is ideal (Fig. **B**). That way, the bristles can reach a little below the gum line as well. There is no need to move the brush around as with a manual brush. Simply guide it slowly over your teeth to have them brushed one by one. Spend a few seconds on each tooth. Make sure to also clean the inner and biting faces. Proper technique helps prevent gingivitis and ensures effective plaque removal.

CLEANING AND MAINTENANCE

Turn off the toothbrush before any cleaning or maintenance. For sanitary reasons, clean the toothbrush after each use.

Warning:

Do not immerse the toothbrush or charging cable in water or other liquids.

Do not use plastic or metal scrapers, harsh cleaners, chemicals, petrol, solvents or other similar substances to clean any parts of the appliance. The surface finish could be seriously damaged.

Daily cleaning

After use, rinse the brush head under running water to clean it and remove paste residue. You may also remove it from the toothbrush. Then let it dry properly before putting it back on the shaft.

Clean the toothbrush with a dry cloth.

Storage

After use and cleaning, place the toothbrush in the storage case. Store that in a dry place away from children and pets.

Charge the toothbrush regularly (at least once every six months), even during prolonged disuse. This helps extend the battery service life.

Removing the Battery

When disposing of the product, the battery must first be removed. The appliance must be switched off during removal. The battery is recyclable.

To remove the battery:

1. Insert a flat head screwdriver through the hole in the cover at the bottom of the handle to release the cover.
2. Push the brush head shaft downwards and at the same time use a flat head screwdriver to release the lower motor buckle. Pull the motor unit out of the handle.
3. Cut the battery connection wire with scissors, then remove the battery using a flat head screwdriver.

**Warning:**

Do not damage or short-circuit the battery or throw it into a fire. Doing so risks leakage of harmful substances, short circuit, explosion, fire and other hazards.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alza.cz/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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© 2025 Alza.cz a.s. All Rights Reserved.