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User Manual

Safety Precautions

Please assemble and use the Air3 IPL Hair Removal Device strictly according to the instructions in this manual. Carefully read the User Manual and keep it for future reference. If you have any questions about using the product, contact Ulike customer service for assistance.

- DANGER: Personal injury due to incorrect use
- WARNING: Danger due to incorrect use
- CAUTION: Personal safety hazard / property loss due to incorrect use

Danger

- Do not use the device directly on your eyes; the strong light may cause eye injury.
- While the device is turned on and in use, avoid touching the side ventilation outlet as it may cause burns.
- Do not insert objects into the device, as this could result in personal harm.
- No data is available for usage by children under 18 years. Keep the device away from children under 18 and do not allow them to use it
- The device should not be used by individuals with physical, sensory, or mental impairments, or by those lacking experience and knowledge, unless supervised or instructed by someone responsible for their safety.

Warning

- Do not apply antiperspirant immediately after treating the armpit area.
- To prevent burns and mechanical issues, regularly clean the light window.
- Keep the light window away from items that may easily cause demagnetization to avoid device failure or damage.
- Do not use the device if the ventilation outlet is blocked, as this can lead to smoke or fire.
- Avoid pointing the light window at black objects, clothes, hair, or other items to prevent blackening, discoloration, or breakage.
- Do not use the product on pets.
- Failure to follow the User Manual may result in eye (blindness) or skin damage—read and follow instructions carefully.
- Do not overuse or reuse the device on one area to avoid skin harm.
- Protect your eyes; do not use the device above the eyelids or near the eyes.
- Test the device on your skin before use and check for abnormal reactions after 2 hours. Only use if the area appears normal.
- Wear the provided safety glasses during use and keep the light window close to the skin.
- After use, clean hair, dirt, and oil residue from the device to prevent malfunction and skin damage.
- Adjust the intensity level according to your skin type and start with Gentle Mode for the first use.
- Slight heat and redness after use is normal.
- Never leave the device unattended while on; turn off when not in use.
- If the device experiences a significant temperature change, turn it off and wait about 2 hours before reuse.
- Do not use near flammable or explosive materials.
- Only use the device at temperatures between 5-30°C; check room temperature before use.
- Avoid UV exposure to the treated area for 48 hours after use.
- Do not expose the device to sunlight or UV light for extended periods.
- If the built-in fan does not operate, do not use the device.
- Do not use for extended periods; turn off for 10 minutes after one hour of use.
- If the casing, lens, or power supply is damaged, stop using the device immediately and do not attempt repairs yourself; contact authorized personnel.

Unsuitable Populations and Situations

This device may cause discomfort or adverse effects for individuals with allergies, skin diseases, allergic dermatitis, or sensitive skin. Discontinue use if adverse effects occur. Do not use on damaged skin, such as burns, blisters, wounds, active skin diseases, or infections. If taking photosensitive medications, seek medical advice before use.

- Not suitable for pregnant, lactating, or menstruating women, or individuals with scars, diabetes, or fragile skin.
- Prohibited for people with severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin conditions, burned skin, or bleeding disorders.
- Do not use if you have abnormal collagen metabolism (including scarring and poor wound healing) or vascular diseases (e.g., varicosity, hemangiectasis) in the treatment area.
- Do not use if your skin is light-sensitive or prone to rashes/allergies.
- Do not use on infected areas, herpes simplex, skin lesions, or hematomas.
- Not suitable for those with immune suppression diseases (e.g., HIV, lupus, porphyria).
- Prohibited for individuals using medications such as isotretinoin, AHA, BHA, azalea acid, aspirin, immunosuppressants, tetracycline, sulfa, or recently injected with fillers/toxins.
- Use with caution if you have had chemical hair removal within 30 days or used skin whitening products/Vitamin A cosmetics recently.
- Do not use with essential oil products; ensure the skin in the treatment area is clean and dry.

Areas and Conditions Where Use Is Prohibited

1. Oral cavity, genitals, mucosa, etc.
2. Acnes
3. Inside the navel
4. Plastic surgery or other surgical procedure
5. Spots and moles
6. Large vein
7. Throat
8. Eyebrow, hair
9. Tattoo, body painting
10. Eyeball, eyelid, around the eyes
11. Ear
12. Head
13. Nipple, areola

(Note: This list is not exhaustive. Consult a doctor if unsure.)


Possible Side Effects After Use

1. Burns
2. Swollen hair follicles and blisters
3. Itchiness and pain
4. Local redness and swelling
5. Erythema and pigmentation
6. Scarring

(Note: This list is not exhaustive. If you experience any discomfort, consult Ulike customer service or a doctor.)

Other Notices

- This product is not intended to treat diseases.
- Do not use to remove white or gray hair.
- Not suitable for naturally dark skin.
- Avoid sun exposure before use.
- On first use, test on a small area in Gentle Mode.
- Ensure the light window is close to the skin; do not use in dark conditions or direct sunlight.
- Do not hit or shake the device violently. Turn off before unplugging.
- When finished, position the ventilation outlet upwards and the treatment window downwards for cooling.

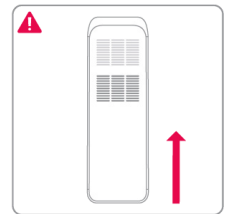Product Introduction

About the Ulike IPL Hair Removal Device

The Air3 IPL Hair Removal Device features clinical sapphire contact cooling technology with Ulike's patented sapphire ice-cooling technology. This system effectively cools the skin, reduces heat, prevents burns, and protects the skin barrier during use. The sapphire flat tip compresses the skin, bringing hair follicles closer to the energy source for improved absorption, resulting in better hair removal. The flat tip is also easy to clean and ensures effective light transmission. With a long lamp lifetime and versatile application (face, lips, underarms, bikini, arms, legs), the device eliminates the need for replacement, saving on accessory costs. It efficiently filters out ineffective wavelengths and harmful light using a top-grade filter.

Product Structure and Accessories






In order of appearance below:

- Host device
- Power adaptor and power cord
- IPL safety glasses
- Razor





Key components include the

1. sapphire light window (1)
2. indicator lights (2) (for auto-glide, cooling system, skin sensor)
3. ice cooling indicator (3)
4. auto-glide indicator (4)
5. treatment button (5)
6. power/intensity adjustment button (6)
7. on/off button (7)
8. power input interface (8)

Air 3 Series

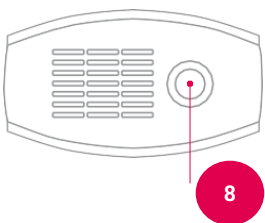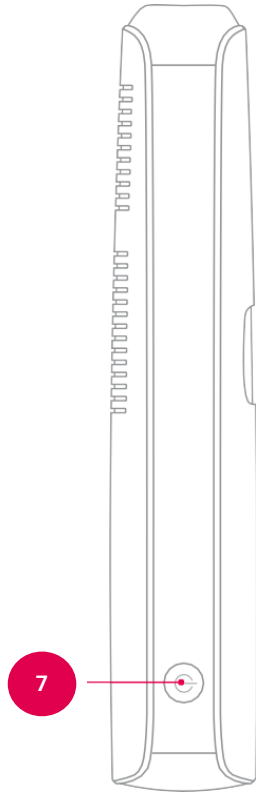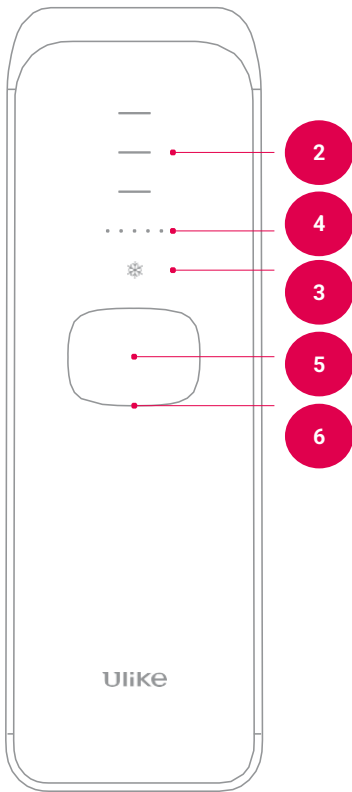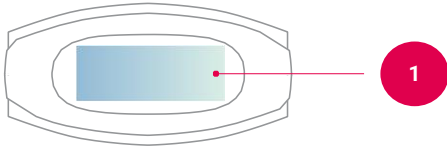
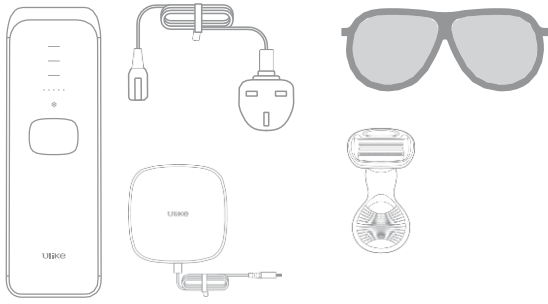
Air 2 Series

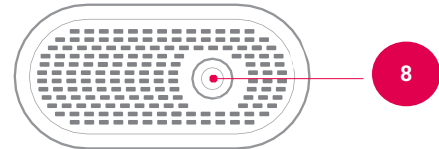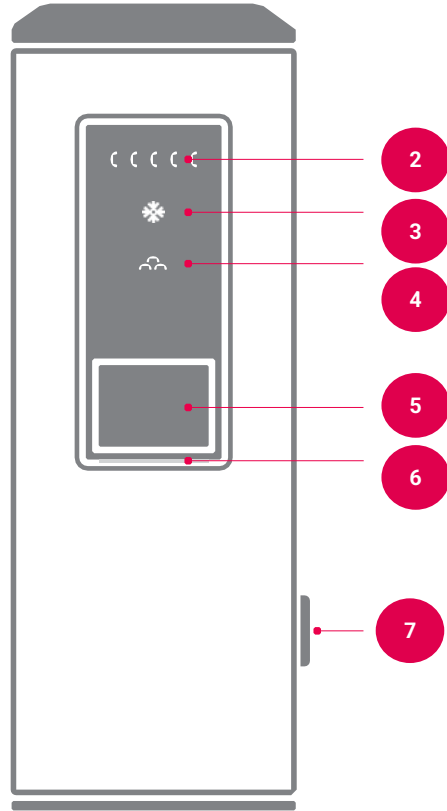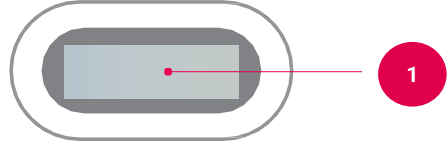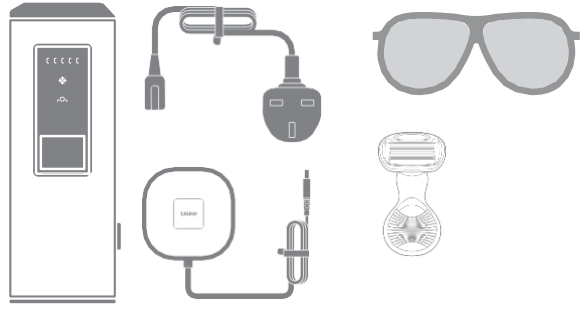Before Use

1. Clean your skin thoroughly. Ensure no hair, oil, or cosmetic residue is left and the skin is completely dry.
2. Make sure the exposure port is clear of dust, hair, and dirt, and is not damaged.
3. Perform a skin test before first use. Check the area after 2 hours; if there are no abnormal signs, proceed with use.

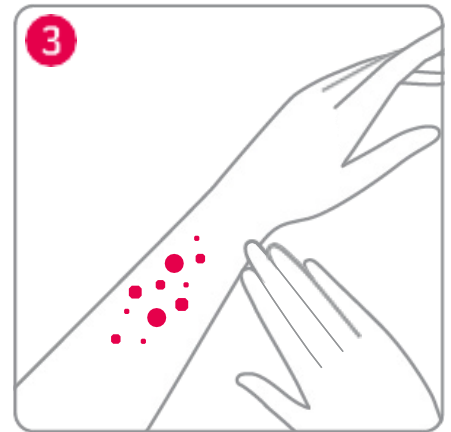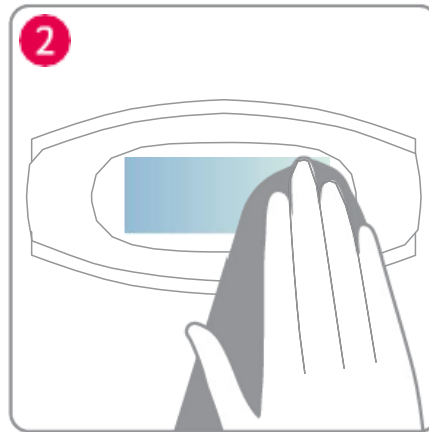Setting Exposure Intensity

In case of Air 3, the device offers three energy modes to accommodate different skin sensitivities and hair types:

- Soft Mode: For sensitive areas (lips), people with dark or sensitive skin.
- Body Mode: For medium or low hardness thick hair (arms, legs, etc.).
- Strong Mode: For thick and hard hair (arms, legs, etc.); may feel hotter but is normal.


In case of Air 2, The product has five gears of energy adjustment, which can meet the needs of different skin sensitivity and different hair colors. The following are the recommended instructions.

In case of skin colour 4, please use the 1-3 setting of Ulike Air 2, and the Soft mode / Body mode of Air 3.

If you feel pain or discomfort during use, please lower the energy gear (if you feel warm or slight burning without pain, it is normal);

People with darker skin are prone to skin allergies, so it is recommended to use it with a lower energy gear;

For people with thick hair, it is recommended to use a higher energy gear (4/5 or strong mode); For people with little hair, it is recommended to use the lower gear (1-3 or soft mode). The specific situation varies from person to person, but for the first time, Please use it from Mode 1 or Soft mode.

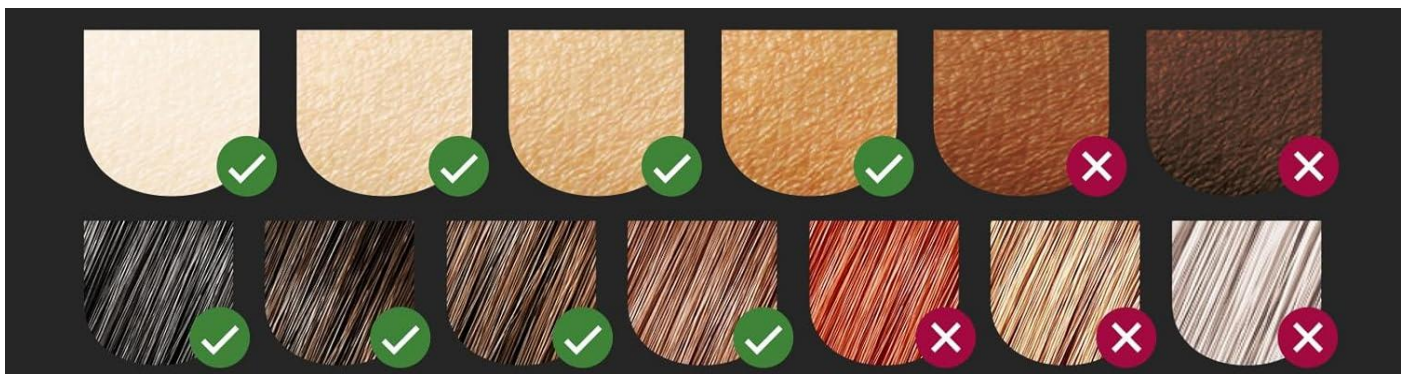
How to Use

1. Shave hair from the treatment area with care, using lather or shaving during a shower to avoid scratches.
2. Connect the power chord to the adapter and then to the device's power input. Plug into a power outlet.
3. Short press the Power On/Off button to turn on the device and display energy level and cooling status.
4. Switch between the modes by short pressing the Power/Mode button.
5. Wear safety glasses before starting hair removal.
6. Place the light outlet on the treatment area, ensuring full contact with the skin so the sensor lamp lights up. Press the flash key to start hair removal. If the lamp is not lit, ensure proper contact.
7. After finishing, long press the Power On/Off button to turn off the device (do not unplug without turning off first).

Treatment Methods


- Auto-Glide Method (for large areas): Place the light outlet on the skin and maintain full contact. The sensor indicator will light up. Long press the flashing button for 1.2 seconds to activate Smart Continuous Play Mode. Move the device back and forth twice over the area for even coverage. Adjust your movement to match the flash frequency.
- Spot Method (for small areas): Place the light outlet on the skin. The sensor indicator lights up. Short press the flash button, then move to the next area and repeat.

Tip: If you feel excessive heat, pause and let the device cool before continuing.

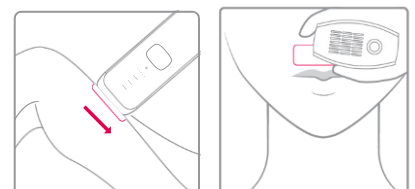Tips for Using Ulike Air Series

- The device heats up due to continuous flashes. For optimal cooling, use the light window on a treated area for about 30 seconds before moving to the next area.
- Press the light window tightly against the skin for better cooling and protection.

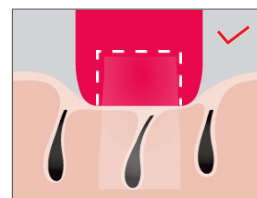Usage Frequency

- Use three times a week for weeks 1–4.
- From week 5 onward, use every half month or month as needed based on hair growth.
- Increase frequency/gear for thicker hair as needed.

Care After Using

- Avoid water contact with the treated area for 6 hours post-treatment; moisturize after that.
- Do not use perfumes, facial masks, or exfoliating cosmetics on the treated area after hair removal.
- Avoid UV exposure for 48 hours following hair removal and apply sunscreen if necessary.
- Artificial tanning is prohibited during treatment.
- Consult a doctor if you plan surgery on a treated area.
- If discomfort occurs, stop using the product and consult a doctor.

Cleaning and Maintenance

- Clean the device before and after each use to prevent hair and dust buildup.
- Always turn off, unplug, and let the device cool before cleaning.
- Use a soft, dry cloth to clean the light window. Store in a dry, dust-free environment between -10 and 55°C.
- Do not wash the device or accessories under water or in the dishwasher.
- Avoid abrasive cleaners or corrosive liquids for cleaning.
- If necessary, use a few drops of high-concentration alcohol (below 95%) for cleaning.

Environmentally Safe Disposal

Dispose of the device at a designated recycling center; do not include with regular household waste.

FAQ

This section addresses common issues that may arise during use. If issues persist, contact Ulike customer service. Keep your invoice and warranty card for service.

Compliance

This device complies with the following standards: EN 55014-1, EN 55014-2, EN 61000-3-2, EN 61000-3-3, EN 62233, EN 60335-1, EN 60335-2-23, IEC 60335-2-113, IEC 62471, BS EN 60335-2-23, BS EN 60335-1, BS EN 62233, BS EN 55014-1, BS EN 55014-2, BS EN IEC 61000-3-2, BS EN 61000-3-3.

Faults and Solutions

Problem	Possible Reason	Solution
The product will not turn on	Not connected to power	Check that it is plugged in
The product does not emit light	Main unit head too hot or not in full contact with skin	Let device cool, ensure full contact and check indicator light
The working indicator light is not on	Product is defective or malfunctioning	Check power supply, repair or return the product
Odor when product flashes	Light filtration mirror contaminated by hair	Clean mirror, shave area, reduce intensity
Unbearable pain during hair removal	Hair not shaven in area	Shave area before use, use lower intensity if needed
Area not suitable for treatment	Refer to guidelines on restricted areas	Do not use on restricted areas; men should not use on face/genitalia
Intensity level too high		Reduce intensity until comfortable
Light window is damaged		Contact Ulike customer service for repair or return
Redness after hair removal	Slight redness is normal	Use aloe vera gel or cool with sapphire lamp head
Unsatisfactory hair removal	Intensity too low	Select higher intensity
Did not press flash button on next area		Ensure to press flash button and overlap treated areas
No effect on certain hair colors	Light blond, gray, red, or white hair	Irradiation not effective for these colors
Device not used frequently enough		Follow recommended regimen for best results
Body's response slower than expected		Use device continuously for at least 5 months
Product becomes hot and humid	Extended use	Let device cool down before reuse

Basic Parameters

Air Series	Air 3 Specification	Air 2 Specification
Product Name	IPL Hair Removal Device	IPL Hair Removal Device
Product Model	UI06	UI04
Product Size	179*58*34mm	170*61*38mm
Rated Input	24V / 2.8A	24V / 2.8A
Power	67.2W	67.2W
Energy Level	3 levels	5 levels
Spot Size	3.3cm ² ±0.25cm ²	3.3cm ² ±0.25cm ²
Spectrum	550-1200nm	550-1200nm
Pulse Duration	1.12ms-7.44ms	0.7ms-1.7ms
Pulse Repetition Rate	0.33Hz-2.5Hz	0.60Hz-1.40Hz
Manufacturer	Shenzhen Ulike Smart Electronics Co., Ltd., 810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen, 518000, Guangdong, P.R. China	Shenzhen Ulike Smart Electronics Co., Ltd., 810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen, 518000, Guangdong, P.R. China

Storage and Operating Conditions

Storage Temperature	-10–55°C	Operating Temperature	5–30°C
Humidity	≤90% RH	Humidity	≤80% RH
Relative Pressure	50–106 kPa	Relative Pressure	80–106 kPa