

SENCOR®

SOI 1010SL



ORAL IRRIGATOR
Translation of the original manual

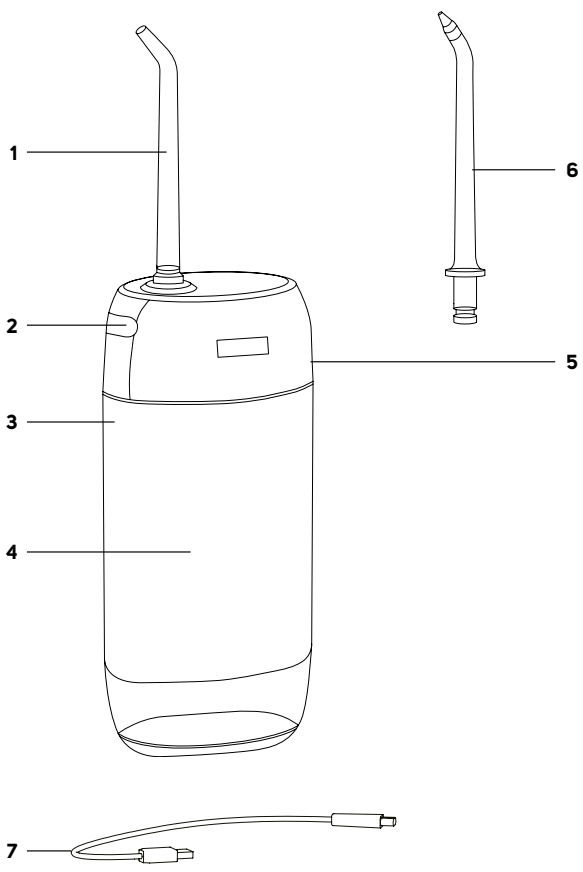


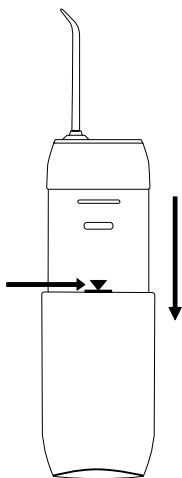
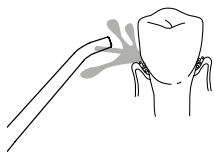
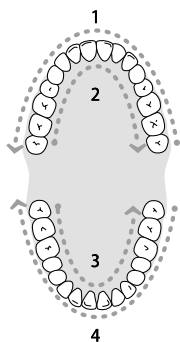
SENCOR®

SOI 1010SL



A



B**C**

IMPORTANT SAFETY INSTRUCTIONS

READ CAREFULLY AND STORE FOR FUTURE USE.

- Oral irrigators may be used by children 8 years of age and older and by persons with physical, sensory or mental impairments or by inexperienced persons, if they are properly supervised.
- Children must not play with the appliance. Cleaning and maintenance must not be performed by unsupervised children.
- If the charging cable is damaged, have it replaced in a professional service centre in order to prevent a dangerous situation. It is forbidden to use the appliance with the charging cable damaged.
- Oral irrigators must only be used with cold or hot water or with special solutions defined in this manual.
- Parts of the appliance that need to be fixed or erected must be fixed or erected in such a way that they cannot fall into the water.



Caution:

To avoid the risk of a sudden resetting of the thermal fuse, this appliance must not be supplied via an external switching device such as a timer or be connected to a circuit with regular on/off switching.

- The appliance may only be powered using a safe small voltage corresponding to the designation on the rating label of the appliance.
 - The power adapter is not supplied.
-

Electrical Safety

- Do not place the power cord over sharp objects. Make sure that the cable does not hang over the edge of a table or touch a hot surface.
- Never use the charging cable supplied with this appliance for other purposes.
- Keep the charging cable dry.
- Do not wind the charging cable around the handle of the appliance.
- Do not connect or disconnect the charging cable with wet hands.
- Before connecting the appliance to a power socket, check that the nameplate rated voltage matches the voltage of the power socket the appliance is to be connected to. Do not use an extension cord.
- Store, charge and use the appliance at an ambient temperature of 15 °C to 35 °C. Do not charge the appliance for longer than necessary to fully charge it. After charging, disconnect it from the power source.

- Do not expose the appliance to temperatures above 50°C. Otherwise damage to the built-in battery may occur.
- The appliance is equipped with a built-in battery with a long service life, which can only be removed by a professional service. Never attempt to remove it from the appliance a switch for a new one by yourself.



Warning:

Never disassemble the appliance in order to remove the battery. Do not short-circuit the appliance or throw it into a fire. There is a risk of leakage of toxic substances, fire or explosion if improperly handled.

- Do not charge, use, store or leave the appliance in a bathroom or any place where it could fall into a bathtub, wash basin or any other container filled with water.
- If the appliance was connected to a mains socket via power adapter and it falls into water, under no circumstances should you touch it, but immediately disconnect the power adapter from the socket.

Safety in Use

- This oral irrigator is designed exclusively for the cleaning of teeth and the oral hygiene of people. The oral irrigator is not intended for commercial purposes or for brushing animals' teeth.
- Use the oral irrigator and its individual attachments in accordance with the instructions in this manual.
- If you have recently undergone intensive surgery such as tooth or gum surgery, have had long-term dental problems, or if you suspect dental problems (fillings, crowns, etc.), consult your dentist before using the oral irrigator.
- If you experience bleeding gums after using the oral irrigator and if the bleeding persists after 2 weeks of use, contact your dentist.
- In case of any health concerns, contact your dentist before using the oral irrigator.
- If any damage to the oral irrigator (handle or mains adapter) occurs, stop using it.
- Use only the original accessories supplied with this oral irrigator. Do not use attachments other than those supplied by the manufacturer.
- The oral irrigator and its power adapter must only be used indoors. Do not use them outdoors.
- Avoid using the oral irrigator in rooms where aerosol dispensers are used or where oxygen is administered.
- Do not expose the oral irrigator and its accessories to the outdoor environment, extreme temperatures, direct sunlight, excessive humidity or excessive dust.
- Do not place the appliance and its accessories near heaters, open fire or other sources of heat.



Warning:

Keep the oral irrigator dry.

- Before each use, make sure that the oral irrigator and its accessories are not damaged. Never use an oral irrigator with a damaged attachment.
- The oral irrigator may heat up when used or recharged. This is completely normal.
- Do not apply too much pressure to the end of the nozzle attachment.
- To avoid injury, do not fill the tank with hot water. Do not insert any objects into the tank.
- Do not point the nozzle of the attachment under the tongue, in the ear, nose, eyes or other sensitive areas, as this could cause injury.

- Do not clean the nozzle attachments, the handle, the mains adapter or the water tank in the dishwasher or microwave.
- Always switch off and unplug the oral irrigator when not in use, when left unattended and before attaching, removing or replacing the attachment, cleaning, maintenance or storage. Store the oral irrigator and its mains adapter in a dry environment.
- Keep the oral irrigator clean. Perform the cleaning and maintenance following the instructions in the Cleaning and Maintenance section.
- Do not perform any maintenance on the oral irrigator other than cleaning its individual components as described in the Cleaning and Maintenance chapter.
- Do not use the oral irrigator if it has been dropped on the floor, does not function properly or is damaged in any way.
- To avoid a hazardous situation, do not repair or modify the oral irrigator or its accessories yourself. Contact an authorised service centre for all repairs. Tampering with the appliance may result in the loss of your legal rights regarding faulty performance or warranty for quality.

Special Safety Instructions for the Li-ion Battery

- The Li-ion battery is built into the oral irrigator and is an integral part of it, which the user is not allowed to remove.
- Follow all the instructions provided in this user's manual.
- In the event of a battery fluid (electrolyte) leak, do not touch the battery with your bare hands. Use protective aids (gloves, a protective mask, goggles, etc.). The liquid inside the battery is toxic and may cause health problems even with slight contact of the bare skin with it.
- We recommend opening windows or doors to ensure adequate ventilation of the space.
- Contact a specialist hazardous waste disposal company or use a hazardous waste collection facility. Do not dispose of the damaged battery or accumulator in ordinary household waste.
- In case of skin or eye splashes, rinse with plenty of water and seek medical attention. If ingestion occurs, seek immediate professional medical help.

EXPLANATION OF SYMBOLS PLACED ON THE PRODUCT OR IN THE ACCOMPANYING DOCUMENTATION



The product meets all the basic requirements of the applicable EU directives.



This symbol on products or accompanying documents means that used electric or electronic products must not be added to ordinary municipal waste.



This symbol on the accumulator or battery means that the used battery must not be added to normal municipal waste.

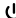
Protection level against harmful water ingress:

IPX7 – The device is protected against harmful ingress of water when submerged to a depth of 1 meter for a period of 30 minutes.

EN Oral irrigator User Manual

- Thank you for purchasing our SENCOR product; we hope it will serve to your satisfaction.
- Prior to using this appliance, please read the user manual thoroughly, even in cases when you are already familiar with the use of similar types of appliances. Use the appliance only as described in this user manual. Keep the manual for future reference. If you pass the appliance on to another person, ensure that this user manual is included.
- Carefully unpack the appliance and make sure not to discard any part of the packaging material until you have removed all its components. It is recommended that you keep the original packaging, packing material, receipt and confirmation of the extent of the seller's liability or warranty certificate at least for the duration of the legal right to file a defect claim or the quality warranty. When transporting the appliance, we recommend repackaging it in the original box provided by the manufacturer.


DESCRIPTION OF THE ORAL IRRIGATOR

- | | |
|---|---|
| A1 Standard attachment with nozzle (1 pc) | A5 Charging socket with cover |
| A2 Nozzle attachment release button | A6 Orthodontic attachment (1 pc) |
| A3  On/Off button | A7 Charging cable |
| A4 Water tank | |

INTENDED USE

- The oral irrigator is designed only for oral hygiene. It helps to remove bacteria from the interdental spaces and the folds of the gum around the teeth, quickly and gently removes food debris from the interdental spaces, while effectively massaging your gums. In combination with careful brushing, it helps you to actively care for the health of your teeth. The oral irrigator is intended for household use and is not designed for commercial or other professional applications.

BEFORE FIRST USE

- Wash the appliance under running water.
- Rinse the water tank above the sink several times and fill it with clean water according to chapter Filling the Water Tank. Insert the attachment into the oral irrigator according to the chapter Attaching and Removing the Attachment. Point the nozzle of the attachment into the sink and turn on the oral irrigator with the  button. Select one of the modes and let the water flow through the nozzle for a while to flush both the suction tube and the nozzle of the attachment.
- Pour the remaining water out of the tank and leave the filling hole cover open to allow the tank to dry out.
- Wipe the oral irrigator with a dry cloth.

CHARGING THE ORAL IRRIGATOR

- The oral irrigator is equipped with a built-in Lithium-ion battery that serves as an internal power source and allows cordless operation of the oral shower. Allow the battery of the oral irrigator to fully charge first before first use or after a long period of non-use. It takes about 3.5 hours to fully charge the battery. Use only the original supplied charging cable to charge.

- Remove the charging socket cover and connect the connector of the charging cable to it. Plug the loose terminal into a suitable mains adapter (not included).
- While the oral irrigator is charging, the charging indicator flashes. When it stops flashing and turns off, the oral irrigator is fully charged.
- Disconnect the cable from the power socket and the charging source. Slide the cover back into the slot.
- If the charging indicator flashes during operation of the oral irrigator, this means that the oral irrigator is almost discharged and needs to be recharged.



Caution:

There must be no moisture around the charging socket or its cover. Therefore, it is recommended to wipe the area with a dry cloth before plugging in the charging cable.



Caution:

Do not leave the oral irrigator connected to an external power source longer than it takes to fully charge it.



Note:

During charging, the use of the oral irrigator is blocked and cannot be started.

- Recommendations on how to achieve the longest possible battery life:
 - To keep the built-in battery working for as long as possible, charge, use and store the oral irrigator at an ambient temperature of 15–35 °C in a dry environment. Lower or higher temperatures adversely affect the life of the battery.
 - Do not wait until the battery is completely discharged. It is better to recharge the battery before it is fully discharged.
 - If you do not use the oral irrigator for a long time, charge the battery at least once every 6 months.

FILLING THE WATER TANK

1. Open the tank filler cap.
2. Pull the water tank downwards until the water level indicator appears – see figure B.
3. Hold the oral irrigator in a horizontal position and pour clean water into the tank.
4. Close the cover properly and then you may return the oral irrigator to a vertical position.



Caution:

The water should be at a comfortable temperature. Do not use water hotter than 40 °C.



Note:


Open the filler cap **BEFORE** pulling the water tank downwards. Otherwise it may be difficult to pull the tank down due to overpressure.

Tip:

- In addition to clean water, you can also use mouthwash or its solution to fill the water tank. Using a mouthwash makes using oral irrigator even more effective. However, do not use mouthwashes containing isopropyl myristate or a high concentration of essential oils. These substances are included in the mouthwashes sold in glass and ceramic containers and could damage the

water tank as well as the internal parts of the oral irrigator. Mouthwashes containing substances harmless to the oral irrigator are sold in plastic containers.

ATTACHING AND REMOVING THE ATTACHMENT

- Insert the attachment perpendicularly into the opening at the top of the oral irrigator until you hear a click. Hold the oral irrigator so that the  button is not pressed.
- To remove the attachment, press and hold button **A2** and pull the attachment straight up out of the opening.






Caution:

Do not remove the attachment while the oral irrigator is in operation.

OVERVIEW OF OPERATING MODES

- The oral irrigator is equipped with three operating modes:
 - **Pulse mode:** water flows out in short intervals, which is ideal for cleaning and massaging teeth and gums;
 - **SOFT mode:** low water flow for gentle cleaning; suitable for people with sensitive teeth and those using an oral irrigator for the first time;
 - **Normal mode:** suitable for everyday use.
-

USE OF THE ORAL IRRIGATOR (FIG. C)

- Check that the oral irrigator is properly assembled, charged and the water tank is full.
 - Brush your teeth and gums and rinse your mouth with clean water.
1. Slide the attachment into the open mouth and lean over the sink so that water and rinsed food remains can flow out.
 2. Turn on the oral irrigator with the  button.
 3. To change the mode, press the  button repeatedly within 5 seconds of switching on the oral irrigator. The indicator light for the selected mode will illuminate.
 4. Direct the stream of water to the place you want to spray, i.e. to the interdental spaces, the edges of your teeth or to your gums.
 5. Move the attachment nozzle slowly from one interdental space to another. In this way, clean all the outer and inner sides.
 6. Thoroughly spray around braces, crowns and bridges.
 7. When massaging the gums, aim the nozzle perpendicularly to the gums and spray the gum line.
 8. The oral irrigator is equipped with a timer that automatically turns it off after 1.5 minute of operation. If you need to switch off the oral irrigator earlier, press the  button.
 9. Remove the nozzle attachment from the mouth when the oral irrigator is switched off.




Caution:

Make sure the attachment is not damaged before each use. Never use an oral irrigator with a damaged attachment.



Note:

The oral irrigator will start in the mode in which it was used for more than 10 seconds before you last turned it off.

- As soon as no water comes out of the nozzle, the water tank is empty. Turn off the oral irrigator with  button and add water following the instructions in the "FILLING THE WATER TANK" section. When the oral irrigator is turned on again, it will automatically start in the last selected mode.
- When used for the first time, it's recommended to use the **Soft** operating mode.

**Caution:**

During the first few uses of the oral irrigator, you may notice slight bleeding from the gums. This is normal and usually disappears after 1 to 2 weeks.

**Caution:**

Consult the use of an oral irrigator with your attending physician if you notice heavy bleeding from the gums, have undergone a medical procedure in the oral cavity, have an injured gum or an open wound in the oral cavity.

**Caution:**

Do not use an oral irrigator if it causes you pain, discomfort, dizziness or vomiting.

CLEANING AND MAINTENANCE

**Caution:**

To clean the oral irrigator, use water up to 50 °C or neutral cleaning agent. Do not use wire cloths, abrasive cleaners, thinners, etc. for cleaning to prevent the cleaned surface from being damaged. When cleaning, do not immerse the oral irrigator in water.

WATER TANK

- Open the tank filler cap and pour out the remaining water from the tank. Close the cap.
- Start the oral irrigator to remove residual water from the tank and hose. When water stops flowing out, switch off the oral irrigator.
- Open the tank filler cap and allow the tank to dry out.

ORAL IRRIGATOR

- Wipe the oral irrigator with a clean, dry cloth. In the event of heavy soiling, you can use a slightly damp cloth. Clean the surface then properly wipe dry.

**Caution:**

The upper part of the oral irrigator in which the battery is stored is not intended for washing under the tap. Do not spray the oral irrigator and the charging cable with water or any other liquid or submerge it in water or any another liquid.

ATTACHMENTS

- Rinse the used attachments with water after each use to avoid the formation of deposits. Then wash the surface of the attachment with water and then always wipe it dry with a soft cloth. For hygiene reasons, we recommend replacing the used attachments with new ones every six months.
-

REMOVING CALCIUM DEPOSITS

- Depending on the mineral content of the used water, hard water deposits may accumulate in the tank, in the suction hose or extension and in the nozzle. These can negatively affect the operation of the oral irrigator and reduce its performance. It is therefore important to remove deposits regularly.
 - To remove mineral deposits, proceed as follows:
 1. Empty the water tank and pour clean warm water into it (max. temperature 85 °C).
 2. Point the attachment with the nozzle into the sink and start the oral irrigator.
 3. Let the contents of the tank be emptied completely.
 4. Repeat the procedure if needed.
 - Depending on the degree of hardness of the water you use, repeat the deposits removal process every 1–3 months to maintain the performance of your oral irrigator.
-

STORAGE

- After use, prepare the oral irrigator as instructed in the CLEANING AND MAINTENANCE chapter. Store the oral irrigator in a dry place away from direct sunlight. If you want to put the oral irrigator away for a long time, we recommend removing the water tank from the handle.

Troubleshooting

Problem	Possible cause	Solution
Water escapes from the space between the end of the attachment and the body of the oral irrigator.	The sleeve gasket is damaged or worn.	Have the gasket replaced or replace the nozzle.
	The attachment is not correctly fitted.	Remove the attachment and reattach it.
Water is leaking from the tank.	The filler cap is not properly closed.	Close the filler cap.
Water pressure is too low.	The nozzle is clogged with dirt.	Please clean or replace the nozzle.
	The nozzle is deformed.	Change the nozzle.
	The oral irrigator is not sufficiently charged.	Allow the oral irrigator to charge.
The oral irrigator does not work when unpacked.	The oral irrigator is not sufficiently charged.	Allow the oral irrigator to charge.
Water is not splashing out of the oral irrigator.	The water tank is empty.	Please refill water.
	The irrigator is too tilted when used.	Keep the irrigator in a vertical position as much as possible.

.....

TECHNICAL SPECIFICATIONS

Oral irrigator

Dimensions incl. attachment (H x D x W).....	135 x 45 x 70 mm
Weight incl. attachment.....	235 g
Tank volume.....	120 ml
Pump frequency.....	1,600 pulses/min
Water pressure.....	30–90 PSI
Charging time.....	3,5 hours
Operating time.....	40 minutes
Protection level against harmful water ingress.....	IPX7

Built-in battery

Type.....	Li-ion 1,100 mAh, 3.7 A
-----------	-------------------------

.....

BATTERY DISPOSAL

- Before disposing of the product or returning it to a designated collection point, ensure that the battery is removed. The battery contains environmentally harmful substances. Proper disposal shall be ensured by the relevant collection point.
-

INSTRUCTIONS AND INFORMATION ON DISPOSAL OF USED PACKAGING MATERIAL

Dispose of used packaging material at the designated municipal waste disposal site.

.....

DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT

Used electrical and electronic products must not be disposed of in ordinary municipal waste. To enable the proper disposal, renewal and recycling of these products, deliver them to the designated collection points. Alternatively, in some European Union states or other European countries the products can be returned to the local retailer when buying an equivalent new product. By disposing of this product correctly, you help to conserve valuable natural resources and prevent any potential adverse impact on the environment and human health that could result from improper waste disposal. For more details, contact your local authority or nearest collection point. Fines may be imposed for the improper disposal of this type of waste as per national regulations.

For business entities in European Union States

If you wish to dispose of electrical and electronic equipment, please request the necessary information from your retailer or supplier.

Disposal in other countries outside the European Union

To dispose of this product elsewhere, request the necessary information on the correct disposal method from local authorities or your retailer.

.....

Changes to text and technical parameters are reserved.