SIGURO

Massage Shawl

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


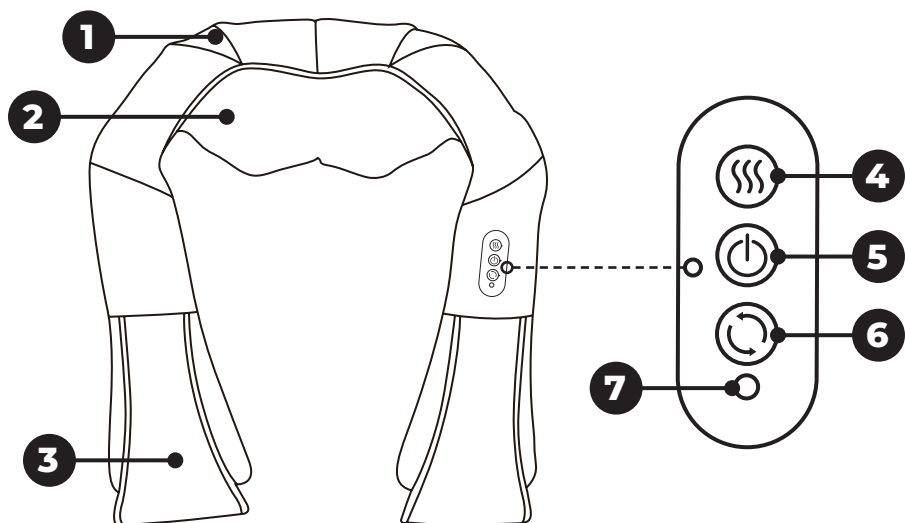
A

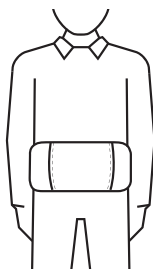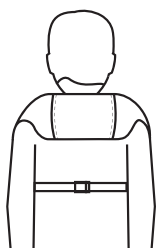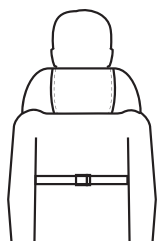Thank you for purchasing a SIGURO appliance. We are grateful for your trust and are happy to introduce on the following pages the device, its functions and uses.

We believe in a fair and responsible society and only work with suppliers who meet our strict criteria on employee protection, exploitation prevention and fair working conditions.

If you need help with extensive maintenance or repair of the product that requires intervention in its internal parts, our authorized service centre is available at siguro@alza.cz. You may also contact the vendor's helpline. For your convenience in solving potential issues with the product, the unified contact points above can also be used for any complaints or post-warranty service.

The current version of this manual is available at www.siguro.net.

SAFETY INFORMATION

Please read this manual carefully before use and keep it for future reference.

1. Children aged at least 8, persons with physical, sensory or mental impairment, and the untrained or inexperienced may use this appliance with instruction or supervision in its safe use, provided they understand the potential dangers.
2. Children are not allowed to play with the product.
3. Cleaning and user maintenance may not be performed by unsupervised children.
4. If the power adapter is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons to avoid danger.
5. This product has a heated surface. Persons with impaired heat sensitivity should take care when using this appliance.
6. Warning: Incorrect use of the product may lead to damage to the product or serious injury.
7. Use only the adapter supplied with the product. Do not use different adapters. There is a risk of electric shock or property damage.
8. Do not use the supplied adapter to power another product.
9. Use the product as described in this manual.
10. This product is intended for body/neck and shoulder massage. Do not use it for any other purpose. It is not intended for use in a commercial environment or for laboratory purposes.
11. The product is designed for indoor use. Do not use it outdoors.
12. Do not place the product in an environment with high temperature or

humidity or in direct sunlight. Do not use the product at temperatures above 40 °C.

13. Avoid using sharp objects.
14. Before use, make sure the appliance is correctly assembled. Do not use the product if it or the woven foot pads are damaged.
15. Use only the accessories supplied with the appliance.
16. Do not use the product when your feet are wet. Make sure your feet are dry before using the product. Do not use massage cream.
17. Turn off and unplug the product after use and before cleaning, storage or moving.
18. Do not cover the product while in use. Otherwise, overheating and damage may occur.
19. Do not disassemble or modify the product. It may become damaged, create a dangerous situation (e.g. fire) or cause personal injury. If necessary, contact an authorized service centre.
20. Stop using the product if it makes you feel unwell, dizzy or uncomfortable.
21. Do not touch the product with wet or damp hands, especially when plugging or unplugging it.
22. Consult your physician before using the product if:
 - your legs are injured, post-surgery, or you have another leg condition;
 - you have high blood pressure, heart disease, cerebral arteriosclerosis, cerebral thrombosis, malignant tumour, or osteoporosis;
 - you are pregnant;
 - you use a pacemaker or other medical aids.
23. Do not use the appliance continuously for more than 30 minutes.
24. Do not use the product while sleeping.
25. Use only dry cloth for cleaning. Follow the instructions in this manual.

Special instructions for rechargeable Li-ion battery

1. The massage collar is equipped with a built-in rechargeable Li-ion battery.
2. Do not expose the massage collar to direct sunlight, frost or high temperatures (e.g. car parked in the sun).
3. Keep the massage collar in a dry place. Avoid moisture and water contact.
4. Check appliance condition regularly. If you notice bulging, odour, unusual heat or other unusual signs, stop using the massage collar immediately.
5. Avoid mechanical damage. Do not puncture, squeeze or deform the battery.
6. Do not leave the massage collar unattended when charging.
7. Do not use damaged or old batteries.
8. Do not disassemble the battery.
9. Do not expose the battery to fire or high pressure. Do not throw batteries in the fire.
10. Charging and discharging outside the recommended limits (e.g. full discharge) may irreversibly damage the battery.


11. Do not use or store the massage collar where flammables, volatiles or explosives are used or stored.
12. Avoid contact with leaked electrolyte. Wear protective gloves, goggles, and possibly a breather mask. Place the massage collar on a non-flammable surface, ideally outdoors or in a well-ventilated area. Do not handle the battery directly with bare hands. Contact a professional battery disposal service or hazardous waste collection facility. Never dispose of damaged batteries in household waste. In case of skin or eye contact, rinse with plenty of water and seek medical attention. In case of ingestion, seek medical attention immediately.

The “interior use only” symbol means the appliance is not designed for outdoor use.

Use only with the adapter marked IVP0900-2000G.

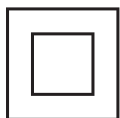Protection class II equipment

Do not use if the adapter plug is damaged.

OVERVIEW

- 1 Massage collar
 - 2 Massage head
 - 3 Straps
 - 4  Heating on / off button
 - 5  Massage collar on / off button
 - 6  Massage head rotation direction selection button
 - 7 Charging indicator light
- Charging adapter (not pictured)

TECHNICAL SPECIFICATIONS

Massage collar

Rated power	18 W
Input voltage	7.4 V ---
Battery	Li-Ion 2200 mAh
Charging time	3.5 hours
Operating time per charge	90 minutes
Charging adapter model	IVP0900-2000G
Dimensions	400 × 160 × 170 mm
Weight	1.63 kg (excluding the adapter)

Adapter

Manufacturer name or trademark, trade registration number and address	Alza.cz a.s. Jankovcova 1522/53 170 00 Prague 7, Czech Republic
Product identifier	Model: IVP0900-2000G
Input voltage	100–240 V~
Input frequency	50/60 Hz

Output voltage	9.0 V
Output current	2.0 A
Output power	18 W
Mean efficiency in active mode	85.31%
Efficiency at low load (10%)	85.14%
Standby power consumption	0.062 W

BEFORE FIRST USE

Remove the massage collar and its accessories from the packaging. Before disposing of packaging, make sure no components are left inside. Check the massage collar and accessories for damage.

Wipe the outer surface of the massage collar with a soft and dry sponge.

Method of use

This product is designed to massage parts of the body (e.g. neck, shoulders, thighs, etc.).

You can find examples of use in Figure A.

Charging

- The massage collar must be charged before first use.
- Insert the charging adapter connector into the charging socket of the massage collar. Plug the charging adaptor into a properly earthed mains socket. The charging indicator will light up red, indicating that the massage collar is charging.
- Charging time is usually 3.5 hours.
- When the massage collar is fully charged, the charging indicator will light up green. Disconnect the charging adapter connector from the charging socket and unplug the adapter from the mains socket.

Warning:

The massage collar cannot be used while it is charging.
Only use the supplied charging adaptor to charge the device.

USE

1. Put on the massage collar.
2. Press the  button to switch on the massage collar.
3. Press the  button repeatedly to select the rotation direction of the massage heads.
4. Press the  button to switch on the heating. Press it again to switch it off.
5. The massage collar switches off automatically after 15 minutes of continuous use. If you need to switch it off sooner, press the  button.

Once you have finished, clean the massage collar as described in the 'Cleaning and Maintenance' section.

Warning:

We recommend not using the massage collar for longer than 15 minutes. Do not use it for longer periods.

CLEANING AND MAINTENANCE

Always switch off and unplug the appliance before cleaning or any maintenance.

Warning:

Do not immerse the appliance, adapter, or power cord in water or other liquids.

Warning:

Do not use plastic or metal scrapers, harsh cleaners, chemicals, petrol, solvents, or other similar substances to clean any parts of the appliance. The surface finish could be seriously damaged.

Outer surface

Wipe the outer surface of the massage collar with a soft and dry sponge.

Storage

Make sure the massage collar is cool and clean before storing.

Store in a dry, airy place out of reach of children and pets. Do not place any objects on top of the product.

TROUBLESHOOTING

Problem	Cause	Solution
The massage collar wouldn't fasten.	The massage collar is dis-charged.	Let it charge.
	The ⏻ button was not pressed.	Press the ⏻ button to switch on the massage collar.
The massage collar suddenly switched off.	Timer has expired.	Press the ⏻ button to switch on the massage collar.
	Operating time has been too long, increasing risk of overheating.	This is normal behaviour. The thermal fuse protects the massage collar from damage caused by overheating. Let it cool down before using it again.

Problem	Cause	Solution
The massage collar is hot.	You have switched on the heating function.	Press the }} button to switch off the heating function.

If the problem is not listed in the table above or if it persists, switch off the appliance, unplug it from the mains and contact an authorized service centre.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alza.cz/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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R20260707

SIGURO

Alza.cz a.s., Jankovcova 1522/53, 170 00 Prague 7, Czech Republic
© 2026 Alza.cz a.s. All Rights Reserved.

www.siguro.net
www.alza.cz/contact