MASSAGE GUN

M2

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1. The instructions for appliances incorporating a battery that contains materials hazardous to the environment shall give details on how to remove the battery and shall state that:
 - the battery must be removed from the appliance before it is scrapped,
 - the appliance must be disconnected from the supply mains when removing the battery,
 - the battery is to be disposed of safely.
2. Instructions for appliances containing non user-replaceable batteries shall state the substance of the following:

This appliance contains batteries that are only replaceable by skilled persons.
3. Instructions for appliances containing non-replaceable batteries shall state the substance of the following:

This appliance contains batteries that are non-replaceable.
4. For appliances intending to be supplied from a detachable supply unit for the purposes of recharging the battery, the type reference of the detachable supply unit shall be stated along with the substance of the following:

WARNING: For the purposes of recharging the battery, only use detachable supply unit provided with this appliance.

The instructions shall include the substance of the following:

This appliance can be used by persons aged from 18 years and above with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children.

Thank you for purchasing this product. Please read this manual carefully before use to ensure proper operation and use. Please keep this manual in a safe place for future reference.

Massage gun application: suitable for physical massage before or after sports and other similar physical activities.

Important Information

Statement

To reduce the risk of malfunction, electric shock, injury, fire, death and equipment damage, please read all instructions and warnings in this manual carefully before using this product.

Prohibited

1. This product is not to be used by minors. If you are injured, do not use this product. If you have muscle damage, consult your doctor before using this device.
2. Do not use on non-dry, non-clean skin surfaces. Press and move on the skin for about 60 seconds on each part.
3. Do not use above the neck and with less muscle coverage. Do not use the product on bones. Use the product only on the soft tissues of the body without causing pain or discomfort.
4. Do not use when bruising, signs of pain or discomfort occur.
5. Do not place any objects inside the ventilation holes of the product.
6. Do not immerse the product in water, and do not allow water to enter the ventilation holes of the product.
7. Do not put the product into fire or near a fire source.
8. Do not intentionally drop this product.
9. Do not use other chargers for charging the product. You must use the charger provided by our company for charging.
10. Do not disassemble this product in any way.
11. Do not place fingers, hair, or other parts of the body near the shaft and back of the massage head, as they may be pinched.

Do not touch this part when starting the machine without the massage head installed.

Do not touch this part when starting the machine with the massage head installed.

Warnings and Reminders

1. Keep away from liquid or heat sources. Keep the vents free of dust.
2. After charging is complete, unplug the charger.
3. The product cannot be used while charging.
4. Please use this product according to the instruction manual.
5. Before each use, carefully check whether the product can be used in the normal way.


Important

This product cannot be used in any of the following situations without the doctor's prior consent:


Pregnancy, complications of diabetes (such as neuropathy or retinal damage), pacemaker wearing, recent surgery, epilepsy or migraine, disc herniation, spine advancement, vertebral detachment, spinal joint stiffness, recent joint replacement surgery or IUD, metal needle or steel plate for any health related issues. These contraindications do not mean that you cannot use this product, we recommend that you consult your doctor before use.

Research currently underway is about the impact of vibration training on specific medical diseases. This is likely to reduce the above-mentioned contraindications. Practical experience shows that in many cases, combining vibration training into a treatment plan is also a treatment modality. This must be done on the advice of or accompanied by a doctor, expert, or physical therapist.

Product Description

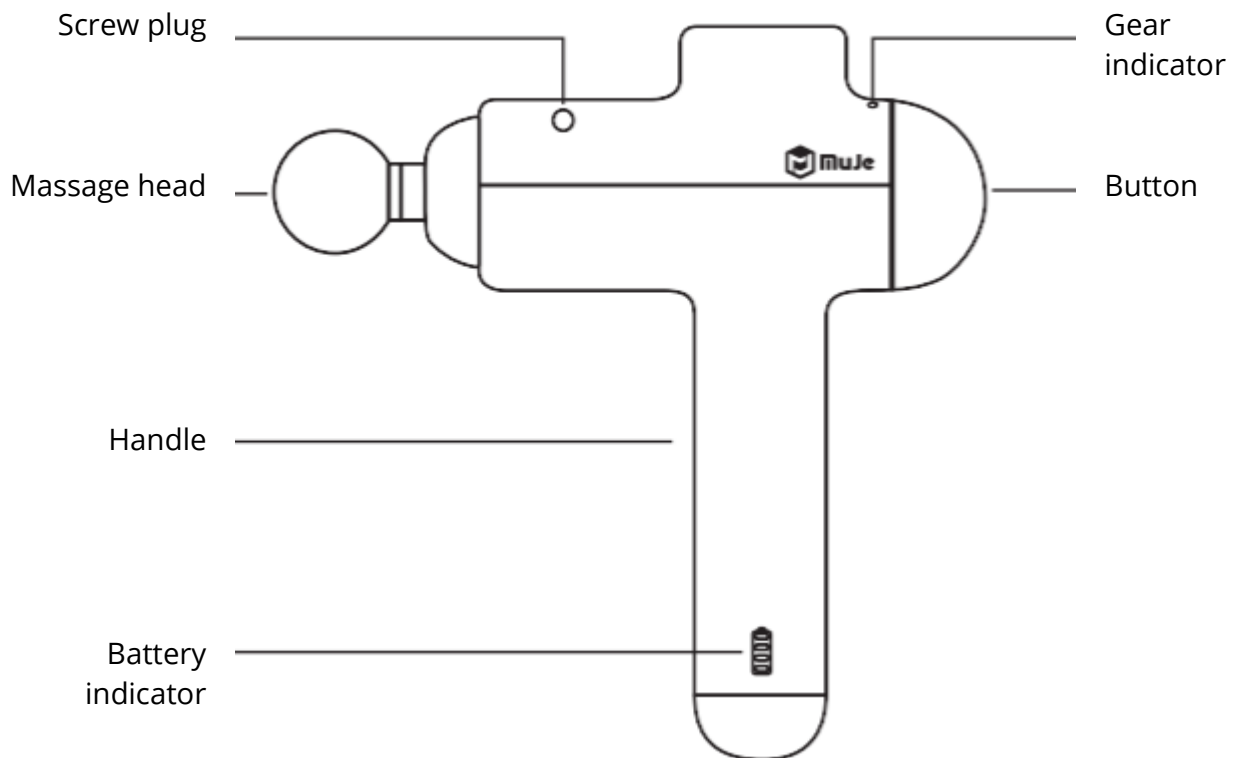
1. High-speed brushless balance motor, adopting dynamic balance and silencing technology.
2. EVA soft massage head provides a comfortable massage experience while reducing damage to bones.
3. Rechargeable lithium-ion battery (approximately 4 hours per full charge).
4. Small handle, more comfortable grip (41mm).
5. Wireless, lightweight and portable (0.8kg).
6. 3 gears can be adjusted at will.
7. Sweat and moisture resistant humanized design.
8. The battery and body are an integrated design, no risk of loose handle, very safe to use.

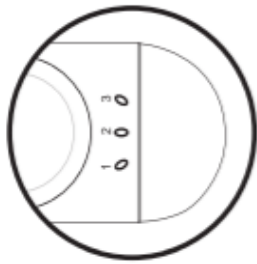Gear indicator

Button

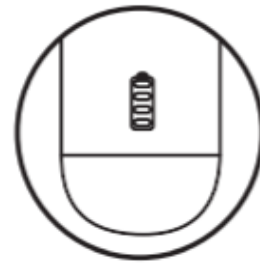Battery indicator

1. Gear indicator: 1-3 gears.
2. Power indicator: after powering on one of the following will happen:
 - a. 4 lights are on: 100% power.
 - b. 3 lights are on: 65% - 80% power.
 - c. 2 lights are on: 40-60% power.
 - d. 1 light is on: power is low, the product needs charging.
3. Keys: press and hold for 2 seconds to turn the machine on/off. Press to switch gears.

The product comes equipped with four interchangeable massage heads

Round head

Bullet head

Flat head

Fork head

Installation method of massage head:

While pulling/installing, rotate the massage head slightly to pull out/install. Do not pull/install directly (except for Fork head). The Fork head has alignment lines.

When installing a Fork head, when the inside of the horizontal alignment line (the red line in the figure) is aligned with the machine as shown in the figure, it indicates that it is in place. If there is a gap between the line and the machine, there is a risk of it falling off during use.

Basic Parameters

Product Name: MUJE Massage Fascia Gun	Product model: M2
Product Net Weight: 0.8kg	Input voltage: 100-240V ~ 50/60Hz
Interface type: DC, 16.8V = 0.8A	Battery capacity: 2000mAh
Gear strength: 3 speeds to adjust	Charging time: about 3 hours
Power supply mode: Rechargeable lithium-ion battery	Operating environment: 0°C ~ 40°C, not waterproof
Product size: 180mm * 62mm * 238mm	Persons suitable for use: adults over 18 years of age. Pregnant women are forbidden from using this product.

Charging

1. Before initial use, charge for four hours to fully charge the product.
2. When charging, connect the DC end of the supplied 26V adapter to the battery charging port, and plug the AC-side cable into a 100 ~ 240V power outlet.
3. The LED light on the device flashes regularly, showing the battery level, indicating that it is charging.
4. The first indicator is on, and the remaining battery charge level is 20% - 35%. While the second indicator is on, the remaining battery charge level is 40% - 60%. While the third indicator is on, the remaining battery charge level is 65% - 80%. While the fourth indicator is on, the remaining battery charge level is 85% - 100%.
5. The battery can be changed at any time and at any power level, and the device cannot be turned on when charging.
6. The battery is exhausted to less than 5% (the first LED is flashing) and needs to be charged.
7. The average running time of the product when fully charged is about 3 hours. There are some differences in gears and strength.

Operation

1. Clean the device with a slightly damp towel and then dry it with a soft cloth.
2. When storing or traveling, put the product in the package provided.
3. The product should be stored in a dry, well-ventilated place without corrosive gases. The distance from the ground should not be less than 250mm, and the number of stacking layers should not exceed 5.
4. After a long period of non-use, the product should be stored in a clean package and recharged every six months.
5. This product is suitable for transportation by car, train, ship and aircraft.
6. This product should be fixed in place during transportation to prevent violent vibration, collision, rolling and falling, It should also be fixed upward and in a rainproof place.

Environmental Protection

At the end of the product's lifespan, please do not discard it with general household waste. Put this product in a government designated recycling bin, or send this product to the manufacturer's after-sales service center so it can be handled in an environmentally friendly way.

Name and Content of Noxious Substances in Product

At the end of the product's lifespan, please do not discard it with general household waste. Put this product in a government designated recycling bin, or send this product to the manufacturer's after-sales service center so it can be handled in an environmentally friendly way.

Part names	Toxic and Noxious Substances or Elements					
	Pb	Hg	Cd	Cr(VI)	PBB	PBDE
Product Shell	O	O	O	O	O	O
Circuit Board Components	×	O	O	O	O	O
Metallic Parts	O	O	O	O	O	O
Rubber Parts	O	O	O	O	O	O
Adapter	×	O	O	O	O	O
Connection Line	O	O	O	O	O	O
Handle Cover	O	O	O	O	O	O
<i>This form was prepared in accordance with SJ /T 11364</i> ※: Circuit board components include printed circuit board machine components, such as resistors, capacitors, integrated circuits, connectors, etc. O: indicates that the content of noxious substances in all homogenous materials of the part is below the limit requirement specified in GB /T 26572. X: indicates that the content of noxious substances in at least one homogenous material of the part exceeds the limit requirement specified in GB /T 26572.						

Explanation for Environmental Protection Period

The use of this electronic and electrical product under the conditions described in the user manual will not cause leakage or escape of harmful substances or elements, and will not cause environmental pollution or serious damage to people and property for a period of 10 years.

Instructions For Use


Ball Head

Mainly suitable for large muscle groups.


Bullet Head

Mainly suitable for deep tissue.

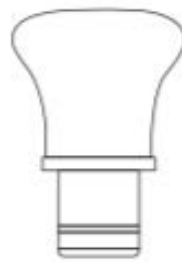
Flat Head

Mainly used to relax and shape various muscle parts.

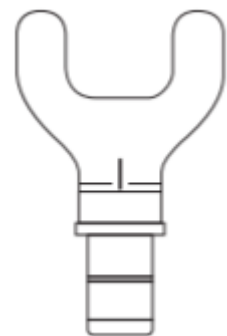
Fork Head

Mainly suitable for the sides of the neck, the spine and achilles tendon.

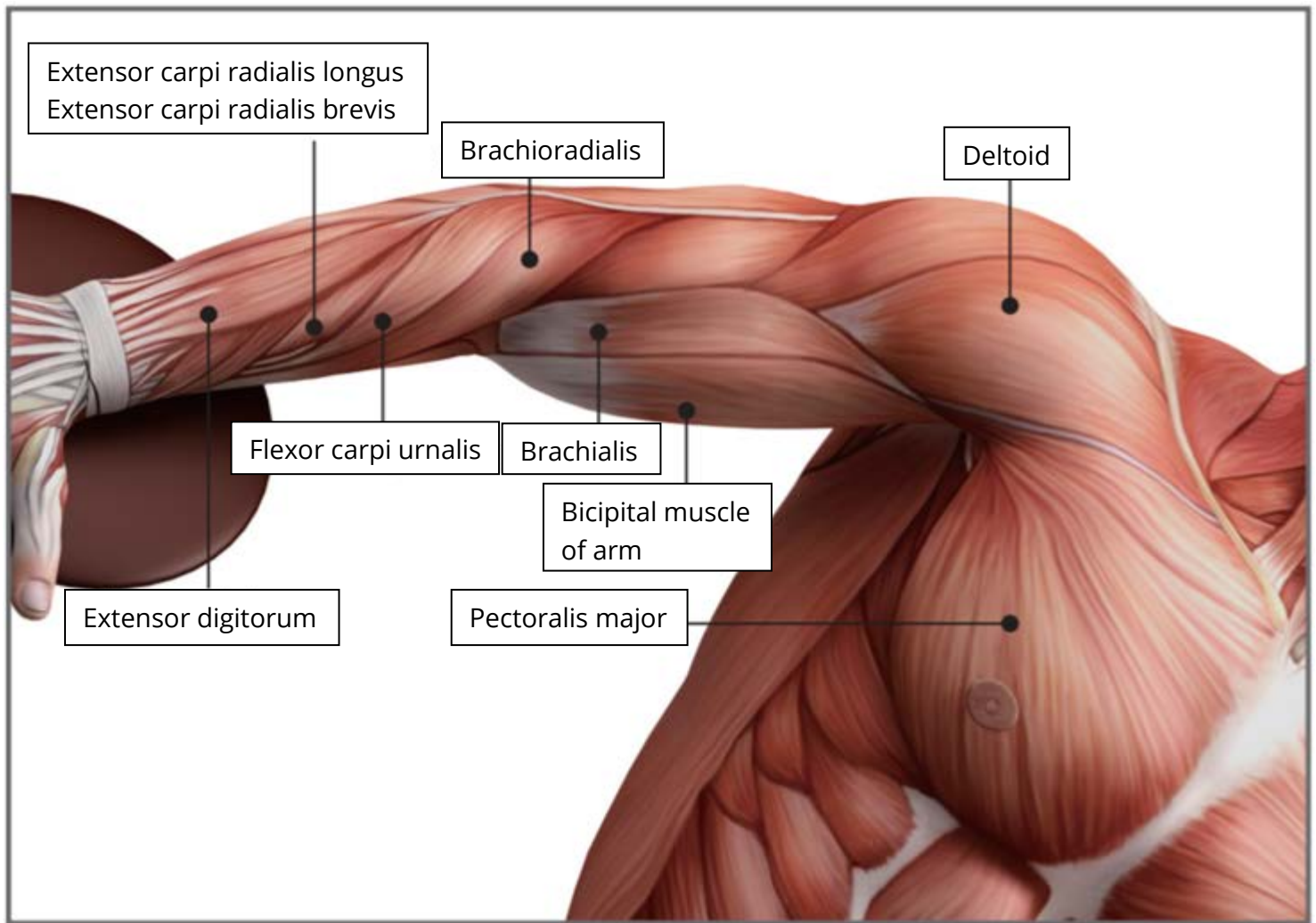Arm


Brachioradialis	10-15s * 3 times
Flexor carpi radialis	10-15s * 3 times
Palmaris longus	12-15s * 3 times
Flexor carpi ulnaris	10-15s * 3 times
Extensor digitorum	5-10s * 3 times
Extensor carpi radialis longus	5-10s * 3 times
Extensor carpi radialis brevis	5-10s * 3 times
Bicipital muscle of arm	12-20s * 3 times
Brachialis	12-20s * 3 times
Deltoid	12-20s * 3 times
Pectoralis major	12-20s * 3 times

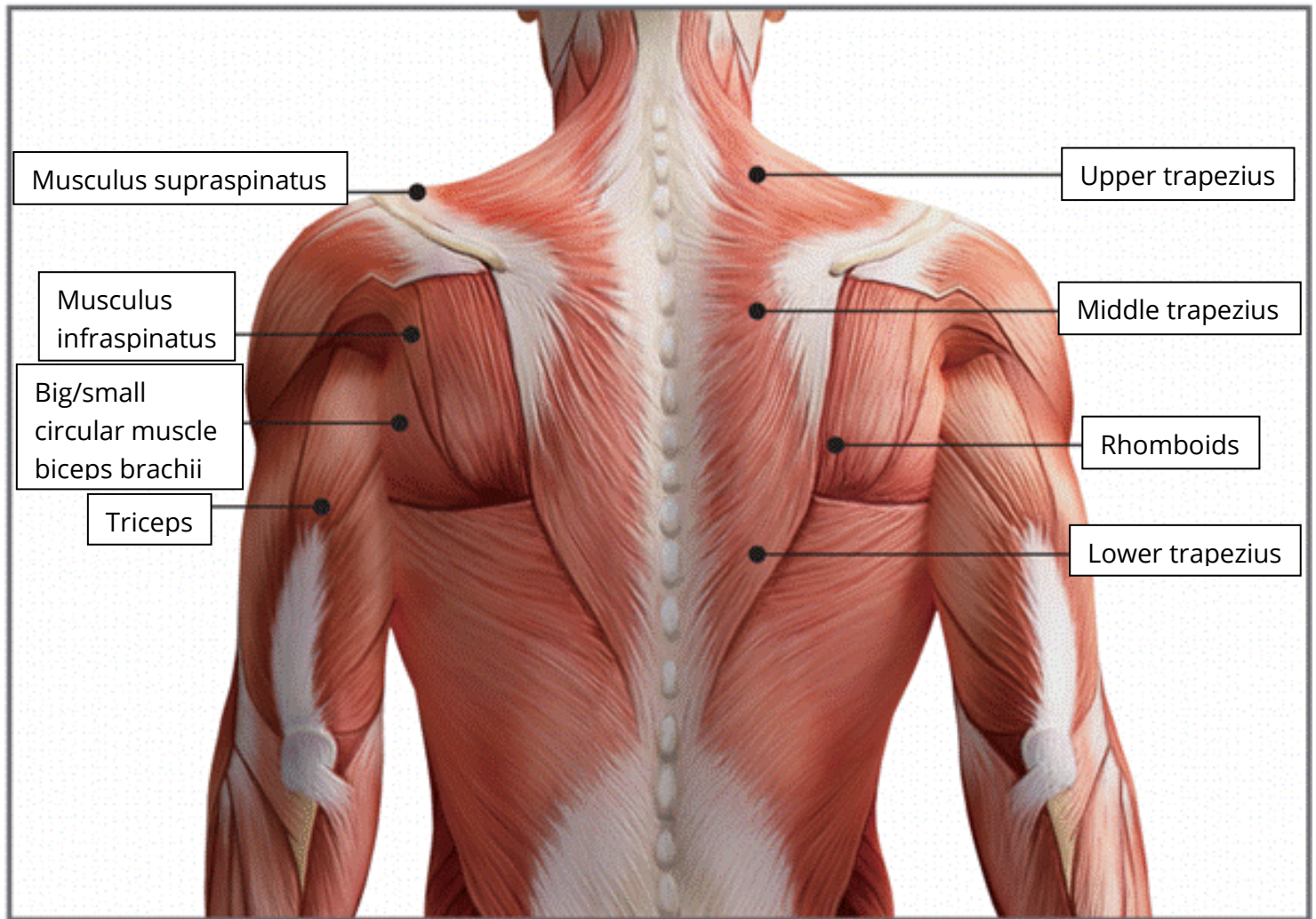Shoulder

Triceps	15-20s * 3 times
Rhomboids	20-30s * 3 times
Upper Trapezius	15-20s * 3 times
Middle Trapezius	15-20s * 3 times
Lower Trapezius	15-20s * 3 times
Musculus supraspinatus	10-15s * 3 times
Levator Scapulae	10-12s * 3 times
Big circular muscle biceps brachii	10-12s * 3 times
Small cilcular muscle biceps brachii	10-12s * 3 times
Musculus infraspinatus	10-12s * 3 times


Advice: Massage the muscles from the middle to the sides, from top to bottom.

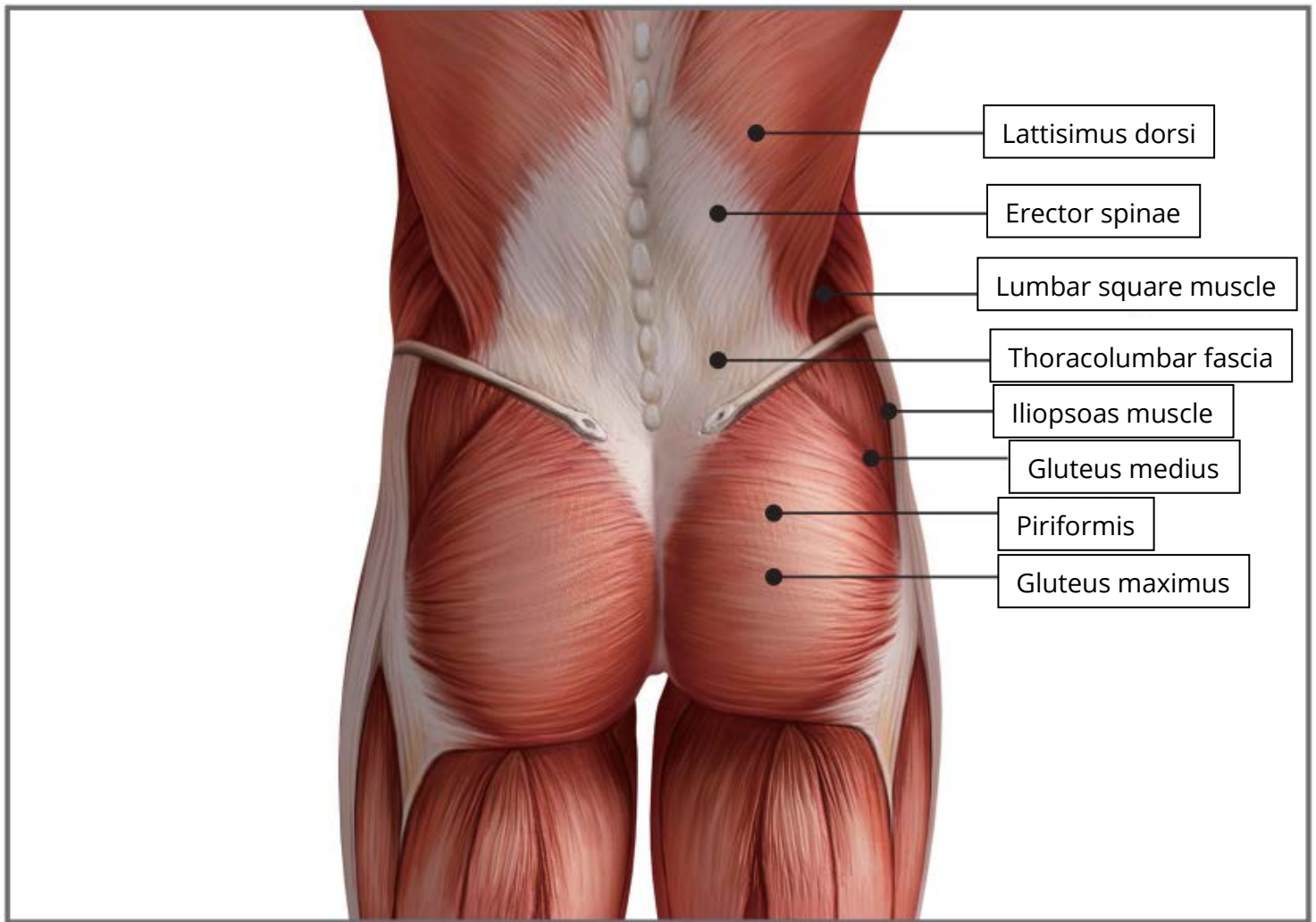Waist

Erector spinae (waist)	30-60s * 4 times
Lattissimus dorsi	30-60s * 4 times
Thoracolumbar fascia	40-60s * 3 times
Psoas major	30-60s * 3 times
Abdominal oblique muscle	10-15s * 2 times
Vertical spinal muscle (thoracic segment)	15-20s * 2 times
Lumbar square muscle	20-30s * 2 times
Iliopsoas muscle	20-30s * 3 times
Gluteus maximus	30-60s * 4 times
Gluteus medius	30-60s * 4 times
Piriformis	60-90s * 1 time


Advice: Massage the muscles from the middle to the sides, from top to bottom. Large muscle groups can run along the muscles.

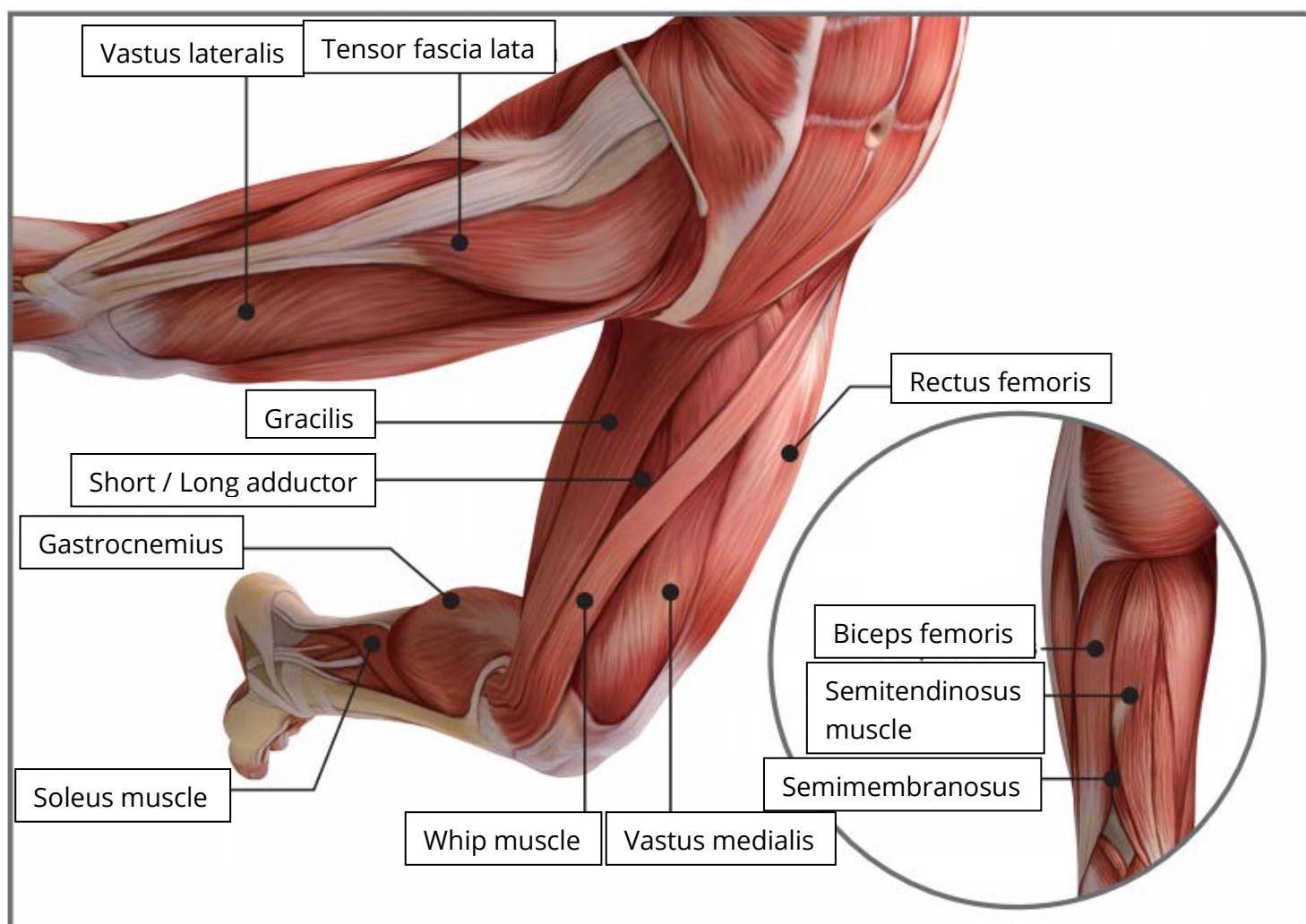Leg

Biceps femoris	20-30s * 3 times
Semitendinosus muscle	20-30s * 3 times
Semimembranosus	20-30s * 2 times
Gracilis	10-15s * 2 times
Short adductor	10-15s * 2 times
Long adductor	10-15s * 2 times
Rectus femoris	30-45s * 2 times
Vastus lateralis	30-45s * 2 times
Vastus medialis	30-45s * 2 times
Whip muscle	30-45s * 1 time
Tensor fascia lata	15-20s * 2 times

Advice: Massage the muscles from top to bottom.

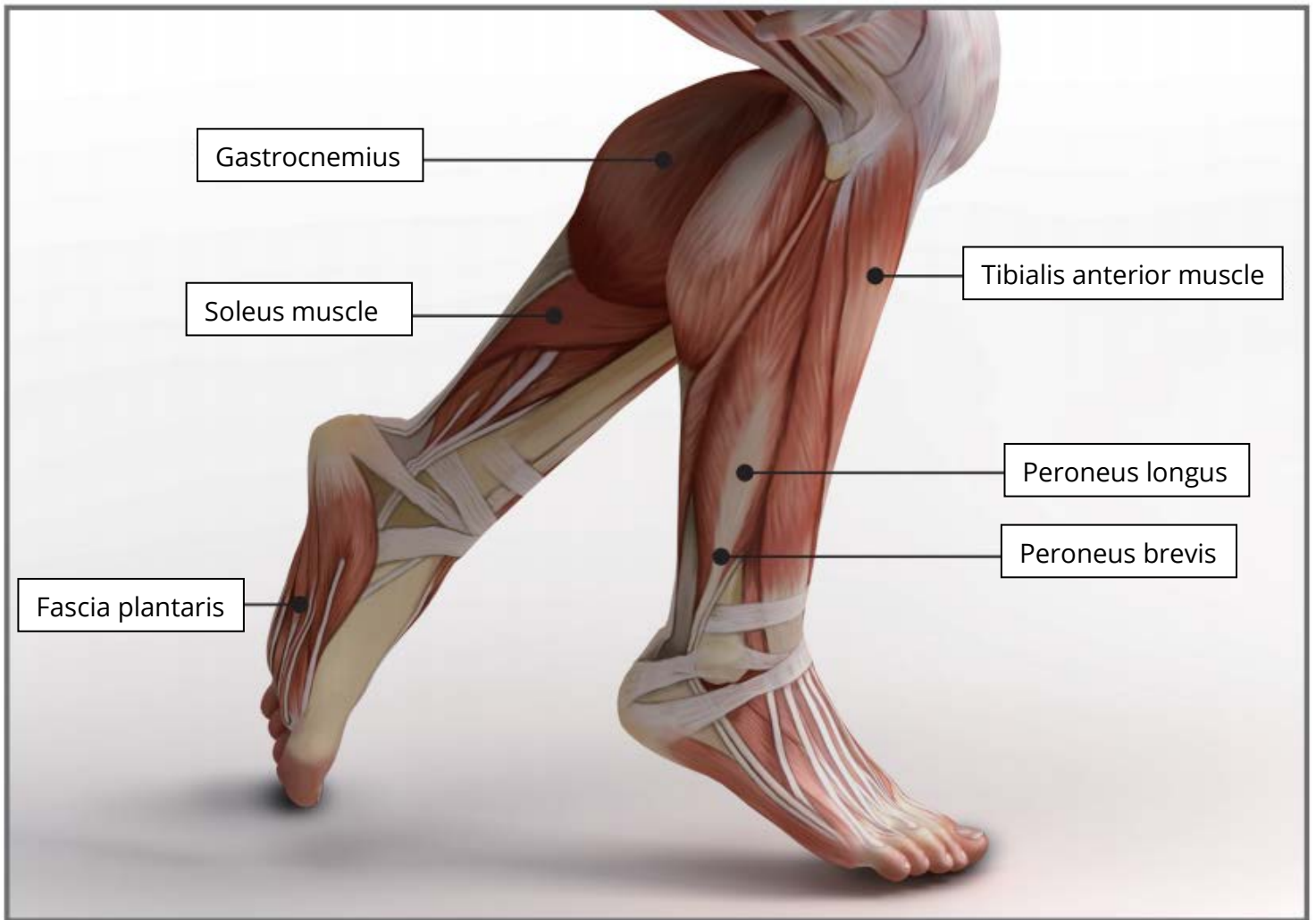Foot

Gastrocnemius	20-30s * 3 times
Soleus muscle	20-30s * 3 times
Tibialis anterior muscle	15-20s * 2 times
Peroneus longus	15-20s * 2 times
Peroneus brevis	15-20s * 2 times
Fascia plantaris	20-30s * 3 times

Advice: Massage the muscles from top to bottom.

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