# **SIGURO**

## **SC110x**

Personal Scale

## **User Manual**

<b>ENGLISH</b>	<b>4</b>
<b>DEUTSCH</b>	<b>9</b>
<b>ČESKY</b>	<b>15</b>
<b>SLOVENSKY</b>	<b>20</b>
<b>MAGYAR</b>	<b>25</b>
WARRANTY & SUPPORT	30
GARANTIE & BETREUUNG	30
ZÁRUKA & PODPORA	30
ZÁRUKA & PODPORA	30
GARANCIA & TÁMOGATÁS	30

# SAFETY INFORMATION

**Please read this manual carefully before use and keep it for future reference. The manufacturer declines any liability for damage caused by improper use or failure to follow the instructions in this manual. Improper use or non-compliance with the instructions will void the warranty.**


1. This appliance may be used by children 8 years of age and older and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, provided they have been given supervision or instruction concerning use of the appliance in a safe manner and understand the potential hazards. Children must not play with the device. Cleaning and user maintenance must not be performed by children without supervision.
2. The digital personal scale is intended for household use only.
3. **Have all repairs carried out by a specialized service!** Do not disassemble the scale or remove any components.
4. Do not use the scale near devices that generate electromagnetic fields. Place the scale at a sufficient distance from such a product to avoid mutual interference.
5. The weighing result has only a guiding value. This product is not a medical device. Please always consult your doctor about medication or diet to achieve your optimal weight.
6. Do not insert any objects into the openings of the device.
7. The scale is designed for weighing people from 10 years of age.
8. Pets must not be weighed on the scale.
9. Protect the scale from falls and impacts and do not jump on it.
10. You or the device may slip on a wet surface. Also, never step on the edge of the scale on a smooth surface and keep the scale in balance when measuring.
11. Do not immerse the scale in water.

## Notes about safety when weighing

1. Always take off your shoes and socks before weighing.
2. The surface of the scale is made of glass. Use extra caution not to slip on the surface. Do not step on the edges of the scale.
3. We recommend weighing at the same time every day.
4. The weighing results can be misleading after intense exercise, a strong diet or at extremely high temperatures.
5. Do not subject the scale to shocks, vibrations or drop it on the floor. Handle it with care.
6. Do not use the scale for commercial purposes.

## Notes regarding batteries

1. Insert the batteries with the correct polarity.
2. Do not throw the batteries in a fire.
3. Always replace the entire set of batteries at once. Do not mix old and new batteries or different types of batteries.
4. Replace empty batteries in time to prevent electrolyte leakage.
5. If you will not be using the appliance for a long time, remove the batteries.

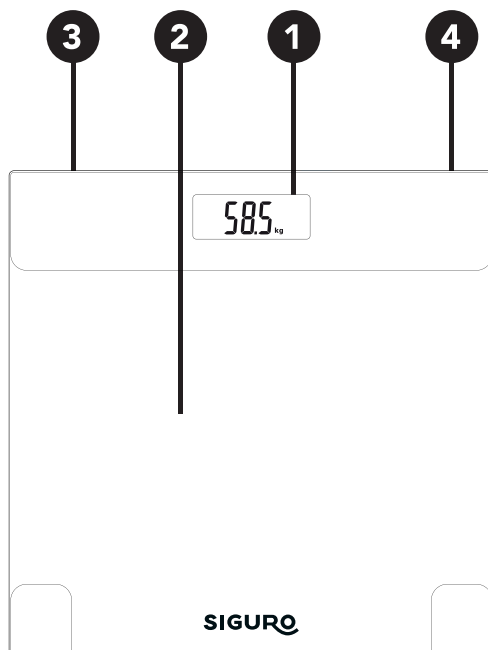- 1 Display
- 2 Weighing area
- 3 Battery cover (on the bottom side)
- 4 Unit setting switch (on the bottom side)

## SPECIFICATIONS

<b>Dimensions</b>	30.2 × 30.2 × 2.2 cm
<b>Weight</b>	1.5 kg
<b>Power</b>	2× 1.5 V batteries type AAA (included)
<b>Maximum load</b>	180 kg
<b>Weight measurement deviation</b>	100 g
<b>Thickness of the weighing area</b>	6 mm

# USING THE SCALE

## Before use


- The device is powered by 2 × 1.5 V AAA batteries. Remove the insulating foil from the battery contacts before use.
- If the scale cannot be switched on, check that the batteries are inserted or that the batteries are not discharged. Insert the batteries or replace the batteries with new ones.

## Inserting / replacing batteries

1. Open the battery compartment cover on the bottom of the scale.
2. Remove discharged batteries; use a suitable object if necessary.
3. Insert new batteries, always with one end first, and then push the other end. Close the battery compartment cover. Make sure the polarity is correct.


## Weight measurement

1. To achieve accurate weighing results, place the scale on a flat, hard and stable surface and avoid placing it on carpets and soft surfaces. Make sure that no foreign object is lying under the scale before weighing.
2. Switch on the scale by briefly pressing with your foot on the middle of the weighing surface, "0.0" appears on the display.
3. Carefully step on the scale. Stand evenly on the scale, do not move and wait until the weight indication on the display stabilizes.
4. The device switches off automatically after you step off the scale. The automatic switch-off occurs if "0.0" appears on the display or the same weight is displayed for about 8 seconds.
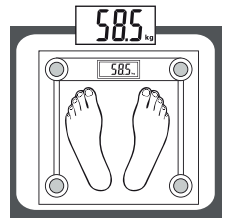
## Switching measuring units

You can switch between the units of measurement (kilograms / pounds / stones) using the button on the bottom of the scale. When the scale switched on, press the button. The LCD display shows the currently used unit of measurement. Press the button again to change the unit.
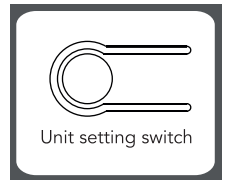
## Error messages on the display

**Overload.** If you load the scale more than the permissible limit, this error message will appear on the display. Reduce the load to prevent damage to the scale.

**Low battery warning.** Insert new batteries according to the instructions above.

**Erroneous measurement.** An error occurred during the weighing. Please, repeat weighing.

# CLEANING AND MAINTENANCE

---

- Use a damp cloth to clean the scale and make sure that no water gets inside the scale. DO NOT use chemical / abrasive cleaning agents. Do not immerse the scale in water!
- If an error appears on the display or if the scale cannot be switched off for a long time, please remove the batteries for approx. 3 seconds and then reinsert them. If this does not solve the problem, contact a service center or the vendor.

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## **ENGLISH**

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on [www.alzashop.com/doc](http://www.alzashop.com/doc).

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

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© 2020 Alza.cz, a. s. All Rights Reserved.