
STORMRED

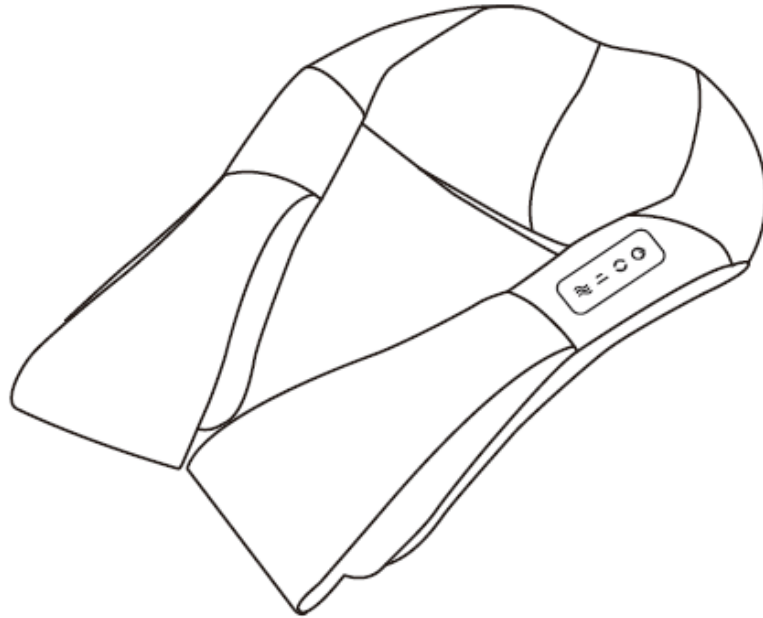NECK & SHOULDER SHIATSU MASSAGER WITH HEAT

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt


☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Instruction

Congratulations on your purchase of the Snailax Massager. It combines lighter heat therapy together with 8 massage nodes which rotate in unison to provide you with a deep-penetrating warmth and rejuvenating massage at home or away. Your body will feel more relaxed and rejuvenated.

Features

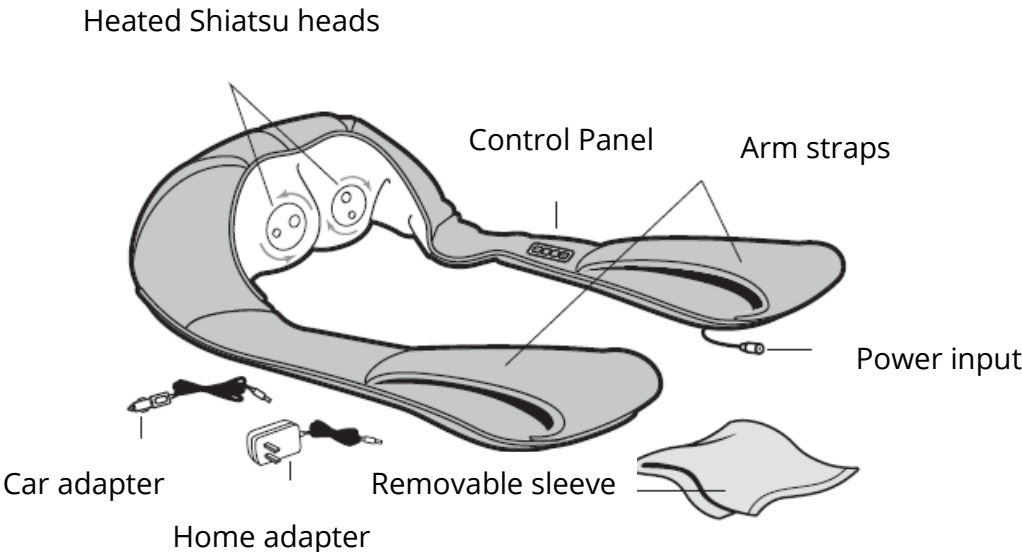Contents


- Neck & Shoulder Shiatsu massager with heat
- Home adapter
- 12V car adapter
- Removable sleeve

Technical data

Measurements:	19.81 x 12.2 inches
Weight:	4.85lbs
Voltage:	Input: AC 100-240V ~50/60Hz Output: 12VDC 1500mA
Nominal Power:	Max. 18 Watt
Automatic runtime:	15 minutes

Setup and Operation

1. Connect the adapter cable to the cable on the cushion.

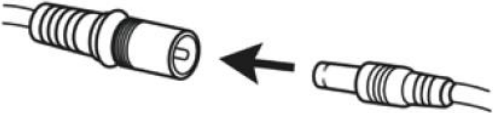Cushion side

Adaptor side


- 2. Plug the home adapter into an outlet or the car adapter into your car's power receptacle.
- 3. Turn on the massager and select your desired settings through below controller.
- 4. Additional dust-proof sleeve is provided to cover the massage heads for softening purpose when needed.

Power Button

Press to turn on the massager and activate shiatsu clockwise rotation at high level, Press again to turn off the massager

Intensity Button

When the massager is turned on, Press to set the shiatsu intensity at high or low level

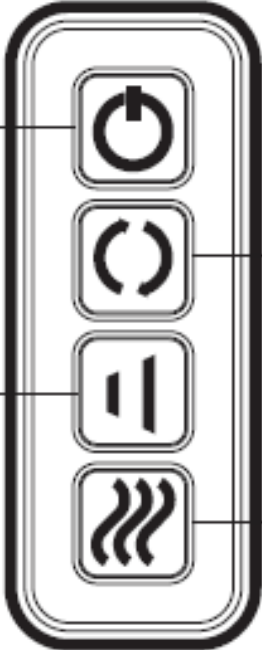Direction Button

When the massager is turned on, Press to shift Shiatsu clockwise rotation or Shiatsu anti-clockwise rotation

Heat Button


When the massager is turned on. Press to turn on or off the heat

Note: Press the power button to turn on the massager and activate the default setting massage mode i.e. the massage heads will rotate clockwise at high level of intensity and with heating, until it auto shuts off or press the corresponding button to change it.

This device will switch off automatically after 15 minutes timer runs out.

How to use your massager


Please use in the manner shown below

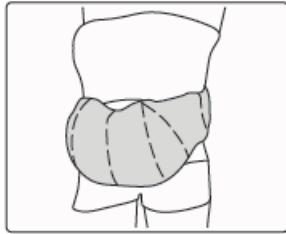
Abdomen


Neck & Shoulder

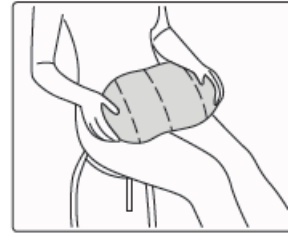
Thigh

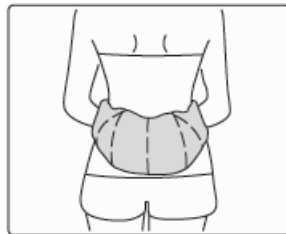Waist

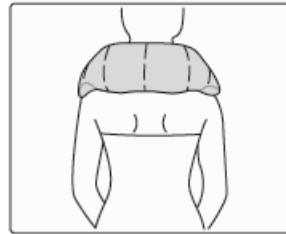Neck

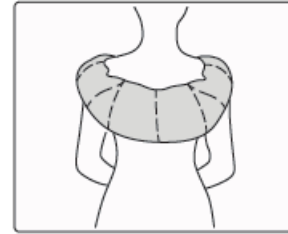Back

- Remove the adapter from the outlet when not in use.
- Do not setup or use the massager in a bathroom or similarly wet / damp areas.
- For your own safety we equip the massage mattress with a protection against over-heating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk and after a runtime of 15 minutes. The massager will switch off automatically and will only allow switching on again after an appropriate downtime.
- Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 minutes.

Safety Instructions

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This Appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- In plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- To avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 30 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kind of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

Care and Cleaning Instructions

CARE AND MAINTENANCE

Always unplug the unit before cleaning it. Never touch the power plug with wet hands. Before cleaning, allow unit to cool.

Wipe the unit with a soft, dry cloth. DO not use cloths containing any kind of chemical or alcohol and solvent liquids.

Never submerge any part of the unit in liquid.

If the unit is particularly dirty, use the following instructions:

1. Dampen a cloth in water or a 3%-5% mild detergent solution.
2. Wipe the dirty areas with the wet cloth.
3. Clean the cloth in water and wring it out, then wipe away any remaining detergent off the surface of the unit.
4. Wipe the areas with soft, dry cloth to dry.
5. Wait for the unit become fully dry before use.

Place massager in a safe, dry, and cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the surface.

Removable Sleeve

- Sleeve wraps around unit, covering massage heads, and fastens with loop strip.
- Wash in cold water with mild detergent. Allow to air dry.

Troubleshooting

PROBLEM	SOLUTIONS
Massage heads rotating at reduced speed.	The load on the massage heads is excessive. Reduce pressure and try again.
Device switches off suddenly.	Device may have been in use for 15 minutes. Auto-shutoff function activated. Continue using by pressing the ON button after resting the unit for 15 minutes.
	Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 mins before restarting.
Device fails to start.	Check if device has been used for more than 15 mins. If so, let device rest for 15 mins before continuing. Device may require longer time to cool in warmer environments.
	Check if adapter is plugged into socket.
	Check if adapter output is connected to device power inlet wire.
Jerky motion of vibrations from device.	Excessive pressure may have applied onto the device. Reduce the load and try again.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g. during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

IČO: 27082440

Subject of the declaration:

Title: Neck & shoulder shiatsu massager with heat

Model / Type: SL-632N

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. (EU) 2014/35/EU

Directive No. (EU) 2014/30/EU


Directive No. (EU) 2011/65/EU as amended 2015/863/EU

Prague, 27.1.2021

WEEE

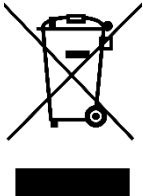This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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