
SIGURO

**Ultrasonic
skin scrubber**

SGR-SK-U530W

| | |
|------------------|-----------|
| ENGLISH | 5 |
| DEUTSCH | 9 |
| ČESKY | 13 |
| SLOVENSKY | 17 |
| MAGYAR | 21 |
| FRANÇAIS | 25 |
| ITALIANO | 29 |
| ROMÂNĂ | 33 |

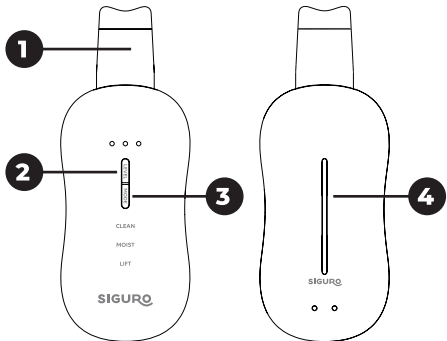SAFETY INFORMATION

1. This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, even if they have been given supervision or instruction concerning use of the appliance in a safe manner and understand the potential hazards. Children should be supervised to ensure that they do not play with the appliance.
2. Never immerse the appliance in water!
3. The appliance is intended for household and similar use only. It is not intended for use in a beauty salon or other commercial use. Use of this device is not intended as a substitute for medical care.
4. Do not use the appliance if you have hypersensitive, damaged or skin with eczema (e.g. fresh scars, various wounds, burns, redness, swollen cheeks or skin diseases). To prevent the spread of micro-inflammation, do not use it even if you have acne on your skin.
5. If you feel pain or discomfort during cleaning, or if you experience any other health problems, switch off the appliance immediately and consult your doctor or dermatologist about its further use.
6. Do not use the appliance on children, legally incompetent persons or people who are unable to express their feelings.
7. Do not use the appliance with exfoliating or otherwise aggressive gels (e.g. with alcohol)!
8. Do not use the appliance close to sensitive body parts, not intended for this device.
9. Do not use the appliance in one time longer than specified in the operating instructions! Prolonged use may cause skin injuries.

- 1 **Ultrasonic skin spatula**
- 2 **Intensity setting button**
- 3 **Power and mode switch**
- 4 **Sensor ring**

Included in package but not shown: charging base and cable and cleaning cloth.

| | |
|------------------------|-----------------------------|
| Dimensions (H x W x D) | 158 × 58 × 28 mm |
| Weight | 190 g |
| Power | 5 V / 0.5 A |
| Power consumption | 2.4 W |
| Battery | Rechargeable 3.7 V, 800 mAh |

USE

Charging

- You can charge the device using the charging base.
- Fully charge the device before using it for the first time, which takes about 4 hours. When fully charged, the battery lasts approximately 60 minutes of continuous operation.
- The charge status will be indicated by the intensity indicators above the **2** button. When charging is complete, all 3 indicators will be lit continuously.

Method of use

The device is used to treat the face using ultrasonic waves. Always remove all makeup and moisten your face with water before using the device. Apply the ultrasonic spatula to the face at an angle of approx. 30 ° and move it slowly along the contours of the face.

The device offers 3 modes of operation:

- **CLEAN** (cleaning mode): serves for deep pore cleansing and cell regeneration.
- **MOIST** (moistening mode): improves the absorption of facial cosmetics.
- **LIFT** (firming and lifting mode): firms the skin and smoothes wrinkles.

Controls

- Turn on the device by pressing the **3** button. The default mode is **CLEAN**.
- Press briefly the **3** button and you can switch between the individual modes. The active mode will be shown by the indicator.
- Use the **2** button to set the required intensity. The set intensity will be shown by indicators.
- You can turn off the device by pressing and holding the **3** button. It also switches off automatically after 5 minutes of use and after a short press of button 3 at the end of the mode switching cycle (that is after the **LIFT** mode).

CLEANING AND MAINTENANCE

- Always make sure that the appliance is switched off before cleaning or storing it.
- Rinse the device under running water. Do not immerse in water or any other liquid. Do not use any cleaning agents, lubricants or polishes.
- Do not dry the device in direct sunlight.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc. This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic

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© 2021 Alza.cz, a. s. All Rights Reserved.