
SIGURO

**Multi-therapy
beauty device**

SGR-SK-R760W

ENGLISH	5
DEUTSCH	10
ČESKY	16
SLOVENSKY	21
MAGYAR	27
FRANÇAIS	33
ITALIANO	39
ROMÂNĂ	45

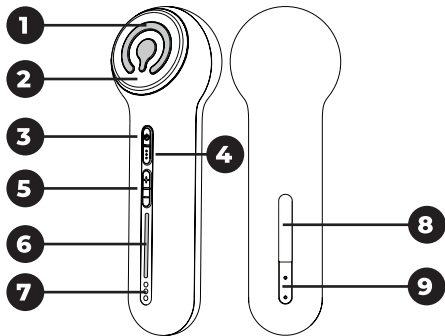SAFETY INFORMATION

1. This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, even if they have been given supervision or instruction concerning use of the appliance in a safe manner and understand the potential hazards. Children should be supervised to ensure that they do not play with the appliance.
2. Never immerse the appliance in water!
3. The appliance is intended for household and similar use only. It is not intended for use in a beauty salon or other commercial use. Use of this device is not intended as a substitute for medical care.
4. Do not use the appliance if you have hypersensitive, damaged or skin with eczema (e.g. fresh scars, various wounds, burns, redness, swollen cheeks or skin diseases). To prevent the spread of micro-inflammation, do not use it even if you have acne on your skin.
5. If you feel pain or discomfort during cleaning, or if you experience any other health problems, switch off the appliance immediately and consult your doctor or dermatologist about its further use.
6. Do not use the appliance on children, legally incompetent persons or people who are unable to express their feelings.
7. Do not use the appliance with exfoliating or otherwise aggressive gels (e.g. with alcohol)!
8. Do not use the appliance close to sensitive body parts, not intended for this device.
9. Do not use the appliance in one time longer than specified in the operating instructions! Prolonged use may cause skin injuries.

- 1 **Electrodes**
- 2 **LED lights**
- 3 **On / Off button**
- 4 **Mode switching button**
- 5 **+/- buttons for intensity control**
- 6 **Intensity indicator**
- 7 **Operating mode indicators**
- 8 **USB charging connector (with cover)**
- 9 **Charging contacts for wall base**

Package parts not shown: charging bas and cable, cleaning cloth and protective case.

Dimensions (H x W x D)	57 × 170 × 52 mm
Weight	212 g
Power	5.0 V / 2.0 A
Power consumption	4.0 W
Battery	Rechargeable Li-ion 7.4 V, 1000 mAh

USE

Charging

- You can charge the device using the included USB cable or a

- charging base that can be glued to the wall.
- Fully charge the device before using it for the first time, which takes about 3 hours. When fully charged, the battery lasts approximately 60 minutes of continuous operation.
- The charge status will be indicated by the intensity indicator (6). It flashes while charging and once fully charged it is continuously lit.
- When the battery is low, the mode indicator flashes and the device turns off.
- You cannot use the device when it is charging.

Method of use


The device is used to treat the face and support collagen fibers by electro stimulation. It offers three modes: radio frequency (RF), micro current (MC) and LED mask.

1. Take off all the jewelry you wear, including earrings, piercing or necklaces.
2. Treat the face in five zones - forehead, left and right face, chin and neck.


3. Thoroughly cover the electrodes (1) with a layer of treatment

gel before using them. Then place the device on your face and move it slowly in a circular motion.

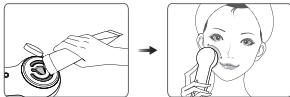4. Each time the treated zone of the face is changed, apply a new layer of treatment gel to the electrodes.
5. Move the device in the way shown in the figure below for lifting effect. For the cheeks, move from the inside towards the ears, while for the forehead, move from the eyebrows to the hairline. Follow the order shown in the figure.

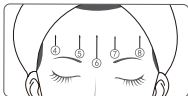Notes on use:

- Use the device in one mode only for one session and not longer than 8 minutes.
- Do not use the device for several days in a row, always let the

face rest for at least 1 day.

- We recommend that you do not use the device more than 2-3 times a week.
- Do not look directly into the LED lamp during use.

Controls

- Turn on the device using **3**.
- Press the **4** button to switch between the required modes of treatment (in the order of RF, MC and LED)
- Use the **5** buttons to set the required intensity. If you feel that the electrodes on your face are too hot, set a lower intensity.
- You can manually turn off the device by pressing the power button again (**3**). The device also switches off automatically depending on the operating mode (after 4 minutes in RF mode, 6 minutes in MC mode and 8 minutes in LED mode).

CLEANING AND MAINTENANCE

- Always make sure that the appliance is switched off before cleaning or storing it.
- Rinse the device under running water. Do not immerse in water or any other liquid. Do not use any cleaning agents, lubricants or polishes.
- Do not dry the device in direct sunlight.

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ENGLISH

The use of the WEEE Symbol indicates that this product should not be treated as household waste. By ensuring this product is disposed of correctly, you will help protect the environment. For more detailed information about the recycling of this product, please contact your local authority, your household waste disposal service provider, or the shop where you purchased the item.

This product meets all of the related basic EU regulation requirements. The EU declaration of conformity is available on www.alzashop.com/doc.

This product sold in the European Union meet the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic

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© 2021 Alza.cz, a. s. All Rights Reserved.